

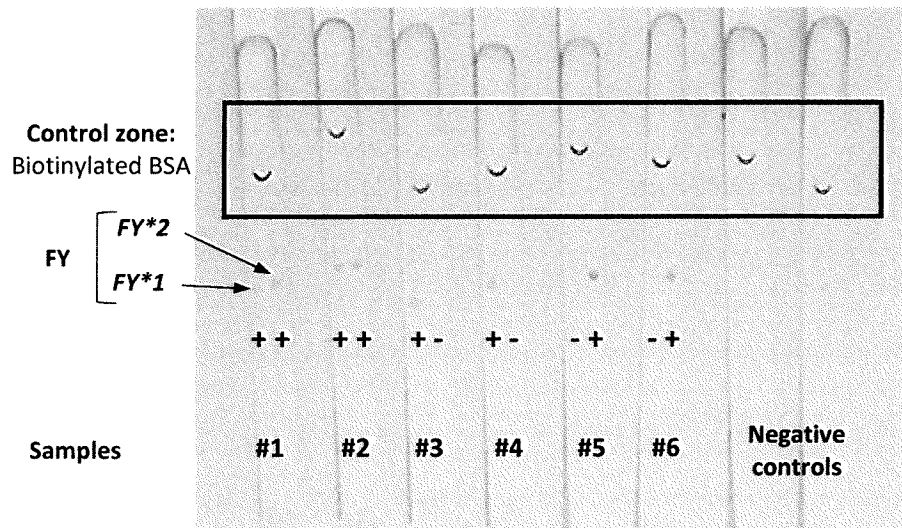


Fig. 1A

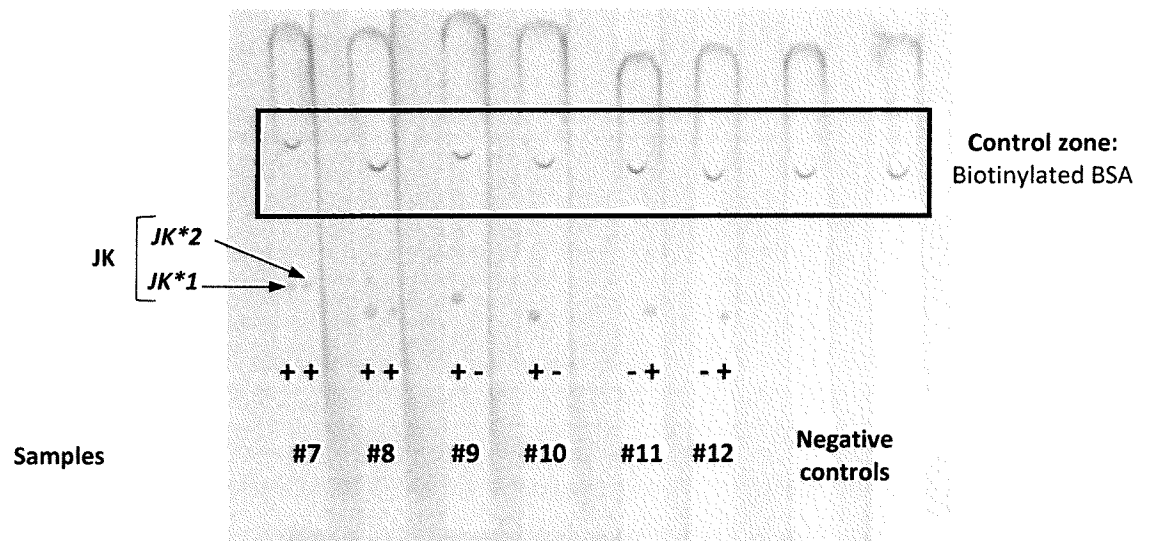


Fig. 1B

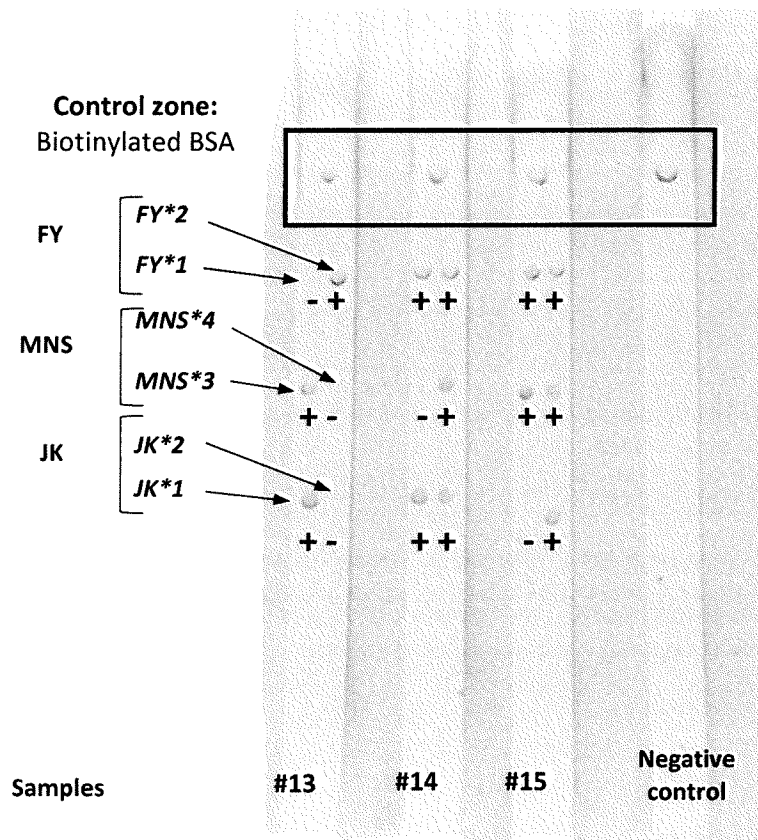


Fig. 1C

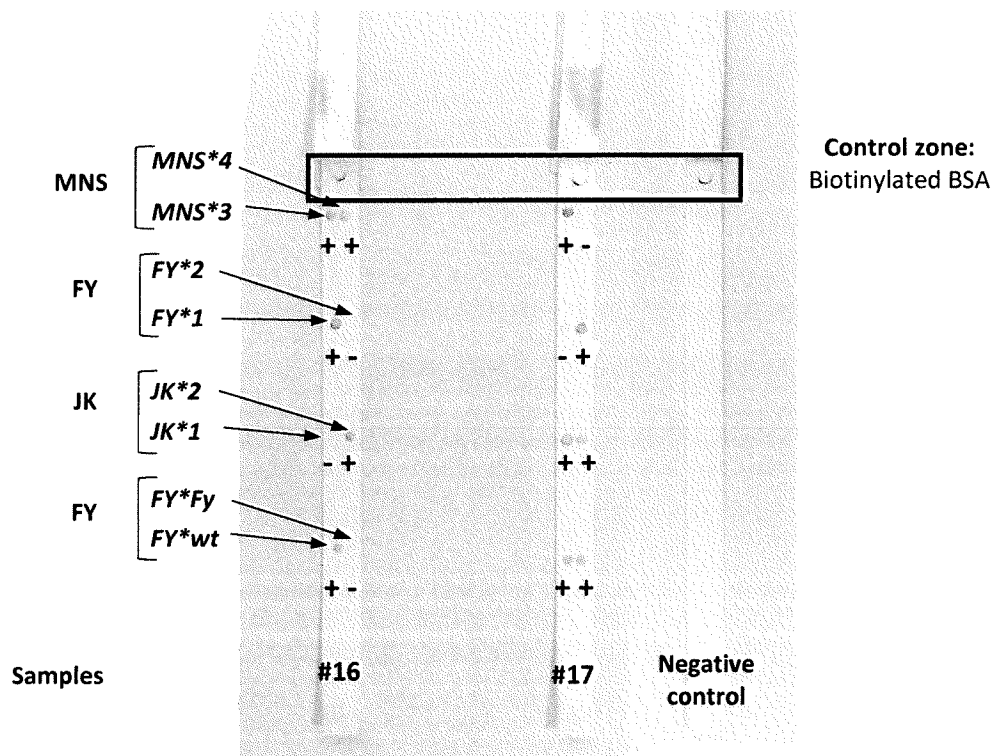


Fig. 1D

		Phenotypes FY system
Sample #1	Serological technique	FY: 1,2
	Reference genotyping technique	FY: 1,2 (FY*1/*2)
	Migration on nitrocellulose membrane technique	FY: 1,2 (FY*1/*2)
Sample #2	Serological technique	FY: 1,2
	Reference genotyping technique	FY: 1,2 (FY*1/*2)
	Migration on nitrocellulose membrane technique	FY: 1,2 (FY*1/*2)
Sample #3	Serological technique	FY: 1,-2
	Reference genotyping technique	FY: 1,-2 (FY*1/*1)
	Migration on nitrocellulose membrane technique	FY: 1,-2 (FY*1/*1)
Sample #4	Serological technique	FY: 1,-2
	Reference genotyping technique	FY: 1,-2 (FY*1/*1)
	Migration on nitrocellulose membrane technique	FY: 1,-2 (FY*1/*1)
Sample #5	Serological technique	FY: -1,2
	Reference genotyping technique	FY: -1,2 (FY*2/*2)
	Migration on nitrocellulose membrane technique	FY: -1,2 (FY*2/*2)
Sample #6	Serological technique	FY: -1,2
	Reference genotyping technique	FY: -1,2 (FY*2/*2)
	Migration on nitrocellulose membrane technique	FY: -1,2 (FY*2/*2)

Figure 2

		Phenotypes JK system
Sample #7	Serological technique	JK: 1, 2
	Reference genotyping technique	JK: 1, 2 (JK*1/*2)
	Migration of amplified targets on nitrocellulose membrane technique	JK: 1, 2 (JK*1/*2)
Sample #8	Serological technique	JK: 1, 2
	Reference genotyping technique	JK: 1, 2 (JK*1/*2)
	Migration of amplified targets on nitrocellulose membrane technique	JK: 1, 2 (JK*1/*2)
Sample #9	Serological technique	JK: 1, -2
	Reference genotyping technique	JK: 1, -2 (JK*1/*1)
	Migration of amplified targets on nitrocellulose membrane technique	JK: 1, -2 (JK*1/*1)
Sample #10	Serological technique	JK: 1, -2
	Reference genotyping technique	JK: 1, -2 (JK*1/*1)
	Migration of amplified targets on nitrocellulose membrane technique	JK: 1, -2 (JK*1/*1)
Sample #11	Serological technique	JK: -1, 2
	Reference genotyping technique	JK: -1, 2 (JK*2/*2)
	Migration of amplified targets on nitrocellulose membrane technique	JK: -1, 2 (JK*2/*2)
Sample #12	Serological technique	JK: -1, 2
	Reference genotyping technique	JK: -1, 2 (JK*2/*2)
	Migration of amplified targets on nitrocellulose membrane technique	JK: -1, 2 (JK*2/*2)

Figure 3

		Phenotypes		
		JK	FY	MNS
Sample #13	Serological technique	JK: 1, -2	FY: -1,2	MNS: 3, -4
	Reference genotyping technique	JK: 1, -2 (JK*1/*1)	FY: -1,2 (FY*2/*2)	MNS: 3, -4 (MNS*3/*3)
	Migration on nitrocellulose membrane technique	JK: 1, -2 (JK*1/*1)	FY: -1,2 (FY*2/*2)	MNS: 3, -4 (MNS*3/*3)
Sample #14	Serological technique	JK: 1, 2	FY: 1,2	MNS: -3, 4
	Reference genotyping technique	JK: 1, 2 (JK*1/*2)	FY: 1,2 (FY*1/*2)	MNS: -3, 4 (MNS*4/*4)
	Migration on nitrocellulose membrane technique	JK: 1, 2 (JK*1/*2)	FY: 1,2 (FY*1/*2)	MNS: -3, 4 (MNS*4/*4)
Sample #15	Serological technique	JK: -1, 2	FY: 1,2	MNS: 3, 4
	Reference genotyping technique	JK: -1, 2 (JK*2/*2)	FY: 1,2 (FY*1/*2)	MNS: 3, 4 (MNS*3/*4)
	Migration on nitrocellulose membrane technique	JK: -1, 2 (JK*2/*2)	FY: 1,2 (FY*1/*2)	MNS: 3, 4 (MNS*3/*4)
Sample #16	Serological technique	JK: -1, 2	FY: 1,-2	MNS: 3, 4
	Reference genotyping technique	JK: -1, 2 (JK*2/*2)	FY: 1,-2 (FY*1/*1)	MNS: 3, 4 (MNS*3/*4)
	Migration of amplified targets on nitrocellulose membrane technique	JK: -1, 2 (JK*2/*2)	FY: 1,-2 (FY*1/*1)	MNS: 3, 4 (MNS*3/*4)
Sample #17	Serological technique	JK: 1, 2	FY: -1,2	MNS: 3, -4
	Reference genotyping technique	JK: 1, 2 (JK*1/*2)	FY: -1,2 (FY*2/*Fy)	MNS: 3, -4 (MNS*3/*3)
	Migration of amplified targets on nitrocellulose membrane technique	JK: 1, 2 (JK*1/*2)	FY: -1,2 (FY*2/*Fy)	MNS: 3, -4 (MNS*3/*3)

Figure 4

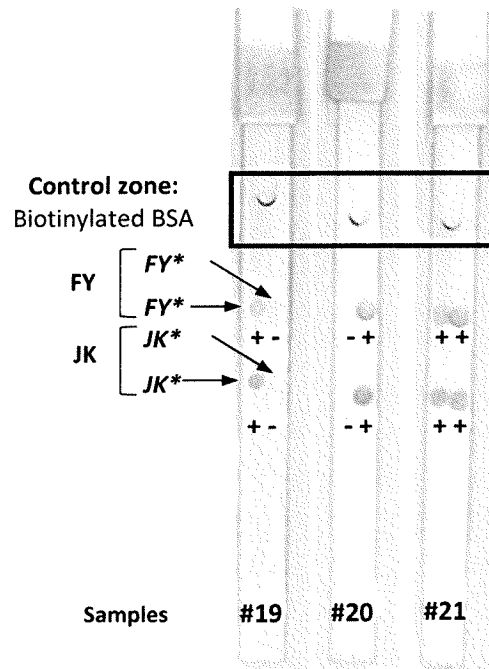


Figure 5

		Phenotypes	
		JK	FY
Sample #19	Serological technique	JK: 1, -2	FY: 1,-2
	Reference genotyping technique	JK: 1, -2 (<i>JK*1/*1</i>)	FY: 1,-2 (<i>FY*1/*1</i>)
	Migration on nitrocellulose membrane technique	JK: 1, -2 (<i>JK*1/*1</i>)	FY: 1,-2 (<i>FY*1/*1</i>)
Sample #20	Serological technique	JK: -1, 2	FY: -1,2
	Reference genotyping technique	JK: -1, 2 (<i>JK*2/*2</i>)	FY: -1,2 (<i>FY*2/*2</i>)
	Migration on nitrocellulose membrane technique	JK: -1, 2 (<i>JK*2/*2</i>)	FY: -1,2 (<i>FY*2/*2</i>)
Sample #21	Serological technique	JK: 1, 2	FY: 1,2
	Reference genotyping technique	JK: 1, 2 (<i>JK*1/*2</i>)	FY: 1,2 (<i>FY*1/*2</i>)
	Migration on nitrocellulose membrane technique	JK: 1, 2 (<i>JK*1/*2</i>)	FY: 1,2 (<i>FY*1/*2</i>)

Figure 6

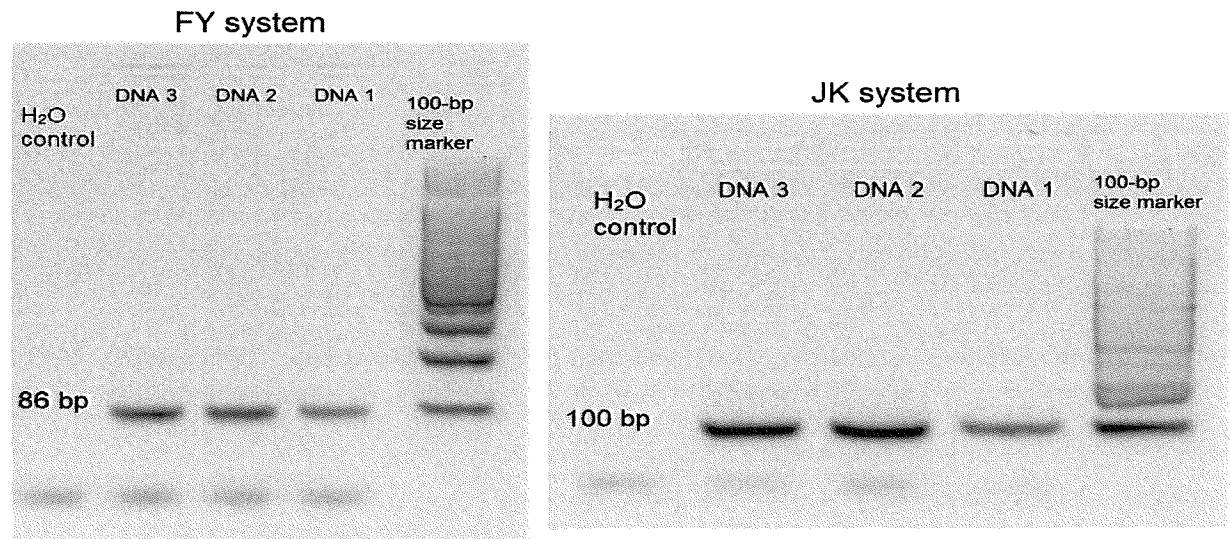


Figure 7

		Phenotypes		
		JK	FY	MNS
Sample #22	Serological technique	JK: 1, 2	FY: -1,2	MNS: 3, -4
	Reference genotyping technique	JK: 1, 2 (JK*1/*2)	FY: -1,2 (FY*2/*Fy)	MNS: 3, -4 (MNS*3/*3)
	Migration on nitrocellulose membrane technique	JK: 1, 2 (JK*1/*2)	FY: -1,2 (FY*2/*Fy)	MNS: 3, -4 (MNS*3/*3)
Sample #23	Serological technique	JK: 1, -2	FY: -1, -2 or FYnull	MNS: -3, 4
	Reference genotyping technique	JK: 1, -2 (JK*1/*1)	FY: -1, -2 (FY*Fy/*Fy)	MNS: -3, 4 (MNS*4/*4)
	Migration on nitrocellulose membrane technique	JK: 1, -2 (JK*1/*1)	FY: -1, -2 (FY*Fy/*Fy)	MNS: -3, 4 (MNS*4/*4)
Sample #24	Serological technique	JK: 1, -2	FY: 1, -2	MNS: -3, 4
	Reference genotyping technique	JK: 1, -2 (JK*1/*1)	FY: 1, -2 (FY*1/*Fy)	MNS: -3, 4 (MNS*4/*4)
	Migration on nitrocellulose membrane technique	JK: 1, -2 (JK*1/*1)	FY: 1, -2 (FY*1/*Fy)	MNS: -3, 4 (MNS*4/*4)
Sample #25	Serological technique	JK: -1, 2	FY: -1,2	MNS: 3, 4
	Reference genotyping technique	JK: -1, 2 (JK*2/*2)	FY: -1,2 (FY*2/*2)	MNS: 3, 4 (MNS*3/*4)
	Migration on nitrocellulose membrane technique	JK: -1, 2 (JK*2/*2)	FY: -1,2 (FY*2/*2)	MNS: 3, 4 (MNS*3/*4)

Figure 8

B2091PC00 Seq Listing_ST25_GB.txt
SEQUENCE LISTING

<110> Etablissement Français du Sang
<120> METHOD AND DEVICE FOR SNP GENOTYPING
<130> B2091PC
<160> 20
<170> PatentIn version 3.5
<210> 1
<211> 22
<212> DNA
<213> Artificial Sequence
<220>
<223> Excess primer for LATE-PCR amplification of the SNP specific
of JK system
<220>
<221> misc_feature
<223> 5'biotinylated primer
<400> 1
cagtcctttca gccccatttg ag 22
<210> 2
<211> 25
<212> DNA
<213> Artificial Sequence
<220>
<223> Limiting primer for LATE-PCR amplification of the SNP specific
of JK system
<400> 2
ggtgagcgcc atgaacattc ctccc 25
<210> 3
<211> 21
<212> DNA
<213> Artificial Sequence
<220>
<223> Excess primer for LATE-PCR amplification of the SNP
specific of MNS system
<220>
<221> misc_feature
<223> 5'biotinylated primer
<400> 3
acctggtaca gtgaaacgat g 21
<210> 4
<211> 27
<212> DNA
<213> Artificial Sequence
<220>
<223> Limiting primer for LATE-PCR amplification of the SNP specific
of MNS system

B2091PC00 Seq Listing_ST25_GB.txt

<400> 4
aggaaacccg cagaacagtt tgattcc 27

<210> 5
<211> 23
<212> DNA
<213> Artificial Sequence

<220>
<223> Excess primer for LATE-PCR amplification of the SNP specific
of FY system

<220>
<221> misc_feature
<223> 5'biotinylated primer

<400> 5
atgattcctt cccagatgga gac 23

<210> 6
<211> 29
<212> DNA
<213> Artificial Sequence

<220>
<223> Limiting primer for LATE-PCR amplification of the SNP specific
of FY system

<400> 6
tgcagagtca tccagcaggt tacaggagt 29

<210> 7
<211> 22
<212> DNA
<213> Artificial Sequence

<220>
<223> Excess primer for LATE-PCR amplification of the SNP specific
of FY (FY*FY) system

<220>
<221> misc_feature
<223> 5'biotinylated primer

<400> 7
ccctcattag tccttggtc tt 22

<210> 8
<211> 24
<212> DNA
<213> Artificial Sequence

<220>
<223> Limiting primer for LATE-PCR amplification of the SNP specific
of FY (FY*FY) system

<400> 8
ctcaccctgt gcagacagtt cccc 24

B2091PC00 Seq Listing_ST25_GB.txt

```

<210> 9
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the JK*01 allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 9
agtagatgtc ctcaaattgg                                     19

<210> 10
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the JK*02 allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 10
aagtagatgt tctcaaattgg g                                  21

<210> 11
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the FY*01 allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 11
aggttggcac catagtctc                                     19

<210> 12
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the FY*02 allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 12
aggttggcat catagtctc                                     19

```

B2091PC00 Seq Listing_ST25_GB.txt

```

<210> 13
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the FY*wt allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 13
gcttccaaga taagagcca 19

<210> 14
<211> 19
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the FY*02M.02 (FY*Fy) allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 14
gcttccaagg taagagcca 19

<210> 15
<211> 20
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the MNS*03 allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 15
aggagaaatg ggacaacttg 20

<210> 16
<211> 21
<212> DNA
<213> Artificial Sequence

<220>
<223> specific probe of the MNS*04 allele

<220>
<221> misc_feature
<223> 5' spacer group C6 NH2-(dT)15

<400> 16
ataggagaaa cgggacaact t 21

```

B2091PC00 Seq Listing_ST25_GB.txt

<210>	17	
<211>	12	
<212>	DNA	
<213>	Artificial Sequence	
<220>		
<223>	specific probe of the JK*01 allele	
<220>		
<221>	misc_feature	
<223>	5' spacer group C6 NH2-(dT)15	
<400>	17	
	agatgtcctc aa	12
<210>	18	
<211>	14	
<212>	DNA	
<213>	Artificial Sequence	
<220>		
<223>	specific probe of the JK*02 allele	
<220>		
<221>	misc_feature	
<223>	5' spacer group C6 NH2-(dT)15	
<400>	18	
	tagatgttct caat	14
<210>	19	
<211>	11	
<212>	DNA	
<213>	Artificial Sequence	
<220>		
<223>	specific probe of the FY*01 allele	
<220>		
<221>	misc_feature	
<223>	5' spacer group C6 NH2-(dT)15	
<400>	19	
	tggcaccata g	11
<210>	20	
<211>	12	
<212>	DNA	
<213>	Artificial Sequence	
<220>		
<223>	specific probe of the FY*02 allele	
<220>		
<221>	misc_feature	
<223>	5' spacer group C6 NH2-(dT)15	
<400>	20	
	tggcatcata gt	12

B2091PC00 Seq Listing_ST25_GB.txt