

TITLE OF THE INVENTION

BACTERIALCIDAL METHODS AND COMPOSITIONS

STATEMENT OF RIGHTS TO INVENTIONS MADE UNDER FEDERALLY5 SPONSORED RESEARCH

[0001] This work was supported by a National Institutes of Health Grant No. R01AI132638. The government may have certain rights to the invention.

RELATED APPLICATIONS/PATENTS

- 10 **[0002]** This application claims priority to U.S. Application Ser. No. 62/799,328, filed on January 31, 2019 as Attorney Docket No. 701586-094400PL01, the contents of which are incorporated herein by reference.

BACKGROUND OF THE INVENTION

- 15 **[0003]** Antibiotic resistance kills an estimated 700,000 people each year worldwide, and studies predict that this number could rise to 10 million by 2050, if efforts are not made to curtail resistance (Willyard, C. J. N. N. The drug-resistant bacteria that pose the greatest health threats. 543, 15 (2017)). Yet, the pace of resistance acquisition from mutation in pathogens is faster than clinical introduction of new antibiotics. There is an urgent need to
20 develop unconventional ways to combat the resistance.

SUMMARY OF THE INVENTION

- [0004]** The lethal effect of certain antibiotics occurs through the generation of Reactive Oxygen Species (ROS). Catalase, the ubiquitous key defense enzyme existing in most of the
25 aerobic pathogens, is utilized to scavenge hydrogen peroxide, thus preventing downstream oxidative damage. It has now been shown that catalase can be optimally photoinactivated by blue light having a wavelength of about 400nm to about 430nm, and specifically, a wavelength of about 410 nm. Photoinactivation of catalase renders broad-SPECTRUM catalase-positive microbial pathogens highly susceptible to ROS-generating antimicrobials
30 and/or immune cell attack. It has now been further determined that the antimicrobial effect of

photoinactivation is significantly and unexpectedly increased upon administration of a low-concentration of H_2O_2 and/or a ROS-generating agent.

[0005] In one aspect, the invention provides a method of treating a tissue of a subject infected with a catalase-positive microbe, said method comprising the steps of: applying light to the
5 tissue of the subject infected with the catalase-positive microbe at a wavelength of about 400nm to about 430nm, wherein the catalase is inactivated, and contacting the tissue with a composition comprising a diluted peroxide solution, thereby treating the tissue of the subject infected with the catalase-positive microbe.

[0006] In one embodiment, the wavelength is about 410nm.

10 [0007] In another embodiment, the dose of the light is about 5 J/cm^2 to about 200 J/cm^2 .

[0008] In yet another embodiment, the dose of the light is about 15 J/cm^2 .

[0009] In yet another embodiment, the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a pulsed nanosecond laser.

15 [0010] In yet another embodiment, the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a continuous wave LED.

[0011] In yet another embodiment, the diluted peroxide solution is a hydrogen peroxide solution.

[0012] In yet another embodiment, the hydrogen peroxide solution is between about 0.03% and about 0.3% hydrogen peroxide.

20 [0013] In yet another embodiment, the method further comprises administering a ROS generating agent to the infected tissue of the subject.

[0014] In yet another embodiment, the ROS generating agent is tobramycin, silver cation, iodine tincture, a gold nanoparticle, methylene blue, a β -lactam antibiotic, an aminoglycoside, a fluoroquinolone, an azole, a membrane-targeting polyene antifungal or a cell-wall targeting
25 antifungal.

[0015] In yet another embodiment, the tissue is skin, scalp or nails.

[0016] In yet another embodiment, the catalase-positive microbe is eradicated.

[0017] In another aspect, the invention provides a method of disinfecting an inanimate surface contaminated with a catalase-positive microbe, said method comprising the steps of:
30 applying light to the inanimate surface at a wavelength of about 400nm to about 430nm,

wherein the catalase is inactivated, and contacting the inanimate surface with a composition comprising a diluted peroxide solution, thereby disinfecting the inanimate surface.

[0018] In one embodiment, the wavelength is about 410nm.

[0019] In another embodiment, the dose of the light is about 5 J/cm² to about 200 J/cm².

5 **[0020]** In yet another embodiment, the dose of the light is about 15 J/cm².

[0021] In yet another embodiment, the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a pulsed nanosecond laser.

[0022] In yet another embodiment, the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a continuous wave LED.

10 **[0023]** In yet another embodiment, the diluted peroxide solution is a hydrogen peroxide solution.

[0024] In yet another embodiment, the hydrogen peroxide solution is between about 0.03% and about 0.3% hydrogen peroxide.

[0025] In yet another embodiment, the method further comprises administering a ROS
15 generating agent to the infected tissue of the subject.

[0026] In yet another embodiment, the ROS generating agent is tobramycin, silver cation, iodine tincture, a gold nanoparticle, methylene blue, a β -lactam antibiotic, an aminoglycoside, a fluoroquinolone, an azole, a membrane-targeting polyene antifungal or a cell-wall targeting antifungal.

20 **[0027]** In yet another embodiment, the inanimate surface is a material comprising metal, plastic, fabric, rubber, stone, composite surfaces or wood.

[0028] In yet another embodiment, the catalase-positive microbe is eradicated.

[0029] In yet another aspect, the invention provides a method of treating a tissue of a subject infected with a catalase-positive microbe, said method comprising the steps of: applying light
25 from a pulsed nanosecond laser to the tissue of the subject infected with the catalase-positive microbe at a wavelength of about 400nm to about 460nm, wherein the catalase is inactivated, and contacting the tissue with a composition comprising a ROS generating agent, thereby treating the tissue of the subject infected with the catalase-positive microbe.

[0030] In one embodiment, the wavelength is about 410nm.

30 **[0031]** In another embodiment, the dose of the light is about 5 J/cm² to about 200 J/cm².

[0032] In yet another embodiment, the dose of the light is about 15 J/cm².

[0033] In yet another embodiment, the catalase-positive microbe is a fungal or bacterial microbe.

5 [0034] In yet another embodiment, the diluted peroxide solution is a hydrogen peroxide solution.

[0035] In yet another embodiment, the hydrogen peroxide solution is between about 0.03% and about 0.3% hydrogen peroxide.

10 [0036] In yet another embodiment, the ROS generating agent is tobramycin, silver cation, iodine tincture, a gold nanoparticle, methylene blue, a β -lactam antibiotic, an aminoglycoside, a fluoroquinolone, an azole, a membrane-targeting polyene antifungal or a cell-wall targeting antifungal.

[0037] In yet another embodiment, the tissue is skin, scalp or nails.

[0038] In yet another embodiment, the catalase-positive microbe is eradicated.

15 [0039] In yet another aspect, the invention provides a method of producing a synergistic antimicrobial effect in a tissue of a subject infected with a catalase-positive microbe, said method comprising the steps of: applying light to the tissue of the subject infected with the catalase-positive microbe at a wavelength of about 400nm to about 460nm, wherein the catalase is inactivated, and contacting the tissue with a composition comprising a diluted peroxide solution, thereby producing the synergistic antimicrobial effect in the tissue of the
20 subject infected with the catalase-positive microbe.

[0040] In one embodiment, the wavelength is about 410nm.

[0041] In another embodiment, the dose of the light is about 5 J/cm² to about 200 J/cm².

[0042] In yet another embodiment, the dose of the light is about 15 J/cm².

25 [0043] In yet another embodiment, the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a pulsed nanosecond laser.

[0044] In yet another embodiment, the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a continuous wave LED.

[0045] In yet another embodiment, the diluted peroxide solution is a hydrogen peroxide solution.

[0046] In yet another embodiment, the hydrogen peroxide solution is between about 0.03% and about 0.3% hydrogen peroxide.

[0047] In yet another embodiment, the method further comprises administering a ROS generating agent to the infected tissue of the subject.

- 5 [0048] In yet another embodiment, the ROS generating agent is tobramycin, silver cation, iodine tincture, a gold nanoparticle, methylene blue, a β -lactam antibiotic, an aminoglycoside, a fluoroquinolone, an azole, a membrane-targeting polyene antifungal or a cell-wall targeting antifungal.

[0049] In yet another embodiment, the tissue is skin, scalp or nails.

- 10 [0050] In yet another embodiment, the catalase-positive microbe is eradicated.

BRIEF DESCRIPTION OF THE DRAWINGS

- [0051] The following Detailed Description, given by way of example, but not intended to limit the invention to specific embodiments described, may be understood in conjunction
15 with the accompanying figures, incorporated herein by reference.

- [0052] Figure 1 depicts the effect of ns-410 nm exposure on pure catalase solution. (a). Absorption spectra of pure catalase solution under ns-410 nm exposure. Catalase solution: 3 mg/ml, filtered with a 0.2 μ m filter. (b). Percent of remaining active catalase after different treatment schemes (different wavelengths under the same dosage). Quantification of catalase
20 was obtained by an Amplex Red Catalase kit. Data: Mean $\pm$ standard deviation (N=3).

- [0053] Figure 2 depicts the effect of ns-410 nm exposure on active catalase percentages from MRSA USA300 and *P. aeruginosa*. (a-b). Percent of active catalase remained inside MRSA USA300 (a) and *P. aeruginosa* (b) after different treatment schemes (different wavelengths under the same dosage). Quantification of catalase was obtained by an Amplex Red Catalase
25 kit. Data: Mean $\pm$ standard deviation (N=3).

[0054] Figure 3 depicts Resonance Raman spectra of bovine liver catalase powder with and without 410 nm exposure. 410 nm dose: 250 mW/cm². Raman spectrum acquisition time: 25 s. 532 nm excitation. Data: Mean $\pm$ SD from five spectra.

- [0055] Figure 4 depicts the comparison between ns-410 nm and CW-410 nm exposure on the
30 catalase photoinactivation effect from pure catalase solution (a), catalase from MRSA

USA300 (b), and catalase from *P. aeruginosa* (c). Quantification of catalase was obtained by an Amplex Red Catalase kit. Data: Mean±standard deviation (N=3). Student unpaired t-test, ***: $p < 0.001$; **: $p < 0.01$.

[0056] Figure 5 depicts CFU ml⁻¹ of stationary-phase MRSA USA300 methicillin-resistant *Staphylococcus aureus* (a), *Pseudomonas aeruginosa* (b), and *Salmonella enterica* (c) under the treatment of 22 mM H₂O₂ with/without the combination with various light exposure. Data: Mean±standard deviation (N=3). Student unpaired t-test, ***: $p < 0.001$; **: $p < 0.01$. 250 CFUs: limit of detection. Figure 5 further depicts the synergistic effect between photoinactivation of catalase and low-concentration hydrogen peroxide to eliminate stationary-phase MRSA USA300 and stationary-phase *Pseudomonas aeruginosa*. Left and right: CFU ml⁻¹ of stationary-phase MRSA and *P. aeruginosa* under different treatment schemes, respectively. N=3. Data: Mean±SD. ***: significant difference. $p < 0.001$. 250 CFUs: detection of limit.

[0057] Figure 6 depicts the killing efficacy comparison between CW-410 nm and ns-410 nm combined with H₂O₂ in both stationary-phase MRSA USA300 and *Pseudomonas aeruginosa*. Left and right: CFU ml⁻¹ of stationary-phase MRSA and *P. aeruginosa* under different treatment schemes, respectively. N=3. Data: Mean±SD. ***: significant difference. $p < 0.001$. 250 CFUs: detection of limit.

[0058] Figure 7 depicts CFU ml⁻¹ of *E. coli* BW25113 under different treatment schemes. Tobramycin: 2 µg/ml, 4-hour incubation. ***: $p < 0.001$, student unpaired t-test.

[0059] Figure 8 depicts CFU ml⁻¹ of *Enterococcus faecalis* NR-31970 under different treatment schemes. Tobramycin: 2 µg/ml, 4-hour incubation.

[0060] Figure 9 depicts confocal laser scanning microscopy of intracellular MRSA. (a-c). Fluorescence images of intracellular live MRSA (a), and dead MRSA (b), along with the transmission images (c) after MRSA infecting RAW 264.7 macrophage cells for 1 hour. (d-f). Fluorescence images of intracellular live MRSA (d), and dead MRSA (e), along with the transmission images (f) after ns-410 exposed MRSA infecting RAW 264.7 macrophage cells for 1 hour. (g-h). Quantitative analysis of live/dead MRSA from the above two groups. Scalar bar=10 µm.

[0061] Figure 10 depicts active catalase percent of various fungal strains with or without 410 nm light exposure. Dose: 410 nm, 150 mW/cm², 5 min. Fungal concentration: 10⁶ cells/ml. *C. albicans* CASC5314: wild-type *Candida albicans*.

[0062] Figure 11 depicts CFU results of *C. albicans* CASC5314 after different treatment schemes. (a). Time-killing assay of CASC5314 after various treatment schemes. (b). Spread plates of CASC5314 after 1-hour incubation at different treatment schemes.

[0063] Figure 12 depicts the synergistic effect between photoinactivation of catalase under various wavelengths and low-concentration hydrogen peroxide to eliminate stationary-phase CASC5314. CFU ml⁻¹ of CASC5314 after treatments under the combination between H₂O₂ and various wavelengths. Dosage: 40 mW/cm², 24 J/cm². H₂O₂: 44 mM, 1.5-hour incubation. Data: Mean±SEM (N=3). ##: detection limit.

[0064] Figure 13 depicts fluorescence signals of PrestoBlue from CASC5314 under various treatment schemes. (a, c). H₂O₂-alone treated stationary-phase CASC5314 and log-phase CASC5314, respectively. (b, d). 410nm plus H₂O₂ treated stationary-phase CASC5314 and log-phase CASC5314, respectively.

[0065] Figure 14 depicts Fluorescence signals of PrestoBlue of three different *C. auris* strains under different treatment schemes. (a, c, e). Amp B alone-treated groups. (b, d, f). 410 nm plus amp B-treated groups.

[0066] Figure 15 depicts confocal laser scanning imaging of live/dead *C. albicans* after infecting RAW264.7 macrophage cells.

[0067] Figure 16 depicts photoinactivation of catalase in combination with silver cation kills MRSA. Shown in the images are spread agar plates of MRSA USA300 under different treatment schemes.

[0068] Figure 17 depicts the comparison between CW-410 and ns-410 to inactivate catalase and eliminate *E. coli* BW25113 by synergizing with silver cation. (a-b). CFU ml⁻¹ of *E. coli* BW25113 after different treatment schemes: 30 min for (a) and 60 min for (b). Dose: 22 J/cm². Silver cation: 0.5 μM. Data: Mean ± SEM (N=3). ***: p<0.001, significant difference. Student unpaired t-test.

DETAILED DESCRIPTION OF THE INVENTION

Definitions

[0069] Unless otherwise defined, all technical and scientific terms used herein have the same meaning as commonly understood by one of ordinary skill in the art to which this invention belongs. In case of conflict, the present application, including definitions will control.

[0070] As used herein, the phrase “treating an infected tissue” refers to curing, alleviating or partially arresting the clinical manifestations of the infection or its complications. Treating an infected tissue achieves a medically desirable result. In some cases, this is a complete eradication of infection. In other cases, it is an improvement in the symptoms of the infection.

[0071] A “ROS-generating agent” is any biological or chemical agent that produces Reactive Oxygen Species (ROS). ROS-generating agents as defined herein, exclude exogenous photosensitizer agents that have been light-activated. A “photosensitizer” is a chemical compound, or a biological precursor thereof, that produces a phototoxic or other biological effect on biomolecules upon photoactivation.

[0072] A “subject” is a vertebrate, including any member of the class mammalia, including humans, domestic and farm animals, and zoo, sports or pet animals, such as mouse, rabbit, pig, sheep, goat, cattle and higher primates.

[0073] A “microbe” is a multi-cellular or single-celled microorganism, including bacteria, protozoa, and some fungi and algae. The term microbe, as used herein, includes pathogenic microorganisms such as bacterium, protozoan, or fungus.

[0074] The term “inanimate surface” refers any non-living surface.

[0075] The term “disinfecting” refers to destroying or eliminating pathogenic microorganisms that cause infections.

[0076] Unless specifically stated or clear from context, as used herein, the term “about” is understood as within a range of normal tolerance in the art, for example within 2 standard deviations of the mean. “About” is understood as within 10%, 9%, 8%, 7%, 6%, 5%, 4%, 3%, 2%, 1%, 0.5%, 0.1%, 0.05%, or 0.01% of the stated value. Unless otherwise clear from context, all numerical values provided herein are modified by the term about.

[0077] Ranges provided herein are understood to be shorthand for all of the values within the range. A range of 1 to 50 is understood to include any number, combination of numbers, or sub-range from the group consisting of 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17,



18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36, 37, 38, 39, 40, 41, 42, 43, 44, 45, 46, 47, 48, 49, or 50 (as well as fractions thereof unless the context clearly dictates otherwise). For example, the wavelengths from about 400nm to about 460nm include the wavelengths 400, 401, 402, 403, 404, 405, 406, 407, 408, 409, 410, 411, 412, 413, 414, 415, 416, 417, 418, 419, 420, 421, 422, 423, 424, 425, 426, 427, 428, 429, 430, 431, 432, 433, 434, 435, 436, 437, 438, 439, 440, 441, 442, 443, 444, 445, 446, 447, 448, 449, 450, 451, 452, 453, 454, 455, 456, 457, 458, 459 and 460nms. The light from about 5 J/cm² – to about 200 J/cm² includes 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36, 37, 38, 39, 40, 41, 42, 43, 44, 45, 46, 47, 48, 49, 50, 51, 52, 53, 54, 55, 56, 57, 58, 59, 60, 61, 62, 63, 64, 65, 66, 67, 68, 69, 70, 71, 72, 73, 74, 75, 76, 77, 78, 79, 80, 81, 82, 83, 84, 85, 86, 87, 88, 89, 90, 91, 92, 93, 94, 95, 96, 97, 98, 99, 100, 101, 102, 103, 104, 105, 106, 107, 108, 109, 110, 111, 112, 113, 114, 115, 116, 117, 118, 119, 120, 121, 122, 123, 124, 125, 126, 127, 128, 129, 130, 131, 132, 133, 134, 135, 136, 137, 138, 139, 140, 141, 142, 143, 144, 145, 146, 147, 148, 149, 150, 151, 152, 153, 154, 155, 156, 157, 158, 159, 160, 161, 162, 163, 164, 165, 166, 167, 168, 169, 170, 171, 172, 173, 174, 175, 176, 177, 178, 179, 180, 181, 182, 183, 184, 185, 186, 187, 188, 189, 190, 191, 192, 193, 194, 195, 196, 197, 198, 199, and 200 J/cm².

[0078] In this disclosure, “comprises,” “comprising,” “containing” and “having” and the like can have the meaning ascribed to them in U.S. Patent law and can mean “includes,” “including,” and the like; “consisting essentially of” or “consists essentially” likewise has the meaning ascribed in U.S. Patent law and the term is open-ended, allowing for the presence of more than that which is recited so long as basic or novel characteristics of that which is recited is not changed by the presence of more than that which is recited, but excludes prior art embodiments.

[0079] Other definitions appear in context throughout this disclosure.

Compositions and Methods of the Invention

[0080] Hydrogen peroxide (H₂O₂) is continuously produced inside microbes from autoxidation of the redox enzyme, and it diffuses quickly into the intracellular environment, causing an acutely detrimental effect (e.g. lipid peroxidation, DNA and protein damage) as a result of the Fenton reaction:

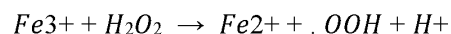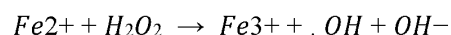

[0081] Photoinactivation of catalase creates potent antimicrobial effects due to a lethal accumulation of ROS. Photoinactivation of catalase further assists immune cells to eliminate intracellular pathogens. Neutrophils and macrophages are highly motile phagocytic cells that serve as the first line of defense of the innate immune system (Segal, A. W., *Annu Rev Immunol* 23, 197-223, doi:10.1146/annurev.immunol.23.021704.115653 (2005)). These cells play an essential role in providing resistance to bacterial and fungal infections through releasing ROS burst (e.g., superoxide, hydroxyl radicals, and singlet oxygen (Hampton, M. B., *Blood* 92, 3007-3017 (1998))). However, pathogens possess an array of elaborate strategies to invade and survive inside neutrophils or macrophages, thus acting as the ‘Trojan horses’ responsible for further dissemination and recurrent infections (Lehar, S. M. et al. *Nature* 527, 323-328 (2015)). Catalase, which is encoded by gene, *katA*, confers indispensable resistance for antimicrobial agents or reactive oxygen species released by immune cells (Flannagan, R., *Pathogens* 4, 826-868 (2015)). Photoinactivation of catalase assists macrophage and neutrophils to reduce the intracellular and extracellular bacterial burden.

[0082] In conducting the methods of the present invention, photoinactivation of catalase is preferably conducted with light having a wavelength of about 400nm to about 430nm, in combination with administration of a low-concentration peroxide solution and/or an ROS generating agent. Methods of the invention exclude the use of exogenous photosensitizing agents that have been activated by light.

[0083] Peroxide solutions include, but are not limited to solutions containing hydrogen peroxides, metal peroxides, and organic peroxides. Hydrogen peroxides include, but are not limited to, peroxy acids, peroxy monosulfuric acid, peracetic acid, peroxydisulfuric acid, peroxy nitric acid, peroxy nitrous acid, perchloric acid, and phthalimidoperoxycaproic acid. Metal peroxides include but are not limited to ammonium periodate, barium peroxide, sodium peroxide, sodium perborate, sodium persulfate, lithium peroxide, magnesium peroxide, magnesium perchlorate and zinc peroxide. Organic peroxides include but are not limited to acetone peroxide, acetozone, alkenyl peroxide arachidonic acid 5-hydroperoxide,

- artelinic acid, artemether, artemisinin, artemotil, arterolane, artesunate, ascaridole, benzoyl peroxide, bis(trimethylsilyl) peroxide, tert-butyl hydroperoxide tert-butyl peroxybenzoate, CSPD ([3-(1-chloro-3'-methoxyspiro[adamantane-4,4'-dioxetane]-3'-yl)phenyl] dihydrogen phosphate), cumene hydroperoxide, di-tert-butyl peroxide, diacetyl peroxide, diethyl ether
- 5 peroxide, dihydroartemisinin, dimethyldioxirane, 1,2-dioxane, 1,2-dioxetane, 1,2-dioxetanedione, dioxirane, dipropyl peroxydicarbonate, ergosterol peroxide, hexamethylene triperoxide diamine, methyl ethyl ketone peroxide, nardosinone, paramenthane hydroperoxide, perfosfamide, peroxyacetyl nitrate, peroxyacyl nitrates, prostaglandin h2, 1,2,4-trioxane, and verruculogen.
- 10 Other peroxides include, potassium peroxydisulfate, bis(trimethylsilyl) peroxide ($\text{Me}_3\text{SiOOSiMe}_3$), phosphorus oxides, ammonium peroxide, copper(II) peroxide, sodium peroxide, cobalt(II) peroxide, mercury(I) peroxide, iron(II) peroxide potassium peroxide, copper(I) peroxide, rubidium peroxide, cesium peroxide, iron(III) peroxide, beryllium peroxide, magnesium peroxide, nickel(II) peroxide, cadmium peroxide, barium
- 15 peroxide, benzoyl peroxide, calcium peroxide, diacetyl peroxide, cesium superoxide, lead(IV) peroxide, lithium peroxide, gallium(II) peroxide, chromium(III) peroxide, mercury(II) peroxide, gold(I) peroxide, strontium peroxide, zinc peroxide, potassium superoxide, and chromium(VI) peroxide.
- [0084]** In other specific embodiments, the diluted peroxide solution is a hydrogen peroxide
- 20 solution formulated with between about 0.030% and about 0.3% hydrogen peroxide (which converts to about 8.8 mM to about 88 mM hydrogen peroxide).
- [0085]** Photoinactivation of catalase and administration of the peroxide solution can also be provided in combination with ROS generating agents including antibiotics, such as tobramycin. Other ROS generating agents include, but are not limited to, silver cation, iodine
- 25 tincture, gold nanoparticles, methylene blue (non-photoactivated), β -lactam antibiotics, aminoglycosides, fluoroquinolones, antifungal azoles, membrane-targeting polyene antifungals, such as amphotericin B, and cell-wall targeting antifungals, such as caspofungin.
- [0086]** Typically following photoinactivation, the peroxide solution can be administered to the site of the infection for a duration of about 10 to about 30 minutes. In alternate
- 30 embodiments, the peroxide solution, the ROS generating agent and/or the photoinactivating

light can be administered concomitantly or sequentially to the site of infection. For example, in specific embodiments, the ROS-generating agent is administered prior to photoinactivation of catalase. In other specific embodiments, the ROS-generating agent is administered after the photoinactivation of catalase. Preferably, the peroxide solution is topically administered
5 (e.g., as a liquid or a spray). Administration of the ROS generating agent can be according to all modes of local or systemic administration known in the art.

[0087] In one embodiment, methods of the invention comprising photoinactivation of catalase are directed to an infected external tissue of a subject, including, but not limited to, skin, hair and nails. In other embodiments, internal tissues, such as gastrointestinal organs or
10 cavities (oral, vaginal or nasal cavities), may be targeted as well.

[0088] Peroxide solutions and/or ROS generating agents can be administered alone or as a component of a pharmaceutical formulation. The compounds may be formulated for administration, in any convenient way for use in human or veterinary medicine. Wetting agents, emulsifiers and lubricants, such as sodium lauryl sulfate and magnesium stearate, as
15 well as coloring agents, release agents, and preservatives can also be present in the compositions.

[0089] Pharmaceutical formulations of the invention include those suitable for intradermal, inhalation, oral/ nasal, topical, parenteral, rectal, and/or intravaginal administration. The formulations may conveniently be presented in unit dosage form and may be prepared by any
20 methods well known in the art of pharmacy. The amount of active ingredient which can be combined with a carrier material to produce a single dosage form will vary depending upon the host being treated, and the particular mode of administration, e.g., intradermal.

[0090] The formulations can include a pharmaceutically acceptable carrier. The amount of active ingredient which can be combined with a carrier material to produce a single dosage
25 form will generally be that amount of the compound which produces a therapeutic effect. Formulations of the invention can be administered parenterally, intraperitoneally, subcutaneously, topically, orally (e.g., the ROS generating agent) or by local administration, such as by aerosol or transdermally. Formulations can be administered in a variety of unit dosage forms depending upon the severity of the infection or the site of the infection and the
30 degree of illness, the general medical condition of each patient, the resulting preferred

method of administration and the like. Details on techniques for formulation and administration of pharmaceuticals are well described in the scientific and patent literature, see, e.g., the latest edition of Remington's Pharmaceutical Sciences, Maack Publishing Co, Easton PA ("Remington's").

5 **[0091]** Pharmaceutical formulations of this invention can be prepared according to any method known to the art for the manufacture of pharmaceuticals. A formulation can be admixed with nontoxic pharmaceutically acceptable excipients which are suitable for manufacture. Formulations may comprise one or more diluents, emulsifiers, preservatives, buffers, excipients, etc. and may be provided in such forms as liquids, powders, emulsions,
10 lyophilized powders, sprays, creams, lotions, controlled release formulations, tablets, pills, gels, on patches, in implants, etc.

[0092] In practicing this invention, the pharmaceutical formulations can be delivered transdermally, by a topical route, formulated as applicator sticks, solutions, suspensions, emulsions, gels, creams, ointments, pastes, jellies, paints, powders, and aerosols. In specific
15 embodiments, delivery can be mediated by a transdermal patch, bandage or dressing impregnated with compositions comprising the peroxide solution and/or ROS generating agent. Sustained release can be provided by transdermal patches, for slow release at the site of infection.

[0093] The amount of pharmaceutical formulation adequate to reduce or eradicate pathogenic
20 microbes is a therapeutically effective dose. The dosage schedule and amounts effective for this use, i.e., the dosing regimen, will depend upon a variety of factors, including the stage of the infection, the severity of the infection, the general state of the patient's health, the patient's physical status, age and the like. In calculating the dosage regimen for a patient, the mode of administration also is taken into consideration.

25 **[0094]** The dosage regimen also takes into consideration pharmacokinetics parameters well known in the art, i.e., the active agents' rate of absorption, bioavailability, metabolism, clearance, and the like (see, e.g., Hidalgo-Aragones (1996) *J. Steroid Biochem. Mol. Biol.* 58:611-617; Groning (1996) *Pharmazie* 51:337-341; Fotherby (1996) *Contraception* 54:59-69; Johnson (1995) *J. Pharm. Sci.* 84:1144-1146; Rohatagi (1995) *Pharmazie* 50:610-613;
30 Brophy (1983) *Eur. J. Clin. Pharmacol.* 24:103-108; the latest Remington's, supra). The state

of the art allows the clinician to determine the dosage regimen for each individual patient, active agent and disease or condition treated. Guidelines provided for similar compositions used as pharmaceuticals can be used as guidance to determine the dosage regimen, i.e., dose schedule and dosage levels, administered practicing the methods of the invention are correct and appropriate.

[0095] Single or multiple administrations of pharmaceutical formulations of the invention can be given depending on for example: the dosage and frequency as required and tolerated by the patient, the persistence of infection, or lack thereof, after each administration, and the like. The formulations should provide a sufficient quantity of peroxide solution to effectively treat, prevent or ameliorate the infection.

[0096] Methods of the invention target catalase positive microbes which are associated with, or may give rise to, infection. Both Gram-negative and Gram-positive bacteria serve as infectious pathogens in vertebrate animals. Such catalase positive Gram-positive bacteria include, but are not limited to, *Staphylococci* species. Catalase positive Gram-negative bacteria include, but are not limited to, *Escherichia coli*, *Pasteurella* species, *Pseudomonas* species (e.g., *P. aeruginosa*), and *Salmonella* species. Specific examples of infectious catalase positive bacteria include but are not limited to, *Helicobacter pylori*, *Borelia burgdorferi*, *Legionella pneumophila*, *Mycobacteria* species (e.g. *M. tuberculosis* complex, *M. avium* complex, *M. gordonae* clade, *M. kansasii* clade, *M. nonchromogenicum/terrae* clade, *Mycobacterium* producing mycobacteria, *M. simiae* clade, *M. abscessus* clade, *M. chelonae* clade, *M. fortuitum* clade, *M. mucogenicum* clade, *M. parafortuitum* clade, *M. vaccae* clade, *M. ulcerans*, *M. vanbaalenii*, *M. gilvum*, *M. bovis*, *M. leprae*, *M. spyr1*, *M. kms*, *M. mcs*, *M. jls*, *M. intracellulare*, and *M. gordonae*), *Acinetobacter baumannii*, *Staphylococcus aureus*, *Neisseria gonorrhoeae*, *Neisseria meningitidis*, *Listeria monocytogenes*, pathogenic *Campylobacter* species, *Haemophilus influenzae*, *Bacillus anthracis*, *Corynebacterium diphtheriae*, *Corynebacterium* species, *Erysipelothrix rhusiopathiae*, *Chlamydia trachomatis*, *Clostridium perfringens*, *Clostridium tetani*, *Klebsiella pneumoniae*, *Pasteurella multocida*, *Bacteroides* species, *Fusobacterium nucleatum*, *Treponema pallidum*, *Treponema pertenue*, *Leptospira*, *Rickettsia*, and *Actinomyces israeli*. *Mycoplasma* and *Chlamydia* species.

[0097] Examples of catalase positive fungi include, but are not limited to, *Aspergillus fumigatus*, *Cryptococcus neoformans*, *Histoplasma capsulatum*, *Coccidioides immitis*, *Blastomyces dermatitidis*, *Candida glabrata*, *Candida tropicalis*, *Candida parapsilosis*, and other catalase-positive *Candida* spp. *Candida auris*, and *Trichophyton rubrum*.

5 [0098] The light for photoactivation of catalase can be produced and delivered to the site of infection by any suitable means known in the art.

[0099] While it has now been determined that the antimicrobial effect of photoinactivation is significantly and unexpectedly increased upon administration of a low-concentration of H₂O₂ and/or a ROS-generating agent, the antimicrobial effectiveness is also significantly improved
10 when the light is provided by a pulsed nanosecond laser compared to continuous wavelength LED. Accordingly, in specific embodiments, the light source is a pulsed nanosecond laser. Pulsed operation of lasers refers to any laser not classified as continuous wave, so that the optical power appears in pulses of some duration at some repetition rate. Nanosecond laser families can range from the UV to the IR with wavelengths up to 1064 nm, repetition
15 rates up to 2 kHz, and pulse energy up to 20 mJ. Photoinactivation of catalase can be conducted with light having a wavelength of about 400nm to about 460nm. In specific embodiments, the wavelength is about 400nm to about 430nm, applied at a dosage of about 5 J/cm² – to about 200 J/cm², and in other specific embodiments, about 14 J/cm² to about 32 J/cm². In other specific embodiments, the pulse duration is about 5 nanoseconds. Light
20 delivered in this range by a pulsed nanosecond laser is clinically advantageous because thermal damage is minimal, temporary or otherwise non-existent. In more specific embodiments, the wavelength is 410nm (delivered using about 15 J/cm²), applied by a pulsed nanosecond laser according to methods known in the art for operation of such lasers.

[00100] Exposure times range from about 5 to about 10 minutes in length, and can be
25 repeated weekly as needed, for example, about twice per week for several months. In clinical applications, patients may receive treatment for between one to 3 months or longer as determined by the practicing physician.

[00101] In other embodiments, photoactivating light can be delivered to the site of infection through various optical waveguides, such as an optical fiber or implant. In some
30 embodiments, the photoinactivating light is delivered by optical fiber devices that directly

illuminate the site of infection. For example, the light can be delivered by optical fibers threaded through small gauge hypodermic needles. In addition, light can be transmitted by percutaneous instrumentation using optical fibers or cannulated waveguides. For open surgical sites, suitable light sources include broadband conventional light sources, broad arrays of LEDs, and defocused laser beams. The light source can be operated in the Continuous Wave (CW) mode. Photoinactivation of catalase is preferably conducted with light having a wavelength of about 400nm to about 430nm and a dosage of about 5 J/cm² – to about 200 J/cm², and in specific embodiments, about 14 J/cm² to about 32 J/cm². In other specific embodiments, the wavelength is 410nm (delivered using about 15 J/cm²), applied by a CW LED according to methods known in the art for operation of such LED sources. Exposure times range from any light source range from about 5 to about 10 minutes in length.

[00102] In other embodiments of the invention, the photoinactivation of catalase is performed on an inanimate surface including but not limited to metal, plastic, fabric, rubber, stone, composite surfaces or wood. In specific embodiments, the inanimate surface comprises objects such as instruments, catheters, medical and military equipment, furniture, handrails, textiles, fixtures such as sinks and plumbing materials, building materials, industrial or electronic equipment, and food product or food processing equipment. Photoinactivation of catalase on inanimate surfaces is preferably conducted with light having a wavelength of about 400nm to about 430nm, and in combination with administration of a solution having a low-concentration of a peroxide. In specific embodiments, the wavelength is 410nm (delivered at 15 J/cm²), applied by a pulsed nanosecond laser. Exposure times range from about 5 to about 10 minutes in length.

[00103] The following examples are put forth for illustrative purposes only and are not intended to limit the scope of what the inventors regard as their invention.

25

EXAMPLES

[00104] The following materials and methods were employed throughout Examples 1.-4.

[00105] Bacterial strains: *Enterococcus faecalis* NR-31970, *Enterococcus faecalis* HM-325, *Escherichia coli* BW 25113, *Escherichia coli* ATCC 25922. *Klebsiella pneumoniae* ATCC BAA 1706. *Klebsiella pneumoniae* ATCC BAA 1705. *Salmonella enterica* ATCC

30

70720. *Salmonella enterica* ATCC 13076. *Acinetobacter baumannii* ATCC BAA 1605. *Acinetobacter baumannii* ATCC BAA-747. *Pseudomonas aeruginosa* ATCC 47085 (PAO-1). *Pseudomonas aeruginosa* 1133. *Pseudomonas aeruginosa* ATCC 15442. *Pseudomonas aeruginosa* ATCC 9027.

- 5 **[00106]** Quantitation of catalase by Amplex red catalase kit: Quantification of catalase both from the pure catalase solution and bacteria was achieved by a fluorescent amplex red catalase kit. 25 μ l of analyte were incubated with 25 μ l (40 μ M of H_2O_2) for 30 min at room temperature. Then 50 μ l of working solution (100 μ M Amplex Red reagent containing 0.4 U/ml horseradish peroxidase) were added to the abovementioned mixture, and the subsequent
- 10 mixture were incubated for another 30–60 min in the dark. After that, the fluorescence was recorded at an emission of 590 nm when excited at 560 nm.

- [00107]** Resonance Raman spectrum of dried catalase film: Catalase was measured by its Raman peaks at around 1300–1700 cm^{-1} measured by resonance Raman spectroscopy (1221, LABRAM HR EVO, Horiba) with a 40 $\times$ objective (Olympus) and an excitation wavelength
- 15 of 532 nm. Samples (dried ‘coffee ring’ were sandwiched between two glass cover slides (48393-230, VWR international) with a spatial distance of $\sim$ 80 μ m. To study the photoinactivation (by a continuous-wave LED), the same samples were measured after each laser irradiation.

- [00108]** CFU experiments to test the potential synergy between photoinactivation of
- 20 catalase and H_2O_2 : Overnight-cultured bacteria was centrifuged, the supernatant was discarded, and the pellet was resuspended with the same amount of PBS. The laser source used in the study is nanosecond (ns) pulsed OPO laser purchased from OPOTEK Inc, model number as Opolette HE355 LD, having the following key specifications: wavelength range, 410–2400 nm; pulse repetition rate, 20 Hz; maximum pulse energy at 460 nm, 8 mJ; pulse
- 25 duration, 5 nanosecond; spectral linewidth, 4–6 cm^{-1} ; and pulse-to-pulse stability, <5%. For each bacterial strain there were four groups: untreated one, ns-410 nm-treated group, H_2O_2 (22 or 44 mM)-treated group, ns- 410 nm plus H_2O_2 (22 or 44 mM)-treated group. Dose for ns-410 nm exposure was 15 J/cm². H_2O_2 was incubated with bacteria for 30 min at 37°C with the shaking speed of 200 rpm. After incubation, bacterial burden from each group was serial
- 30 diluted, inoculated onto TSA plates, then counted by enumeration of these plates.

[00109] CFU experiments to test the potential synergy between photoinactivation of catalase and ROS-generating antibiotics: Overnight-cultured bacteria was centrifuged, and then the supernatant was discarded and re-suspended with the same amount of fresh TSB. Then prior to any treatments, the above solution was incubated with antibiotics (10 µg/ml) for 1 hour. For each bacterial strain, four groups were tested: untreated one, ns-410 nm-treated group, antibiotic (2 µg/ml)-treated group, ns-410 nm plus antibiotic (2 µg/ml)-treated group. Dose for ns-410 nm exposure was 15 J/cm². Antibiotic was incubated with bacteria for up to 6 hours at 37°C with the shaking speed of 200 rpm. At each time interval, bacterial burden from each group was serially diluted, inoculated onto TSA plates, then counted by enumeration of these plates.

[00110] Confocal imaging of intracellular bacteria assay: As described elsewhere (Yang, X., et al. International journal of nanomedicine 13, 8095 (2018)), murine macrophage cells (RAW 264.7) were cultured in DMEM supplemented with 10% FBS at 37 degrees C with CO₂ (5%). Cells were exposed to MRSA USA300 or Salmonella enterica (with/without ns-410 nm exposure) at a multiplicity of infection (MOI) of approximately 100:1 at serum-free DMEM medium. 1 or 2-hour post-infection, RAW 264.7 cells were washed with gentamicin (50 µg/mL, for one hour) to kill extracellular bacteria in DMEM + 10% FBS. After that, RAW 264.7 cells were washed with gentamicin (50 µg/mL) and subsequently lysed using 0.1% Triton-X 100 for 3 min. After membrane permeabilization, infected RAW 264.7 cells were stained with Live/Dead stain for 15 min, then samples were fixed in 10% formalin for 10 min prior to confocal imaging.

Example 1: Pulsed Blue Laser Effectively Inactivates Pure Catalase and Catalase from Bacteria

[00111] Pure catalase solution (bovine liver catalase, 3 mg/ml in the PBS) was prepared in PBS using a protocol previously published to examine the effect of 410-nm exposure on the absorption spectrum of catalase solution (Cheng, L., Photochemistry and Photobiology 34, 125-129 (1981)). Catalase shows a pronounced absorption at around 410 nm, and its absorption at this wavelength gradually decreases as the 410-nm exposure elongates (Figure 1a). This suggests that the secondary structure of catalase might be changed, especially in the active heme-containing domain. In addition, this photoinactivation effect was examined by an

Amplex Red Catalase kit at different wavelengths (Figure 1b). The photoinactivation trend is similar to the absorption spectrum of catalase, with the 410 nm being the most effective, where 5-min exposure depleted ~70% active catalase.

[00112] Since most of the aerobic bacteria and facultative anaerobes express catalase (Mishra, S. & Imlay, J. Arch Biochem Biophys 525, 145-160, doi:10.1016/j.abb.2012.04.014 (2012)), whether one could photoinactivate catalase in situ from the catalase-positive bacteria was examined. MRSA USA300 and *P. aeruginosa* (PAO-1) were selected as the representative for Gram- positive and Gram-negative bacteria, respectively. Noteworthy, catalase from both MRSA USA300 (Figure 2a) and *P. aeruginosa* (Figure 2b) were photoinactivated by blue light exposure region, especially 410-nm exposure. The dose utilized was about 15 J/cm², well below the ANSI safety limit of 200 J/cm², and the specimens were stationary-phase cultured bacteria (~10⁸ cells/ml). ANSI is the American National Standard for Safe Use of Lasers, see ANSI Z136.1, Laser Institute of America 2014.

[00113] To further understand how 410 nm exposure could cause the structural change of catalase, Resonant Raman spectroscopy was performed to capture the Raman signature of dried catalase film (Figure 3). Apparently, 410 nm exposure significantly drops the Raman intensity at 750 cm⁻¹, and the Raman bands ranging from 1300 cm⁻¹ to 1700 cm⁻¹, which are typical vibrational bands of heme ring from catalase (Chuang, W.-J., Heldt, J. & Van Wart, H. J. J. o. B. C. Resonance Raman spectra of bovine liver catalase compound II. Similarity of the heme environment to horseradish peroxidase compound II. 264, 14209-14215 (1989)). These data further consolidate the fact that 410 nm exposure could cause structural change of catalase.

[00114] In addition, the efficacy between ns-410 nm and CW-410 nm to inactivate catalase was compared. ns-410 nm is significantly more effective both in the pure solution form (Figure 4a), or from MRSA USA300 (Figure 4b) and *P. aeruginosa* (Figure 4c) compared to CW-410 nm. Moreover, ns-410 exposure eliminates the necessity of heating tissue during future clinical study.

Example 2: Photo-inactivation of Catalase Sensitizes a Wide Range of Bacteria to Low-Concentration H₂O₂

[00115] Catalase is an essential detoxifying enzyme in bacteria encountering various endogenous or exogenous stress (Nakamura, K. et al. Microbiology and immunology 56, 48-55 (2012)). When the gene encoding the expression of catalase is mutant, pathogens are more susceptible to the environmental stress (Mandell, G. L., J Clin Invest 55, 561-566, doi:10.1172/jci107963 (1975)). Whether exogenous addition of low-concentration H₂O₂ could eliminate those ‘traumatized’ pathogens was investigated. As shown in Figure 5, photoinactivation of catalase (15 J/cm²) alone didn’t reduce the MRSA burden (Figure 5a), *P. aeruginosa* burden (Figure 5b), and *Salmonella enterica* burden (Figure 5c) significantly. Moreover, low-concentration H₂O₂ (22 mM) didn’t exert any significant antimicrobial effect against both MRSA and *P. aeruginosa* (Figure 5). However, subsequent administration of low-concentration H₂O₂ after photoinactivation of catalase significantly reduced the MRSA and *P. aeruginosa* burden (≥ 3 -log₁₀ reduction, Figure 5). Interestingly, the bacterial killing trend versus irradiance wavelength is similar to that of photoinactivation of catalase versus irradiance wavelength. Noteworthy, low-concentration H₂O₂ combined with 410 nm exposure (15 J/cm²) achieved total eradication of *P. aeruginosa* (Figure 5b).

Example 3. Photoinactivation of Catalase and Low-concentration Hydrogen Peroxide Create a Synergistic Effect

[00116] There is a synergistic effect between photoinactivation of catalase and low-concentration hydrogen peroxide to eliminate stationary-phase MRSA USA300 and stationary-phase *Pseudomonas aeruginosa*. Figure 5a depicts the synergistic results in a bar-graph. CFU ml⁻¹ (colony-forming unit) designates the bacterial burden. ‘Untreated’ refers to the original stationary-phase MRSA without any exogenous treatment. ‘H₂O₂ (22 mM, 0.075%)’ and ‘ns-light’ refer to stationary-phase MRSA with H₂O₂ and ns light alone, respectively. As shown in the graph, H₂O₂ alone and ns-light alone do not exert any significant killing effect on MRSA, however, ns-410 nm in combination with H₂O₂ reduces approximately four orders of magnitude of bacterial burden. The same phenomenon happens with other wavelengths as well. Noteworthy, ns-430 nm or ns-430 nm combined with H₂O₂ reduces around 99% of the bacterial burden under the same conditions. ns-450 or ns-460 nm combined with H₂O₂ together reduces around 90% of the bacterial burden. ns-470 nm combined with H₂O₂ together reduces around 50% of the bacterial burden. ns-480 nm

combined with H₂O₂ barely exerts an antimicrobial effect. Altogether, the killing effect of H₂O₂ is significantly enhanced by blue light photoinactivation of catalase, especially when applied using ns-410 nm. A similar phenomenon occurred with stationary-phase *Pseudomonas aeruginosa*, which is a representative of Gram-negative bacteria (Figure 5b) and *Salmonella enterica*. By employing ns-410 nm combined with H₂O₂ to *Salmonella enterica*, an enhanced killing effect of around five orders of magnitude was observed (Figure 5c).

[00117] In addition, ns-410 nm combined with H₂O₂ is significantly more effective in eliminating microbes compared to CW-410 nm combined with H₂O₂ (Figure 6).

10 Example 4: Photoinactivation of Catalase Revives Conventional Antibiotics Against a Wide Range of Bacteria

[00118] Besides H₂O₂, whether photoinactivation of catalase could synergize with conventional antibiotics was investigated, especially for antibiotics that can generate the downstream intracellular ROS. Tobramycin, a representative of aminoglycoside, is an example. Tobramycin can induce downstream ROS burst (Dwyer, D. J. et al. Proceedings of the National Academy of Sciences 111, E2100-E2109, doi:10.1073/pnas.1401876111 (2014)), thus the combination of photoinactivation of catalase and tobramycin administration, together, was tested to see whether an enhanced effect was observed.

[00119] Interestingly, enhanced killing effect was observed in the combination-treated group (Figure 7). More than 100 times enhancement suggests that photoinactivation of catalase indeed accelerates the antimicrobial effect of ROS-generating antibiotics. As a control, the same treatment schemes were tested on a catalase-negative *Enterococcus* strain, *Enterococcus faecalis* NR-31970, which did not produce any enhanced killing effect (Figure 7). Altogether, this indicates that photoinactivation of catalase helps to revive traditional antibiotics against catalase-positive pathogens.

25 Example 5: Photoinactivation of Catalase Assists Macrophage Cells Against Intracellular Pathogens

[00120] Neutrophils and macrophage cells are highly essential phagocytic cells that serve as the first line of defense of the innate immune system (Segal, A. W., Annu Rev Immunol 23, 197-223, doi:10.1146/annurev.immunol.23.021704.115653 (2005)). Catalase, which is

encoded by gene, *katA*, confers indispensable resistance to the antimicrobial agents released by immune cells (Flannagan, R., Heit, B. & Heinrichs, D., *Pathogens* 4, 826-868 (2015)).

Based on these facts, it was hypothesized that photoinactivation of catalase could assist immune cells to eliminate extracellular and intracellular pathogens. To test the potential

5 assistance effect, a fluorescent Live/Dead assay was used to visualize the intracellular live/dead bacteria after ns-410 nm exposure. A higher percent of dead MRSA was observed intracellularly (Figure 9).

[00121] In conclusion, photoinactivation of catalase significantly boosts the efficacy of low-concentration H₂O₂, ROS-generating antibiotics, and immune cells against broad-

10 spectrum bacteria, including the notorious drug-resistant gram-negative bacteria.

[00122] The following materials and methods were employed throughout Examples 5.-9.

[00123] Chemicals and fungal strains: DMSO (W387520, Sigma Aldrich), amphotericin B (A9528-100 MG, Sigma Aldrich), ergosterol (AC1178100050, 98%, ACROS Organics).

YPD broth (Y1375, Sigma Aldrich). YPD agar (Y1500, Sigma Aldrich). PrestoBlue cell
15 viability assay (A13262, Thermo Fisher Scientific). *Candida albicans* SC5314, the test of fungal strains used see Table 1.

Table 1. Fungal strains utilized for amp B imaging experiments.

Pathogen	Strain	Number
<i>C. albicans</i>	SC5314	SC5314
<i>C. glabrata</i>	ATCC2001	ATCC2001
<i>C. tropicalis</i>	C22	H3222861
<i>C. parapsilosis</i>	C23	F825987
<i>C. lusitanae</i>	C30	S1591976
<i>Candida auris</i>	Lung	From MGH-- commonly used in MKM lab
<i>Candida haemulonii</i>	CAU-13	AR-0393
<i>Candida duobushaemulonii</i>	CAU-14	AR-0394
<i>Candida haemulonii</i>	CAU-15	AR-0395
<i>Kodameae ohmeri</i>	CAU-16	AR-0396
<i>C. albicans</i>	Ca C13	
<i>C. albicans</i>	Ca C14	
<i>C. glabrata</i>	Cg C1	
<i>C. glabrata</i>	Cg C2	
<i>Candida krusei</i>	CAU-17	AR-0397
<i>C. lusitanae</i>	CAU-18	AR-0398
<i>Saccharomyces cerevisiae</i>	CAU-19	AR-0399
<i>C. albicans</i>	Ca C15	
<i>C. albicans</i>	Ca C16	
<i>C. albicans</i>	Ca C17	
<i>Candida krusei</i>	CAU-17	AR-0397
<i>C. lusitanae</i>	CAU-18	AR-0398
<i>Saccharomyces cerevisiae</i>	CAU-19	AR-0399
<i>C. albicans</i>	Ca C15	
<i>C. albicans</i>	Ca C16	
<i>C. albicans</i>	Ca C17	
<i>Candida auris</i>	CAU4	AR-0384
	CAU5	AR-0385
	CAU6	AR-0386
	CAU7	AR-0387
	CAU8	AR-0388
	CAU9	AR-0389

[00124] Quantification of catalase from fungus before and after 410 nm exposure:

Quantification of catalase both from the pure catalase solution and fungal solution were achieved by a fluorescent amplex red catalase kit. Basically, 25 µl of analyte were incubated with 25 µl (40 µM of H₂O₂) for 30 min at room temperature. Then 50 µl of working solution (100 µM Amplex Red reagent containing 0.4 U/ml horseradish peroxidase) were added to the abovementioned mixture, and the subsequent mixture was incubated for another 30-60 min in the dark. After that, the fluorescence was recorded at an emission of 590 nm when excited at 560 nm.

[00125] CFU test to quantify the treatment efficacy: Quantification of antifungal treatment

schemes were achieved as following: overnight cultured fungal specimen was washed by sterile PBS. And log-phase fungal pathogens were prepared by dilution into fresh YPD broth at a ratio of 1:50 and cultured for another 2-3 hours at 30 °C with the shaking speed of 200 rpm. After that, the fungal concentration was adjusted to be around 1×10⁸ cells/ml by centrifuging or further dilution with PBS. 10 µl of the above fungal solution was exposed to 410 nm for 5 min (150 mW/cm²). After that, the exposed sample was collected into 990 µl of sterile PBS, then supplemented with treatment agents. Later, CFU of fungal cells was enumerated by serial dilution and cultured in YPD agar plates for 48 hours.

[00126] PrestoBlue viability assay: First log-phase fungal pathogens were prepared by diluting overnight-cultured fungal pathogens into fresh YPD broth at a ratio of 1:50 and

cultured for another 2-3 hours at 30 C with the shaking speed of 200 rpm. After that, the fungal concentration was adjusted to be around 1×10⁸ cells/ml by centrifuging or further dilution with PBS. 10 µl of the above fungal solution was exposed to 410 nm for 5 min (150 mW/cm²). After that, the exposed sample was collected into 990 µl of sterile PBS, then supplemented with treatment agents. Aliquots were made from the above sample into a 96-well plate, with each well containing 100 µl. Then 100 µl sterile YPD broth and 23 µl of PrestoBlue were added into the same well. Fluorescence signal at 590 nm from each well was recorded in a time-course (up to 18 hours with the interval of 30 min) manner at an excitation of 560 nm. For each strain, in order to know the exact number of fungal pathogens in each well, the corresponding fluorescence signals were recorded from fungal pathogens with

known numbers, however no external treatments.

[00127] Macrophage-Candida albicans interaction unveiled by confocal laser scanning microscopy: As described elsewhere (Yang, X., et al. International journal of nanomedicine 13, 8095 (2018)), murine macrophage cells (RAW 264.7) were cultured in DMEM supplemented with 10% FBS plus penicillin and streptomycin at 37 C with CO₂ (5%). Cells were exposed to Candida albicans SC5314 (with/without 410 nm exposure) at a multiplicity of infection (MOI) of approximately 10:1 at serum-free DMEM medium. 1 or 2-hours post-infection, RAW 264.7 cells were washed with gentamicin (50 µg/mL, for one hour) to kill extracellular pathogens in DMEM + 10% FBS. After that, RAW 264.7 cells were washed with gentamicin (50 µg/mL) and subsequently lysed using 0.1% Triton- X 100 for 3 min. After membrane permeabilization, infected RAW 264.7 cells were stained with Live/Dead stain for 15 min, then samples were fixed in 10% formalin for 10 min. Formalin was washed away prior confocal imaging.

Example 5: 410 nm Exposure Reduces Intracellular Catalase Amount

[00128] It is known that most fungal pathogens are catalase positive (Hansberg, W., et al. Arch Biochem Biophys 525, 170-180 (2012)). To test whether 410 nm exposure could cause the loss of catalase activity, the same approach to quantify the intracellular catalase amount by the amplex red catalase kit was utilized before and after 410 nm exposure. Catalase from various fungal pathogens, either log- phase or stationery-phase could be significantly inactivated by 410 nm exposure (Figure 10). Noteworthy, catalase from notorious Candida auris strain reduced by 60% after only 5- min 410 nm exposure.

Example 7: Photoinactivation of Catalase in Combination with H₂O₂ Achieved Total Eradication of C. albicans SC5314 by CFU Assay

[00129] Since catalase was effectively inactivated among various fungal strains, whether photoinactivation of catalase could sensitize fungal strains to external H₂O₂ attack was investigated. With further administration of low-concentration H₂O₂ after 410 nm exposure, eradication was achieved after combinational treatments (Figure 11). Noteworthy, there was more than five orders of magnitude enhancement of the function of H₂O₂ after photoinactivation of catalase (Figure 11).

[00130] Synergism between photoinactivation of catalase and H₂O₂ to eliminate Candida albicans SC5314 was also observed. The result is shown by a scatter plot in Figure 12. CFU

ml-1 (colony-forming unit) refers to the number of bacterial burden. 'Untreated' means the original stationary-phase SC5314 without any exogenous treatment. 'H₂O₂ (44 mM, 0.15%) and 'ns-light' means stationary-phase SC5314 with H₂O₂ and ns-light alone, respectively.

As shown in Figure 12, H₂O₂ alone and ns-light alone doesn't exert significant killing effect

- 5 on CASC5314, however, ns-light in combination with H₂O₂ reduces around four orders of magnitude of bacterial burden. Especially, ns-410 or ns-420, ns-430 combined with H₂O₂ achieved total eradication. ns-450 or ns-480 nm combined with H₂O₂ reduced a similar amount of fungal burden as H₂O₂-alone. Altogether, the killing effect of H₂O₂ is significantly enhanced by photoinactivation of catalase by blue light, especially by ns-410-ns-430 nm.
- 10 Therefore, an effective synergy exists between photoinactivation of catalase under the blue light range and H₂O₂ to eliminate CASC5314.

Example 8. Photoinactivation of Catalase in Combination with H₂O₂ Achieved Efficient Eradication of Broad-spectrum Fungal Species by PrestoBlue Assay

[00131] To further confirm that this combinational therapy works as well for other fungal strains, more clinical fungal strains were tested for feasibility of this synergistic therapy.

Unlike bacteria, fungal cells growth is slower, with each colony forming after around 48 hours. Thus, a high-throughput method, PrestoBlue viability assay, was used to measure the treatment efficacy. As shown in Figure 13, the utilization of PrestoBlue could achieve the same killing effect as the CFU assay. Interestingly, log-phase and stationary-phase

- 20 CASC5314 demonstrate different behavior towards the combinational killing, presumably because of the difference in metabolic activity between these two states. However, either log-phase or stationary-phase, photoinactivation of catalase always boosts the killing effect of low-concentration H₂O₂. This synergistic therapy was tested among more than twenty different clinical fungal isolates, and significant killing was consistently found among them.

25 Example 9. *Candida auris* Strains Are Sensitive to 410 nm Light Exposure

[00132] Apart from H₂O₂, whether photoinactivation of catalase was capable of synergizing with conventional antifungal agents, such as azoles or amphotericin B (amp B) was investigated. Similar to some classes of antibiotics, amp B kills fungi partly due to the oxidative damage (Belenky, P. et al. Fungicidal drugs induce a common oxidative-damage cellular death pathway. Cell Rep 3, 350-358, doi:10.1016/j.celrep.2012.12.021 (2013).

30

Therefore, to test our hypothesis, the PrestoBlue assay was conducted after the treatments of photoinactivation of catalase and subsequent addition of amp B against various clinical fungal isolates, including *C. auris* strains.

[00133] Interestingly, without the assistance of photoinactivation of catalase, some *C.*

- 5 *auris* strains were resilient to amp B (Figure 14). Nonetheless, photoinactivation of catalase achieved total eradication of *C. auris* strains regardless of the addition of amp B. Ten *C. auris* strains were tested and they demonstrated the same behavior. This means *C. auris* strains are exceptionally sensitive to blue light exposure.

Example 10. Photoinactivation of Catalase Inhibits the Formation of Hyphae of *C. albicans*,

- 10 and Assists Macrophage Cells to Phagocytose

[00134] Host immune cells play important roles against external evasive pathogens.

Catalase holds an essential role during the battle between *C. albicans* and neutrophils or macrophage cells (Pradhan, A. et al. Elevated catalase expression in a fungal pathogen is a double-edged sword of iron. Plos Pathog 13, e1006405 (2017). Thus, whether

- 15 photoinactivation of catalase could assist macrophage cells against *C. albicans* was examined. To visualize this effect, RAW 264.7 cells were infected with *C. albicans* and 410 nm-exposed *C. albicans* at a MOI of 10 and labeled with live/dead fluorescence stains.

[00135] As shown in Figure 15, untreated *C. albicans* stay as hyphae form and pierced through macrophage cells. Whereas 410 nm-exposed *C. albicans* remained as dead 'yeast'

- 20 form intracellularly.

[00136] In summary, photoinactivation of catalase in combination with low-concentration H_2O_2 presents an effective and novel approach to eliminate broad-spectrum fungus and fungal infections.

Example 11. Photoinactivation of Catalase in Combination With ROS Activating Agent

- 25 Silver Cation Synergistically Kills Microbes

[00137] Electromagnetic energy having a wavelength of ns-410 nm combined with 10 μ M of silver cation eliminated about 90% of MRSA one hour after treatment, whereas ns-410 nm alone or silver cation alone does not exert any significant antimicrobial effect (Figure 16).

- [00138]** Photoinactivation of catalase and low-concentration silver cation synergistically
30 eliminate *E. coli* BW25113 as well. The result is shown by scatter plots in Figure 17. CFU

- ml-1 (colony-forming unit) is designated as the amount of bacterial burden. 'Untreated' refers to the original E. coli BW25113 without any exogenous treatment. '0.5 μM Ag⁺ and 'CW-410' or 'ns-light' refers to E. coli BW25113 with 0.5 μM Ag⁺ and ns-410 alone, respectively. 0.5 μM Ag⁺ alone and CW-410 alone or ns-410 alone doesn't exert any significant killing effect on E. coli, however, ns-410 in combination with 0.5 μM Ag⁺ reduces around 99% of bacterial burden (Figure 17). The same phenomenon happens at other wavelengths as well. Noteworthy, CW-410 combined with 0.5 μM Ag⁺ didn't significantly reduce bacterial burden under the same conditions. Our results are consistent for both 30 and 60 minutes after treatments.
- 10 **[00139]** From the foregoing description, it will be apparent that variations and modifications may be made to the invention described herein to adopt it to various usages and conditions. Such embodiments are also within the scope of the following claims. The recitation of a listing of elements in any definition of a variable herein includes definitions of that variable as any single element or combination (or subcombination) of listed elements.
- 15 The recitation of an embodiment herein includes that embodiment as any single embodiment or in combination with any other embodiments or portions thereof.

REFERENCES

- 20 **[00140]** All patents, patent applications and publications mentioned in this specification are herein incorporated by reference to the same extent as if each independent patent and publication was specifically and individually indicated to be incorporated by reference.

CLAIMS

1. A method of treating a tissue of a subject infected with a catalase-positive microbe, said method comprising the steps of: applying light to the tissue of the subject infected with the catalase-positive microbe at a wavelength of about 400nm to about 430nm, wherein the catalase is inactivated, and contacting the tissue with a composition comprising a diluted peroxide solution, thereby treating the tissue of the subject infected with the catalase-positive microbe.
2. The method of claim 1, wherein the wavelength is about 410nm.
3. The method of claim 1, wherein the dose of the light is about 5 J/cm² to about 200 J/cm².
4. The method of claim 1, wherein the dose of the light is about 15 J/cm².
5. The method of claim 1, wherein the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a pulsed nanosecond laser.
6. The method of claim 1, wherein the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a continuous wave LED.
7. The method of claim 1, wherein the diluted peroxide solution is a hydrogen peroxide solution.
8. The method of claim 7, wherein the hydrogen peroxide solution is between about 0.03% and about 0.3% hydrogen peroxide.
9. The method of claim 8, wherein the hydrogen peroxide solution is 0.3% hydrogen peroxide.
10. The method of claim 1, further comprising administering a ROS generating agent to the infected tissue of the subject.
11. The method of claim 1, wherein the tissue is skin, scalp or nails.
12. A method of disinfecting an inanimate surface contaminated with a catalase-positive microbe, said method comprising the steps of: applying light to the inanimate surface at a wavelength of about 400nm to about 430nm, wherein the catalase is inactivated, and contacting the inanimate surface with a composition comprising a diluted peroxide solution, thereby disinfecting the inanimate surface.

13. The method of claim 12, wherein the wavelength is about 410nm.
14. The method of claim 13, wherein the dose of the light is about 5 J/cm^2 – to about 200 J/cm^2 .
15. The method of claim 12, wherein the dose of the light is about 15 J/cm^2 .
16. The method of claim 12, wherein the catalase-positive microbe is a fungus or bacteria and the light is provided by a pulsed nanosecond laser.
17. The method of claim 12, wherein the catalase-positive microbe is a fungus or bacteria and the light is provided by a continuous wave LED.
18. The method of claim 12, wherein the diluted peroxide solution is a hydrogen peroxide solution.
19. The method of claim 18, wherein the hydrogen peroxide solution is between about 0.03% and about 0.3% hydrogen peroxide.
20. The method of claim 19, wherein the wherein the hydrogen peroxide solution is 0.3% hydrogen peroxide.
21. The method of claim 12, further comprising administering a ROS generating agent to the infected tissue of the subject.
22. The method of claim 12, wherein the inanimate surface is a material comprising metal, plastic, fabric, rubber, stone, composite surfaces or wood.
23. A method of treating a tissue of a subject infected with a catalase-positive microbe, said method comprising the steps of: applying light from a pulsed nanosecond laser to the tissue of the subject infected with the catalase-positive microbe at a wavelength of about 400nm to about 460nm, wherein the catalase is inactivated, and contacting the tissue with a composition comprising a ROS generating agent, thereby treating the tissue of the subject infected with the catalase-positive microbe.
24. The method of claim 23, wherein the catalase-positive microbe is a fungus or bacteria.
25. The method of claim 23, wherein the dose of the light is about 5 J/cm^2 – to about 200 J/cm^2 .
26. The method of claim 25, wherein the dose of the light is about 15 J/cm^2 .

27. The method of claim 23, wherein the pulse duration of the light is about 5 nanoseconds.

28. A method of producing a synergistic antimicrobial effect in a tissue of a subject infected with a catalase-positive microbe, said method comprising the steps of: applying light to the tissue of the subject infected with the catalase-positive microbe at a wavelength of about 400nm to about 460nm, wherein the catalase is inactivated, and contacting the tissue with a composition comprising a diluted peroxide solution, thereby producing the synergistic antimicrobial effect in the tissue of the subject infected with the catalase-positive microbe.

29. The method of claim 28, wherein the wavelength is about 410nm.

30. The method of claim 28, wherein the dose of the light is about 5 J/cm² to about 200 J/cm².

31. The method of claim 28, wherein the dose of the light is about 15 J/cm².

32. The method of claim 28, wherein the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a pulsed nanosecond laser.

33. The method of claim 28, wherein the catalase-positive microbe is a fungal or bacterial microbe and the light is provided by a continuous wave LED.

34. The method of claim 28, wherein the diluted peroxide solution is a hydrogen peroxide solution.

35. The method of claim 34, wherein the hydrogen peroxide solution is between about 0.03% and about 0.3% hydrogen peroxide.

36. The method of claim 35, wherein the hydrogen peroxide solution is 0.3% hydrogen peroxide.

37. The method of claim 28, further comprising administering a ROS generating agent to the infected tissue of the subject.

38. The method of claim 28, wherein the tissue is skin, scalp or nails.

39. The method of claim 1, 12, 23 or claim 28, wherein the catalase-positive microbe is eradicated.