

1/8


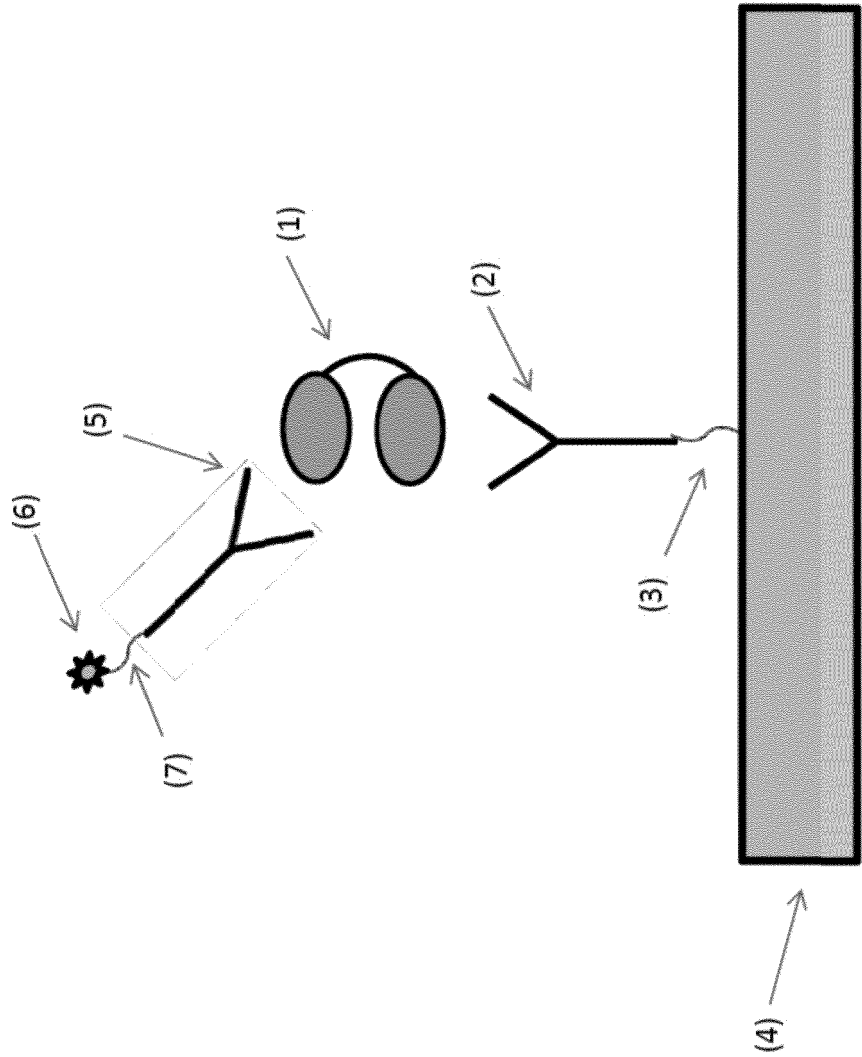


Figure 1

2/8

Figure 2

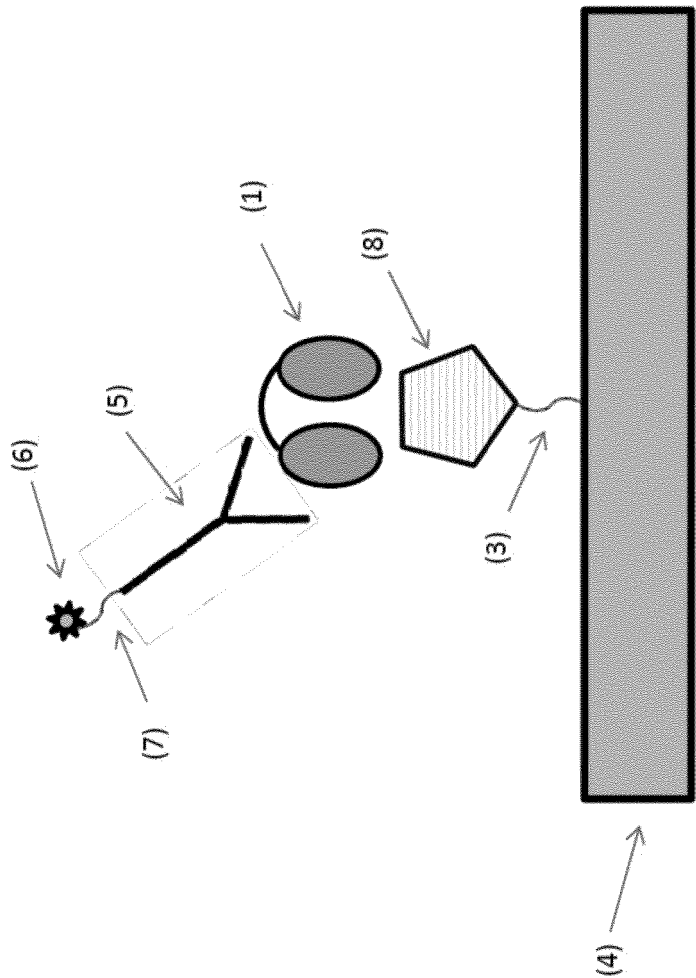

3/8

Figure 3


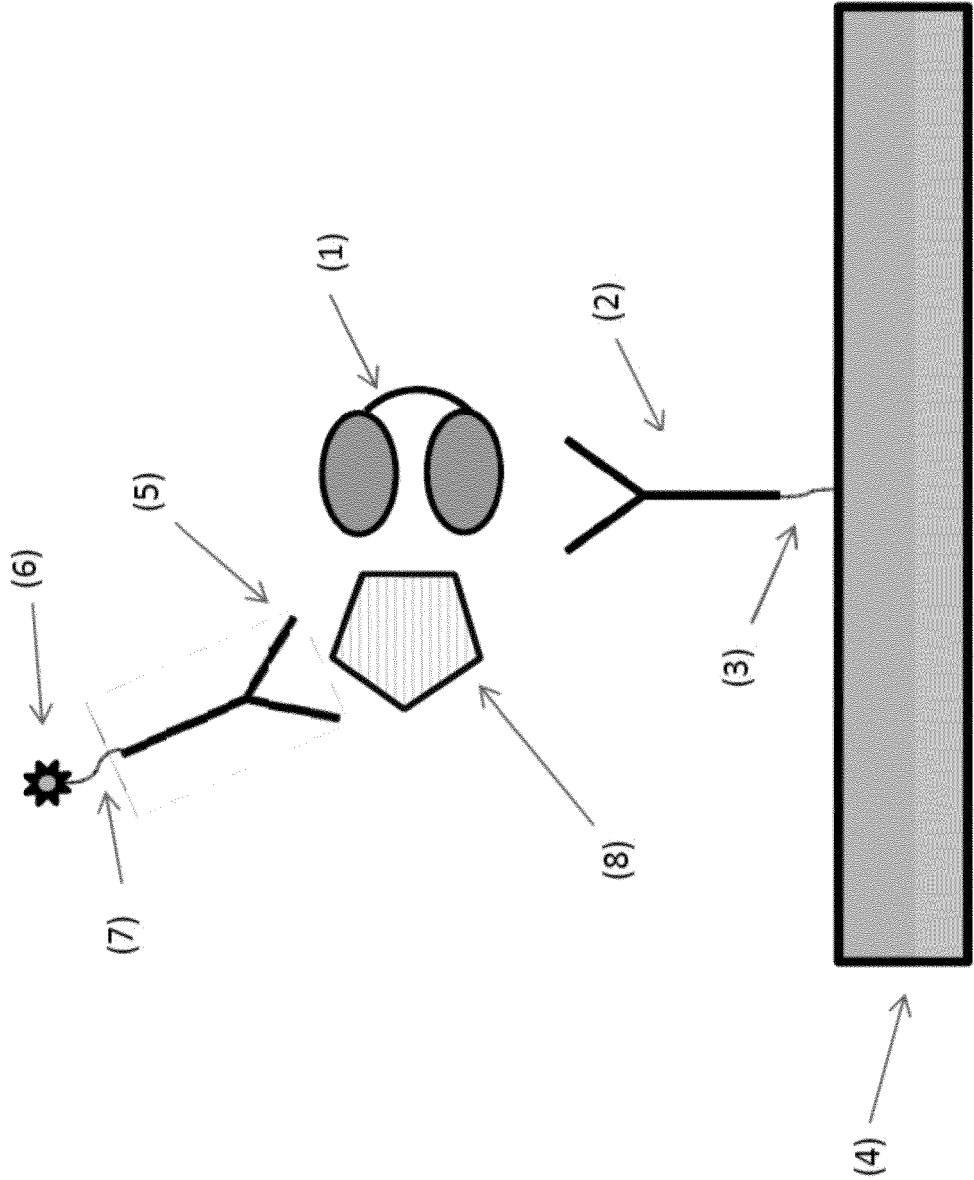


Figure 4

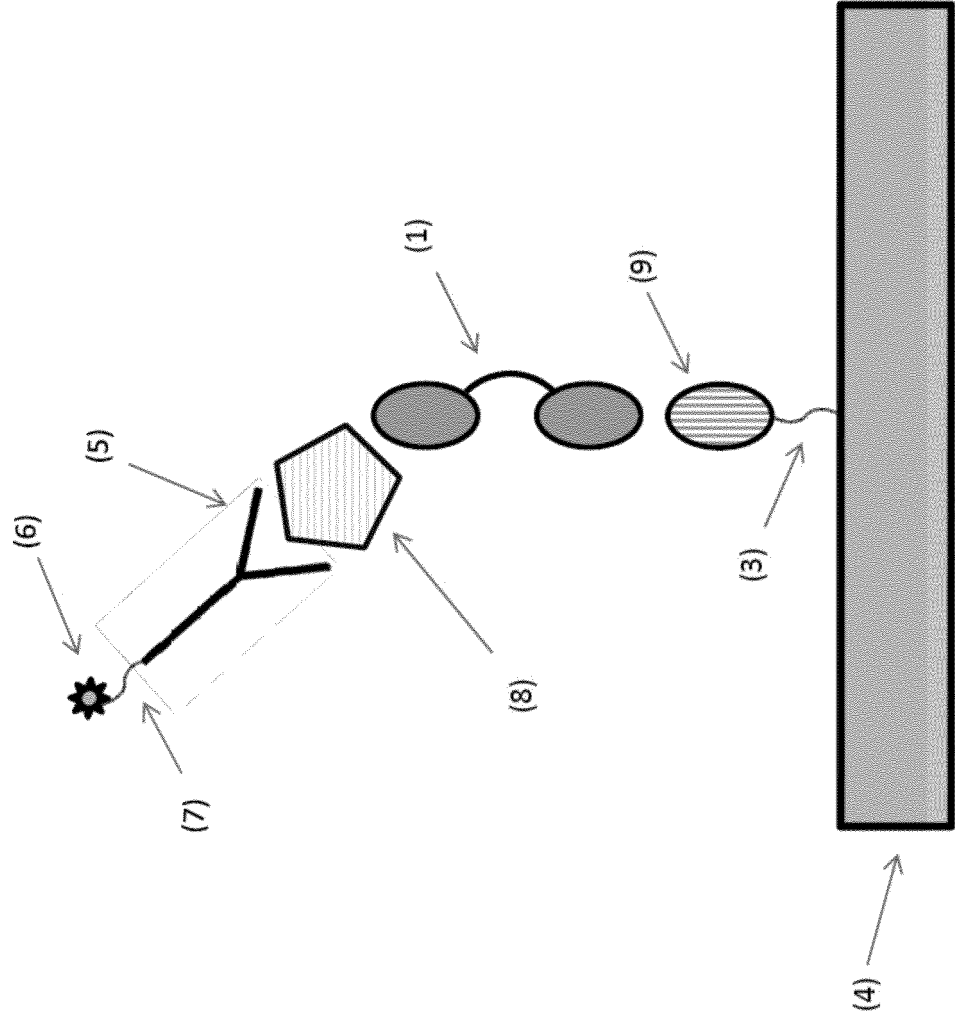

5/8

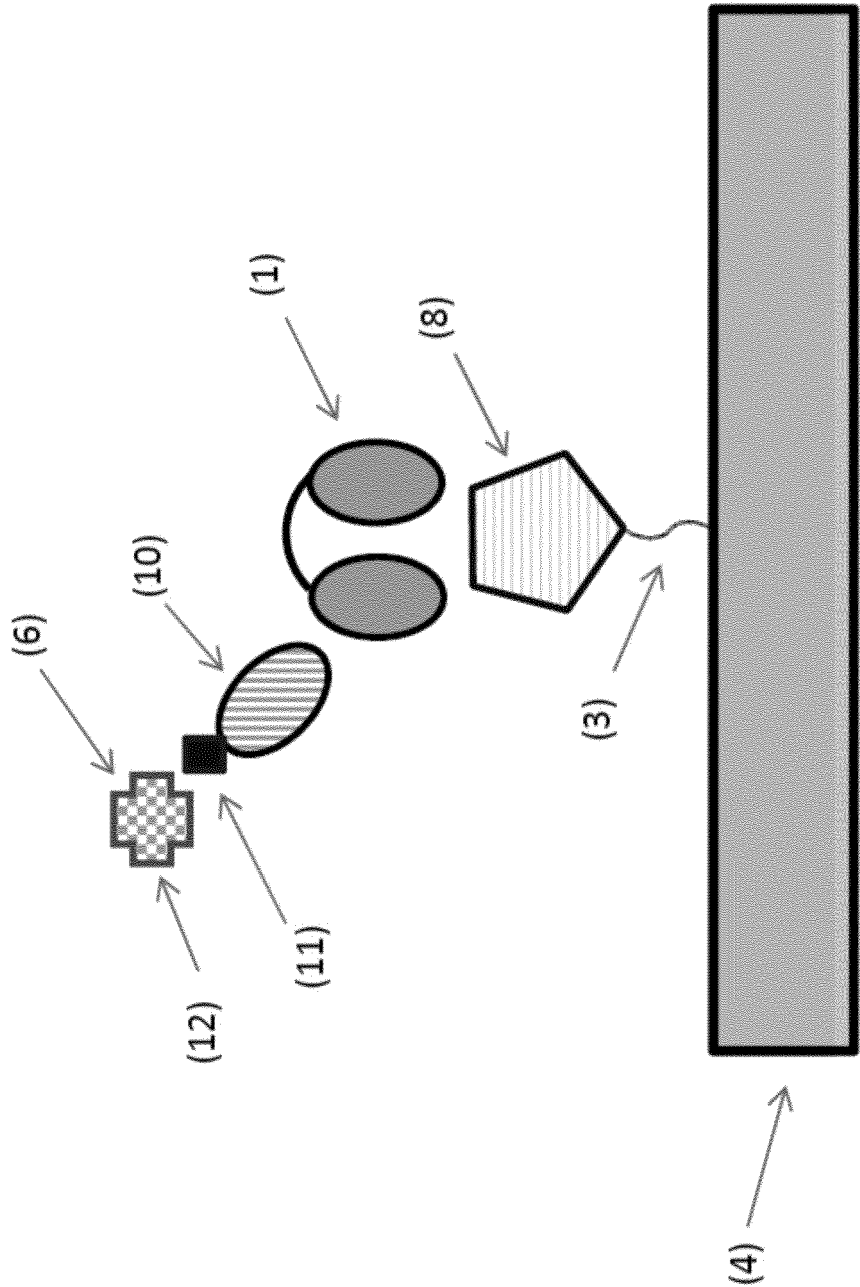


Figure 5

6/8

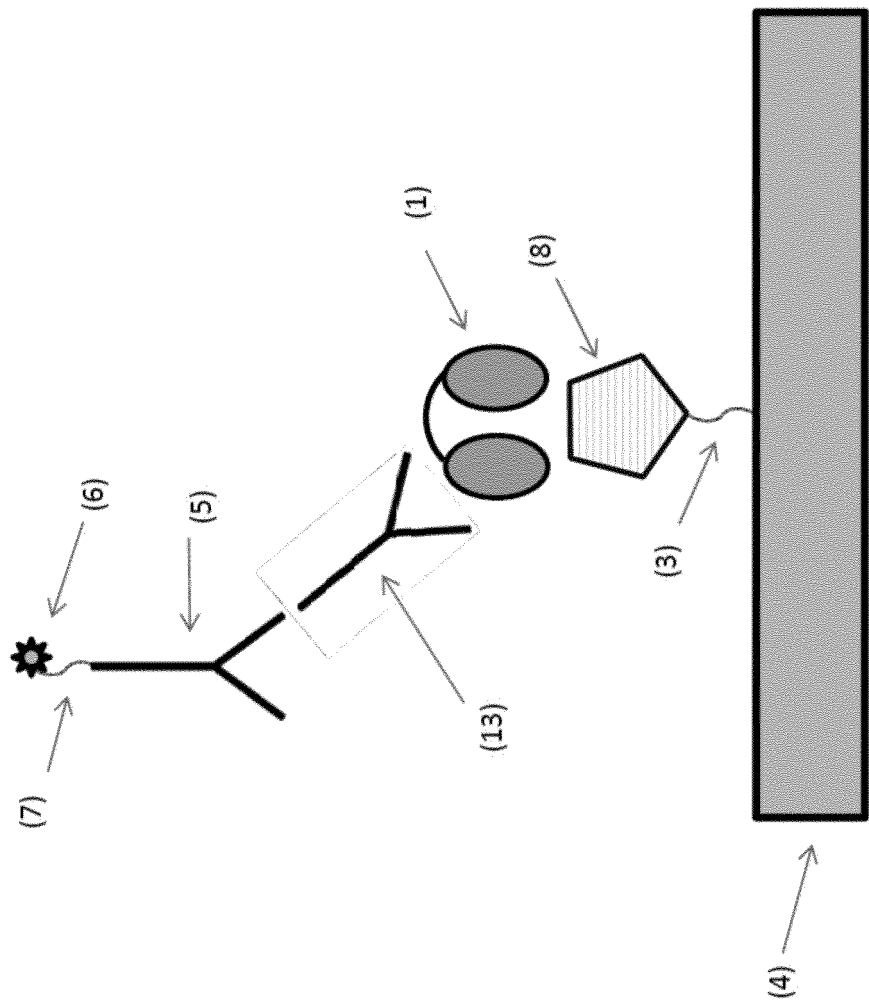


Figure 6

Figure 7

SEQ ID NO:	Sequence
1	VTVSS
2	HHHHHHDVQLVESGGGLVQPGGSLRLSCAASGRTFNYYAMGWF RQAPGKEREFVAAITRSGVRSVSAIYGDSVKDRFTISRDNKNT LYLQMNSLRPEDTAVYYCAASAIGSGALRRFEYDYSQGQTLTV SSGGGGSGGGGSGGGGSGGGGSGGGGSGGGGSGGGGSEVQLVE SGGGLVQPGGSLRLSCAASGRTFNYYAMGWFRQAPGKEREFV AITRSGVRSVSAIYGDSVKDRFTISRDNKNTLYLQMNSLRPED TAVYYCAASAIGSGALRRFEYDYSQGQTLTVSSGGGGSGGGG GGGGSGGGGSGGGGSGGGGSGGGGSEVQLVESGGGLVQPGNSL RLSCAASGFTFSSFGMSWVRQAPGKGLEWVSSISGSGSDTLYADS VKGRFTISRDNKNTLYLQMNSLRPEDTAVYYCTIGGSLSRSSQG TLTVSS
3	HHHHHHEVQLVESGGDLVQPGNSLRLSCAASGFTFSSFGMSWVR QAPGKGLEWVSSISGSGSDTLYADSVKGRFTISRDNKNTLYLQM NSLRPEDTAVYYCTIGGSLSRSSQGTLTVSS
4	EVQLVESGGGLVQPGNSLRLSCAASGFTFSSFGMSWVRQAPGKG LEWVSSISGSGSDTLYADSVKGRFTISRDNKNTLYLQMNSLRPE DTAVYYCTIGGSLSRSSQGTLTVSS
5	EVQLVESGGGLVQPGGSLRLSCAASGFTFSSYAMSWVRQAPGKG LEWVSGIKSSGDSTRYAGSVKGRFTISRDNKNTLYLQMNSLRPE DTAVYYCAKSRVSRTGLYTYDNRGQGTTLTVSSGGGGSGGGG GGGGSGGGGSEVQLVESGGGLVQPGGSLRLSCAASGRTFNYYA MGWFRQAPGKEREFVAAITRSGVRSVSAIYGDSVKDRFTISRDN AKNTLYLQMNSLRPEDTAVYYCAASAIGSGALRRFEYDYSQGQ TLTVSS

Figure 8

Capture	
Flow path	2
Flow rate (µl/min)	5
Contact time (s)	180
Concentration (µg/ml)	Suitably chosen from a range between 1 and 100 µg/ml
Binding and dissociation	
Flow path	1,2
Flow rate (µl/min)	45
Sample contact time (s)	120
Sample concentration (nM)	Ranging from 0 to 6 micromolar
Dissociation time (s)	1200
Regeneration1	
Flow path	1,2
Flow rate (µl/min)	10
Regeneration contact time (s)	180
Regeneration buffer	10mM Glycine-HCl pH1.7
Stabilization time (s)	60
If ...Then...Else	If after regeneration1 >15RU on Fc2
	Else exit cycle
Regeneration2	
Flow path	1,2
Flow rate (µl/min)	10
Regeneration contact time (s)	240
Regeneration buffer	10mM Glycine-HCl pH1.7