

FIELD OF INVENTION

[0001] The present invention relates to small molecules or metabolites that are found
5 to have significantly different abundances or intensities between clinically diagnosed
MULTIPLE SCLEROSIS or other neurological disorders, and normal patients. The present
invention also relates to methods for diagnosing MULTIPLE SCLEROSIS and other
neurological disorders, or individuals at risk of getting MULTIPLE SCLEROSIS or other
neurological disorders.

10 BACKGROUND OF THE INVENTION

[0002] MULTIPLE SCLEROSIS is the most common neurological disorder effecting
people under the age of 30, and is second only to epilepsy as the most common
disease of the central nervous system (CNS) [1]. It is generally accepted that MULTIPLE
SCLEROSIS is an autoimmune disorder that results in focal and discrete areas of
15 inflammation and demyelination throughout the white matter of the CNS.

[0003] The prevalence rate of MULTIPLE SCLEROSIS throughout North America ranges
from 1 per 500 to 1 per 1000, affecting an estimated 50,000 Canadians and 400,000
Americans; there are approximately 2 million people affected world-wide.
Epidemiological studies have revealed females are twice as likely to develop the
20 disease, the age of onset is relatively early (peak age of 30), and there is a greater
susceptibility in people of northern European descent [2]. Although differing theories
have implicated the involvement of various environmental factors [3-6], immune
dysfunction [3,4], and genetic anomalies [3,4] in the development of this disorder, the
etiology is still unknown. It is reasonable to assume that any factor that results in an
25 autoimmune reaction against myelin proteins results in MULTIPLE SCLEROSIS. The most
accepted theory involving its etiology takes into account several factors and suggests
genetically susceptible individuals are exposed to a foreign entity, such as a virus or a
toxin, and through some type of molecular mimicry, an autoimmune reaction against
myelin proteins is initiated. Approximately five to fifteen years later, the first clinical
30 symptoms become apparent/evident [7].

[0004] The pathological hallmark of MULTIPLE SCLEROSIS is discrete and focal areas of myelin loss, known as plaques or lesions. These plaques can consist of varying amounts of demyelination, gliosis, inflammation, edema and axonal degradation [8]. Although the exact locations of the plaques vary among patients, a general anatomical pattern is evident. Plaques within the human brain are located periventricular, within the temporal lobe, corpus callosum, optic nerves, brain stem, and/or cerebellum and tend to surround one or more blood vessels [7,9]. More than half of MULTIPLE SCLEROSIS patients have plaques within the cervical portion of the spinal cord [10]. The physiological consequence of the plaques is the slowing or blocked transmission of nerve impulses which manifests itself as sensory and/or motor impairment. In 2000, Lucchinetti et al [11] described four distinct patterns of MULTIPLE SCLEROSIS plaques in terms of their histological features. Two of these patterns suggest that demyelination results from the destruction of the myelin-producing cells within the CNS, oligodendrocytes, whereas the other two patterns indicate that myelin destruction results from T-cell or T-cell plus antibody targeting of the myelin sheath. The two patterns where oligodendrocytes are destroyed differ from one another by the selective destruction of specific myelin proteins in one pattern. The demyelinated lesions that contain T cells differ due to immunoglobulin-containing deposition and activated complement characteristic of one pattern. The discovery of the four patterns of MULTIPLE SCLEROSIS plaques was important since it indicates that the process of demyelination within this disorder can be achieved in several ways, and, hence, supports the notion that any process which triggers the formation of these plaques results in the clinical manifestation of MULTIPLE SCLEROSIS.

[0005] However, the pathological examination of MULTIPLE SCLEROSIS plaques is problematic in that it is derived primarily from post-mortem tissue, which represents only a snapshot of the disease at a given time. The majority of this tissue is acquired from individuals who had MULTIPLE SCLEROSIS for several years, and therefore represent tissue from the chronic stage of the disease. While post-mortem tissue may provide some information about the pathology of the disease, but it cannot elucidate how the disease progresses or where the lesions began. Magnetic resonance imaging (MRI) is commonly used to visualize MULTIPLE SCLEROSIS lesions *in vivo*. The use of

MRI to study MULTIPLE SCLEROSIS lesions is limited, however, because it cannot provide information about the pathological composition of the lesions.

[0006] The initial diagnosis of MULTIPLE SCLEROSIS is typically either relapsing-remitting (RR-MULTIPLE SCLEROSIS) or primary-progressive (PP-MULTIPLE SCLEROSIS). PP-MULTIPLE SCLEROSIS is the initial diagnosis in 10-15% of patients and is defined as a gradual worsening of symptoms throughout the course of the disease without any clinical remissions [4,12]. RR-MULTIPLE SCLEROSIS is the most common form as it is the initial diagnosis in 80% of patients, and is defined by clinical attacks (relapses) that last at least 24 hours followed by partial or complete recovery (remission). Within 20 years of initial diagnosis, 90% of RR-MULTIPLE SCLEROSIS patients will proceed to the secondary-progressive form of MULTIPLE SCLEROSIS (SP-MULTIPLE SCLEROSIS), where the symptoms worsen and remission periods eventually disappear. Some RR-MULTIPLE SCLEROSIS patients within a 15-year time period experience few relapses with no worsening of symptoms and long remission periods; these patients would have developed benign-MULTIPLE SCLEROSIS (BN-MULTIPLE SCLEROSIS). Currently, there is no evidence that indicates why a patient would initially manifest either PP-MULTIPLE SCLEROSIS or RR-MULTIPLE SCLEROSIS.

[0007] In 2001, the McDonald Criteria [13] was published to standardize the diagnosis of MULTIPLE SCLEROSIS. The fundamental feature of the criteria involves the objective evidence of lesions disseminated in both time and space. Clinical evidence alone can be adequate to secure a diagnosis if: 1) the individual has experienced two attacks/relapses and 2) there is clinical evidence of two or more lesions separated by time and space. If the individual does not reach this clinical criterion, additional paraclinical tests from MRI, cerebrospinal fluid (CSF) analysis and/or visual evoked potentials (VEP) are performed. MRI is the most sensitive and specific paraclinical test as it can provide objective evidence for dissemination of lesions in both time and space. CSF analysis can provide evidence of immune or inflammatory reactions of lesions and can aid in diagnosis when the clinical presentation and MRI criteria are not met, but it cannot provide information about dissemination of lesions or events in time or space. VEP in MULTIPLE SCLEROSIS are delayed, but exhibit a well-preserved waveform and can be used to provide evidence of a second lesion if the first lesion does not affect the visual pathway. The supplemental evidence provided by the

paraclinical tests might result in a diagnosis of either: a) having MULTIPLE SCLEROSIS, b) not having MULTIPLE SCLEROSIS, or c) having possible MULTIPLE SCLEROSIS. The majority of individuals diagnosed with having MULTIPLE SCLEROSIS exhibit the RR-MULTIPLE SCLEROSIS form, so the dissemination of lesions in time and space is often
5 evident. However, since there are no remission periods in PP-MULTIPLE SCLEROSIS, paraclinical tests are particularly important to secure a diagnosis. CSF analysis and either MRI or VEP must be obtained to provide objective evidence about space, whereas the use of MRI and continued progression of clinical symptoms for one year could provide evidence about dissemination over time.

10 [0008] Prior to the utilization of these paraclinical tests, it took an average of seven years before a physician could secure a diagnosis. Today, the use of these tests can secure a diagnosis of RR-MULTIPLE SCLEROSIS within months. The McDonald Criteria decreased the time required for diagnosis substantially, but for those individuals who are diagnosed with possible MULTIPLE SCLEROSIS, or will eventually receive a
15 diagnosis of PP-MULTIPLE SCLEROSIS, it has fallen short.

[0009] While the paraclinical tests may aid in the diagnosis of multiple sclerosis and provide information regarding the dissemination of lesions, no specific information regarding the pathological composition of the lesions is obtained. In addition, the interpretation of paraclinical test results is subjective and requires the expertise of
20 trained personnel. Furthermore, tools such as the pathological examination of multiple sclerosis plaques and the paraclinical test do not provide any information on susceptibility to the disease, but rather are used once symptoms become apparent.

SUMMARY OF THE INVENTION

25 [0010] The present invention relates to small molecules or metabolites that are found to have significantly different abundances or intensities between persons with MULTIPLE SCLEROSIS or other neurological disorders, and normal patients. The present invention also relates to small molecules or metabolites that have significantly different abundances or intensities between persons with neuropathology associated
30 with MULTIPLE SCLEROSIS and persons absent of such pathology such that these small

molecules or metabolites may be indicative of a pre-clinical pathological state. The present invention also relates to methods for diagnosing MULTIPLE SCLEROSIS and other neurological disorders.

5 [0011] The present invention provides novel methods for discovering, validating, and implementing a diagnostic method for one or more diseases or particular health-states. In particular, the present invention provides a method for the diagnosis and differential diagnosis of MULTIPLE SCLEROSIS in humans by measuring the levels of specific small molecules present in a sample and comparing them to "normal" reference levels.

10 [0012] The type of neurological disorder diagnosed by the above method may be MULTIPLE SCLEROSIS, or other type of demyelinating disorder. The sample obtained from the human may be a blood sample.

15 [0013] A method is provided for the diagnosis of subjects afflicted with MULTIPLE SCLEROSIS (relapsing-remitting or primary-progressive) and/or for the differential diagnosis of subjects transitioning from relapsing-remitting to secondary progressive MULTIPLE SCLEROSIS.

[0014] The methods of the present invention, including high throughput screening (HTS) assays, can be used for the following, wherein the specific "health-state" in this application may refer to, but is not limited to, MULTIPLE SCLEROSIS:

20 [0015] 1. identifying small-molecule metabolite biomarkers that can discriminate between multiple health-states using any biological sample taken from an individual;

[0016] 2. specifically diagnosing a health-state using metabolites identified in serum, plasma, whole blood, CSF, and/or other tissue biopsy as described in this application;

25 [0017] 3. selecting the minimal number of metabolite features required for optimal diagnostic assay performance statistics using supervised statistical methods such as those mentioned in this application;

[0018] 4. identifying structural characteristics of biomarker metabolites selected from non-targeted metabolomic analysis using LC-MS/MS, MSⁿ and NMR;

- [0019] 5. developing a high-throughput LC-MS/MS method for assaying selected metabolite levels in serum, plasma, whole blood, CSF, saliva, urine, hair, and/or other tissue biopsy; and
- [0020] 6. diagnosing a given health-state, or risk for development of a health-state by determining the levels of any combination of metabolite features disclosed from the Fourier Transform Mass Spectrometry (FTMS) analysis patient serum or other biological fluid or tissue, using any method including, but not limited to, mass spectrometry, NMR, UV detection, ELISA (enzyme-linked immunosorbant assay), chemical reaction, image analysis, or other.
- 10 [0021] The present invention provides for the longitudinal monitoring or screening of the general population for one or more health-states using any single or combination of features disclosed in the method, described above.
- [0022] The present invention also provides several hundred metabolite masses that have statistically significant differential abundances between clinically diagnosed RR-
15 MULTIPLE SCLEROSIS, clinically diagnosed PP-MULTIPLE SCLEROSIS, clinically diagnosed SP-MULTIPLE SCLEROSIS, and normal samples, also referred to herein as a reference sample. Of the metabolite masses identified, an optimal panel of between four and 45 metabolite masses can be used, or any number there between; for example, an optimal panel of 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20,
20 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36, 37, 38, 39, 40, 41, 42, 43, 44, or 45 metabolite masses can be used to differentiate between clinically diagnosed RR-MULTIPLE SCLEROSIS, clinically diagnosed PP-MULTIPLE SCLEROSIS, clinically diagnosed SP-MULTIPLE SCLEROSIS and normal states. In a specific, non-limiting example, an optimal panel of 36 metabolite masses can be used.
- 25 [0023] The present invention also provides a panel of about 257 metabolite masses that can be used as a diagnostic indicator of RR-MULTIPLE SCLEROSIS disease course in serum samples compared to normal samples (see Table 1); in a further example, the panel may contain about 240 metabolite masses. In a more specific example, an optimal panel of nine metabolite masses can be extracted and used as a diagnostic
30 indicator of RR-MULTIPLE SCLEROSIS disease course in serum samples compared to

normal samples; for example, the panel of nine metabolites can include those with masses (measured in Daltons) 452.3868, 496.4157, 524.4448, 540.4387, 578.4923, 580.5089, 594.4848, 596.5012, 597.5062 where a +/- 5 ppm difference would indicate the same metabolite.

5 [0024] Also, the invention provides a panel of about 100 metabolite masses that can be used as a diagnostic indicator of PP-MULTIPLE SCLEROSIS disease course in serum samples compared to normal samples (see Table 2); in a further example, the panel may contain about 60 metabolite masses. In a more specific example, an optimal panel of five metabolite masses can be extracted and used as a diagnostic indicator of
10 PP-MULTIPLE SCLEROSIS disease course in serum samples compared to normal samples; for example, the optimal panel of five metabolites can include those with masses (measured in Daltons) 202.0453, 216.04, 243.0719, 244.0559, 857.7516, where a +/- 5 ppm difference would indicate the same metabolite.

[0025] In addition, the invention provides a panel of about 226 metabolite masses that
15 can be used as a diagnostic indicator of SP-MULTIPLE SCLEROSIS disease course in serum samples compared to normal samples (see Table 3); in a further example, the panel may contain about 129 metabolite masses. In a more specific example, an optimal panel of eighteen metabolite masses can be extracted and used as a diagnostic indicator of SP-MULTIPLE SCLEROSIS disease course in serum samples compared to
20 normal samples; for example, the optimal panel of eighteen metabolites can include those with masses (measured in Daltons) 194.0803, 428.3653, 493.385, 541.3415, 565.3391, 576.4757, 578.4923, 590.4964, 594.4848, 495.4883, 596.5012, 596.5053, 597.5062, 597.5068, 805.5609, 806.5643, 827.5446, 886.5582, where a +/- 5 ppm difference would indicate the same metabolite.

25 [0026] Furthermore, the invention provides a panel of about 142 metabolite masses that can be used as a diagnostic indicator of RR-MULTIPLE SCLEROSIS disease course in serum samples compared to SP - MULTIPLE SCLEROSIS samples (see Table 4); in a further example, the panel may contain about 135 metabolite masses. In a more specific example, an optimal panel of six metabolite masses that can be extracted and
30 used as an indicator of RR-MULTIPLE SCLEROSIS disease course in serum samples compared to SP- MULTIPLE SCLEROSIS samples, also referred to herein as a reference

sample; for example, the optimal panel of six metabolites can include those with masses (measured in Daltons) 540.4387, 576.4757, 594.4848, 595.4883, 596.5012, 597.5062, where a +/- 5 ppm difference would indicate the same metabolite.

[0027] The present invention further provides a panel of about 148 metabolite masses that can be used as a diagnostic indicator of the transition from RR-MULTIPLE SCLEROSIS patients to SP-MULTIPLE SCLEROSIS compared to RR-MULTIPLE SCLEROSIS, also referred to herein as a reference sample (see Table 5); in a more specific example, an optimal panel of 5 metabolites masses that can be extracted and used as an indicator of early neuropathology changes within the transition from RR-MULTIPLE SCLEROSIS patients to SP-MULTIPLE SCLEROSIS compared to RR-MULTIPLE SCLEROSIS; for example, the optimal panel of five metabolites can include those with masses (measured in Daltons) 576.4757, 578.4923, 594.4848, 596.5012, 597.5062, where a +/- 5 ppm difference would indicate the same metabolite.

[0028] Moreover, the invention provides a panel of about 309 metabolite masses that can be used as a diagnostic indicator of the transition from RR-MULTIPLE SCLEROSIS to SP-MULTIPLE SCLEROSIS compared to SP-MULTIPLE SCLEROSIS (see Table 6), also referred to herein as a reference sample; in a further example, the panel may contain about 42 metabolite masses. In a more specific example, an optimal panel of eight metabolite masses that can be extracted and used as an indicator of early neuropathology changes within the transition from RR-MULTIPLE SCLEROSIS to SP-MULTIPLE SCLEROSIS compared to SP-MULTIPLE SCLEROSIS; for example, the optimal panel of eight metabolites can include those with masses (measured in Daltons) 617.0921, 746.5118, 760.5231, 770.5108, 772.5265, 784.5238, 786.5408, and 787.5452, where a +/- 5 ppm difference would indicate the same metabolite.

[0029] The present invention further provides a method for diagnosing RR-MULTIPLE SCLEROSIS, PP-MULTIPLE SCLEROSIS, and SP-MULTIPLE SCLEROSIS, comprising the steps of: introducing one or more samples from one or more patients with clinically diagnosed RR-MULTIPLE SCLEROSIS, clinically diagnosed PP-MULTIPLE SCLEROSIS or clinically diagnosed SP-MULTIPLE SCLEROSIS, introducing said sample containing a plurality of metabolites into a high resolution mass spectrometer, for example, a Fourier Transform Ion Cyclotron Resonance Mass Spectrometer (FTICR-MS);

obtaining, identifying and quantifying data for the metabolites; creating a database of said identifying and quantifying data; comparing, identifying and quantifying data from the sample with corresponding data from a sample from normal subject (one who does not have MULTIPLE SCLEROSIS); identifying one or more metabolites that differ; and selecting the minimal number of metabolite markers needed for optimal diagnosis.

[0030] In a further embodiment of the present invention there is provided a method for identifying specific biomarkers for RR-MULTIPLE SCLEROSIS, PP-MULTIPLE SCLEROSIS, and SP-MULTIPLE SCLEROSIS, comprising the steps of: introducing one or more samples from one or more patients with clinically diagnosed RR-MULTIPLE SCLEROSIS, clinically diagnosed PP-MULTIPLE SCLEROSIS, or clinically diagnosed SP-MULTIPLE SCLEROSIS, said sample containing a plurality of metabolites into an FT/IC-MS; obtaining, identifying, and quantifying data for the metabolites; creating a database of said identifying and quantifying data; comparing the identifying and quantifying data from the sample with corresponding data from a sample from a normal subject (one who does not have MULTIPLE SCLEROSIS) identifying one or more metabolites that differ; and selecting the minimal number of metabolite markers needed for optimal diagnosis. The metabolite markers needed for optimal diagnosis of RR-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 452.3868, 496.4157, 524.4448, 540.4387, 578.4923, 580.5089, 594.4848, 596.5012, 597.5062, where a +/- 5 ppm difference would indicate the same metabolite. The metabolite markers needed for optimal diagnosis of PP-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 202.0453, 216.04, 243.0719, 244.0559, 857.7516, where a +/- 5 ppm difference would indicate the same metabolite. The metabolite markers needed for optimal diagnosis of SP-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 194.0803, 428.3653, 493.385, 541.3415, 565.3391, 576.4757, 578.4923, 590.4964, 594.4848, 495.4883, 596.5012, 596.5053, 597.5062, 597.5068, 805.5609, 806.5643, 827.5446, 886.5582, where a +/- 5 ppm difference would indicate the same metabolite. The metabolite markers needed for optimal differentiation of RR-MULTIPLE SCLEROSIS patients from

SP-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 540.4387, 576.4757, 594.4848, 595.4883, 596.5012, 597.5062, where a ± 5 ppm difference would indicate the same metabolite. The metabolite markers needed for optimal

5 differentiation of RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS (RR-SP) as compared to SP-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 617.0921, 746.5118, 760.5231, 770.5108, 772.5265, 784.5238, 786.5408, and 787.5452, where a ± 5 ppm difference would indicate the same metabolite. The

10 metabolite markers needed for optimal differentiation of RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS (RR-SP) as compared to RR-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 576.4757, 578.4923, 594.4848, 596.5012, 597.5062, where a ± 5 ppm difference would indicate the same

15 metabolite.

[0031] In a further embodiment of the present invention there is provided a method for diagnosing a patient for RR-MULTIPLE SCLEROSIS, PP-MULTIPLE SCLEROSIS, and SP-MULTIPLE SCLEROSIS, comprising the steps of: screening a sample from said

20 patient for quantification of one or more metabolic markers and comparing the amounts of metabolite markers to corresponding data from a sample from a normal subject (one who does not have MULTIPLE SCLEROSIS). The metabolite markers for diagnosis of RR-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 452.3868, 496.4157, 524.4448, 540.4387, 578.4923, 580.5089, 594.4848, 596.5012, 597.5062,

25 where a ± 5 ppm difference would indicate the same metabolite. The metabolite markers for diagnosis of PP-MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 202.0453, 216.04, 243.0719, 244.0559, 857.7516, where a ± 5 ppm difference would indicate the same metabolite. The metabolite markers for diagnosis of SP-MULTIPLE

30 SCLEROSIS from healthy controls in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 194.0803, 428.3653, 493.385, 541.3415, 565.3391, 576.4757, 578.4923, 590.4964, 594.4848,

495.4883, 596.5012, 596.5053, 597.5062, 597.5068, 805.5609, 806.5643, 827.5446, 886.5582,, where a +/- 5 ppm difference would indicate the same metabolite. The metabolite markers for diagnosis of RR-MULTIPLE SCLEROSIS from SP- MULTIPLE SCLEROSIS in a serum sample may be selected from the group consisting of

5 metabolites with accurate masses (measured in Daltons) 540.4387, 576.4757, 594.4848, 595.4883, 596.5012, 597.5062, where a +/- 5 ppm difference would indicate the same metabolite. The metabolite markers needed for optimal differentiation of RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS (RR-SP) as compared to SP-MULTIPLE SCLEROSIS in a serum sample may be selected

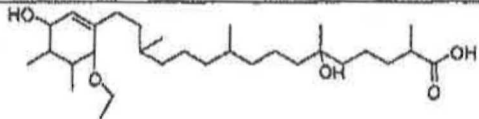
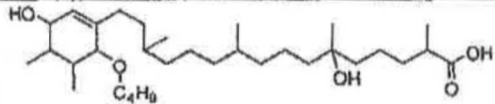
10 from the group consisting of metabolites with accurate masses (measured in Daltons) 617.0921, 746.5118, 760.5231, 770.5108, 772.5265, 784.5238, 786.5408, and 787.5452, where a +/- 5 ppm difference would indicate the same metabolite. The metabolite markers needed for optimal differentiation of RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS (RR-SP) as compared to RR-MULTIPLE

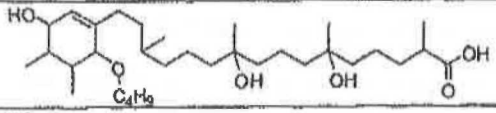
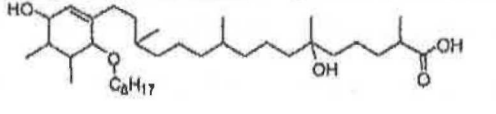
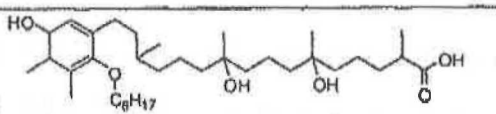
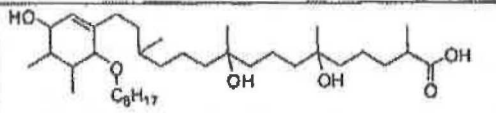
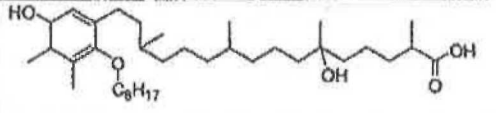
15 SCLEROSIS in a serum sample may be selected from the group consisting of metabolites with accurate masses (measured in Daltons) 576.4757, 578.4923, 594.4848, 596.5012, 597.5062, where a +/- 5 ppm difference would indicate the same metabolite.

[0032] The molecular formulae and proposed structure for some of the MULTIPLE

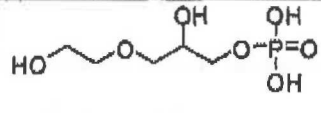
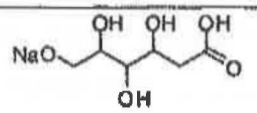
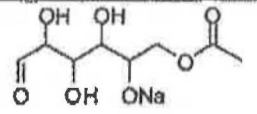
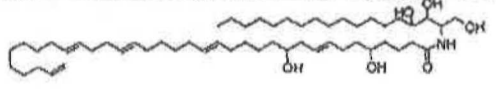
20 SCLEROSIS biomarkers referred to above were determined in one embodiment of the present invention. These are summarized below. According to results the biomarkers are thoughts to be derivatives of sugars, phospholipids and tocopherols.

[0033] RR-MULTIPLE SCLEROSIS as compared to a Normal patient

Mass	Formula	Structure
496.4157	$C_{30}H_{56}O_5$	
524.4448	$C_{32}H_{60}O_5$	

540.4387	$C_{32}H_{60}O_6$	
580.5089	$C_{36}H_{68}O_5$	
594.4848	$C_{36}H_{66}O_6$	
596.5012	$C_{36}H_{68}O_6$	
578.4923	$C_{36}H_{66}O_5$	

[0034] PP-MULTIPLE SCLEROSIS as compared to a Normal patient

Mass	Formulae	Structure
216.04	$C_5H_{13}O_7P$	
202.0453	$C_6H_{11}O_6Na$	
244.0559	$C_8H_{13}O_7Na$	
857.7516	$C_{54}H_{99}NO_6$	

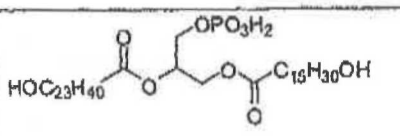
[0035] SP-MULTIPLE SCLEROSIS as compared to a Normal patient

Mass	Formulae	Structure
541.3415	$C_{25}H_{52}NO_9P$	
565.3391	$C_{27}H_{52}NO_9P$	
428.3653	$C_{29}H_{48}O_2$	
805.5609	$C_{48}H_{80}NO_8P$	
194.0803	$C_7H_{14}O_6$	
578.4923	$C_{36}H_{66}O_5$	

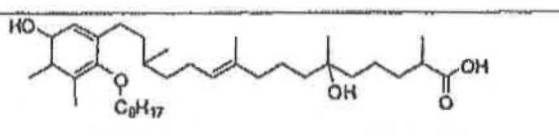
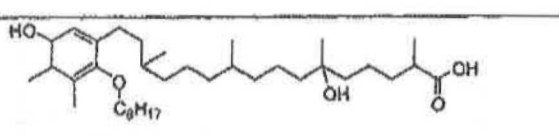
[0036] RR-MULTIPLE SCLEROSIS as compared to a SP -MULTIPLE SCLEROSIS patient

Mass	Formulae	Structure
540.4387	$C_{32}H_{60}O_6$	
576.4757	$C_{36}H_{64}O_5$	

[0037] RR-MULTIPLE SCLEROSIS transitioning to SP -MULTIPLE SCLEROSIS as compared to a SP -MULTIPLE SCLEROSIS patient

Mass	Formulae	Structure
786.5408	$C_{43}H_{79}O_{10}P$	

5 [0038] RR-MULTIPLE SCLEROSIS transitioning to SP -MULTIPLE SCLEROSIS as compared to a RR -MULTIPLE SCLEROSIS patient

Mass	Formulae	Structure
576.4757	$C_{36}H_{64}O_5$	
578.4923	$C_{36}H_{66}O_5$	

[0039] The identification of MULTIPLE SCLEROSIS-specific biomarkers in human serum is extremely useful since it is minimally invasive, and can be used to detect the presence of MULTIPLE SCLEROSIS pathology prior to the manifestation of clinical

10 symptoms. A serum test is minimally invasive and would be accepted by the general population. The metabolite masses presently identified were found to have statistically significantly differential abundances between RR-MULTIPLE SCLEROSIS, PP-MULTIPLE SCLEROSIS, SP-MULTIPLE SCLEROSIS and normal serum, of which an optimal panels can be extracted and used as a diagnostic indicator of disease presence. A diagnostic

15 assay based on small molecules or metabolites in serum can be developed into a relatively simple and cost-effective assay that is capable of detecting specific metabolites. Translation of the method into a clinical assay, compatible with current clinical chemistry laboratory hardware, is commercially acceptable and effective, and

could result in a rapid deployment worldwide. Furthermore, the requirement for highly trained personnel to perform and interpret the test would be eliminated.

[0040] Since the present invention relates to panels of molecules that are increased in individuals with RR-MULTIPLE SCLEROSIS and PP-MULTIPLE SCLEROSIS as compared to
5 healthy individuals, there the test can be used as an indicator of susceptibility to the specific type of MULTIPLE SCLEROSIS or, alternatively, an indicator of very early disease onset. The possibility of a highly accurate MULTIPLE SCLEROSIS predisposition assay in serum would be the first of its kind.

[0041] The impact of the present invention on the diagnosis of MULTIPLE SCLEROSIS
10 would be tremendous, as literally everyone could be screened longitudinally throughout their lifetime to assess risk. Given that the performance characteristics of the test of the present invention are representative for the general population, this test alone may be superior to any other currently available screening method, as it may have the potential to detect disease progression prior to the emergence of clinical
15 symptoms.

[0042] This summary of the invention does not necessarily describe all features of the invention.

BRIEF DESCRIPTION OF THE DRAWINGS

20 [0043] These and other features of the invention will become more apparent from the following description in which reference is made to the appended drawings, wherein:

[0044] FIGURE 1A shows a Prediction Analysis of Microarray (PAM) training error plot and FIGURE 1B shows a cross validated misclassification error plot, in accordance with an embodiment of the present invention.

25 [0045] FIGURE 2 shows cross-validated diagnostic probabilities for clinically diagnosed RR-MULTIPLE SCLEROSIS patients and controls, in accordance with an embodiment of the present invention.

[0046] FIGURE 3 shows a receiver-operator characteristic (ROC) curve based on cross-validated probabilities, in accordance with a further embodiment of the present invention.

5 [0047] FIGURE 4 shows diagnostic predictions for blinded test set, in accordance with a further embodiment of the present invention.

[0048] FIGURE 5 shows a ROC curve based on predicted test set of clinically diagnosed RR-MULTIPLE SCLEROSIS patients and controls, in accordance with a further embodiment of the present invention.

10 [0049] FIGURE 6 shows a ROC curve based on clinically diagnosed PP-MULTIPLE SCLEROSIS and controls, in accordance with a further embodiment of the present invention.

[0050] FIGURE 7 shows a ROC curve based on clinically diagnosed SP-MULTIPLE SCLEROSIS and controls, in accordance with a further embodiment of the present invention.

15 [0051] FIGURE 8 shows a ROC curve based on clinically diagnosed RR-MULTIPLE SCLEROSIS and SP-MULTIPLE SCLEROSIS in accordance with a further embodiment of the present invention.

20 [0052] FIGURE 9 shows a ROC curve based on clinically diagnosed RR-MULTIPLE SCLEROSIS patients and RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS, in accordance with a further embodiment of the present invention.

[0053] FIGURE 10 shows a ROC curve based on clinically diagnosed SP-MULTIPLE SCLEROSIS and RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS, in accordance with a further embodiment of the present invention.

25 [0054] FIGURE 11 shows a mean signal-to-noise +/-SEM of the RR-MULTIPLE SCLEROSIS 9 serum biomarker panel relative to controls, in accordance with a further embodiment of the present invention.

- [0055] FIGURE 12 shows a mean signal-to-noise $\pm$ SEM of the PP-MULTIPLE SCLEROSIS 5 serum biomarker panel relative to controls, in accordance with a further embodiment of the present invention.
- [0056] FIGURE 13 shows a mean signal-to-noise $\pm$ SEM of the SP-MULTIPLE SCLEROSIS 18 serum biomarker panel relative to controls, in accordance with a further embodiment of the present invention.
- [0057] FIGURE 14 shows a mean signal-to-noise $\pm$ SEM of the RR-MULTIPLE SCLEROSIS 6 serum biomarker panel relative to SP-MULTIPLE SCLEROSIS, in accordance with a further embodiment of the present invention.
- [0058] FIGURE 15 shows a mean signal-to-noise $\pm$ SEM of the RR-MULTIPLE SCLEROSIS 5 serum biomarker panel relative to RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS, in accordance with a further embodiment of the present invention.
- [0059] FIGURE 16 shows a mean signal-to-noise $\pm$ SEM of the RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS 8 serum biomarker panel relative to SP-MULTIPLE SCLEROSIS, in accordance with a further embodiment of the present invention.

DETAILED DESCRIPTION

- [0060] The present invention relates to small molecules or metabolites that are found to have significantly different abundances or intensities between clinically diagnosed MULTIPLE SCLEROSIS or other neurological disorders, and normal patients. The present invention also relates to methods for diagnosing MULTIPLE SCLEROSIS and other neurological disorders.
- [0061] The present invention provides novel methods for discovering, validating, and implementing a diagnosis method for one or more diseases or particular health-states. In particular, the present invention provides a method for the diagnosis and differential diagnosis of MULTIPLE SCLEROSIS in humans by measuring the levels of

specific small molecules present in a sample and comparing them to "normal" reference levels. A reference sample can be a normal sample or a sample from a patient with other forms of MULTIPLE SCLEROSIS. The sample may be any biological sample, including, but not exclusive to blood, urine, saliva, hair, cerebrospinal fluid (CSF), biopsy or autopsy samples. The methods measure the intensities of specific small molecules, also referred to as metabolites, in the sample from patients with MULTIPLE SCLEROSIS and compare these intensities to the intensities observed in a population of healthy (non-MULTIPLE SCLEROSIS) individuals.

[0062] The small molecules measured in a sample may also be referred to herein as "markers", "biomarkers", or "metabolites". The metabolites may be characterized in any manner known in the art, for example but not limited to, by mass (also referred to as "metabolite mass" or "accurate mass"), molecular formula, polarity, acid/base properties, NMR spectra, MS/MS or MSⁿ spectra, molecular structure, or any combination thereof. The term "metabolite feature" refers to a metabolite, a fragment thereof, an analogue thereof, or a chemical equivalent thereof.

[0063] The diagnosis or the exclusion of any type(s) of neurological disorders is contemplated by the present invention, using all or a subset of the metabolites disclosed herein. The types of neurological disorders include, but are not limited to: Alzheimer's disease (AD), dementia with Lewy bodies (DLB), frontotemporal lobe dementia (FTD), vascular induced dementia (e.g. multi-infarct dementia), anoxic event induced dementia (e.g. cardiac arrest), trauma to the brain induced dementia (e.g. dementia pugilistica [boxer's dementia]), dementia resulting from exposure to an infectious (e.g. Creutzfeldt-Jakob Disease) or toxic agent (e.g. alcohol-induced dementia), Acute Disseminated Encephalomyelitis, Guillain-Barré Syndrome, Adrenoleukodystrophy, Adrenomyeloneuropathy, Leber's Hereditary Optic Neuropathy, HTLV-associated Myelopathy, Krabbe's Disease, phenylketonuria, Canavan Disease, Pelizaeus-Merzbacher Disease, Alexander's Disease, Neuromyelitis Optica, Central Pontine Myelinolysis, Metachromatic Leukodystrophy, Schilder's Disease, Autism, Multiple Sclerosis, Parkinson's Disease, Bipolar Disorder, Ischemia, Huntington's Chorea, Major Depressive Disorder, Closed Head Injury, Hydrocephalus, Amnesia, Anxiety Disorder, Traumatic Brain Injury, Obsessive Compulsive Disorder, Schizophrenia, Mental Retardation, Epilepsy and/or any other

condition that is associated with an immune response, demyelination, myelitis or encephalomyelitis.

[0064] The present invention provides a method of diagnosing MULTIPLE SCLEROSIS and its subtypes by measuring the levels of specific small molecules present in a sample obtained from a human and comparing them to "normal" reference levels.

[0065] In order to determine whether there are biochemical markers of a given health-state in particular population, a group of patients representative of the health state (i.e. a particular disease) and a group of "normal" counterparts are required. Biological samples taken from the patients in a particular health-state category are then compared to the same samples taken from the normal population as well as to patients in similar health-state categories to identify biochemical differences between the two groups, by analyzing the biochemicals present in the samples using FTMS and/or LC-MS. The biological samples could originate from anywhere within the body, including, but not limited to, blood (serum/plasma), cerebrospinal fluid (CSF), urine, stool, saliva, or biopsy of any solid tissue including tumor, adjacent normal, smooth and skeletal muscle, adipose tissue, liver, skin, hair, brain, kidney, pancreas, lung, colon, stomach, or other. Of particular interest are samples that are serum. While the term "serum" is used herein, those skilled in the art will recognize that plasma, whole blood, or a sub-fraction of whole blood may be used.

[0066] The method of the present invention, based on small molecules or metabolites in a sample, makes an ideal screening test as the development of assays capable of detecting specific metabolites is relatively simple and cost effective. The test is minimally invasive and is indicative of MULTIPLE SCLEROSIS pathology, and may be useful to differentiate MULTIPLE SCLEROSIS subtypes from each other. Translation of the method into a clinical assay compatible with current clinical chemistry laboratory hardware is commercially acceptable and effective. Furthermore, the method of the present invention does not require highly trained personnel to perform and/or interpret the test.

[0067] The present invention also provides several hundred metabolite masses that were found to have statistically significantly differential abundances between

clinically diagnosed RR-MULTIPLE SCLEROSIS, clinically diagnosed PP-MULTIPLE SCLEROSIS, clinically diagnosed SP-MULTIPLE SCLEROSIS and normal serum.

[0068] *Non-Targeted Metabolomic Strategies.* Multiple non-targeted metabolomics strategies have been described in the scientific literature including NMR [14], GC-MS
5 [15-17], LC-MS, and FTMS strategies [14,18-20]. The metabolic profiling strategy employed for the discovery of differentially expressed metabolites in this application was the non-targeted FTMS strategy developed by Phenomenome Discoveries [17,20-23; see also US Published Application No. 2004-0029120 A1, Canadian Application No. 2,298,181, and WO 01/57518]. Non-targeted analysis involves the measurement
10 of as many molecules in a sample as possible, without any prior knowledge or selection of components prior to the analysis. Therefore, the potential for non-targeted analysis to discover novel metabolite biomarkers is high versus targeted methods, which detect a predefined list of molecules. The present invention uses a non-targeted method to identify metabolite components in serum samples that differ between:

15 [0069] 1. Clinically diagnosed RR-MULTIPLE SCLEROSIS patients and healthy controls;

[0070] 2. Clinically diagnosed PP-MULTIPLE SCLEROSIS patients and healthy controls;

[0071] 3. Clinically diagnosed SP-MULTIPLE SCLEROSIS patients and healthy controls;

[0072] 4. Clinically diagnosed RR-MULTIPLE SCLEROSIS patients and clinically diagnosed SP-MULTIPLE SCLEROSIS patients;

20 [0073] 5. Clinically diagnosed RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS patients and clinically diagnosed SP-MULTIPLE SCLEROSIS patients; and

[0074] 5. Clinically diagnosed RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS patients and clinically diagnosed RR-MULTIPLE SCLEROSIS patients.

[0075] *Sample Processing.* When a blood sample is drawn from a patient there are
25 several ways in which the sample can be processed. The range of processing can be as little as none (i.e. frozen whole blood) or as complex as the isolation of a particular cell type. The most common and routine procedures involve the preparation of either serum or plasma from whole blood. All blood sample processing methods, including

spotting of blood samples onto solid-phase supports, such as filter paper or other immobile materials, are also contemplated by the present invention.

[0076] *Sample Extraction.* The processed blood sample described above is then further processed to make it compatible with the methodical analysis technique to be employed in the detection and measurement of the biochemicals contained within the processed serum sample. The types of processing can range from as little as no further processing to as complex as differential extraction and chemical derivatization. Extraction methods may include sonication, soxhlet extraction, microwave assisted extraction (MAE), supercritical fluid extraction (SFE), accelerated solvent extraction (ASE), pressurized liquid extraction (PLE), pressurized hot water extraction (PHWE), and/or surfactant assisted extraction (PHWE) in common solvents such as methanol, ethanol, mixtures of alcohols and water, or organic solvents such as ethyl acetate or hexane. The preferred method of extracting metabolites for FTMS non-targeted analysis is to perform a liquid/liquid extraction whereby non-polar metabolites dissolve in an organic solvent and polar metabolites dissolve in an aqueous solvent.

[0077] *Mass spectrometry analysis of extracts.* Extracts of biological samples are amenable to analysis on essentially any mass spectrometry platform, either by direct injection or following chromatographic separation. Typical mass spectrometers are comprised of a source, which ionizes molecules within the sample, and a detector for detecting the ionized molecules or fragments of molecules. Examples of common sources include electron impact, electrospray ionization (ESI), atmospheric pressure chemical ionization, atmospheric pressure photo ionization (APPI), matrix assisted laser desorption ionization (MALDI), surface enhanced laser desorption ionization (SELDI), and derivations thereof. Common mass separation and detection systems can include quadrupole, quadrupole ion trap, linear ion trap, time-of-flight (TOF), magnetic sector, ion cyclotron (FTMS), Orbitrap, and derivations and combinations thereof. The advantage of FTMS over other MS-based platforms is its high resolving capability that allows for the separation of metabolites differing by only hundredths of a Dalton, many of which would be missed by lower resolution instruments.

[0078] *Training classifier.* Cross-validated training classifier was created using the Prediction Analysis of Microarrays (PAM) (<http://www-stat.stanford.edu/~tibs/PAM/>)

algorithm [24]. The method involves training a classifier algorithm using samples with known diagnosis that can then be applied to blinded diagnosed samples (i.e. a test set). Several supervised methods exist, of which any could have been used to identify the best feature set, including artificial neural networks (ANNs), support
5 vector machines (SVMs), partial least squares discriminative analysis (PLSDA), sub-linear association methods, Bayesian inference methods, supervised principle component analysis, shrunken centroids, or others (see [25] for review).

[0079] With reference to Examples 1 to 4, and based on the similarity of molecular formula, MS/MS fragmentation patterns, and NMR data, the metabolites identified in
10 serum, or subsets thereof, comprising the diagnostic feature set may be chemically related. In addition, there are many other related compounds present in the FTMS dataset that also show increased abundance in the MULTIPLE SCLEROSIS population, and which share similar molecular formulas to the subset identified. Therefore, the results suggest that an entire family of metabolites sharing common structural
15 properties is abnormal in MULTIPLE SCLEROSIS patients. Without wishing to be bound by theory, the biochemical pathway responsible for regulating the levels of these metabolites may be perturbed in MULTIPLE SCLEROSIS patients, and consequently may be a putative interventional target for treatment. Possible types of intervention include the development of agonists or antagonists for proteins involved in the implicated
20 pathways and/or the development of nutritional supplements that would decrease the concentration of the implicated metabolites or the development of pro-drugs or pro-nutrients to decrease the concentration of these metabolites.

[0080] The present invention also provides the structural characteristics of the metabolites used for the differential diagnosis of RR-MULTIPLE SCLEROSIS, PP-
25 MULTIPLE SCLEROSIS, and SP-MULTIPLE SCLEROSIS, which may include accurate mass and molecular formula determination, polarity, acid/base properties, NMR spectra, and MS/MS or MSⁿ spectra. Techniques used to determine these characteristics include, but are not limited to, reverse phase LC-MS using a C18 column followed by analysis by MS, MS/MS fragmentation using collision induced dissociation (CID),
30 nuclear magnetic resonance (NMR), and extraction. The characteristics of the metabolites obtained by various methods are then used to determine the structure of the metabolites.

[0081] The present invention also provides high throughput methods for differential diagnosis of MULTIPLE SCLEROSIS and normal states. The method involves fragmentation of the parent molecule; in a non-limiting example, this may be accomplished by a Q-Trap™ system. Detection of the metabolites may be performed using one of various assay platforms, including colorimetric chemical assays (UV, or other wavelength), antibody-based enzyme-linked immunosorbant assays (ELISAs), chip- and PCR-based assays for nucleic acid detection, bead-based nucleic-acid detection methods, dipstick chemical assays or other chemical reaction, image analysis such as magnetic resonance imaging (MRI), positron emission tomography (PET) scan, computerized tomography (CT) scan, nuclear magnetic resonance (NMR), and various mass spectrometry-based systems.

[0082] The metabolites and the methods of the present invention may also be combined with the current diagnostic tools for MULTIPLE SCLEROSIS, which include clinical history, neuroimaging analysis, evoked potentials, and cerebrospinal fluid analysis of proteinaceous and inflammatory components within the cerebrospinal fluid. Imaging techniques include, but are not limited to, structural magnetic resonance imaging (MRI), contrast-enhanced MRI, positron emission tomography (PET), computerized tomography (CT), functional magnetic resonance imaging (fMRI), electroencephalography (EEG), single positron emission tomography (SPECT), event related potentials, magnetoencephalography, and/or multi-modal imaging. The clinical assessment may include, but is not limited to, the Kurtzke's extended disability status scale (EDSS), multiple sclerosis impact scale (MSIS), Scripps neurologic rating scale (NRS), ambulation index (AI), MS-related symptoms scale, 15-item activities of daily living self-care scale for MS Persons, Incapacity status scale, functional independent measure, and/or internuclear ophthalmoplegia. A person skilled in the art would recognize that the combination of metabolites and methods as described herein with current techniques has the potential to diagnosis or differentiate any form of multiple sclerosis and/or its pathology.

[0083] The present invention will be further illustrated in the following examples.

30 Example 1: Identification of Differentially Expressed Metabolites

[0084] Differentially expressed metabolites are identified in individuals with clinically diagnosed RR-MULTIPLE SCLEROSIS, clinically diagnosed PP-MULTIPLE SCLEROSIS, clinically diagnosed SP-MULTIPLE SCLEROSIS, as well as healthy controls.

[0085] *Clinical Samples.* For the MULTIPLE SCLEROSIS serum diagnostic assay described, samples were obtained from representative populations of healthy individuals and those with clinically diagnosed RR-MULTIPLE SCLEROSIS, clinically diagnosed PP-MULTIPLE SCLEROSIS, and clinically diagnosed SP-MULTIPLE SCLEROSIS patients. The biochemical markers of RR-MULTIPLE SCLEROSIS described in the invention were derived from the analysis of 93 serum samples from patients clinically diagnosed with RR-MULTIPLE SCLEROSIS, serum samples from 18 patients with clinically diagnosed PP-MULTIPLE SCLEROSIS, serum samples from 22 patients with clinically diagnosed SP-MULTIPLE SCLEROSIS, and 51 serum samples from controls. The 93 patients with RR-MULTIPLE SCLEROSIS were further divided into one of two groups: those still exhibiting a relapsing-remitting disease course (mean disease duration 5.9 years, n = 46) and those transitioning into the chronic secondary-progressive disease course (mean disease duration 11.4 years, n = 47). Samples in the four groups were from a diverse population of individuals, ranging in age, demographic, weight, occupation, and displaying varying non-MULTIPLE SCLEROSIS-related health-states. All samples were single time-point collections

[0086] The metabolites contained within the 184 serum samples used in this example were separated into polar and non-polar extracts through sonication and vigorous mixing (vortex mixing).

[0087] Analysis of serum extracts collected from 184 individuals (93 clinically diagnosed RR-MULTIPLE SCLEROSIS, 18 clinically diagnosed PP-MULTIPLE SCLEROSIS, 22 clinically diagnosed SP-MULTIPLE SCLEROSIS, and 51 healthy controls) was performed by direct injection into a FTMS and ionization by either ESI or atmospheric pressure chemical ionization (APCI) in both positive and negative modes. Sample extracts were diluted either three or six-fold in methanol:0.1%(v/v) ammonium hydroxide (50:50, v/v) for negative ionization modes, or in methanol:0.1% (v/v) formic acid (50:50, v/v) for positive ionization modes. For APCI, sample extracts were directly injected without diluting. All analyses were

performed on a Bruker Daltonics APEX III Fourier transform ion cyclotron resonance mass spectrometer equipped with a 7.0 T actively shielded superconducting magnet (Bruker Daltonics, Billerica, MA). Samples were directly injected using electrospray ionization (ESI) and APCI at a flow rate of 1200 μ L per hour. Ion transfer/detection
5 parameters were optimized using a standard mix of serine, tetra-alanine, reserpine, Hewlett-Packard tuning mix and the adrenocorticotrophic hormone fragment 4-10. In addition, the instrument conditions were tuned to optimize ion intensity and broadband accumulation over the mass range of 100–1000 amu according to the instrument manufacturer's recommendations. A mixture of the abovementioned standards was
10 used to internally calibrate each sample spectrum for mass accuracy over the acquisition range of 100–1000 amu.

[0088] In total six separate analyses comprising combinations of extracts and ionization modes were obtained for each sample:

Aqueous Extract

- 15 1. Positive ESI (analysis mode 1101)
 2. Negative ESI (analysis mode 1102)

Organic Extract

3. Positive ESI (analysis mode 1201)
 4. Negative ESI (analysis mode 1202)
20 5. Positive APCI (analysis mode 1203)
 6. Negative APCI (analysis mode 1204)

[0089] *Mass Spectrometry Data Processing.* Using a linear least-squares regression line, mass axis values were calibrated such that each internal standard mass peak had a mass error of <1 ppm compared with its theoretical mass. Using XMASS software
25 from Bruker Daltonics Inc., data file sizes of 1 megaword were acquired and zero-filled to 2 megawords. A sinm data transformation was performed prior to Fourier transform and magnitude calculations. The mass spectra from each analysis were integrated, creating a peak list that contained the accurate mass and absolute intensity of each peak. Compounds in the range of 100–2000 m/z were analyzed. In order to
30 compare and summarize data across different ionization modes and polarities, all detected mass peaks were converted to their corresponding neutral masses assuming hydrogen adduct formation. A self-generated two-dimensional (mass vs. sample

intensity) array was then created using DISCOVAmetrics™ software (Phenomenome Discoveries Inc., Saskatoon, SK, Canada). The data from multiple files were integrated, and this combined file was then processed to determine all of the unique masses. The average of each unique mass was determined, representing the y axis.

- 5 This value represents the average of all of the detected accurate masses that were statistically determined to be equivalent. Considering that the mass accuracy of the instrument for the calibration standards is approximately 1 ppm, a person skilled in the art will recognize that these average masses may include individual masses that fall within +/- 5 ppm of this average mass. A column was created for each file that
- 10 was originally selected to be analyzed, representing the x axis. The intensity for each mass found in each of the files selected was then filled into its representative x,y coordinate. Coordinates that did not contain an intensity value were left blank. Once in the array, the data were further processed, visualized and interpreted, and putative chemical identities were assigned. Each of the spectra were then peak picked to obtain
- 15 the mass and intensity of all metabolites detected. These data from all of the modes were then merged to create one data file per sample. The data from all 184 samples were then merged and aligned to create a two-dimensional metabolite array in which each sample is represented by a column and each unique metabolite is represented by a single row. In the cell corresponding to a given metabolite sample combination, the
- 20 intensity of the metabolite in that sample is displayed. When the data is represented in this format, metabolites showing differences between groups of samples can be determined.

[0090] *Advanced Data Interpretation - Serum Biomarkers.* A student's T-test was used to select for metabolites which differed significantly between the following

- 25 different clinical groups in serum:

1. clinically diagnosed RR-MULTIPLE SCLEROSIS patients (n=46) and controls (n=51), [240 metabolites, see Table 1];
2. clinically diagnosed PP-MULTIPLE SCLEROSIS patients (n=18) and controls (n=51), [60 metabolites, see Table 2];
- 30 3. clinically diagnosed SP-MULTIPLE SCLEROSIS patients (n=22) and controls (n=51), [129 metabolites, see Table 3];

4. clinically diagnosed RR-MULTIPLE SCLEROSIS patients (n=46) and clinically diagnosed SP-MULTIPLE SCLEROSIS (n=22), [135 metabolites, see Table 4];
5. clinically diagnosed RR-MULTIPLE SCLEROSIS patients (n=46) and RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] (n=47), [148 metabolites, see Table 5];
6. RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] (n=47) and SP-MULTIPLE SCLEROSIS patients (n=22), [42 metabolites, see Table 6].
- 10 [0091] Metabolites that were less than $p < 0.05$ were considered significant.
- [0092] Tables 1-6 show metabolite features whose concentrations or amounts in serum are significantly different ($p < 0.05$) between the tested populations and therefore have potential diagnostic utility for identifying each of the aforesaid populations. The features are described by their accurate mass and analysis mode, which together are
- 15 sufficient to provide the putative molecular formulas and chemical characteristics (such as polarity and putative functional groups) of each metabolite.
- [0093] For each clinical pairing, a cross-validated training classifier was created using the PAM algorithm, previously described. The classifier algorithm was trained using samples with known diagnosis and then applied to blinded sample (i.e. a test set).
- 20 [0094] The lowest training classifier obtained with the fewest number of metabolites was selected for each clinical pairing. The graph in Figure 1A shows the number of metabolites required to achieve a given training error at various threshold values (a user-definable PAM parameter). The plot shows that a training classifier with less than 22% error rate (0.22 training error) is possible with five metabolite features
- 25 (threshold value of approximately 3.59, see arrow). The graph in Figure 1B is conceptually similar to that in 1A, however, the graph in 1B shows the misclassification error of the trained classifier for clinically diagnosed RR-MULTIPLE SCLEROSIS patients and control patients following the cross-validation procedure integral to the PAM program. The line connected by the diamonds mirrors the
- 30 previous result, showing that minimal cross-validated misclassification error for controls were achieved using as few as five metabolites. It also shows that clinically

diagnosed RR-MULTIPLE SCLEROSIS patients, depicted by the squares, were 93% accurately diagnosed as having RR-MULTIPLE SCLEROSIS using only three metabolite feature, but at this threshold, the misclassification for the controls was 66% (see arrows). The individual cross-validated diagnostic probabilities for each of the RR-MULTIPLE SCLEROSIS patients and controls are shown in Figure 2. All of clinically diagnosed RR-MULTIPLE SCLEROSIS patients are listed on right side of the graph, and the controls are on the left. Each sample contains two points on the graph, one showing the probability of having RR-MULTIPLE SCLEROSIS (squares), and one showing the probability of not having RR-MULTIPLE SCLEROSIS (i.e. normal, diamonds). From the graph, six RR-MULTIPLE SCLEROSIS samples were classified as non-MULTIPLE SCLEROSIS and five control samples were classified as RR-MULTIPLE SCLEROSIS. The five metabolites are listed in Table 7. The predicted probabilities were then used to create the receiver-operating characteristic (ROC) curve in Figure 3 using JROCFIT (<http://www.rad.jhmi.edu/jeng/javarad/roc/JROCFITi.html>), which shows the true positive fraction (those with RR-MULTIPLE SCLEROSIS being predicted to have RR-MULTIPLE SCLEROSIS) versus the false positive fraction (control individuals predicted as having RR-MULTIPLE SCLEROSIS). The area under the curve is 81.4%, with a sensitivity of 94.3%, and a specificity of 72.5%. Overall, the diagnostic accuracy is 81.4% based on the cross-validated design.

[0095] The above first principle component analysis allowed the initial identification of the optimal metabolites for each clinical pairing. In order to confirm these findings, a second PAM analysis was performed. For each clinical pairing, the second analysis (discussed below) generally provided a larger number of metabolites than the first principle component analysis. From this expanded set of metabolites, the best candidates for differentiation between clinical health states, which generally correspond to the initially identified metabolites, were identified.

[0096] In the second PAM method, the samples for each clinical pairing were randomly split in half, using one half to generate a classifier and other half as a blinded “test set” for diagnosis. Since the first method creates the classifier using more samples, its predictive accuracy would be expected to be higher than the second approach, and consequently requires a fewer number of metabolites for high diagnostic accuracy. Following the previous example of all clinically diagnosed RR-

MULTIPLE SCLEROSIS patients and controls, the training set was comprised of 30 clinically diagnosed RR-MULTIPLE SCLEROSIS patients and 26 controls. The predicted probabilities of the blinded test samples as either being RR-MULTIPLE SCLEROSIS-specific or controls are plotted in Figure 4. The results show four of the clinically-
5 diagnosed RR-MULTIPLE SCLEROSIS samples were given a higher probability of being controls and four of the controls were given a higher probability of being RR-MULTIPLE SCLEROSIS. The optimal number of metabolites required for the lowest misclassification error using these samples was 16, listed in Table 8. The classifier was next used to predict the diagnosis of the remaining samples (blinded; 17 clinically
10 diagnosed RR-MULTIPLE SCLEROSIS patients and 25 controls). Table 9 contains the patients that were used in the test set and their actual and predicted diagnosis. The probabilities from Figure 4 were then translated into a ROC curve (Figure 5). The performance characteristics based on classification of the blinded test set were sensitivity of 76.5%, specificity of 84.0%, and overall diagnostic accuracy of 81.0%.

15 [0097] The PAM analysis was repeated for each of the clinical pairings. The sample numbers used in each training set as well as the optimal number of metabolites required for the lowest misclassification error are listed in Table 10. The classifiers for the training sets were next used to predict the diagnosis of the remaining samples for each clinical pairing.

20 [0098] *i) Clinically diagnosed PP-MULTIPLE SCLEROSIS patients and controls.* Table 11 contains the expanded set of metabolites and the actual and predicted diagnosis of the patients that were used in the test set. The probabilities from Table 11 were translated into a ROC curve (Figure 6). The performance characteristics based on the classification of the blinded test set were: sensitivity of 44.4%, specificity of 92%, and
25 overall diagnostic accuracy of 79.4%.

[0099] *ii) Clinically diagnosed SP-MULTIPLE SCLEROSIS patients and controls.* Table 12 contains the expanded set of metabolites and the actual and predicted diagnosis of the patients that were used in the test set. The probabilities from Table 12 were translated into a ROC curve (Figure 7). The performance characteristics based
30 on the classification of the blinded test set were: sensitivity of 63.6%, specificity of 100%, and overall diagnostic accuracy of 88.9%.

[00100] **iii) Clinically diagnosed RR-MULTIPLE SCLEROSIS patients and SP-MULTIPLE SCLEROSIS patients.** Table 13 contains the expanded set of metabolites and the actual and predicted diagnosis of the patients that were used in the test set. The probabilities from Table 13 were translated into a ROC curve (Figure 8). The performance characteristics based on the classification of the blinded test set were: sensitivity of 88.9%, specificity of 100%, and overall diagnostic accuracy of 97.1%.

[00101] **iv) Clinically diagnosed RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] and RR-MULTIPLE SCLEROSIS patients.** Table 14 contains the expanded set of metabolites and the actual and predicted diagnosis of the patients that were used in the test set. The probabilities from Table 14 were translated into a ROC curve (Figure 9). The performance characteristics based on the classification of the blinded test set were: sensitivity of 100%, specificity of 92.3%, and overall diagnostic accuracy of 95.7%.

[00102] **v) Clinically diagnosed RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] and SP-MULTIPLE SCLEROSIS patients.** Table 15 contains the expanded set of metabolites and the actual and predicted diagnosis of the patients that were used in the test set. The probabilities from Table 15 were translated into a ROC curve (Figure 10). The performance characteristics based on the classification of the blinded test set were: sensitivity of 72.7%, specificity of 95.5%, and overall diagnostic accuracy of 87.9%.

[00103] Using an initial panel of about 240 metabolites, and an expanded set of about 16 metabolites, it was determined that a combination of nine metabolites fulfills the criteria for a serum diagnostic test of RR-MULTIPLE SCLEROSIS compared to normal samples. The best combination of nine metabolites includes the metabolites with masses (measured in Daltons) 452.3868, 496.4157, 524.4448, 540.4387, 578.4923, 580.5089, 594.4848, 596.5012, 597.5062. Although these are the actual masses, a person skilled in the art of this technology would recognize that +/- 5 ppm difference would indicate the same metabolite.

[00104] Using an initial panel of about 60 metabolites, and an expanded set of about 7 metabolites, it was determined that a combination of five metabolites fulfills

the criteria for a serum diagnostic test of PP-MULTIPLE SCLEROSIS compared to normal samples. The best combination of five metabolites includes the metabolites with masses (measured in Daltons) 202.0453, 216.04, 243.0719, 244.0559, 857.7516, where a +/- 5 ppm difference would indicate the same metabolite.

- 5 [00105] Using an initial panel of about 129 metabolites, and an expanded set of about 16 metabolites, it was determined that a combination of eighteen metabolites fulfills the criteria for a serum diagnostic test of SP-MULTIPLE SCLEROSIS compared to normal samples. The best combination of eighteen metabolites includes the metabolites with masses (measured in Daltons) 194.0803, 428.3653, 493.385,
10 541.3415, 565.3391, 576.4757, 578.4923, 590.4964, 594.4848, 495.4883, 596.5012, 596.5053, 597.5062, 597.5068, 805.5609, 806.5643, 827.5446, 886.5582, where a +/- 5 ppm difference would indicate the same metabolite.

- [00106] Using an initial panel of about 135 metabolites, and an expanded set of about 16 metabolites, it was determined that a combination of six metabolites fulfills
15 the criteria for a serum indicator of RR-MULTIPLE SCLEROSIS compared to SP-MULTIPLE SCLEROSIS. The best combination of six metabolites includes the metabolites with masses (measured in Daltons) 540.4387, 576.4757, 594.4848, 595.4883, 596.5012, 597.5062, where a +/- 5 ppm difference would indicate the same metabolite.

- 20 [00107] Using an initial panel of about 148 metabolites, and an expanded set of about 9 metabolites, it was determined that a combination of 5 metabolites fulfills the criteria for a serum indicator of RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] compared to RR-MULTIPLE SCLEROSIS patients. The best combination of five metabolites includes the metabolites with masses (measured in
25 Daltons) 576.4757, 578.4923, 594.4848, 596.5012, 597.5062, where a +/- 5 ppm difference would indicate the same metabolite.

- [00108] Using an initial panel of about 42 metabolites, and an expanded set of about 17 metabolites, it was determined that a combination of 8 metabolites fulfills the criteria for a serum indicator of RR-MULTIPLE SCLEROSIS patients transitioning to
30 SP-MULTIPLE SCLEROSIS [RR-SP] compared to SP-MULTIPLE SCLEROSIS patients. The

best combination of eight metabolites includes the metabolites with masses (measured in Daltons) 617.0921, 746.5118, 760.5231, 770.5108, 772.5265, 784.5238, 786.5408, 787.5452, where a ± 5 ppm difference would indicate the same metabolite.

[00109] Bar graphs representing the mean $\pm$ SEM of the biomarkers for the
5 different clinical groups are shown in Figures 11-16. Relative to control individuals, the three non-control states can be described as follows:

[00110] 1. RR-MULTIPLE SCLEROSIS vs. control:

- a. Biomarker 452.3868 – increased
- b. Biomarker 496.4157 – increased
- 10 c. Biomarker 524.4448 – increased
- d. Biomarker 540.4387 – increased
- e. Biomarker 578.4923 – increased
- f. Biomarker 580.5089 – increased
- g. Biomarker 594.4848 – increased
- 15 i. Biomarker 596.5012 – increased
- h. Biomarker 597.5062 – increased

[00111] 2. PP-MULTIPLE SCLEROSIS vs. control:

- a. Biomarker 202.0453 – increased
- b. Biomarker 216.0400 – increased
- 20 c. Biomarker 243.0719 – increased
- d. Biomarker 244.0559 – increased
- e. Biomarker 857.7516 – increased

[00112] 3. SP-MULTIPLE SCLEROSIS vs. control:

- a. Biomarker 194.0803 – decreased
- 25 b. Biomarker 428.3653 – increased
- c. Biomarker 493.3850 – decreased
- d. Biomarker 541.3415 – decreased
- e. Biomarker 565.3391 – decreased
- f. Biomarker 576.4757 – decreased
- 30 g. Biomarker 578.4923 – decreased

- h. Biomarker 590.4964 – decreased
- i. Biomarker 594.4848 – decreased
- j. Biomarker 595.4883 – decreased
- k. Biomarker 596.5012 – decreased
- 5 l. Biomarker 596.5053 – decreased
- m. Biomarker 597.5062 – decreased
- n. Biomarker 597.5068 – decreased
- o. Biomarker 805.5609 – increased
- p. Biomarker 806.5643 – increased
- 10 q. Biomarker 827.5446 – increased
- r. Biomarker 886.5582 – decreased

[00113] Relative to RR-MULTIPLE SCLEROSIS patients, the two chronic clinical groups can be described as follows:

- [00114] 1. SP-MULTIPLE SCLEROSIS vs. RR-MULTIPLE SCLEROSIS:
- 15 a. Biomarker 540.4387 – decreased
 - b. Biomarker 576.4757 – decreased
 - c. Biomarker 594.4848 – decreased
 - d. Biomarker 595.4883 – decreased
 - e. Biomarker 596.5012 – decreased
 - 20 f. Biomarker 597.5062 – decreased

[00115] 2. RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] vs. RR-MULTIPLE SCLEROSIS:

- a. Biomarker 576.4757 – decreased
- b. Biomarker 578.4923 – decreased
- 25 c. Biomarker 594.4848 – decreased
- d. Biomarker 596.5012 – decreased
- e. Biomarker 597.5062 – decreased

[00116] Relative to SP-MULTIPLE SCLEROSIS patients, the RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] can be described as follows:

30

[00117] 1. RR-MULTIPLE SCLEROSIS transitioning to SP-MULTIPLE SCLEROSIS

[RR-SP] vs. SP-MULTIPLE SCLEROSIS:

- a. Biomarker 617.0921 -- increased
- b. Biomarker 746.5118 -- increased
- 5 c. Biomarker 760.5231 -- increased
- d. Biomarker 770.5108 -- increased
- e. Biomarker 772.5265 -- increased
- f. Biomarker 784.5238 -- increased
- g. Biomarker 786.5408 -- increased
- 10 e. Biomarker 787.5452 -- increased

[00118] The biomarker panels were then applied to the various clinical groups and the ten patients for each clinical group that showed the best separation were selected. A student's T-test was performed on all the serum metabolites using only ten patients per clinical group.

- 15 1. Clinically diagnosed RR-MULTIPLE SCLEROSIS patients (n=10) and controls (n=10), [257 metabolites, see Table 16];
- 2. Clinically diagnosed PP-MULTIPLE SCLEROSIS patients (n=10) and controls (n=10), [100 metabolites, see Table 17];
- 3. Clinically diagnosed SP-MULTIPLE SCLEROSIS patients (n=10) and controls
- 20 (n=10), [226 metabolites, see Table 18];
- 4. Clinically diagnosed RR-MULTIPLE SCLEROSIS patients (n=10) and clinically diagnosed SP-MULTIPLE SCLEROSIS (n=10), [142 metabolites, see Table 19];
- 5. RR- MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] (n=10) and clinically diagnosed RR-MULTIPLE SCLEROSIS patients
- 25 (n=10), [148 metabolites, see Table 20];
- 6. RR-MULTIPLE SCLEROSIS patients transitioning to SP-MULTIPLE SCLEROSIS [RR-SP] (n=10) and clinically diagnosed SP-MULTIPLE SCLEROSIS patients (n=10), [309 metabolites, see Table 19].

30 [00119] The sample set (184 individuals) used for this example was comprised of individuals of various geographical backgrounds, and of varying age and health

status. Therefore, it is expected that the findings are representative of the general MULTIPLE SCLEROSIS population.

Example 2: Independent Method Confirmation of Discovered Metabolites

5 [00120] The metabolites and their associations with the clinical variables described in this invention are further confirmed using an independent mass spectrometry system. Representative sample extracts from each variable group are re-analyzed by LC-MS using an HP 1050 high-performance liquid chromatography (HPLC), or equivalent, interfaced to an ABI Q-Star, or equivalent, mass spectrometer
10 to obtain mass and intensity information for the purpose of identifying metabolites that differ in intensity between the clinical variables under investigation.

[00121] By determining the levels of the identified metabolites in a person's blood and comparing these levels to levels in a normal "reference" population, a prediction is made whether the person has RR-MULTIPLE SCLEROSIS, PP-MULTIPLE
15 SCLEROSIS, or early stages of SP-MULTIPLE SCLEROSIS. This is carried out in one of several ways: 1) Using a prediction algorithm to classify the test sample, as previously described, which outputs a percentage probability for having a form of MULTIPLE SCLEROSIS. A predictive approach would work independently of the assay method, as long as the intensities of the metabolites are measured. 2) Applying a method based on
20 setting a threshold intensity level from the mass spectrometer, and determining whether a person's profile is above or below the threshold, which indicates their disease status. 3) Using a quantitative assay to determine the molar concentration of the 36 serum metabolites in the normal and disease population. An absolute threshold concentration is then determined for MULTIPLE SCLEROSIS-positivity versus non-
25 MULTIPLE SCLEROSIS-positivity. In a clinical setting, this means that if the measured levels of the metabolites, or combinations of the metabolites, are above a certain concentration, there would be an associated probability that the individual is positive for a type of MULTIPLE SCLEROSIS.

Example 3: Structure Elucidation of the Primary Metabolite Biomarkers

[00122] Characteristics that can be used for structure elucidation of metabolites include accurate mass and molecular formula, polarity, acid/base properties, NMR spectra, and MS/MS or MSⁿ spectra. These data can be used as fingerprints of a particular metabolite and are unique identifiers of a particular metabolite regardless of whether the complete structure has been determined. The data include:

- [00123] 1. *LC retention time.* The extracts containing the metabolites of interest are subjected to reverse phase LC-MS using a C18 column and analysis by MS to determine their retention time under standardized conditions.
- 10 [00124] 2. *MS/MS spectra.* Metabolites of interest are further characterized by performing MS/MS fragmentation using collision induced dissociation (CID). This MS/MS analysis is performed in real time (i.e. during the chromatographic elution process) or off-line on fractions collected from the chromatographic separation process. The structure of a given molecule dictates a specific fragmentation pattern under defined conditions and is specific for that molecule (equivalent to a person's fingerprint). Even slight changes to the molecule's structure can result in a different fragmentation pattern. In addition to providing a fingerprint of the molecule's identity, the fragments generated by CID are used to gain insights about the structure of a molecule, and for generating a very specific high-throughput quantitative detection method (see [26-29] for examples).

[00125] 3. *NMR spectra.* The MS/MS fragmentation provides highly specific descriptive information about a metabolite. However, NMR can solve and confirm the structures of the molecules. As NMR analysis techniques are typically less sensitive than mass spectrometry techniques, multiple injections are performed on the HPLC and the retention time window corresponding to the metabolites of interest collected and combined. The combined extract is then evaporated to dryness and reconstituted in the appropriate solvent for NMR analysis.

[00126] Multiple NMR techniques and instruments are available, for example, NMR spectral data are recorded on Bruker Avance 600 MHz spectrometer with cryogenic probe after the chromatographic separation and purification of the

metabolites of interest. ¹H NMR, ¹³C NMR, no-difference spec, as well as 2-D NMR techniques like heteronuclear multiple quantum correlation (HMQC), and heteronuclear multiple bond correlation (HMBC) are used for structure elucidation work on the biomarkers.

- 5 [00127] 4. *Extraction conditions.* The conditions of extraction also provide insights about the chemical properties of the biomarkers. All nine metabolites in the serum (from Example 1) were ionized in negative mode (APCI), which is indicative of a molecule containing an acidic moiety such as a carboxylic acid or phosphate. Any moiety capable of losing a hydrogen atom can be detected in negative ionization
10 mode. The metabolite markers were extracted into an organic ethyl acetate fraction, indicating that these metabolites are non-polar under acidic condition.

- [00128] All chemicals and media were purchased from Sigma-Aldrich Canada Ltd., Oakville, ON., Canada. All solvents were HPLC grade. HPLC analysis were carried out with a high performance liquid chromatograph equipped with quaternary
15 pump, automatic injector, degasser, and a Hypersil ODS column (5 µm particle size silica, 4.6 i.d × 200 mm) with an inline filter. Mobile phase: linear gradient H₂O-MeOH to 100% MeOH in a 52 min period at a flow rate of 1.0 ml/min. High resolution (HR) mass spectra (MS) were recorded on Bruker apex 7T Fourier transform ion cyclotron resonance (FT-ICR) spectrometer and MS/MS data collected
20 using QStar XL TOF mass spectrometer with atmospheric pressure chemical ionization (APCI) and electro spray ionization (ESI) sources in both positive and negative modes.

[00129] Metabolite Characterization Data

[00130] Biomarker 1

- 25 [00131] HRAPCI-MS *m/z*: [M - H]⁻, C₂₈H₅₁O₄⁻ measured; 451.3795, calcd. 451.3793. MS/MS *m/z* (relative intensity): 451 ([M - H]⁻, 20%), 433 (100%), 407(30%), 389 (90%), 281 (10%), 279 (25%), 183 (20%), 169 (10%), 153 (10%), 125 (20%), 111(25%), 97 (25%).

[00132] Biomarker 2

[00133] HRAPCI-MS m/z : $[M - H]^-$, $C_{30}H_{35}O_5^-$ measured; 495.4054, calcd. 495.4055 MS/MS m/z (relative intensity): 495 ($[M - H]^-$, 5%), 451 (5%), 477 (15%), (433 (15%), 415(5%), 307 (5%), 297 (45%), 279 (100%), 235 (5%), 223 (20%), 215 (70%), 197 (90%), 179(50%), 181 (10%), 169 (100%), 157 (25%), 155 (10%), 153 (5%), 141 (10%), 139 (5%), 127 (10%), 125 (10%), 113 (5%).

[00134] Biomarker 3

[00135] HRAPCI-MS m/z : $[M - H]^-$, $C_{32}H_{39}O_5^-$ measured; 523.4375, calcd; 523.4368. MS/MS m/z (relative intensity): 523 ($[M - H]^-$, 30%), 505 (100%), 487 (25%), 479 (40%), 463 (40%), 461 (45%), 443 (40%), 365 (30%), 337 (20%), 299 (25%), 297 (25%), 281 (25%), 279 (40%), 271 (65%), 269 (20%), 253 (35%), 251 (55%), 243 (30%), 225 (65%), 197 (55%), 171 (20%), 169(25%), 157 (20%), 155 (10%), 143 (10%), 141 (20%), 139 (20%).

[00136] Biomarker 4

[00137] HRAPCI-MS m/z : $[M - H]^-$, $C_{32}H_{39}O_6^-$ measured; 539.4312, calcd; 539.4317. MS/MS m/z (relative intensity): 539 ($[M - H]^-$, 20%), 521 (100%), 503 (50%), 495 (40%), 477 (40%), 461 (30%), 459 (40%), 419 (30%), 335 (70%), 315 (40%), 313 (40%), 297 (60%), 279 (90%), 259 (40%), 255 (40%), 253 (20%), 243 (20%), 241 (30%), 225 (20%), 223 (30%), 213 (30%), 179 (20%), 171 (40%), 155 (30%), 141 (50%), 127 (40%).

20 [00138] Biomarker 5

[00139] HRAPCI-MS m/z : $[M - H]^-$, $C_{36}H_{63}O_5^-$ measured; 575.4678, calcd; 575.4681. MS/MS m/z (relative intensity): 575 ($[M - H]^-$, 45%), 557 (75%), 539 (70%), 531 (30%), 513 (60%), 495 (100%), 417 (50%), 403 (60%), 371 (25%), 297 (15%), 279 (40%).

25 [00140] Biomarker 6

[00141] HRAPCI-MS m/z : $[M - H]^-$, $C_{36}H_{65}O_5^-$ measured; 577.4850, calcd; 577.4837. MS/MS m/z (relative intensity): 577 ($[M - H]^-$, 45%), 559 (75%), 541

(70%), 533 (30%), 515 (60%), 497 (100%), 419 (50%), 405 (60%), 387 (25%), 373 (25%), 297 (15%), 281 (25%), 279 (40%).

[00142] Biomarker 7

[00143] HRAPCI-MS m/z : $[M - H]^-$, $C_{36}H_{67}O_5^-$ measured; 579.5016, calcd; 579.4994. MS/MS m/z (relative intensity): 579 ($[M - H]^-$, 45%), 561 (90%), 543 (40%), 535 (25%), 517 (60%), 499 (100%), 421 (20%), 407 (20%), 389 (20%), 375 (20%), 299 (25%), 281 (30%), 279 (40%), 263 (10%), 253 (15%), 185 (10%), 171 (25%).

[00144] Biomarker 8

10 [00145] HRAPCI-MS m/z : $[M - H]^-$, $C_{36}H_{65}O_6^-$ measured; 593.4775, calcd; 593.4787. MS/MS m/z (relative intensity): 593 ($[M - H]^-$, 50%), 575 (55%), 557 (30%), 549 (15%), 531 (20%), 513 (25%), 495 (10%), 421 (15%), 371 (30%), 315 (50%), 297 (100%), 279 (90%). 201 (30%), 171 (60%), 141 (25%), 127 (25%).

[00146] Biomarker 9

15 [00147] HRAPCI-MS m/z : $[M - H]^-$, $C_{36}H_{67}O_6^-$ measured; 595.4939, calcd; 595.4943. MS/MS m/z (relative intensity): 595 ($[M - H]^-$, 20%), 577 (20%), 559 (15%), 551 (5%), 515 (15%), 497 (5%), 423 (5%), 373 (15%), 315 (75%), 297 (70%), 281 (40%), 279 (100%), 269 (5%), 251 (5%), 171 (25%), 155 (15%), 153 (10%), 141 (15%), 139 (10%), 127 (15%).

20 [00148] Biomarker 10

[00149] HRAPCI-MS m/z : $[M - H]^-$, $C_{43}H_{78}O_{10}P^-$ measured; 785.5329, calcd; 785.5338. MS/MS m/z (relative intensity): 758 ($[M - H]^-$, 100%), 529 (10%), 425 (20%), 273 (73%), 169 (5%), 125 (100%), 97 (5%).

[00150] Biomarker 11

25 [00151] HRAPCI-MS m/z : $[M - H]^-$, $C_5H_{12}O_7P^-$ measured; 215.0322, calcd; 215.0326. MS/MS m/z (relative intensity): 215 ($[M - H]^-$, 100%), 197 (30%), 171 (40%), 153 (90%), 135 (20%).

- [00152] Biomarker 12
- [00153] HRAPCI-MS m/z : $[M - H]^-$, $C_{25}H_{51}NO_9P^-$ measured; 540.3337, calcd; 540.3301. MS/MS m/z (relative intensity): 540 ($[M - H]^-$, >1%), 480 (17%), 255 (100%), 242 (>1%), 224 (5%), 168 (>1%), 153 (>1%), 78 (>1%).
- 5 [00154] Biomarker 13
- [00155] HRAPCI-MS m/z : $[M - H]^-$, $C_{27}H_{51}NO_9P^-$ measured; 564.3313, calcd; 564.3307. MS/MS m/z (relative intensity): 564 ($[M - H]^-$, >1%), 504 (10%), 279 (100%), 242 (>1%), 224 (5%), 168 (>1%), 153 (>1%), 78 (>1%).
- [00156] Biomarker 14
- 10 [00157] HRAPCI-MS m/z : $[M + H]^+$, $C_6H_{12}O_6Na^+$ measured; 203.0531, calcd; 205.0526. MS/MS m/z (relative intensity): 203 ($[M + H]^+$, 100%), 159 (15%), 115 (23%), 89 (38%), 97 (5%).
- [00158] Biomarker 15
- [00159] HRAPCI-MS m/z : $[M + H]^+$, $C_8H_{13}O_7Na^+$ measured; 245.0637, calcd; 15 245.0631. MS/MS m/z (relative intensity): 245 ($[M + H]^+$, 100%), 227 (5%), 209 (5%), 155 (10%), 125 (15%), 83 (5%).
- [00160] Biomarker 16
- [00161] HRAPCI-MS m/z : $[M + H]^+$, $C_{29}H_{49}O_2^+$ measured; 429.3732, calcd; 429.3727. MS/MS m/z (relative intensity): 429 ($[M + H]^+$, 1%), 205 (5%), 165 20 (100%).
- [00162] Biomarker 17
- [00163] HRAPCI-MS m/z : $[M + H]^+$, $C_{46}H_{81}NO_8P^+$ measured; 806.5687, calcd; 806.5694. MS/MS m/z (relative intensity): 806 ($[M + H]^+$, 21%), 478 (>1%), 237 (>1%), 184 (100%).
- 25 [00164] Biomarker 18

[00165] HRAPCI-MS m/z : $[M + H]^+$, $C_7H_{15}O_6^+$ measured; 195.0881, calcd; 195.0863. MS/MS m/z (relative intensity): 195 ($[M + H]^+$, 2%), 177 (>1%), 165 (>1%), 163 (>1%), 138 (100%), 123 (6%).

[00166] Biomarker 19

5 [00167] HRAPCI-MS m/z : $[M + H]^+$, $C_{54}H_{100}NO_6^+$ measured; 858.7594, calcd; 858.7545. MS/MS m/z (relative intensity): 858 ($[M + H]^+$, 100%), 576 (10%), 314 (12%), 165 (7%), 151 (10%), 95 (2%).

[00168] The accurate masses of the biomarkers were used to deduce the molecular formulae. Tandem mass spectrometry on the biomarkers were used to
10 propose the structures that are summarized in Table 22. The biomarkers were thought to be derivatives of sugars, phospholipids and tocopherols.

[00169] The MS/MS spectral data obtained for each of the multiple sclerosis biomarkers was used to deduce their structures. Upon comparing the MS/MS fragmentation patterns of MS biomarkers 1 - 9 against that of the CRC panel (see
15 applicant's co-pending application PCT/CA 2006/001502; published as WO/CA2007/030928 on March 22, 2007) a number of similarities were observed. In addition to the common ionization modes for both CRC and these MS biomarkers, their MS/MS spectra also showed signals due to fragment ions corresponding to phytol chain type fatty acid entities, C18:1 or C18:2 (m/z 281, 279) for all of the
20 detected biomarkers as well as fragment losses due to $[M-H-CO_2]$, $[M-H-H_2O]$ and $[M-H-CO_2-H_2O]^+$. Another similarity is that, the MS/MS spectra of MS biomarkers 1 - 9 showed fragment ions deduced as loss of chroman type ring system after cleavage of phytol side chain [(153 (1), 197 (2), 225 (3, 4), 279 (5, 6, 7) and 281 (8, 9), Tables 23 - 31)]. These observations led to the assignment of tocopherol type structures for
25 biomarkers 1 - 9. The loss(es) of water and carbon dioxide suggest the presence of free hydroxyl and carboxylic acid groups. The main differences between MS biomarkers and the CRC's as observed in the MS/MS spectra are the open chroman ring system and chain elongation proposed at position 1.

[00170] The molecular formula of 1 was determined as $C_{28}H_{52}O_4$ by HRAPCI-
30 MS, with three degrees of unsaturation. As indicated above, MS/MS spectra of 1

showed fragment ions due to loss of water (m/z 433), carbon dioxide (m/z 407) and presence of phytyl side chain (m/z 279). Fragment ion observed at m/z 153 was assigned as a cyclohexenyl ring system generated after the loss of the phytyl side chain. Based on these deductions the structure of metabolite 1 was assigned as shown in Table 22.

[00171] As indicated above, metabolites 2 - 9 have all the structural similarities to 1 and additional hydroxylations and chain elongations via ether linkages with the oxygen atom at position 1. The cyclohexenyl ring unit left after the cleavage of the phytyl side chains of these biomarkers gave unique fragment ions having some variation in the degrees of unsaturation and the number of hydroxylations. These ions observed at m/z 197 for 2, m/z 225 for 3 and 4, m/z 279 for 5, 6 and 8 and m/z 281 for 7 and 9 (See Tables 23 – 31) were used to assign the different alkyl chain elongations; ethyl, butyl and octyl respectively with the appropriate hydroxylations. In some detail, these fragmentation patterns clearly show the differences between each cyclohexenyl ring system. For 1 where there is no chain elongation at position 1, the cyclohexenyl ring fragment resulted when cleaved at C2-C3, generating the formula $C_{10}H_{17}O$ (m/z 153). In 2 where the ethylation is thought to occur at position 1, and with an additional hydroxy group on the ring, the formula of the cyclohexenyl ring fragment showed an increase by C_2H_4O entity compared to 1, thus the fragment having $C_{12}H_{21}O_2$ (m/z 197) as formula. These predictions complied with the observation in the MS/MS spectra of 2 thus validating the structural assignments. In 3 and 4, the chain elongation was thought to occur with a butyl unit (C_4H_9), thus an increase by C_2H_4 entity with formula $C_{14}H_{25}O_2$ (m/z 225) observed when compared to 2. For biomarkers 5 - 9 the alkoxy chain elongation at position 1 was by C_8H_{17} entity. Upon comparison of their formulae and MS/MS spectra, 7 and 9 ($C_{36}H_{68}O_5$ and $C_{36}H_{68}O_6$) showed similar features except for an additional oxygen atom in 9. This was consequently assigned on the phytyl chain. Therefore for 7 and 9 the cyclohexenyl ring component fragment ion was observed at m/z 281 ($C_{18}H_{33}O_2$). In the same vane, biomarkers 6 ($C_{36}H_{66}O_5$) and 8 ($C_{36}H_{66}O_6$) showed similarity like 7 and 9, the only difference being an added unsaturation, thus their cyclohexenyl fragment was at m/z 279 ($C_{18}H_{31}O_2$). An additional degree of unsaturation in 5 ($C_{36}H_{64}O_5$) compared to 6 and 8 but with ring fragment m/z 279 ($C_{18}H_{31}O_2$) suggested the additional unsaturation was on the phytyl

chain. Based on these deductions, the structures of metabolites 2 to 9 were assigned as shown in Table 22.

[00172] Biomarker 10 which was detected in the same mode as 1 -9 suggested a different class of metabolite based on the molecular formula FT-ICRMS data. The
5 obtained formula, $C_{43}H_{79}O_{10}P$ suggests a hydroxylated diacylglycerol-phospholipid type structure. The proposed structure and the MS/MS fragments are given in Table 22 and 32 respectively.

[00173] MS/MS data obtained on aqueous extracts of serum in the negative mode with electro spray ionization for biomarkers 11 - 13 were individually analyzed
10 to deduce their structures. The biomarkers identified in this panel were with the formulae of $C_5H_{13}O_7P$, $C_{25}H_{52}NO_9P$, and $C_{27}H_{52}NO_9P$. MS/MS data of 11 ($C_5H_{13}O_7P$) shows the fragments due to loss of two water molecules as well as a HPO_3 group (Table 33), which can be assigned using the proposed structure. Biomarkers 12 and
13, (m/z 541.3415, $C_{25}H_{52}NO_9P$ and m/z 565.3391, $C_{27}H_{52}NO_9P$) were found to be the
15 same as two Prostrate cancer biomarkers (see applicant's co-pending application PCT/CA 2007/000469, filed on March 23, 2007). In the negative mode with electro spray ionization, (ESI), the most commonly observed ions are the acidic phospholipids such as glycerophosphoinositol, glycerophosphoserine, glycerophosphatidic acid and glycerophosphoethanolamine. But under certain
20 circumstances it is possible that the phosphocholines can be detected as an adduct of $[M + Cl]^-$ or $[M + \text{acetate/formate}]^-$ as ion species in the negative ESI mode. Since the laboratory procedure of ESI aqueous extractions involves the use of formic acid there is a good probability that these ions could be the formate adduct of phosphocholines. As a result of the addition of the formate group forms a neutral cluster of
25 glycerophosphocholine which forms the corresponding molecular ion ($[M-H]^+$) upon subjected to negative ESI now that the ionization site is the phosphatidic group. This suggests the de-protonation of the phosphate group leaving the negatively charged phosphate ion as the parent ion. The fragmentation analysis of biomarkers 12 and 13 are given in Tables 34 and 35.

30 [00174] MS/MS data was obtained on organic and aqueous extracts of serum in positive mode with ESI and APCI for biomarkers 14 – 19. Biomarkers 14 and 15

(aqueous extract) were identified as sodium adducts of small monosaccharide related metabolites using their MS/MS fragment fingerprint (Tables 36 and 37). Biomarker 16 (Table 38), (m/z 428.3653, $C_{29}H_{48}O_2$) from organic extracts was assigned as a derivative of α tocopherol since its MS/MS spectra was quite similar to that of α tocopherol standard except for an additional degree of unsaturation. Biomarker 17 (Table 39), (m/z 805.5609, $C_{46}H_{80}NO_8P$) also from organic extracts of serum was proposed as Oleyl, eicosapentenoic(EPA), N-methyl phosphoethanolamine since the MS/MS data showed fragment ions for the presence of EPA and oleyl groups as well as the N-methyl substituted phosphoethanolamine back bone. The MS/MS spectral data of metabolites 18 (Table 40) and 19 (Table 41) using APCI source, were putatively assigned as monosaccharide and sphingolipid derived biomarkers respectively.

15 Example 4: High Throughput Commercial Method Development

[00175] For routine analysis of a subset of the metabolites described, a high throughput analysis method is developed. There are multiple types of cost-effective assay platform options currently available depending on the molecules being detected. These include colorimetric chemical assays (UV, or other wavelength), antibody-based enzyme-linked immunosorbant assays (ELISAs), chip- and PCR-based assays for nucleic acid detection, bead-based nucleic-acid detection methods, dipstick chemical assays, image analysis such as magnetic resonance imaging (MRI), positron emission tomography (PET) scan, computerized tomography (CT) scan, and various mass spectrometry-based systems.

25 [00176] The method involves the development of a high-throughput MS/MS method that is compatible with current laboratory instrumentation and triple-quadrupole mass spectrometers which are readily in place in many labs around the world. A Q-Trap™ system is used to isolate the parent molecule, fragment it; and then the fragments are measured.

[00177] All citations are hereby incorporated by reference.

[00178] The present invention has been described with regard to one or more
embodiments. However, it will be apparent to persons skilled in the art that a number
of variations and modifications can be made without departing from the scope of the
5 invention as defined in the claims.

[00179] Table 1: Accurate mass features differing between clinically diagnosed RR-MULTIPLE SCLEROSIS patients and controls ($p < 0.05$).

581.5126	1204	5.573	0.279	3.307	0.250	1.685	6.02E-09
452.3868	1204	3.933	0.137	2.758	0.161	1.426	9.63E-09
496.4157	1204	10.848	0.581	6.751	0.455	1.607	3.72E-08
524.4448	1204	5.474	0.257	3.608	0.241	1.517	7.83E-08
469.3863	1204	5.090	0.204	3.536	0.206	1.439	1.25E-07
580.5089	1204	13.697	0.678	8.725	0.641	1.570	1.26E-07
534.4645	1204	3.935	0.175	2.771	0.144	1.420	1.81E-07
510.3937	1204	6.354	0.265	4.445	0.265	1.429	1.99E-07
552.4784	1204	8.370	0.465	5.140	0.379	1.628	2.64E-07
468.384	1204	18.514	0.778	12.977	0.754	1.427	2.69E-07
506.2853	1201	4.224	0.345	3.261	0.273	1.295	2.95E-07
541.4422	1204	12.488	0.710	7.651	0.668	1.632	4.16E-07
484.3788	1204	8.292	0.379	5.832	0.359	1.422	5.37E-07
450.3729	1204	9.249	0.323	6.881	0.364	1.344	5.51E-07
494.3968	1204	14.347	0.613	10.018	0.645	1.432	5.77E-07
540.4387	1204	36.603	2.086	22.346	1.749	1.638	5.92E-07
522.4313	1204	15.891	0.597	11.438	0.681	1.389	9.00E-07
508.3782	1204	5.343	0.248	3.819	0.212	1.399	1.27E-06
578.4923	1204	41.017	2.169	26.583	2.068	1.544	1.37E-06
466.3656	1204	13.731	0.553	10.077	0.552	1.363	1.46E-06
610.482	1204	9.001	0.477	6.136	0.359	1.467	1.46E-06
536.41	1204	11.063	0.461	7.896	0.469	1.401	1.75E-06
566.454	1204	8.805	0.369	5.967	0.419	1.475	1.81E-06
440.3526	1204	4.439	0.187	3.273	0.180	1.356	1.84E-06
579.4958	1204	16.171	0.852	10.661	0.805	1.517	2.19E-06
480.3473	1204	3.955	0.153	3.030	0.142	1.305	2.27E-06
562.4989	1204	8.091	0.419	7.162	0.349	1.130	2.47E-06
482.3604	1204	5.195	0.241	3.772	0.196	1.377	2.55E-06
512.4079	1204	16.119	0.815	10.772	0.864	1.496	2.64E-06
568.4723	1204	13.768	0.729	8.786	0.743	1.567	3.70E-06
448.3562	1204	10.559	0.375	8.119	0.373	1.301	4.90E-06
569.4769	1204	5.509	0.287	3.508	0.309	1.570	6.73E-06
523.4337	1204	5.325	0.209	3.934	0.237	1.354	8.23E-06
495.4018	1204	4.656	0.199	3.441	0.203	1.353	8.62E-06
550.4602	1204	10.983	0.490	7.791	0.513	1.410	8.79E-06
538.4257	1204	30.014	1.397	21.271	1.294	1.411	9.34E-06
327.0307	1204	7.570	0.224	6.421	0.164	1.179	1.08E-05
513.4116	1204	5.468	0.275	3.821	0.296	1.431	1.18E-05
521.4188	1204	5.877	0.205	4.437	0.265	1.324	1.49E-05
564.513	1204	4.550	0.230	3.320	0.217	1.371	1.52E-05
493.385	1204	3.460	0.129	2.549	0.165	1.357	1.94E-05
467.3711	1204	4.522	0.194	3.431	0.188	1.318	2.01E-05
598.5107	1204	14.432	1.001	8.823	0.842	1.636	2.07E-05
520.4131	1204	16.144	0.600	12.295	0.696	1.313	2.10E-05
590.4585	1204	11.382	0.604	8.041	0.547	1.415	2.17E-05
548.4438	1204	7.502	0.289	5.581	0.347	1.344	2.51E-05
537.4142	1204	4.383	0.178	3.289	0.207	1.329	2.53E-05
596.5053	1202	15.513	0.919	9.864	0.894	1.573	2.78E-05
438.3354	1204	3.474	0.150	2.674	0.147	1.299	3.05E-05
597.5062	1204	64.543	4.659	39.473	3.816	1.636	3.06E-05
564.4396	1204	3.613	0.169	2.625	0.183	1.377	3.32E-05
596.5012	1204	181.033	13.876	108.540	10.921	1.668	3.44E-05
378.9906	1204	4.143	0.107	3.556	0.099	1.165	4.07E-05
492.3832	1204	9.413	0.388	7.222	0.415	1.303	4.33E-05
618.4834	1201	6.333	0.434	3.880	0.406	1.632	4.41E-05

570.4903	1204	4.537	0.344	2.752	0.258	1.649	4.42E-05
597.5068	1202	5.874	0.358	3.833	0.321	1.532	4.66E-05
188.0143	1102	4.309	0.492	2.424	0.354	1.778	0.0001
253.8165	1101	13.068	0.529	11.088	0.490	1.179	0.0001
462.3348	1204	3.673	0.137	2.982	0.148	1.232	0.0001
464.3524	1204	9.234	0.360	7.293	0.403	1.266	0.0001
478.4044	1204	3.745	0.161	2.814	0.173	1.331	0.0001
539.4274	1204	9.897	0.627	6.743	0.574	1.468	0.0001
551.4646	1204	4.024	0.184	2.946	0.193	1.366	0.0001
563.5013	1204	4.904	0.196	3.762	0.208	1.303	0.0001
576.4757	1204	45.791	2.161	33.088	2.440	1.384	0.0001
577.4795	1204	16.958	0.784	12.368	0.894	1.371	0.0001
594.4848	1204	116.663	6.054	80.027	6.363	1.458	0.0001
595.4883	1204	46.584	2.416	32.020	2.556	1.455	0.0001
462.3716	1204	3.016	0.089	2.549	0.120	1.183	0.0002
534.3912	1204	4.631	0.213	3.585	0.213	1.292	0.0002
546.3413	1204	3.789	0.217	2.844	0.150	1.332	0.0002
576.4765	1202	4.474	0.255	3.174	0.267	1.410	0.0002
594.4875	1202	10.080	0.505	7.204	0.579	1.396	0.0002
612.4994	1204	6.141	0.393	4.493	0.288	1.367	0.0002
616.4675	1201	4.908	0.294	3.410	0.291	1.439	0.0003
255.8135	1101	16.977	0.699	14.569	0.623	1.165	0.0005
384.3389	1203	69.859	1.997	62.624	1.789	1.116	0.0005
595.4928	1202	4.276	0.204	3.142	0.244	1.381	0.0006
366.3284	1203	26.455	0.868	23.895	0.747	1.107	0.0007
519.3998	1204	4.316	0.197	3.389	0.235	1.273	0.0007
572.4455	1204	4.831	0.251	3.693	0.223	1.308	0.0007
592.4717	1204	38.481	2.350	27.621	2.399	1.393	0.0007
769.5638	1204	124.642	6.236	106.685	4.871	1.168	0.0007
518.3969	1204	12.360	0.553	9.780	0.625	1.264	0.0008
593.4736	1204	16.057	1.005	11.506	1.011	1.396	0.0008
763.5153	1204	21.461	2.719	12.389	1.529	1.732	0.0008
770.569	1204	56.365	2.587	49.027	2.067	1.150	0.0008
591.4614	1204	3.999	0.218	3.096	0.219	1.292	0.001
476.3869	1204	4.746	0.185	3.852	0.239	1.232	0.0011
502.4054	1204	6.672	0.292	5.301	0.330	1.258	0.0011
716.4323	1204	14.005	0.923	10.954	0.494	1.279	0.0011
381.311	1203	68.023	2.436	59.906	2.515	1.136	0.0014
385.3428	1203	22.370	0.617	20.336	0.543	1.100	0.0017
446.341	1204	13.375	0.593	10.949	0.621	1.222	0.0017
271.8051	1102	6.853	0.380	5.799	0.449	1.182	0.0018
1018.9399	1203	14.299	1.077	10.464	0.789	1.367	0.0018
211.8495	1102	5.945	0.287	5.125	0.303	1.160	0.0019
367.3325	1203	11.717	0.255	10.717	0.271	1.093	0.0019
1254.1311	1203	6.597	0.466	4.773	0.406	1.382	0.002
713.5097	1204	12.498	0.743	10.692	0.559	1.169	0.0021
765.5316	1204	32.360	2.588	24.079	1.770	1.344	0.0022
1016.9279	1203	23.892	2.182	17.413	1.667	1.372	0.0023
257.8106	1101	8.432	0.356	7.356	0.334	1.146	0.0026
532.4503	1204	4.833	0.212	3.972	0.232	1.217	0.0027
546.4298	1204	4.437	0.233	3.508	0.247	1.265	0.0027
1253.1236	1203	8.820	0.595	6.551	0.561	1.346	0.0028
345.8738	1101	5.102	0.275	5.977	0.235	0.854	0.0033
855.6798	1204	5.567	0.503	3.533	0.454	1.576	0.0033
474.3731	1204	4.936	0.232	4.069	0.229	1.213	0.0034
886.5582	1102	7.913	0.483	9.581	0.608	0.826	0.0035
646.5702	1203	7.748	0.465	6.225	0.383	1.245	0.0038
488.2996	1204	5.809	0.309	4.794	0.190	1.212	0.0041
793.5663	1204	55.191	2.714	46.988	2.378	1.175	0.0043
380.3079	1203	231.313	8.275	207.315	8.508	1.116	0.0046
792.555	1204	29.746	2.033	24.379	1.428	1.220	0.0047
202.0453	1101	27.711	2.242	21.457	0.892	1.291	0.0048
448.3194	1204	3.729	0.121	3.241	0.133	1.151	0.0048
791.5488	1204	57.558	4.280	46.668	2.924	1.233	0.0048
490.3676	1204	6.262	0.288	5.209	0.308	1.202	0.0049

382.1084	1101	3.227	0.518	1.965	0.158	1.642	0.0051
468.3577	1201	3.557	0.279	2.792	0.208	1.274	0.0053
504.4188	1204	6.954	0.301	5.648	0.349	1.231	0.0053
376.2759	1203	18.398	0.666	17.053	0.629	1.079	0.0057
378.2921	1203	40.341	1.816	37.340	1.732	1.080	0.0057
634.3951	1204	8.858	0.797	6.233	0.437	1.421	0.0057
712.5074	1204	28.299	1.742	24.760	1.335	1.143	0.0061
702.4175	1204	10.640	0.718	8.121	0.375	1.310	0.0063
781.6001	1204	9.460	0.461	8.087	0.423	1.170	0.0063
745.5643	1204	120.519	6.555	105.213	5.381	1.145	0.0068
741.5302	1204	30.672	2.256	25.435	1.820	1.206	0.007
832.6022	1102	18.856	1.613	22.020	1.676	0.856	0.0071
306.2568	1204	9.978	0.463	8.702	0.360	1.147	0.0072
736.5031	1204	14.415	0.870	12.305	0.672	1.171	0.0075
556.4497	1204	5.617	0.316	5.615	0.346	1.001	0.0078
460.2681	1204	10.495	0.519	8.750	0.386	1.200	0.0078
610.3691	1201	12.407	0.915	9.705	0.850	1.278	0.0081
530.4379	1204	5.260	0.269	4.259	0.279	1.235	0.0085
559.4688	1204	4.738	0.324	4.029	0.291	1.176	0.0085
575.4628	1204	13.125	0.926	9.966	0.854	1.317	0.0085
766.5372	1204	15.044	1.269	12.085	0.826	1.247	0.009
746.5701	1204	52.385	2.683	46.291	2.263	1.132	0.0092
364.3123	1203	5.078	0.151	4.629	0.136	1.097	0.0094
447.3433	1204	4.155	0.228	3.361	0.217	1.236	0.0096
574.4594	1204	33.564	2.435	25.417	2.231	1.321	0.0096
432.3252	1204	6.136	0.165	5.573	0.150	1.101	0.0098
831.5992	1102	43.218	3.987	50.976	4.257	0.848	0.0098
708.4632	1201	2.976	0.180	2.474	0.195	1.203	0.0101
739.5146	1204	22.411	2.380	17.043	1.472	1.315	0.0103
311.7754	1101	5.487	0.325	4.484	0.299	1.224	0.0107
611.3724	1201	4.409	0.327	3.525	0.273	1.251	0.0108
312.231	1204	5.225	0.185	4.654	0.155	1.123	0.0112
794.6718	1204	28.572	1.318	25.149	1.154	1.136	0.0114
446.2525	1204	5.762	0.262	4.871	0.192	1.183	0.0122
737.5045	1204	6.771	0.410	5.827	0.333	1.162	0.0122
558.4649	1204	6.289	0.344	5.208	0.244	1.208	0.0123
243.0719	1101	33.095	2.902	26.746	1.469	1.237	0.0125
296.2357	1204	10.315	0.333	9.220	0.275	1.119	0.0128
218.0192	1101	8.365	0.788	6.159	0.509	1.358	0.0134
574.4635	1202	3.433	0.250	2.672	0.242	1.285	0.0134
743.5461	1204	365.902	20.760	321.870	20.695	1.137	0.0135
379.2967	1203	10.549	0.454	9.864	0.447	1.069	0.0143
273.8743	1101	9.119	0.350	8.331	0.298	1.095	0.0146
747.5761	1204	16.251	0.779	14.555	0.671	1.117	0.0149
283.8453	1101	8.079	0.306	7.317	0.276	1.104	0.0154
474.2846	1204	9.713	0.515	7.933	0.332	1.224	0.0158
290.1737	1204	3.105	0.157	2.605	0.118	1.192	0.0159
377.2801	1203	6.240	0.228	5.860	0.228	1.065	0.0162
495.3322	1201	3.653	0.262	4.284	0.287	0.853	0.0163
730.4535	1204	27.705	1.725	22.138	0.909	1.251	0.0165
244.0559	1101	9.936	0.826	8.100	0.241	1.227	0.0169
267.811	1102	6.583	0.399	5.420	0.365	1.215	0.0169
775.5514	1204	26.705	1.913	22.066	1.818	1.210	0.0169
833.7541	1203	5.373	0.673	8.021	1.021	0.670	0.0169
557.4527	1204	10.976	0.676	8.715	0.561	1.259	0.0175
744.5516	1204	151.209	7.754	135.492	7.998	1.116	0.0176
734.488	1204	15.280	1.078	12.822	0.881	1.192	0.018
551.4976	1203	38.060	4.291	54.579	6.668	0.697	0.0182
689.5083	1204	11.213	0.726	9.921	0.544	1.130	0.0183
314.2481	1204	5.729	0.255	5.145	0.220	1.114	0.0189
743.5475	1203	17.087	1.399	13.178	0.986	1.297	0.019
205.8867	1101	8.505	0.247	7.734	0.228	1.100	0.0192
260.004	1101	5.043	0.421	3.917	0.255	1.288	0.0197
209.8525	1102	4.762	0.278	3.920	0.272	1.215	0.0198
428.3653	1201	5.297	0.338	4.617	0.315	1.147	0.0203

544.3636	1204	4.414	0.247	3.576	0.180	1.234	0.0207
1017.9316	1203	21.209	2.096	15.726	1.093	1.349	0.0209
855.6009	1102	21.610	1.823	26.526	2.217	0.815	0.0215
552.5008	1203	7.495	0.783	10.398	1.173	0.721	0.022
282.2572	1204	168.937	10.620	141.140	7.379	1.197	0.0227
333.9639	1102	4.018	0.340	3.076	0.264	1.306	0.0229
744.55	1203	7.721	0.649	5.952	0.454	1.297	0.0242
502.3165	1204	33.700	1.944	27.434	1.147	1.228	0.0243
693.631	1204	16.868	1.766	12.232	1.297	1.379	0.0248
550.4954	1203	105.271	12.069	147.523	18.043	0.714	0.0249
503.3194	1204	9.409	0.559	7.602	0.307	1.238	0.025
318.1421	1201	8.664	0.503	7.285	0.537	1.189	0.0252
758.4785	1204	68.915	3.462	58.805	2.587	1.172	0.0256
524.296	1201	4.652	0.284	3.922	0.281	1.186	0.0258
269.8081	1102	10.759	0.685	8.876	0.682	1.212	0.026
268.1287	1201	46.799	2.580	40.474	2.525	1.156	0.0265
277.8861	1101	11.485	0.501	12.370	0.526	0.928	0.0268
688.5048	1204	26.803	1.813	23.685	1.393	1.132	0.0269
304.2398	1202	9.267	0.531	8.084	0.527	1.146	0.0272
694.6323	1204	10.428	0.979	7.889	0.727	1.322	0.0272
632.5038	1204	3.332	0.305	2.653	0.213	1.256	0.0277
283.2602	1204	32.186	1.992	27.097	1.398	1.188	0.0278
648.5861	1203	27.564	1.377	24.101	1.204	1.144	0.0279
374.2613	1203	7.381	0.254	6.920	0.250	1.067	0.028
781.5619	1204	10.474	0.737	8.753	0.595	1.197	0.0281
558.3761	1204	4.507	0.282	3.634	0.242	1.240	0.0295
274.1778	1202	44.020	1.999	37.866	1.523	1.163	0.0305
275.1811	1202	6.901	0.302	5.925	0.241	1.165	0.0306
687.4916	1204	30.334	1.910	27.231	1.812	1.114	0.0307
558.4663	1202	40.924	2.515	32.452	2.372	1.261	0.0308
766.5051	1201	3.411	0.211	2.896	0.207	1.178	0.0309
207.8836	1101	6.751	0.210	6.223	0.230	1.085	0.0316
789.5658	1204	10.254	0.434	9.339	0.370	1.098	0.0321
886.4879	1204	72.574	4.856	64.772	4.573	1.120	0.0329
649.5895	1203	13.365	0.671	11.763	0.593	1.138	0.0331
856.6045	1102	10.911	0.876	13.159	1.067	0.829	0.0341
542.3447	1102	6.632	0.354	7.305	0.314	0.908	0.0347
280.2413	1204	143.389	6.365	127.689	5.733	1.123	0.0351
715.5228	1204	22.831	1.767	19.876	1.283	1.149	0.0353
767.5473	1204	204.597	14.465	176.048	11.501	1.162	0.0355
722.479	1201	3.450	0.199	3.002	0.200	1.149	0.0356
265.8424	1101	7.574	0.284	6.878	0.249	1.102	0.0365
296.1801	1201	95.136	6.065	83.281	5.046	1.142	0.0383
328.2393	1202	4.172	0.377	3.300	0.355	1.264	0.0394
249.9677	1102	6.420	0.393	5.792	0.353	1.109	0.0397
768.5525	1204	93.141	5.972	81.499	4.915	1.143	0.0407
281.2447	1204	28.015	1.237	25.049	1.101	1.118	0.0408
580.478	1203	15.556	0.914	12.279	0.877	1.267	0.0412
1251.1042	1203	7.618	0.586	6.046	0.642	1.280	0.0418
333.8302	1101	7.201	0.306	6.569	0.258	1.096	0.0441
742.5366	1204	15.072	1.088	13.383	0.905	1.126	0.0443
256.24	1202	4.238	0.364	3.529	0.330	1.201	0.0449
246.1467	1202	16.897	0.885	14.462	0.580	1.168	0.045
392.294	1204	5.413	0.407	4.474	0.399	1.210	0.0485
552.3273	1201	6.177	0.323	5.319	0.366	1.161	0.0488

[00180] Table 2: Accurate mass features differing between clinically diagnosed PP- MULTIPLE SCLEROSIS patients and controls ($p < 0.05$).

188.0143	1102	5.494	0.988	2.424	0.354	2.267	7.38E-08
244.0559	1101	11.378	0.934	8.100	0.241	1.405	1.73E-06
202.0453	1101	29.842	2.721	21.457	0.892	1.391	2.59E-05
218.0371	1102	7.872	0.586	6.067	0.256	1.297	3.72E-05
218.04	1102	23.392	1.656	18.040	0.754	1.297	4.97E-05
243.0719	1101	33.520	3.834	26.746	1.469	1.253	0.0003
273.9985	1102	4.594	0.455	3.181	0.267	1.444	0.0003
218.0192	1101	9.506	1.366	6.159	0.509	1.543	0.0004
226.0688	1102	12.009	1.155	10.686	0.540	1.124	0.0004
290.1737	1204	3.446	0.253	2.605	0.118	1.323	0.0006
278.1494	1201	11.298	1.626	6.249	0.625	1.808	0.0008
260.004	1101	5.873	0.706	3.917	0.255	1.499	0.0014
326.1708	1201	6.913	0.832	4.477	0.371	1.544	0.0017
613.3404	1202	6.307	0.559	5.304	0.302	1.189	0.0045
827.5445	1101	5.066	0.511	4.010	0.235	1.263	0.005
546.3413	1204	3.704	0.305	2.844	0.150	1.303	0.0052
246.1467	1202	17.624	1.405	14.462	0.560	1.219	0.0054
269.132	1201	7.977	0.617	5.987	0.374	1.332	0.006
634.3951	1204	8.465	0.758	6.233	0.437	1.358	0.007
506.4338	1204	3.033	0.356	2.406	0.161	1.260	0.0082
268.1287	1201	53.182	4.401	40.474	2.525	1.314	0.0085
273.8743	1101	8.819	0.350	8.331	0.298	1.059	0.01
1228.1101	1203	12.407	1.879	9.994	0.940	1.241	0.0104
257.8106	1101	8.396	0.392	7.358	0.334	1.141	0.0119
474.2846	1204	9.786	0.712	7.933	0.332	1.234	0.0133
432.2385	1204	3.775	0.254	3.251	0.139	1.161	0.0136
623.5003	1203	8.814	0.720	6.953	0.309	1.268	0.0139
333.9539	1102	4.124	0.515	3.076	0.264	1.341	0.0148
611.3724	1201	5.093	0.575	3.525	0.273	1.445	0.0149
828.5479	1101	2.853	0.300	2.322	0.123	1.229	0.015
282.1444	1201	10.169	0.739	7.891	0.482	1.289	0.0162
622.4973	1203	18.911	1.515	14.995	0.686	1.261	0.0174
296.1601	1201	105.884	8.348	83.281	5.046	1.271	0.0175
488.2996	1204	5.767	0.402	4.794	0.190	1.203	0.021
203.1157	1101	5.032	0.697	4.122	0.337	1.221	0.0222
283.8453	1101	7.812	0.304	7.317	0.276	1.068	0.0228
246.1472	1204	12.531	0.803	10.559	0.477	1.187	0.0248
253.8165	1101	12.096	0.565	11.088	0.490	1.091	0.0257
792.555	1204	27.305	2.818	24.379	1.428	1.120	0.0282
161.1051	1101	4.381	0.482	3.749	0.267	1.163	0.0289
793.4936	1204	38.528	3.782	33.573	1.651	1.148	0.0292
791.5488	1204	52.376	6.009	46.668	2.924	1.122	0.03
517.3141	1201	2.684	0.302	2.994	0.162	0.897	0.0315
610.3691	1201	14.143	1.697	9.705	0.850	1.457	0.0322
310.1758	1201	7.318	0.503	5.609	0.387	1.305	0.0323
217.9124	1101	11.644	0.463	11.018	0.336	1.057	0.0347
446.2525	1204	5.623	0.324	4.871	0.192	1.154	0.0356
328.2393	1202	5.686	1.013	3.300	0.355	1.693	0.0365
318.1421	1201	9.677	0.763	7.285	0.537	1.328	0.0383
274.1778	1202	44.628	3.137	37.866	1.523	1.179	0.0389
297.1634	1201	16.418	1.267	13.359	0.793	1.229	0.0391
831.5992	1102	33.924	6.553	50.976	4.257	0.666	0.0391
275.8713	1101	5.707	0.219	5.459	0.187	1.046	0.0393
819.5831	1204	15.703	1.548	14.148	0.734	1.110	0.0423
460.2681	1204	10.050	0.633	8.750	0.386	1.149	0.0429
506.2853	1201	5.202	0.690	3.261	0.273	1.595	0.0431
832.6022	1102	15.107	2.641	22.020	1.676	0.686	0.0439

899.5871	1102	7.493	1.294	10.990	0.852	0.682	0.0462
328.2415	1204	4.683	0.491	3.818	0.259	1.227	0.0465
503.3194	1204	8.986	0.646	7.602	0.307	1.182	0.0475

[00181] Table 3: Accurate mass features differing between clinically diagnosed SP- MULTIPLE SCLEROSIS patients and controls ($p < 0.05$).

428.3653	1201	9.177	0.839	4.617	0.315	1.988	2.84E-05
590.4964	1204	3.690	0.441	5.275	0.443	0.700	0.0003
597.5068	1202	2.135	0.223	3.833	0.321	0.557	0.0003
596.5053	1202	5.345	0.487	9.864	0.894	0.542	0.0005
493.385	1204	1.599	0.143	2.549	0.165	0.627	0.001
594.4875	1202	4.484	0.610	7.204	0.579	0.622	0.0014
763.5153	1204	19.590	3.550	12.389	1.529	1.581	0.0016
764.5196	1204	8.495	1.885	4.784	0.817	1.776	0.0017
194.0803	1203	3.200	0.530	10.851	1.415	0.295	0.0019
872.6715	1204	4.860	0.498	2.578	0.305	1.885	0.0019
597.5062	1204	18.785	1.588	39.473	3.816	0.476	0.0022
618.4675	1201	2.340	0.315	3.410	0.291	0.686	0.0023
495.4018	1204	2.381	0.174	3.441	0.203	0.692	0.0025
595.4928	1202	2.005	0.277	3.142	0.244	0.638	0.0025
523.4337	1204	2.645	0.225	3.934	0.237	0.672	0.0026
598.5107	1204	4.542	0.376	8.823	0.842	0.515	0.0028
596.5012	1204	50.842	4.367	108.540	10.921	0.468	0.003
618.4834	1201	2.583	0.254	3.880	0.406	0.666	0.0032
610.5204	1204	7.767	0.828	12.405	1.447	0.626	0.0033
539.4274	1204	4.835	0.723	6.743	0.574	0.717	0.0037
791.5488	1204	52.523	5.677	46.668	2.924	1.125	0.0038
577.4795	1204	7.187	0.718	12.368	0.894	0.581	0.0041
578.4923	1204	15.340	1.107	26.563	2.068	0.578	0.0042
821.5288	1204	17.973	1.220	15.743	0.732	1.142	0.0043
792.555	1204	27.350	2.825	24.379	1.428	1.122	0.0046
576.4757	1204	19.011	1.969	33.088	2.440	0.575	0.0047
490.3676	1204	3.541	0.400	5.209	0.308	0.680	0.0048
594.4848	1204	43.087	5.319	80.027	6.363	0.538	0.0048
579.4958	1204	6.320	0.458	10.661	0.805	0.593	0.0049
793.4936	1204	37.612	3.233	33.573	1.651	1.120	0.0049
595.4883	1204	17.987	2.221	32.020	2.556	0.561	0.0051
492.3832	1204	4.994	0.446	7.222	0.415	0.692	0.0054
851.5686	1102	6.306	0.461	9.813	0.818	0.643	0.0056
541.4422	1204	4.920	0.418	7.651	0.568	0.643	0.0068
466.3656	1204	7.059	0.607	10.077	0.552	0.700	0.0069
550.4602	1204	5.085	0.560	7.791	0.513	0.653	0.007
606.4872	1204	4.658	0.508	6.560	0.561	0.710	0.0072
806.5643	1201	27.948	1.837	18.717	0.860	1.493	0.0075
522.4313	1204	7.993	0.637	11.438	0.681	0.699	0.0076
551.4646	1204	1.920	0.203	2.946	0.193	0.652	0.0078
495.3321	1101	10.510	0.594	11.115	0.423	0.946	0.0081
440.3526	1204	2.306	0.202	3.273	0.180	0.705	0.0083
558.4663	1202	3.020	0.486	4.029	0.291	0.749	0.0084
467.3711	1204	2.462	0.223	3.431	0.188	0.718	0.009
519.3322	1101	5.123	0.498	6.008	0.350	0.853	0.009
520.4131	1204	8.374	0.748	12.295	0.696	0.681	0.009
548.4438	1204	3.745	0.372	5.581	0.347	0.671	0.0091
805.5609	1201	55.027	3.212	36.921	1.704	1.490	0.0093
468.3577	1201	5.153	0.520	2.792	0.208	1.846	0.0094
538.4257	1204	14.296	1.804	21.271	1.294	0.672	0.0094
464.3524	1204	5.162	0.501	7.293	0.403	0.708	0.0097
542.3447	1102	4.282	0.293	7.305	0.314	0.586	0.0098
446.341	1204	7.886	0.859	10.949	0.621	0.720	0.01
513.4116	1204	2.379	0.211	3.821	0.296	0.623	0.0107
540.4387	1204	14.521	1.258	22.346	1.749	0.650	0.0108
202.0453	1101	24.111	1.347	21.457	0.892	1.124	0.0109
328.2415	1204	4.988	0.862	3.818	0.259	1.308	0.0109
819.5831	1204	15.250	1.244	14.148	0.734	1.078	0.0112

569.3687	1102	4.492	0.482	7.792	0.345	0.576	0.0117
568.4723	1204	5.475	0.537	8.786	0.743	0.623	0.0123
518.3969	1204	6.830	0.808	9.780	0.625	0.698	0.0125
828.5477	1201	7.801	0.619	5.137	0.244	1.519	0.0126
494.3968	1204	8.860	0.516	10.018	0.645	0.685	0.0129
576.4765	1202	2.269	0.250	3.174	0.267	0.715	0.0129
249.9677	1102	6.069	0.455	5.792	0.353	0.876	0.0143
468.384	1204	9.354	0.797	12.977	0.754	0.721	0.0147
382.1084	1101	2.573	0.251	1.965	0.158	1.309	0.015
566.454	1204	3.969	0.446	5.967	0.419	0.665	0.0154
484.3788	1204	4.070	0.305	5.832	0.359	0.698	0.0157
512.4079	1204	6.864	0.551	10.772	0.864	0.637	0.0157
610.482	1204	4.242	0.241	6.136	0.359	0.691	0.0159
537.4142	1204	2.263	0.275	3.299	0.207	0.686	0.0167
720.4696	1204	6.131	0.514	4.968	0.257	1.234	0.0167
580.5089	1204	5.730	0.402	8.725	0.641	0.657	0.017
855.6798	1204	5.942	0.996	3.533	0.454	1.682	0.017
448.3194	1204	4.225	0.221	3.241	0.133	1.304	0.0177
508.3782	1204	2.805	0.263	3.819	0.212	0.735	0.0178
438.3354	1204	2.071	0.176	2.674	0.147	0.774	0.0181
574.4594	1204	15.435	2.101	25.417	2.231	0.607	0.0187
613.3404	1202	5.155	0.432	5.304	0.302	0.972	0.0189
482.3604	1204	2.926	0.227	3.772	0.196	0.776	0.019
827.5446	1201	15.396	1.216	9.932	0.505	1.550	0.0191
564.4396	1204	1.877	0.175	2.625	0.183	0.715	0.0192
448.3562	1204	6.275	0.562	8.119	0.373	0.773	0.0194
541.3415	1102	15.111	1.031	25.470	1.129	0.593	0.0203
622.4973	1203	20.247	1.598	14.995	0.688	1.350	0.0219
311.7754	1101	4.688	0.288	4.484	0.299	1.045	0.022
385.3428	1203	21.477	0.879	20.338	0.543	1.056	0.022
574.4635	1202	1.995	0.279	2.672	0.242	0.747	0.0231
566.3431	1102	4.439	0.519	7.334	0.380	0.605	0.024
521.3478	1101	3.649	0.287	4.014	0.207	0.909	0.0244
328.2393	1202	6.800	1.193	3.300	0.355	2.061	0.0248
480.3473	1204	2.492	0.216	3.030	0.142	0.823	0.0249
253.8165	1101	10.790	0.555	11.088	0.490	0.973	0.025
510.3937	1204	3.222	0.237	4.445	0.265	0.725	0.0251
1228.1101	1203	8.624	1.735	9.984	0.940	0.863	0.0253
565.3391	1102	14.619	1.846	24.344	1.318	0.601	0.0256
593.4736	1204	7.260	0.952	11.506	1.011	0.631	0.0256
519.3998	1204	2.331	0.340	3.389	0.235	0.688	0.0265
886.5582	1102	5.104	0.302	9.581	0.608	0.533	0.0267
694.6323	1204	12.155	1.212	7.889	0.727	1.541	0.0283
820.589	1204	9.217	0.740	8.675	0.401	1.062	0.0291
384.3399	1203	65.198	3.300	62.624	1.789	1.041	0.0303
546.4298	1204	2.524	0.327	3.508	0.247	0.719	0.0305
766.5372	1204	13.388	1.664	12.065	0.926	1.110	0.0308
469.3863	1204	2.523	0.248	3.536	0.206	0.714	0.0312
312.231	1204	5.656	0.280	4.654	0.155	1.215	0.0313
592.4717	1204	17.304	2.411	27.621	2.399	0.626	0.0313
541.3141	1201	2.734	0.234	2.835	0.227	0.964	0.0315
474.3731	1204	3.158	0.377	4.069	0.229	0.776	0.0326
575.4628	1204	6.499	0.860	9.966	0.854	0.652	0.0329
723.6395	1204	9.645	0.668	6.930	0.437	1.392	0.0333
244.0559	1101	9.067	0.491	8.100	0.241	1.119	0.0335
246.1468	1201	6.273	0.755	4.636	0.356	1.353	0.0339
765.5316	1204	26.162	3.135	24.079	1.770	1.087	0.0342
521.4188	1204	3.362	0.267	4.437	0.265	0.758	0.0343
534.3912	1204	2.755	0.256	3.585	0.213	0.769	0.0346
569.4769	1204	2.292	0.255	3.508	0.309	0.653	0.0349
523.3637	1101	2.894	0.313	3.422	0.124	0.846	0.0357
243.0719	1101	24.190	1.697	26.746	1.469	0.904	0.0366
255.8135	1101	14.393	0.624	14.569	0.623	0.988	0.0374
536.41	1204	6.261	0.609	7.896	0.469	0.793	0.0387
541.3141	1101	5.759	0.624	6.840	0.414	0.842	0.0414

768.5468	1102	1.679	0.306	3.219	0.321	0.522	0.0415
590.4585	1204	5.278	0.564	2.834	0.368	1.862	0.0421
684.6037	1203	4.261	0.453	3.097	0.392	1.376	0.0431
852.5724	1102	3.428	0.268	5.434	0.416	0.631	0.0436
552.4784	1204	16.345	1.998	22.583	1.626	0.724	0.0454
560.4821	1204	5.721	0.595	8.041	0.547	0.712	0.049

[00182] Table 4: Accurate mass features differing between clinically diagnosed SP- MULTIPLE SCLEROSIS patients and RR- MULTIPLE SCLEROSIS patients ($p < 0.05$).

452.3868	1204	2.163	0.154	3.906	0.137	0.554	4.29E-11
580.5089	1204	5.730	0.402	13.528	0.711	0.424	1.98E-10
578.4923	1204	15.340	1.107	40.496	2.263	0.379	2.21E-10
493.385	1204	1.599	0.143	3.412	0.131	0.469	2.85E-10
523.4337	1204	2.645	0.225	5.249	0.219	0.504	3.26E-10
522.4313	1204	7.993	0.637	15.695	0.623	0.509	3.42E-10
512.4079	1204	6.864	0.551	15.906	0.838	0.432	4.59E-10
579.4958	1204	6.320	0.458	15.978	0.888	0.396	5.47E-10
494.3968	1204	6.860	0.516	14.131	0.626	0.485	6.51E-10
495.4018	1204	2.381	0.174	4.575	0.202	0.521	8.69E-10
484.3788	1204	4.070	0.305	8.180	0.389	0.497	1.00E-09
513.4116	1204	2.379	0.211	5.401	0.281	0.440	1.14E-09
596.5053	1202	5.345	0.487	15.297	0.951	0.349	1.54E-09
581.5126	1204	2.422	0.183	5.487	0.296	0.441	1.69E-09
466.3656	1204	7.059	0.607	13.494	0.566	0.523	1.72E-09
550.4602	1204	5.085	0.560	10.852	0.521	0.469	1.85E-09
510.3937	1204	3.222	0.237	6.276	0.278	0.513	2.63E-09
468.384	1204	9.354	0.797	18.294	0.792	0.511	2.96E-09
469.3863	1204	2.523	0.248	5.002	0.209	0.504	4.28E-09
597.5068	1202	2.135	0.223	5.788	0.369	0.369	4.59E-09
440.3526	1204	2.306	0.202	4.382	0.192	0.526	5.93E-09
568.4723	1204	5.475	0.537	13.607	0.761	0.402	6.09E-09
618.4834	1201	2.583	0.254	6.236	0.449	0.414	9.48E-09
577.4795	1204	7.187	0.718	16.759	0.829	0.429	1.26E-08
524.4448	1204	2.878	0.241	5.382	0.265	0.535	1.34E-08
450.3729	1204	5.537	0.435	9.090	0.331	0.609	1.36E-08
594.4875	1202	4.484	0.610	9.946	0.528	0.451	1.40E-08
551.4646	1204	1.920	0.203	3.993	0.192	0.481	1.71E-08
568.454	1204	3.969	0.446	8.691	0.393	0.457	1.85E-08
552.4784	1204	3.804	0.378	8.282	0.478	0.459	1.89E-08
598.5107	1204	4.542	0.376	14.230	1.030	0.319	1.97E-08
576.4757	1204	19.011	1.969	45.354	2.261	0.419	2.02E-08
548.4438	1204	3.745	0.372	7.407	0.309	0.506	2.90E-08
448.3582	1204	6.275	0.562	10.433	0.390	0.601	3.15E-08
569.4769	1204	2.292	0.255	5.444	0.298	0.421	3.44E-08
520.4131	1204	8.374	0.748	15.943	0.634	0.525	3.59E-08
467.3711	1204	2.462	0.223	4.437	0.193	0.555	3.62E-08
597.5062	1204	18.785	1.588	63.633	4.798	0.295	4.01E-08
508.3782	1204	2.805	0.263	5.241	0.265	0.535	4.46E-08
564.4396	1204	1.877	0.175	3.584	0.173	0.524	4.64E-08
521.4188	1204	3.362	0.267	5.796	0.222	0.580	4.75E-08
541.4422	1204	4.920	0.418	12.302	0.734	0.400	5.43E-08
496.4157	1204	5.182	0.354	10.685	0.591	0.485	5.52E-08
492.3832	1204	4.994	0.446	9.311	0.399	0.536	5.78E-08
594.4848	1204	43.087	5.319	115.274	6.359	0.374	5.81E-08
537.4142	1204	2.263	0.275	4.317	0.185	0.524	5.85E-08
536.41	1204	6.261	0.609	10.884	0.483	0.575	6.25E-08
540.4387	1204	14.521	1.256	36.081	2.150	0.402	6.32E-08
595.4883	1204	17.967	2.221	46.028	2.536	0.390	6.43E-08
595.4928	1202	2.005	0.277	4.218	0.215	0.475	6.44E-08
616.4675	1201	2.340	0.315	4.851	0.305	0.482	7.05E-08
482.3604	1204	2.926	0.227	5.143	0.245	0.569	7.54E-08
596.5012	1204	50.842	4.367	178.485	14.249	0.285	7.94E-08
610.482	1204	4.242	0.241	8.903	0.494	0.476	1.17E-07
464.3524	1204	5.162	0.501	9.142	0.365	0.565	1.72E-07
480.3473	1204	2.492	0.216	3.911	0.157	0.637	1.78E-07
438.3354	1204	2.071	0.176	3.426	0.150	0.604	2.41E-07

539.4274	1204	4.835	0.723	9.733	0.636	0.497	2.43E-07
576.4765	1202	2.269	0.250	4.415	0.265	0.514	3.05E-07
538.4257	1204	14.296	1.804	29.559	1.464	0.484	3.34E-07
562.4989	1204	7.595	0.589	12.437	0.512	0.611	3.83E-07
590.4585	1204	5.721	0.595	11.249	0.631	0.509	4.39E-07
563.5013	1204	2.892	0.258	4.830	0.203	0.599	5.17E-07
478.4044	1204	2.256	0.210	3.682	0.161	0.613	5.30E-07
518.3969	1204	6.830	0.808	12.262	0.569	0.557	1.07E-06
462.3716	1204	1.914	0.187	3.007	0.090	0.636	1.33E-06
446.341	1204	7.886	0.859	13.291	0.599	0.593	2.02E-06
476.3869	1204	3.097	0.261	4.685	0.186	0.661	2.08E-06
519.3998	1204	2.331	0.340	4.275	0.202	0.545	2.45E-06
593.4736	1204	7.260	0.952	15.915	1.032	0.456	2.57E-06
592.4717	1204	17.304	2.411	38.114	2.424	0.454	2.82E-06
570.4903	1204	1.853	0.317	4.480	0.352	0.414	3.57E-06
534.3912	1204	2.755	0.256	4.595	0.215	0.600	3.63E-06
534.4645	1204	2.436	0.202	3.904	0.180	0.624	4.64E-06
532.4503	1204	3.116	0.284	4.752	0.228	0.656	6.03E-06
490.3676	1204	3.541	0.400	6.198	0.296	0.571	7.12E-06
462.3346	1204	2.328	0.276	3.632	0.138	0.641	1.43E-05
502.4054	1204	4.125	0.456	6.631	0.298	0.622	1.47E-05
591.4614	1204	2.320	0.217	3.990	0.220	0.582	1.63E-05
574.4594	1204	15.435	2.101	32.679	2.206	0.472	2.00E-05
546.4298	1204	2.524	0.327	4.374	0.245	0.577	2.51E-05
504.4188	1204	4.556	0.423	6.780	0.309	0.672	3.45E-05
575.4628	1204	6.499	0.860	12.808	0.831	0.507	3.95E-05
572.4455	1204	2.656	0.324	4.775	0.263	0.556	4.64E-05
574.4635	1202	1.995	0.279	3.453	0.242	0.578	4.77E-05
327.0307	1204	6.113	0.237	7.511	0.224	0.814	0.0001
447.3433	1204	2.679	0.300	4.068	0.211	0.658	0.0001
474.3731	1204	3.158	0.377	4.916	0.233	0.642	0.0001
530.4379	1204	3.465	0.344	5.085	0.240	0.681	0.0001
558.4649	1204	22.470	2.918	39.813	2.242	0.564	0.0001
558.4663	1202	3.020	0.486	4.836	0.312	0.624	0.0001
559.4688	1204	8.680	1.100	15.071	0.830	0.576	0.0001
561.4863	1204	6.581	0.759	9.904	0.461	0.664	0.0001
560.4821	1204	16.345	1.998	24.793	1.209	0.659	0.0002
564.513	1204	2.848	0.261	4.456	0.242	0.639	0.0002
612.4994	1204	3.494	0.237	6.110	0.396	0.572	0.0002
532.1851	1204	1.780	0.243	1.088	0.050	1.636	0.0003
506.4338	1204	2.504	0.192	3.502	0.168	0.715	0.0004
610.5204	1204	7.767	0.828	13.914	1.696	0.558	0.0005
556.4497	1204	6.304	0.714	10.754	0.704	0.586	0.0013
590.4964	1204	4.097	0.371	3.690	0.441	1.110	0.0023
821.5288	1204	17.973	1.220	15.703	0.731	1.145	0.0026
557.4527	1204	2.788	0.293	4.409	0.298	0.632	0.004
606.4872	1204	4.658	0.508	6.362	0.467	0.732	0.0041
340.2407	1204	6.645	0.655	5.152	0.169	1.290	0.0042
851.5666	1102	6.306	0.461	9.698	1.102	0.650	0.0046
886.7896	1203	5.463	1.422	4.340	0.882	1.259	0.0075
378.9906	1204	3.602	0.183	4.124	0.106	0.873	0.0121
852.5724	1102	3.428	0.268	5.055	0.553	0.678	0.0135
194.0803	1203	3.200	0.530	9.406	1.860	0.340	0.0144
834.5963	1201	7.969	0.614	5.250	0.377	1.518	0.0154
264.9759	1204	6.124	0.234	6.909	0.130	0.886	0.0158
872.6715	1204	4.860	0.498	2.907	0.376	1.671	0.0173
477.3218	1201	6.130	0.392	3.865	0.260	1.586	0.0183
551.4991	1201	2.961	0.530	1.953	0.210	1.516	0.0199
833.5931	1201	14.601	1.100	9.580	0.695	1.524	0.0221
539.4286	1204	1.364	0.195	1.714	0.365	0.795	0.0222
428.3653	1201	9.177	0.839	5.455	0.349	1.682	0.024
634.3951	1204	5.491	0.664	8.401	0.739	0.654	0.0305
662.4267	1204	4.861	0.451	6.334	0.360	0.768	0.0305

274.1777	1203	1.377	0.177	2.279	0.299	0.604	0.0311
835.6094	1201	5.680	0.531	3.714	0.244	1.529	0.0314
780.5303	1204	7.194	0.535	9.693	0.351	0.742	0.033
793.4936	1204	37.612	3.233	35.440	1.943	1.061	0.0353
368.1656	1102	1.679	0.306	3.982	0.611	0.422	0.0354
646.5702	1203	6.529	0.770	7.797	0.471	0.837	0.0375
632.5038	1204	2.056	0.470	3.021	0.265	0.680	0.0379
729.5727	1204	6.518	0.512	9.419	0.497	0.692	0.0397
808.5643	1201	27.948	1.637	19.725	1.161	1.417	0.0444
786.51	1204	31.030	2.459	41.726	1.808	0.744	0.0451
805.5609	1201	55.027	3.212	38.808	2.258	1.418	0.0463
541.3141	1201	2.734	0.234	2.773	0.153	0.986	0.0464
856.6045	1102	7.004	0.557	11.014	0.891	0.636	0.0464
366.3284	1203	21.097	1.078	26.496	0.868	0.796	0.0474
501.3217	1201	5.313	0.419	3.495	0.214	1.520	0.0498

[00183] Table 5: Accurate mass features differing between clinically diagnosed RR- MULTIPLE SCLEROSIS patients transitioning to SP- MULTIPLE SCLEROSIS and RR- MULTIPLE SCLEROSIS patients ($p < 0.05$).

580.5089	1204	7.612	0.633	13.528	0.711	0.563	2.52E-15
452.3868	1204	2.747	0.171	3.906	0.137	0.703	6.69E-15
522.4313	1204	10.338	0.808	15.695	0.623	0.659	8.43E-15
578.4923	1204	22.148	2.145	40.496	2.263	0.547	9.92E-15
450.3729	1204	6.629	0.358	9.090	0.331	0.729	1.01E-14
579.4958	1204	8.772	0.838	15.978	0.888	0.549	1.08E-14
581.6126	1204	3.136	0.254	5.487	0.296	0.571	1.80E-14
484.3788	1204	5.246	0.379	8.180	0.389	0.641	1.89E-14
466.3656	1204	9.116	0.702	13.494	0.566	0.676	2.23E-14
494.3968	1204	9.307	0.747	14.131	0.626	0.659	3.02E-14
550.4602	1204	6.431	0.504	10.852	0.521	0.593	3.96E-14
523.4337	1204	3.432	0.263	5.249	0.219	0.654	4.34E-14
610.3937	1204	4.078	0.340	6.276	0.278	0.660	5.04E-14
495.4018	1204	2.979	0.244	4.575	0.202	0.661	8.95E-14
512.4079	1204	9.451	0.830	15.906	0.838	0.594	2.38E-13
448.3562	1204	7.324	0.439	10.433	0.390	0.702	3.03E-13
552.4784	1204	4.480	0.334	8.282	0.478	0.541	4.73E-13
524.4448	1204	3.500	0.239	5.382	0.265	0.650	5.83E-13
536.41	1204	7.097	0.572	10.884	0.483	0.652	6.16E-13
568.4723	1204	7.427	0.835	13.607	0.761	0.546	9.45E-13
577.4795	1204	10.116	1.038	16.759	0.829	0.604	9.63E-13
576.4757	1204	26.727	2.897	45.354	2.261	0.589	1.30E-12
467.3711	1204	3.093	0.231	4.437	0.193	0.697	1.48E-12
468.384	1204	12.450	0.935	18.294	0.792	0.681	1.56E-12
493.385	1204	2.195	0.216	3.412	0.131	0.643	1.71E-12
513.4116	1204	3.325	0.286	5.401	0.281	0.616	1.91E-12
521.4188	1204	3.983	0.338	5.796	0.222	0.687	2.15E-12
596.5053	1202	7.401	1.039	15.297	0.951	0.484	2.15E-12
594.4875	1202	5.256	0.642	9.946	0.526	0.528	3.39E-12
537.4142	1204	2.863	0.247	4.317	0.185	0.663	4.81E-12
548.4438	1204	4.947	0.496	7.407	0.309	0.668	5.93E-12
469.3863	1204	3.607	0.253	5.002	0.209	0.721	7.00E-12
440.3526	1204	3.040	0.235	4.382	0.192	0.694	7.60E-12
551.4648	1204	2.450	0.191	3.993	0.192	0.814	8.25E-12
597.5068	1202	2.933	0.380	5.788	0.369	0.507	9.90E-12
569.4769	1204	2.998	0.330	5.444	0.298	0.551	1.02E-11
520.4131	1204	11.310	0.945	15.943	0.634	0.709	1.03E-11
566.454	1204	5.597	0.588	8.691	0.393	0.644	1.22E-11
492.3832	1204	6.447	0.555	9.311	0.399	0.692	1.26E-11
598.5107	1204	7.070	0.979	14.230	1.030	0.497	1.80E-11
595.4883	1204	25.969	3.279	46.028	2.536	0.564	1.93E-11
595.4928	1202	2.356	0.274	4.218	0.215	0.558	2.23E-11
594.4848	1204	65.843	8.436	115.274	6.358	0.571	2.35E-11
591.4614	1204	2.460	0.213	3.990	0.220	0.617	2.90E-11
482.3604	1204	3.414	0.245	5.143	0.245	0.664	3.32E-11
576.4765	1202	2.471	0.241	4.415	0.265	0.560	3.37E-11
476.3869	1204	3.426	0.183	4.685	0.186	0.731	3.85E-11
496.4157	1204	6.968	0.586	10.685	0.591	0.652	4.23E-11
508.3782	1204	3.587	0.251	5.241	0.265	0.684	4.43E-11
464.3524	1204	6.563	0.496	9.142	0.365	0.718	5.30E-11
618.4834	1201	2.975	0.365	6.236	0.449	0.477	7.22E-11
590.4585	1204	6.938	0.556	11.249	0.631	0.617	7.43E-11
597.5062	1204	32.711	4.840	63.633	4.798	0.514	8.79E-11
438.3354	1204	2.443	0.185	3.426	0.150	0.713	1.03E-10
541.4422	1204	7.441	0.963	12.302	0.734	0.605	1.26E-10
596.5012	1204	89.261	13.368	178.485	14.249	0.500	1.62E-10

540.4387	1204	21.969	2.791	36.081	2.150	0.609	1.82E-10
564.4396	1204	2.322	0.202	3.584	0.173	0.648	2.11E-10
538.4257	1204	19.272	2.025	29.559	1.464	0.652	2.93E-10
592.4717	1204	21.832	2.186	38.114	2.424	0.573	3.20E-10
593.4736	1204	8.928	0.908	15.915	1.032	0.561	3.21E-10
539.4274	1204	5.471	0.803	9.733	0.636	0.562	4.48E-10
534.3912	1204	3.038	0.259	4.595	0.215	0.661	5.76E-10
518.3969	1204	8.400	0.712	12.262	0.568	0.685	6.53E-10
532.4503	1204	3.489	0.190	4.752	0.228	0.734	7.60E-10
610.482	1204	5.842	0.532	8.903	0.494	0.656	8.00E-10
616.4675	1201	2.688	0.263	4.851	0.305	0.554	8.68E-10
462.3346	1204	2.692	0.199	3.632	0.138	0.741	9.35E-10
480.3473	1204	2.852	0.191	3.911	0.157	0.729	1.26E-09
446.341	1204	9.664	0.647	13.291	0.599	0.727	2.74E-09
504.4188	1204	4.955	0.270	6.780	0.309	0.731	3.65E-09
478.4044	1204	2.730	0.154	3.682	0.161	0.742	5.27E-09
570.4903	1204	2.348	0.203	4.480	0.352	0.524	6.68E-09
560.4821	1204	17.670	1.056	24.793	1.209	0.713	1.21E-08
502.4054	1204	4.726	0.317	6.631	0.298	0.713	1.51E-08
561.4863	1204	7.248	0.403	9.904	0.461	0.732	1.82E-08
490.3676	1204	4.358	0.414	6.198	0.296	0.703	2.12E-08
574.4594	1204	19.303	2.070	32.679	2.206	0.591	2.38E-08
575.4628	1204	7.733	0.793	12.808	0.831	0.604	3.58E-08
546.4298	1204	3.032	0.284	4.374	0.245	0.693	4.07E-08
574.4635	1202	1.927	0.171	3.453	0.242	0.558	4.32E-08
519.3998	1204	3.134	0.264	4.275	0.202	0.733	6.07E-08
572.4455	1204	3.261	0.227	4.775	0.263	0.683	8.21E-08
506.4338	1204	2.578	0.156	3.502	0.168	0.736	1.86E-07
474.3731	1204	3.602	0.287	4.916	0.233	0.733	6.45E-07
558.4663	1202	2.842	0.237	4.836	0.312	0.588	8.38E-07
559.4688	1204	10.587	0.932	16.071	0.830	0.703	9.35E-07
447.3433	1204	2.967	0.213	4.068	0.211	0.729	1.19E-06
562.4989	1204	9.973	0.554	12.437	0.512	0.802	1.25E-06
558.4649	1204	27.632	2.484	39.813	2.242	0.694	1.66E-06
534.4645	1204	3.046	0.155	3.904	0.180	0.780	1.89E-06
556.4497	1204	7.612	0.401	10.976	0.668	0.694	2.99E-06
557.4627	1204	3.388	0.213	4.409	0.298	0.769	6.92E-06
563.5013	1204	3.903	0.201	4.830	0.203	0.808	1.74E-05
530.4379	1204	4.240	0.273	5.085	0.240	0.834	0.0001
590.4964	1204	4.097	0.371	4.955	0.456	0.825	0.0003
784.6228	1204	19.037	1.831	10.803	1.101	1.762	0.0004
612.4994	1204	4.706	0.316	6.110	0.396	0.770	0.0005
327.0307	1204	7.030	0.236	7.511	0.224	0.936	0.0012
462.3716	1204	2.793	0.143	3.007	0.090	0.929	0.0012
783.6174	1204	30.598	3.318	16.139	1.805	1.896	0.0014
816.5159	1204	8.856	0.270	8.243	0.279	1.074	0.0027
560.478	1203	6.582	0.404	8.188	0.375	0.804	0.0031
244.2189	1203	6.901	0.266	7.322	0.161	0.943	0.0034
333.9539	1102	3.572	0.276	3.748	0.334	0.953	0.0047
744.55	1203	6.217	0.552	7.484	0.600	0.831	0.0079
747.5121	1204	54.149	2.246	45.567	2.452	1.168	0.0096
564.513	1204	4.140	0.278	4.456	0.242	0.929	0.0099
779.5828	1204	20.561	1.093	23.925	1.363	0.859	0.01
832.5211	1204	6.584	0.407	5.016	0.421	1.313	0.0101
743.5475	1203	14.077	1.093	16.578	1.335	0.849	0.0102
260.2137	1203	5.339	0.154	5.622	0.119	0.950	0.0106
828.5477	1201	5.936	0.361	5.534	0.275	1.073	0.0132
246.2345	1203	7.833	0.456	8.467	0.280	0.925	0.0143
584.2641	1202	4.241	0.433	5.003	0.562	0.848	0.015
821.5288	1204	18.114	1.077	15.703	0.731	1.154	0.0151
216.1877	1203	7.189	0.385	7.496	0.290	0.959	0.0166
831.5758	1201	16.234	0.935	16.394	0.762	0.990	0.017
239.939	1102	4.791	0.275	4.879	0.341	0.982	0.0174
830.5634	1201	5.708	0.313	5.574	0.235	1.024	0.0198
726.5438	1204	5.481	0.552	7.236	0.574	0.757	0.0201

214.1721	1203	9.114	0.523	9.867	0.358	0.924	0.0206
823.5427	1204	9.742	0.414	8.824	0.338	1.104	0.0242
200.1566	1203	7.319	0.271	7.980	0.158	0.917	0.0251
610.5204	1204	12.409	1.532	13.914	1.696	0.892	0.0258
839.6019	1202	7.634	0.269	8.499	0.363	0.898	0.0264
277.8861	1101	13.450	0.660	11.133	0.491	1.208	0.0276
303.108	1202	33.776	3.719	37.528	3.650	0.900	0.0286
731.5464	1201	4.127	0.502	3.918	0.335	1.053	0.0292
181.9806	1102	5.275	0.323	5.215	0.335	1.011	0.0298
188.1567	1203	9.264	0.414	9.796	0.299	0.946	0.0302
834.5372	1204	16.761	0.942	13.473	0.712	1.244	0.0311
781.6001	1204	8.456	0.370	9.429	0.471	0.897	0.0316
835.5417	1204	9.435	0.518	7.667	0.387	1.231	0.0327
202.1721	1203	9.681	0.507	10.182	0.375	0.951	0.0335
345.8738	1101	6.088	0.293	5.198	0.260	1.171	0.0344
331.957	1102	3.787	0.310	3.746	0.376	1.011	0.0373
546.3413	1204	3.436	0.266	3.807	0.217	0.903	0.0374
813.5871	1202	5.325	0.192	6.026	0.297	0.884	0.038
378.9906	1204	3.925	0.109	4.124	0.106	0.952	0.0398
718.4736	1204	6.848	0.489	5.319	0.469	1.287	0.0401
384.3399	1203	68.518	2.885	70.221	1.973	0.976	0.0403
804.5476	1201	16.738	1.052	16.471	0.725	1.016	0.043
174.1411	1203	6.924	0.220	7.701	0.168	0.889	0.0441
780.5872	1204	10.743	0.554	12.203	0.690	0.880	0.0443
793.4936	1204	40.783	2.537	35.440	1.943	1.151	0.0443
834.5963	1201	5.412	0.351	4.956	0.368	1.092	0.0476
541.3141	1201	2.684	0.174	2.773	0.153	0.968	0.048

[00184] Table 6: Accurate mass features differing between clinically diagnosed RR- MULTIPLE SCLEROSIS patients transitioning to SP- MULTIPLE SCLEROSIS and SP- MULTIPLE SCLEROSIS patients ($p < 0.05$).

541.3141	1201	2.876	0.170	2.734	0.234	1.052	0.0007
567.3547	1102	12.556	0.589	7.173	0.794	1.750	0.0022
239.939	1102	4.477	0.259	4.397	0.495	1.018	0.0034
872.6715	1204	2.361	0.259	4.860	0.498	0.486	0.0052
555.3102	1102	6.389	0.294	3.908	0.545	1.635	0.0061
760.5231	1204	95.686	4.233	54.548	5.832	1.754	0.0067
761.529	1204	40.261	1.759	24.002	2.250	1.677	0.0071
788.5549	1204	20.553	0.816	13.393	0.984	1.535	0.0075
566.3431	1102	7.397	0.292	4.439	0.519	1.666	0.0081
786.5408	1204	114.078	4.957	69.266	6.251	1.647	0.0081
565.3391	1102	24.577	1.046	14.619	1.846	1.681	0.0084
784.5238	1204	91.976	3.932	55.254	6.193	1.665	0.0099
783.6174	1204	32.380	3.488	8.590	2.388	3.769	0.0106
746.5118	1204	74.992	3.148	46.823	4.708	1.602	0.0107
303.1081	1102	4.422	0.334	2.767	0.507	1.598	0.0108
249.9677	1102	6.051	0.305	5.069	0.455	1.194	0.0115
787.5452	1204	52.453	2.226	32.632	2.794	1.607	0.0124
305.8792	1102	6.257	0.388	3.671	0.627	1.704	0.0125
784.6228	1204	20.109	1.910	6.221	1.486	3.232	0.0125
684.6037	1204	2.308	0.286	5.278	0.564	0.437	0.0145
785.5287	1204	46.467	1.990	27.795	3.186	1.672	0.0157
718.4738	1204	6.929	0.499	4.007	0.540	1.729	0.0175
770.5108	1204	90.825	3.423	61.142	5.177	1.485	0.0185
331.957	1102	3.358	0.274	3.566	0.398	0.942	0.0212
808.5225	1204	43.754	2.334	26.220	2.895	1.669	0.0215
633.3232	1102	4.455	0.197	2.881	0.303	1.546	0.025
809.5264	1204	22.485	1.215	13.576	1.484	1.656	0.0258
333.9539	1102	3.182	0.238	3.120	0.333	1.020	0.0276
772.5255	1204	117.162	4.503	81.159	5.960	1.444	0.028
733.501	1204	31.675	1.880	16.611	2.461	1.907	0.0297
747.5121	1204	54.603	2.277	41.134	3.281	1.327	0.0317
246.1468	1201	4.033	0.358	6.273	0.755	0.643	0.0334
828.7213	1201	2.551	0.267	3.744	0.802	0.681	0.0346
856.7527	1201	9.061	1.308	14.037	3.456	0.645	0.0369
617.0921	1204	277.341	10.025	201.035	13.034	1.380	0.0378
574.4958	1201	6.712	0.792	12.295	2.617	0.546	0.0382
742.4745	1204	10.479	0.392	7.859	0.576	1.333	0.0398
716.4987	1204	25.168	1.344	20.579	1.653	1.223	0.0403
757.5008	1204	47.638	2.832	25.711	4.126	1.853	0.0403
854.737	1201	6.333	0.933	10.183	2.684	0.622	0.0403
379.2536	1204	3.069	0.171	1.706	0.224	1.799	0.0454
734.508	1204	14.722	1.041	6.522	1.436	2.257	0.0475

[00185] Table 7: Metabolites identified in first principle component analysis for RR-multiple sclerosis.

540.4387	1204	36.603	2.086	22.346	1.749	1.638	5.92E-07
578.4923	1204	41.017	2.169	26.563	2.068	1.544	1.37E-06
596.5012	1204	181.033	13.876	108.540	10.921	1.668	3.44E-05
597.5062	1204	64.543	4.659	39.473	3.816	1.635	3.06E-05
594.4848	1204	116.663	6.054	80.027	6.363	1.458	0.0001

[00186] Table 8: Expanded set of metabolites identified in second PAM analysis for RR-multiple sclerosis.

540.4387	1204	36.603	2.086	22.346	1.749	1.638	5.92E-07
538.4257	1204	30.014	1.397	21.271	1.294	1.411	9.34E-06
594.4848	1204	116.663	6.054	80.027	6.363	1.458	0.0001
578.4923	1204	41.017	2.169	26.563	2.068	1.544	1.37E-06
596.5012	1204	181.033	13.876	108.540	10.921	1.668	3.44E-05
468.384	1204	18.514	0.778	12.977	0.754	1.427	2.69E-07
595.4883	1204	46.584	2.416	32.020	2.556	1.455	0.0001
597.5062	1204	64.543	4.659	39.473	3.816	1.635	3.06E-05
384.3399	1203	69.859	1.997	62.624	1.789	1.116	0.0005
576.4757	1204	45.791	2.161	33.088	2.440	1.384	0.0001
763.5153	1204	21.461	2.719	12.389	1.529	1.732	0.0008
541.4422	1204	12.488	0.710	7.651	0.568	1.632	4.16E-07
522.4313	1204	15.891	0.597	11.438	0.681	1.389	9.00E-07
496.4157	1204	10.848	0.581	6.751	0.455	1.607	3.72E-08
765.5316	1204	32.360	2.588	24.079	1.770	1.344	0.0022
745.5643	1204	120.519	6.555	105.213	5.381	1.145	0.0068

[00187] Table 9: Clinically diagnosed RR- MULTIPLE SCLEROSIS patients and controls used in the test set and their actual and predicted diagnosis.

BB000636	RR-MS	RR-MS
BB000761	RR-MS	RR-MS
BB000775	RR-MS	control
BB000792	RR-MS	RR-MS
BB000796	RR-MS	control
BB000852	RR-MS	RR-MS
BB000855	RR-MS	control
BB000866	RR-MS	RR-MS
BB000870	RR-MS	RR-MS
BB000712	RR-MS	RR-MS
BB000241	RR-MS	RR-MS
BB000246	RR-MS	RR-MS
BB000249	RR-MS	RR-MS
BB000251	RR-MS	RR-MS
BB000633	RR-MS	RR-MS
BB000235	RR-MS	control
BB000259	RR-MS	RR-MS
BB003037	control	control
BB002858	control	RR-MS
BB002859	control	control
BB002862	control	control
BB002865	control	RR-MS
BB003011	control	control
BB003012	control	control
BB003013	control	control
BB003016	control	control
BB003017	control	control
BB002866	control	control
BB002867	control	control
BB002861	control	control
BB002870	control	control
BB002874	control	control
BB003006	control	control
BB003009	control	control
BB003014	control	control
BB003021	control	control
BB003023	control	control
BB002852	control	control
BB002854	control	RR-MS
BB002855	control	control
BB002863	control	RR-MS
BB002864	control	control

[00188] Table 10: Sample numbers and optimal number of metabolites used in training sets for each clinical pairing.

Clinically diagnosed RR-MIS	17	16	16.1%
Controls	25		
Clinically diagnosed RR-MIS	9	7	11.4%
Controls	18		
Clinically diagnosed RR-MIS	11	16	16.6%
Controls	23		
Clinically diagnosed RR-MIS	18	17	5%
Controls	18		
Clinically diagnosed RR-MIS	18	9	6.3%
Controls	15		
Clinically diagnosed RR-MIS	18	17	14.2%
Controls	15		

[00189] Table 11: Optimal Number of Metabolites and Prediction Results for clinically diagnosed PP- MULTIPLE SCLEROSIS and controls.

218.04	1102	23.392	1.656	18.040	0.754	1.297	4.97E-05
202.0453	1101	29.842	2.721	21.457	0.892	1.391	2.59E-05
244.0559	1101	11.378	0.934	8.100	0.241	1.405	1.73E-06
218.0371	1102	7.872	0.566	6.067	0.256	1.297	3.72E-05
831.5992	1102	33.924	6.553	50.976	4.257	0.666	0.0391
243.0719	1101	33.520	3.834	26.746	1.469	1.253	0.0003
832.6022	1102	15.107	2.641	22.020	1.676	0.666	0.0439

BB000816		control
BB000879		
BB000929		control
BB001827		
BB000840		control
BB001432		control
BB001924		
BB001925		control
BB002927		control
BB003021	control	control
BB003023	control	control
BB003026	control	control
BB003027	control	control
BB003028	control	control
BB003030	control	control
BB003032	control	control
BB003034	control	control
BB003037	control	control
BB002858	control	control
BB002856	control	control
BB002857	control	control
BB002861	control	control
BB002870	control	control
BB002874	control	control
BB003013	control	control
BB003016	control	control
BB003017	control	control
BB003018	control	control
BB003019	control	control
BB003022	control	
BB003031	control	control
BB003033	control	control
BB003035	control	
BB002851	control	control

[00190] Table 12: Optimal Number of Metabolites and Prediction Results for clinically diagnosed SP- MULTIPLE SCLEROSIS and controls.

805.5609	1201	55.027	3.212	36.921	1.704	1.490	0.0093
806.5643	1201	27.948	1.637	18.717	0.860	1.493	0.0075
541.3415	1102	15.111	1.031	25.470	1.129	0.593	0.0203
594.4848	1204	43.087	5.319	80.027	6.363	0.538	0.0048
596.5012	1204	50.842	4.367	108.540	10.921	0.468	0.003
597.5062	1204	18.785	1.588	39.473	3.816	0.476	0.0022
827.5446	1201	15.396	1.216	9.932	0.505	1.550	0.0191
538.4257	1204	14.296	1.804	21.271	1.294	0.672	0.0094
576.4757	1204	19.011	1.969	33.088	2.440	0.575	0.0047
595.4883	1204	17.967	2.221	32.020	2.556	0.561	0.0051
886.5582	1102	5.104	0.302	9.581	0.608	0.533	0.0267
578.4923	1204	15.340	1.107	26.563	2.068	0.578	0.0042
540.4387	1204	14.521	1.256	22.346	1.749	0.650	0.0108
428.3653	1201	9.177	0.839	4.617	0.315	1.988	2.84E-05
622.4973	1203	20.247	1.598	14.995	0.688	1.350	0.0219
694.6323	1204	12.155	1.212	7.889	0.727	1.541	0.0283

BB000786
BB000787
BB000847
BB000829
BB000906
BB001744
BB001826
BB001928
BB001942
BB002759
BB002878
BB003014
BB003021
BB003023
BB003026
BB003027
BB002858
BB002859

control
control
control
control
control
control
control
control
control
control
control
control
control
control
control
control
control
control

BB002862
BB002865
BB002866
BB002856
BB002857
BB002861
BB002870
BB002874
BB003007
BB003011
BB003012
BB003013
BB003016
BB003004
BB003015
BB003022
BB003031
BB003033

control
control
control
control
control
control
control
control
control
control
control
control
control
control
control
control
control
control
control

[00191] Table 13: Optimal Number of Metabolites and Prediction Results for clinically diagnosed RR- MULTIPLE SCLEROSIS and SP- MULTIPLE SCLEROSIS patients.

578.4923	1204	15.340	1.107	40.496	2.263	0.379	2.21E-10
594.4848	1204	43.087	5.319	115.274	6.359	0.374	5.81E-08
596.5012	1204	50.842	4.367	178.485	14.248	0.285	7.94E-08
576.4757	1204	19.011	1.969	45.354	2.261	0.419	2.02E-08
595.4883	1204	17.967	2.221	46.028	2.536	0.390	6.43E-08
597.5062	1204	18.785	1.588	63.633	4.798	0.295	4.01E-08
805.5809	1201	55.027	3.212	38.808	2.258	1.418	0.0463
592.4717	1204	17.304	2.411	38.114	2.424	0.464	2.82E-06
512.4079	1204	6.864	0.551	15.906	0.838	0.432	4.59E-10
579.4958	1204	6.320	0.458	15.978	0.888	0.396	5.47E-10
580.5089	1204	5.730	0.402	13.528	0.711	0.424	1.98E-10
460.384	1204	9.354	0.797	18.294	0.792	0.511	2.96E-09
538.4257	1204	14.296	1.804	29.559	1.464	0.484	3.34E-07
577.4795	1204	7.187	0.718	16.759	0.829	0.429	1.26E-08
806.5643	1201	27.948	1.637	19.726	1.161	1.417	0.0444
540.4387	1204	14.521	1.256	36.081	2.150	0.402	6.32E-08

BB000636	RR-MS	RR-MS
BB000761	RR-MS	RR-MS
BB000775	RR-MS	RR-MS
BB000792	RR-MS	RR-MS
BB000796	RR-MS	RR-MS
BB000736	RR-MS	RR-MS
BB000758	RR-MS	RR-MS
BB000763	RR-MS	RR-MS
BB000766	RR-MS	RR-MS
BB000771	RR-MS	RR-MS
BB000246	RR-MS	RR-MS
BB000249	RR-MS	RR-MS
BB000251	RR-MS	RR-MS
BB000633	RR-MS	RR-MS
BB000734	RR-MS	RR-MS
BB000777	RR-MS	RR-MS
BB000780	RR-MS	RR-MS
BB000781	RR-MS	RR-MS
BB000782	RR-MS	RR-MS
BB000793	RR-MS	RR-MS
BB000841	RR-MS	RR-MS
BB000848	RR-MS	RR-MS
BB000857	RR-MS	RR-MS
BB000858	RR-MS	RR-MS
BB000863	RR-MS	RR-MS
BB000867	RR-MS	RR-MS
BB000829		
BB000906		
BB000921		
BB001124		
BB001126		
BB001928		
BB001942		
BB002759		
BB002878		RR-MS

[00192] Table 14: Optimal Number of Metabolites and Prediction Results for RR- MULTIPLE SCLEROSIS patients transitioning to SP- MULTIPLE SCLEROSIS and clinically diagnosed RR- MULTIPLE SCLEROSIS patients.

578.4923	1204	22.148	2.145	40.496	2.263	0.547	9.92E-15
594.4848	1204	65.843	8.436	115.274	6.359	0.571	2.35E-11
576.4757	1204	26.727	2.897	45.354	2.261	0.589	1.30E-12
596.5012	1204	89.261	13.368	178.485	14.249	0.500	1.62E-10
595.4883	1204	25.969	3.279	46.028	2.536	0.564	1.93E-11
597.5062	1204	32.711	4.840	83.633	4.798	0.514	8.79E-11
640.4387	1204	21.969	2.791	36.081	2.150	0.609	1.82E-10
592.4717	1204	21.832	2.186	38.114	2.424	0.573	3.20E-10
579.4958	1204	8.772	0.838	15.978	0.888	0.549	1.08E-14

BB000775	RR-MS	
BB000792	RR-MS	RR-MS
BB000796	RR-MS	RR-MS
BB000799	RR-MS	RR-MS
BB000814	RR-MS	RR-MS
BB000771	RR-MS	RR-MS
BB000773	RR-MS	RR-MS
BB000777	RR-MS	RR-MS
BB000780	RR-MS	RR-MS
BB000781	RR-MS	RR-MS
BB000863	RR-MS	RR-MS
BB000867	RR-MS	RR-MS
BB000223	RR-MS	RR-MS
BB000230	RR-MS	RR-MS
BB000234	RR-MS	RR-MS
BB000793	RR-MS	RR-MS
BB000800	RR-MS	RR-MS
BB000815	RR-MS	RR-MS
BB000832	RR-MS	RR-MS
BB000856	RR-MS	RR-MS
BB000235	RR-MS	
BB000259	RR-MS	RR-MS
BB000636	RR-MS	RR-MS

BB000761	RR-MS	RR-MS
BB000822	RR-MS	RR-MS
BB000841	RR-MS	RR-MS
BB000801		
BB000807		
BB000817		
BB000826		
BB000827		
BB000717		
BB000754		
BB000759		
BB000764		
BB000794		
BB000224		
BB000227		
BB000232		
BB000238		
BB000240		
BB000859		
BB000221		
BB000225		
BB000236		
BB000252		

[00193] Table 15: Optimal Number of Metabolites and Prediction Results for RR- MULTIPLE SCLEROSIS patients transitioning to SP- MULTIPLE SCLEROSIS and clinically diagnosed SP- MULTIPLE SCLEROSIS patients.

760.5231	1204	96.686	4.233	54.548	5.832	1.754	0.0067
746.5118	1204	74.992	3.148	46.823	4.708	1.602	0.0107
786.5408	1204	114.078	4.957	69.256	8.251	1.647	0.0081
565.3391	1102	24.577	1.046	14.619	1.846	1.681	0.0084
808.5225	1204	43.754	2.334	26.220	2.895	1.669	0.0215
761.529	1204	40.261	1.759	24.002	2.250	1.677	0.0071
772.5265	1204	117.162	4.503	81.159	5.960	1.444	0.028
784.5238	1204	91.976	3.932	55.254	6.193	1.665	0.0099
617.0921	1204	277.341	10.025	201.035	13.034	1.380	0.0378
787.5452	1204	52.453	2.226	32.632	2.794	1.607	0.0124
733.501	1204	31.675	1.880	16.611	2.461	1.907	0.0297
785.5287	1204	46.467	1.990	27.795	3.186	1.672	0.0157
770.5108	1204	90.825	3.423	61.142	5.177	1.485	0.0185
809.5264	1204	22.485	1.215	13.576	1.484	1.656	0.0258
783.6174	1204	32.380	3.488	8.590	2.388	3.769	0.0106
734.508	1204	14.722	1.041	6.522	1.436	2.257	0.0475
757.5008	1204	47.638	2.832	25.711	4.126	1.853	0.0403

BB000232
BB000236
BB000238
BB000240
BB000247
BB000834
BB000836
BB000842
BB000846
BB000849
BB000749
BB000752
BB000754
BB000759
BB000764
BB000222
BB000224

BB000225
BB000227
BB000248
BB000801
BB000807
BB000809
BB001124
BB001125
BB001153
BB001386
BB001744
BB000755
BB000784
BB000786
BB000787
BB001942
BB002759

[00194] Table 16: Accurate mass features differing between 10 clinically diagnosed RR- MULTIPLE SCLEROSIS patients and 10 controls ($p < 0.05$).

450.3729	1204	11.764	0.454	4.823	0.364	2.439	1.44E-10
512.3347	1201	3.928	0.689	2.788	0.548	1.409	1.44E-10
580.5089	1204	20.678	1.357	4.648	0.571	4.449	5.92E-10
513.4116	1204	7.981	0.438	2.392	0.342	3.337	1.89E-09
578.4923	1204	63.796	4.758	13.250	1.859	4.815	2.05E-09
679.4958	1204	24.980	1.842	5.626	0.692	4.440	3.44E-09
452.3868	1204	5.142	0.235	1.726	0.172	2.979	4.16E-09
581.5126	1204	8.401	0.625	1.791	0.257	4.691	5.10E-09
541.4422	1204	16.895	0.945	4.594	0.763	3.678	7.55E-09
596.5053	1202	22.164	1.738	4.557	0.776	4.864	1.19E-08
540.4387	1204	49.312	3.081	12.594	2.011	3.916	1.42E-08
448.3562	1204	12.704	0.509	6.139	0.508	2.069	1.69E-08
523.3637	1101	2.792	0.154	3.567	0.155	0.783	2.14E-08
494.3968	1204	19.629	1.092	6.605	0.867	2.972	2.25E-08
522.4313	1204	19.942	1.072	7.469	0.944	2.670	2.85E-08
594.4848	1204	165.993	10.796	45.401	7.979	3.656	4.81E-08
595.4883	1204	65.756	4.493	18.270	3.207	3.599	5.87E-08
597.5068	1202	8.430	0.707	2.001	0.308	4.213	6.69E-08
484.3788	1204	11.357	0.673	3.696	0.563	3.073	7.21E-08
588.4723	1204	19.140	1.666	4.476	0.627	4.276	7.36E-08
510.3937	1204	8.643	0.463	2.755	0.495	3.137	1.05E-07
610.482	1204	13.307	1.141	4.336	0.505	3.069	1.09E-07
552.3273	1201	6.531	1.036	4.365	0.492	1.496	1.11E-07
576.4757	1204	64.886	4.627	19.420	3.428	3.341	1.12E-07
495.4018	1204	6.142	0.341	2.325	0.300	2.642	1.21E-07
506.4338	1204	4.255	0.199	1.532	0.236	2.777	1.42E-07
478.4044	1204	4.567	0.254	2.095	0.208	2.180	1.54E-07
536.41	1204	13.137	0.707	5.601	0.743	2.345	1.66E-07
521.4188	1204	6.886	0.318	3.052	0.382	2.256	1.74E-07
468.384	1204	23.300	1.274	8.827	1.121	2.640	1.97E-07
569.4769	1204	13.342	0.782	16.192	0.834	0.824	2.06E-07
577.4795	1204	23.556	1.747	7.266	1.323	3.242	2.11E-07
508.3782	1204	6.701	0.397	2.628	0.357	2.650	2.43E-07
598.5107	1204	23.222	2.389	3.965	0.530	5.857	2.84E-07
550.4602	1204	14.919	1.141	4.800	0.796	3.108	2.93E-07
469.3863	1204	6.296	0.306	2.276	0.365	2.766	3.36E-07
466.3656	1204	17.573	0.871	7.234	0.974	2.429	3.53E-07
566.454	1204	11.369	0.761	4.079	0.698	2.787	5.16E-07
496.4157	1204	14.715	1.316	4.097	0.438	3.592	6.10E-07
597.5062	1204	105.816	11.948	17.041	2.549	6.209	6.44E-07
596.5012	1204	305.004	35.345	46.149	7.337	6.609	8.06E-07
548.4438	1204	9.290	0.523	3.896	0.575	2.384	8.46E-07
524.4448	1204	7.031	0.602	2.530	0.347	2.779	9.55E-07
467.3711	1204	5.622	0.304	2.318	0.352	2.425	1.13E-06
537.4142	1204	5.123	0.277	2.320	0.330	2.208	1.24E-06
590.4585	1204	15.750	1.532	5.639	0.766	2.793	1.25E-06
440.3526	1204	5.976	0.328	2.323	0.361	2.573	1.56E-06
520.4131	1204	19.275	0.868	8.963	1.213	2.151	1.63E-06
327.0307	1204	8.911	0.403	5.729	0.200	1.555	1.66E-06
562.4989	1204	15.568	0.969	6.545	0.819	2.379	1.91E-06
482.3604	1204	6.472	0.517	2.744	0.359	2.359	2.29E-06
538.4257	1204	35.521	2.265	14.248	2.317	2.493	2.63E-06
492.3832	1204	11.353	0.647	5.057	0.656	2.245	2.78E-06
570.4903	1204	7.316	0.836	1.596	0.209	4.584	2.84E-06
564.4396	1204	4.654	0.336	1.888	0.348	2.465	3.44E-06
551.4646	1204	5.429	0.483	1.857	0.313	2.924	3.83E-06
534.4645	1204	4.848	0.327	1.969	0.266	2.462	4.10E-06
534.3912	1204	5.803	0.287	2.348	0.490	2.471	5.38E-06

563.5013	1204	5.800	0.356	2.661	0.306	2.180	6.24E-06
564.513	1204	5.778	0.469	2.465	0.135	2.344	6.52E-06
594.4875	1202	12.561	1.194	4.445	0.981	2.826	6.62E-06
493.385	1204	3.804	0.230	1.686	0.272	2.256	6.96E-06
595.4928	1202	5.270	0.461	2.080	0.364	2.534	1.05E-05
576.4765	1202	5.732	0.787	1.685	0.289	3.402	1.18E-05
438.3354	1204	4.500	0.240	1.999	0.320	2.251	1.41E-05
518.3969	1204	14.324	0.833	7.105	1.015	2.016	1.55E-05
378.2921	1203	39.401	4.829	35.156	2.294	1.121	1.84E-05
464.3524	1204	10.992	0.512	5.327	0.783	2.063	1.83E-05
476.3869	1204	5.546	0.259	2.924	0.410	1.897	1.94E-05
519.3998	1204	5.121	0.373	2.423	0.396	2.113	2.59E-05
618.4834	1201	8.419	1.292	1.816	0.294	4.636	3.02E-05
480.3473	1204	4.616	0.238	2.369	0.377	1.949	4.26E-05
384.3399	1203	79.789	4.637	54.933	2.431	1.452	4.36E-05
593.4736	1204	21.988	2.537	7.306	1.312	3.010	4.64E-05
253.8165	1101	14.128	1.242	9.713	0.417	1.455	0.0001
284.9759	1204	7.668	0.276	5.890	0.250	1.302	0.0001
504.4188	1204	7.881	0.427	4.454	0.498	1.769	0.0001
591.4614	1204	5.417	0.569	2.303	0.323	2.352	0.0001
592.4717	1204	52.044	8.345	17.230	3.267	3.021	0.0001
612.4994	1204	8.902	0.755	3.979	0.613	2.237	0.0001
255.8136	1101	18.945	1.883	12.863	0.495	1.473	0.0002
385.3428	1203	25.610	1.669	17.850	0.822	1.435	0.0003
569.3687	1102	5.808	0.687	9.114	0.798	0.637	0.0003
616.4675	1201	5.743	0.847	2.166	0.403	2.651	0.0003
769.5638	1204	143.765	14.976	95.450	10.239	1.506	0.0003
770.569	1204	63.988	8.122	44.106	4.234	1.451	0.0003
474.3731	1204	5.273	0.308	2.976	0.401	1.772	0.0004
572.4455	1204	5.865	0.678	2.979	0.269	1.969	0.0004
446.341	1204	14.884	1.020	8.471	1.005	1.757	0.0005
447.3433	1204	4.850	0.350	2.489	0.366	1.949	0.0005
574.4594	1204	43.223	6.002	15.738	3.104	2.746	0.0005
462.3346	1204	4.094	0.171	2.281	0.372	1.795	0.0007
490.3676	1204	6.960	0.460	3.926	0.650	1.773	0.0007
502.4054	1204	7.152	0.530	4.183	0.460	1.710	0.0007
546.4298	1204	4.944	0.458	2.619	0.430	1.888	0.0008
575.4628	1204	16.464	2.289	6.352	1.189	2.592	0.0008
712.5074	1204	33.680	4.456	22.449	2.683	1.500	0.0011
1018.9399	1203	16.341	1.538	8.303	1.195	1.968	0.0011
558.3761	1204	7.224	0.825	4.791	0.249	1.508	0.0013
532.4503	1204	5.319	0.396	2.825	0.426	1.883	0.0015
716.4323	1204	17.596	2.776	9.512	0.571	1.850	0.0015
561.4863	1204	11.883	0.972	6.295	0.992	1.888	0.0016
713.5097	1204	14.360	1.899	9.557	1.160	1.603	0.0017
314.2461	1204	6.654	0.541	4.749	0.383	1.401	0.0018
160.1256	1203	7.861	0.351	6.088	0.622	1.291	0.0021
558.4649	1204	47.704	4.927	23.917	4.147	1.995	0.0021
781.6001	1204	11.335	1.338	7.350	0.701	1.542	0.0022
747.5761	1204	19.020	1.510	13.858	1.410	1.372	0.0024
539.4274	1204	10.318	1.781	3.468	0.856	2.975	0.0025
686.4879	1204	91.719	12.953	62.831	9.366	1.460	0.0025
546.3413	1204	4.730	0.643	2.304	0.359	2.053	0.003
688.5048	1204	31.508	3.980	22.044	3.110	1.429	0.0033
700.4371	1204	8.421	1.098	5.617	0.564	1.499	0.0033
1016.9279	1203	29.801	2.895	14.233	2.907	2.094	0.0037
560.478	1203	10.183	0.676	6.147	1.045	1.657	0.0039
559.4688	1204	17.743	1.887	9.237	1.570	1.921	0.004
367.3325	1203	12.234	0.329	10.201	0.434	1.199	0.0041
687.4916	1204	37.542	5.128	26.691	3.782	1.407	0.0044
523.3637	1101	2.792	0.184	3.567	0.155	0.783	0.0045
381.311	1203	65.856	6.771	55.622	3.160	1.184	0.005
574.4635	1202	4.109	0.785	1.733	0.354	2.371	0.005
376.2759	1203	18.196	1.641	16.364	0.695	1.112	0.0052
793.5663	1204	62.347	7.352	40.458	5.240	1.541	0.0056

746.5701	1204	59.736	5.308	43.321	5.482	1.379	0.0058
249.9677	1102	7.037	0.812	5.282	0.619	1.332	0.0059
644.3636	1204	4.822	0.603	3.228	0.162	1.494	0.0059
737.5045	1204	8.325	1.154	4.953	0.718	1.681	0.006
745.5643	1204	139.059	14.204	98.300	12.989	1.415	0.006
257.8106	1101	9.000	0.921	6.494	0.503	1.386	0.0063
794.5716	1204	32.136	3.430	21.759	2.588	1.477	0.0064
556.4497	1204	12.580	1.583	6.930	0.673	1.815	0.0067
689.5083	1204	13.040	1.560	9.335	1.306	1.397	0.0071
306.2568	1204	11.187	1.061	7.856	0.676	1.424	0.0072
370.351	1203	83.555	5.638	58.717	5.196	1.423	0.0073
205.8867	1101	8.244	0.345	7.034	0.240	1.172	0.0075
378.9906	1204	4.794	0.192	3.175	0.185	1.510	0.0075
557.4527	1204	5.130	0.622	2.610	0.443	1.966	0.008
369.3475	1203	641.745	44.780	441.804	42.240	1.453	0.0081
371.3542	1203	8.222	0.598	5.626	0.546	1.461	0.0088
702.4175	1204	12.357	2.087	7.404	0.542	1.669	0.0091
736.5031	1204	16.701	2.447	10.717	1.478	1.558	0.0092
743.5461	1204	452.234	52.550	321.639	51.754	1.406	0.0092
832.6022	1102	19.129	4.193	27.747	5.113	0.689	0.0092
744.55	1203	10.268	1.731	4.670	0.589	2.199	0.0095
722.5244	1204	11.878	1.868	7.942	0.684	1.496	0.0096
244.2189	1203	7.865	0.432	6.421	0.300	1.225	0.0103
263.8453	1101	8.105	0.644	6.599	0.299	1.228	0.0104
154.0035	1204	28.896	1.974	23.716	0.905	1.218	0.0105
530.3474	1204	54.770	6.770	41.447	1.894	1.321	0.0106
698.4885	1204	16.687	1.264	13.465	0.804	1.239	0.0106
776.556	1204	16.853	2.541	9.051	2.083	1.862	0.0107
779.5828	1204	29.410	4.417	18.711	2.200	1.572	0.0112
778.571	1204	13.161	2.066	7.907	1.215	1.664	0.0113
856.6009	1102	20.936	3.842	30.073	5.796	0.696	0.0113
743.5475	1203	22.097	3.492	10.450	1.319	2.115	0.0114
340.2407	1204	5.090	0.309	5.985	0.256	0.850	0.0115
831.5992	1102	44.121	10.631	65.931	13.152	0.869	0.0116
460.2681	1204	11.360	1.336	8.419	0.492	1.349	0.0117
624.5133	1203	24.990	1.800	16.782	1.739	1.489	0.0117
720.5081	1204	7.819	0.953	5.605	0.546	1.395	0.0117
730.4535	1204	31.096	4.976	20.165	0.999	1.542	0.012
432.2365	1204	3.845	0.328	2.860	0.316	1.344	0.0122
789.5658	1204	11.503	0.837	9.419	0.562	1.221	0.0127
446.2525	1204	6.209	0.643	4.533	0.209	1.370	0.0129
646.5702	1203	8.272	0.830	5.522	0.802	1.498	0.013
758.4785	1204	77.135	9.484	54.605	2.972	1.413	0.013
740.4966	1204	23.136	2.101	18.075	1.188	1.280	0.0131
744.5516	1204	183.026	19.581	136.722	20.173	1.339	0.0135
780.5872	1204	14.936	2.228	9.362	1.249	1.595	0.0135
907.7722	1203	26.436	2.735	17.366	2.098	1.522	0.0147
625.5161	1203	11.098	0.830	7.474	0.759	1.485	0.0148
623.5003	1203	8.971	1.018	5.816	0.536	1.542	0.0156
885.7866	1203	1.000	0.000	8.100	2.845	0.123	0.0158
906.7669	1203	45.385	5.161	29.567	3.715	1.535	0.0167
488.2996	1204	6.572	0.927	4.279	0.292	1.536	0.0168
558.4663	1202	5.415	1.033	3.113	0.430	1.739	0.0168
775.5514	1204	35.635	5.243	20.650	4.220	1.726	0.0168
239.939	1102	5.794	0.741	4.295	0.571	1.349	0.0171
462.3716	1204	3.522	0.290	2.681	0.165	1.314	0.0171
530.3474	1204	54.770	6.770	41.447	1.894	1.321	0.0181
856.6045	1102	10.732	1.955	15.064	2.836	0.712	0.0182
541.3415	1102	22.495	2.750	27.652	2.246	0.814	0.0189
648.5661	1203	28.975	1.870	20.476	3.455	1.415	0.019
211.8495	1102	6.647	0.650	4.835	0.395	1.375	0.0195
516.3324	1204	8.573	1.124	6.171	0.354	1.389	0.0201
729.5727	1204	11.687	1.327	8.167	1.211	1.431	0.0201
380.3079	1203	219.668	23.133	187.016	10.229	1.115	0.0202
232.2189	1203	9.603	0.999	6.990	0.603	1.374	0.0208

502.3165	1204	37.199	5.399	25.495	1.609	1.459	0.0213
570.3766	1201	2.370	0.366	1.292	0.159	1.834	0.0214
726.5438	1204	10.420	1.505	6.298	1.240	1.654	0.0219
146.11	1203	5.606	0.311	4.764	0.525	1.177	0.0221
503.3194	1204	10.490	1.478	6.992	0.517	1.500	0.0222
524.296	1201	5.049	0.769	3.162	0.424	1.597	0.0227
742.5366	1204	18.967	2.099	12.996	2.718	1.459	0.0231
777.5678	1204	26.326	4.120	16.489	2.472	1.597	0.0233
727.5554	1204	27.866	3.724	18.243	3.941	1.527	0.0234
286.2656	1203	10.402	1.248	7.166	0.711	1.452	0.0247
728.5605	1204	13.981	1.722	9.663	1.803	1.447	0.0254
260.2507	1204	21.769	3.306	14.449	1.741	1.507	0.026
265.8424	1101	7.482	0.580	6.008	0.368	1.245	0.026
753.5883	1204	26.430	4.780	16.198	2.044	1.632	0.0263
242.2032	1203	18.369	2.262	12.872	0.982	1.427	0.0272
545.3455	1101	2.905	0.196	3.734	0.258	0.778	0.0272
377.2801	1203	6.047	0.635	5.892	0.316	1.062	0.0275
649.5895	1203	13.892	1.020	9.866	1.723	1.408	0.0281
531.3504	1204	16.841	2.194	12.776	0.643	1.318	0.0285
763.5153	1204	29.196	5.814	13.222	5.672	2.208	0.0285
569.369	1202	13.342	0.782	16.192	0.834	0.824	0.0292
909.7867	1203	16.698	1.762	11.334	1.721	1.473	0.0301
311.7754	1101	5.876	0.894	3.902	0.528	1.506	0.0308
272.2501	1203	8.545	1.006	6.078	0.598	1.406	0.0309
622.4973	1203	19.859	2.358	13.139	1.241	1.511	0.0316
552.3273	1201	6.531	1.036	4.365	0.492	1.496	0.0319
672.586	1203	11.817	1.553	7.815	1.024	1.512	0.0324
340.2621	1204	5.577	0.638	4.305	0.400	1.295	0.0326
271.8051	1102	8.179	0.967	5.695	0.859	1.436	0.0328
858.6798	1204	5.556	1.040	2.729	0.579	2.036	0.0334
715.4864	1204	21.704	1.962	18.368	1.657	1.182	0.0338
899.5871	1102	10.315	2.129	13.414	2.488	0.769	0.0344
244.0559	1101	10.439	1.409	7.511	0.240	1.390	0.0346
512.4079	1204	24.363	1.277	6.085	0.936	4.004	0.0346
181.9806	1102	6.360	0.914	4.833	0.724	1.314	0.0357
754.5724	1204	11.650	1.883	7.945	0.856	1.466	0.0361
783.6174	1204	14.897	2.745	25.452	6.501	0.585	0.0368
379.2957	1203	10.169	1.195	9.312	0.567	1.092	0.0369
725.5376	1204	20.287	3.674	11.990	2.404	1.692	0.0369
764.5196	1204	13.325	3.070	5.424	2.939	2.457	0.037
345.8738	1101	5.195	0.436	6.130	0.367	0.847	0.0372
797.5973	1204	32.079	2.947	25.825	2.485	1.242	0.0381
330.2569	1204	3.288	0.416	2.317	0.210	1.419	0.0385
626.5271	1203	31.047	2.292	22.120	2.012	1.404	0.0385
202.0453	1101	27.142	4.189	18.281	0.929	1.485	0.0386
542.3447	1102	6.520	0.775	7.672	0.576	0.850	0.0395
738.5185	1204	27.693	3.848	20.578	2.843	1.346	0.04
144.0944	1203	6.264	0.355	5.320	0.523	1.177	0.0412
699.4908	1204	7.701	0.616	6.202	0.488	1.242	0.0422
584.2641	1202	5.916	1.460	3.159	0.549	1.873	0.0431
806.413	1204	4.879	1.664	2.072	0.870	2.355	0.0433
305.2439	1204	8.424	0.931	6.262	0.570	1.345	0.0435
207.8836	1101	6.837	0.380	5.851	0.405	1.169	0.044
780.5303	1204	10.924	0.738	8.993	0.851	1.215	0.044
773.5954	1204	32.726	3.281	27.668	2.679	1.183	0.0441
304.241	1204	39.983	4.615	29.496	2.692	1.356	0.0445
634.3951	1204	9.972	2.107	5.214	0.746	1.913	0.0446
792.555	1204	37.287	4.795	25.974	4.202	1.436	0.0446
688.4658	1204	12.687	1.417	9.572	0.751	1.325	0.0447
788.4794	1204	12.925	0.869	10.844	0.654	1.192	0.0447
627.5285	1203	13.706	1.127	9.720	0.884	1.410	0.0451
716.4987	1204	26.851	2.554	22.322	1.404	1.203	0.0452
765.5316	1204	37.075	5.186	23.443	5.659	1.581	0.0463
628.5393	1203	15.926	1.688	10.413	1.515	1.529	0.0466
791.5488	1204	72.427	10.198	49.500	8.854	1.463	0.0468

461.2707	1204	3.218	0.476	2.560	0.170	1.257	0.047
741.5302	1204	38.041	5.335	26.680	5.610	1.426	0.0472
781.5619	1204	12.601	1.666	7.960	1.653	1.570	0.0481
711.4947	1204	21.098	3.889	14.530	2.238	1.452	0.049

[00195] Table 17: Accurate mass features differing between 10 clinically diagnosed PP- MULTIPLE SCLEROSIS patients and 10 controls (p<0.05).

218.0371	1102	9.533	0.524	4.929	0.294	1.934	1.13E-08
244.0559	1101	13.840	1.110	6.736	0.509	2.055	3.93E-08
216.04	1102	28.201	1.565	14.157	0.586	1.992	7.45E-08
202.0453	1101	36.309	3.700	16.259	0.940	2.233	9.69E-07
226.0888	1102	14.690	1.479	8.398	0.768	1.749	2.99E-06
243.0719	1101	41.426	5.686	19.420	1.732	2.133	8.10E-06
273.9985	1102	5.556	0.606	1.951	0.368	2.848	2.31E-05
382.1084	1101	4.696	1.074	1.216	0.115	3.862	4.34E-05
253.8165	1101	12.798	0.768	8.954	0.705	1.429	0.0001
218.0192	1101	11.301	2.271	3.211	0.680	3.519	0.0002
188.0143	1102	6.273	1.597	1.293	0.196	4.852	0.0005
257.8106	1101	8.886	0.670	6.118	0.463	1.452	0.0005
260.004	1101	6.981	1.052	2.623	0.456	2.661	0.0005
333.9539	1102	5.156	0.719	1.864	0.382	2.768	0.0008
806.5643	1201	22.470	2.372	16.863	2.290	1.348	0.001
833.5931	1201	12.180	1.662	7.722	0.947	1.577	0.001
805.5609	1201	44.055	4.902	33.136	4.448	1.330	0.0013
263.8453	1101	8.384	0.323	6.535	0.524	1.283	0.0014
834.5963	1201	6.556	0.888	4.267	0.506	1.536	0.0014
506.2853	1201	5.751	0.873	2.233	0.235	2.575	0.0016
570.3766	1201	2.849	0.382	1.659	0.285	1.717	0.0017
311.7754	1101	5.404	0.343	3.232	0.556	1.672	0.0019
331.957	1102	5.296	0.724	2.066	0.378	2.563	0.0024
205.8867	1101	8.743	0.471	6.644	0.433	1.316	0.003
255.8135	1101	16.177	1.027	12.092	0.980	1.338	0.003
611.3724	1201	5.038	0.732	3.082	0.523	1.635	0.0031
271.8051	1102	8.408	1.095	4.389	0.612	1.916	0.0032
209.8525	1102	5.656	0.606	2.987	0.496	1.894	0.0038
275.8713	1101	5.952	0.369	4.788	0.389	1.243	0.0038
269.8081	1102	12.049	1.596	6.580	0.794	1.831	0.0042
610.3691	1201	14.008	2.142	8.756	1.600	1.600	0.0047
943.7452	1204	5.801	0.947	3.476	0.533	1.669	0.0055
882.7648	1203	131.949	15.074	87.892	8.198	1.501	0.0063
203.1157	1101	6.030	1.049	3.217	0.557	1.874	0.0064
428.295	1204	3.809	0.428	4.489	0.510	0.849	0.0066
828.5477	1201	6.444	0.675	4.716	0.655	1.366	0.0087
758.5655	1201	64.803	8.320	48.518	8.035	1.336	0.0089
267.811	1102	7.363	0.960	4.236	0.420	1.738	0.009
757.5622	1201	128.225	16.744	97.646	16.254	1.313	0.0098
884.7764	1203	65.646	8.766	39.419	4.888	1.665	0.0101
150.1413	1203	4.696	0.369	3.874	0.409	1.212	0.0102
766.5051	1201	3.970	0.401	2.771	0.518	1.433	0.0111
857.7516	1203	132.550	15.354	80.972	9.566	1.637	0.0111
452.244	1201	7.063	0.810	5.156	0.648	1.370	0.012
613.3404	1202	6.275	0.752	4.929	0.666	1.273	0.013
273.8743	1101	9.490	0.470	7.689	0.685	1.234	0.0131
265.8424	1101	7.220	0.470	6.038	0.317	1.196	0.0133
856.7475	1203	246.398	29.810	152.319	19.797	1.618	0.0133
337.2697	1203	5.612	0.449	4.461	0.237	1.258	0.0149
1253.124	1203	8.679	1.209	5.137	0.834	1.690	0.0161
813.5871	1202	4.817	0.487	6.593	0.671	0.731	0.0198
827.5445	1101	5.520	0.780	3.455	0.349	1.598	0.0201
861.5265	1102	6.955	0.882	7.014	1.082	0.992	0.0207
601.5163	1203	137.952	15.032	109.993	11.197	1.254	0.0215
1228.11	1203	14.619	3.088	7.593	1.801	1.925	0.0215
835.6094	1201	4.387	0.625	3.337	0.569	1.315	0.0223
858.7607	1203	104.669	15.454	57.677	8.571	1.815	0.0228
602.5287	1203	435.124	61.111	296.579	40.379	1.467	0.023

134.11	1203	11.334	1.148	10.055	1.025	1.127	0.0235
785.5934	1201	66.177	9.947	48.757	9.132	1.357	0.0238
1254.131	1203	6.626	0.892	3.523	0.537	1.881	0.024
339.2851	1203	12.818	1.706	7.989	0.924	1.604	0.0249
603.532	1203	182.681	25.875	124.378	17.252	1.469	0.0259
827.5446	1201	12.164	1.506	9.243	1.482	1.318	0.0261
600.513	1203	329.443	38.573	257.344	27.668	1.280	0.0267
810.5967	1201	23.957	3.670	16.403	2.568	1.461	0.0268
136.1258	1203	5.285	0.522	4.442	0.479	1.190	0.0274
885.778	1203	29.109	5.738	17.223	3.290	1.690	0.0276
789.5163	1204	20.925	4.537	13.491	2.853	1.551	0.0278
285.1366	1201	3.390	0.833	1.221	0.221	2.776	0.028
859.7662	1203	48.036	6.665	26.489	3.984	1.813	0.0295
162.1412	1203	5.843	0.525	5.260	0.399	1.111	0.0296
211.8495	1102	6.585	0.809	4.119	0.402	1.599	0.0309
628.5393	1203	14.633	1.243	10.951	1.097	1.336	0.0309
828.5479	1101	3.221	0.413	2.144	0.126	1.502	0.031
336.266	1203	22.012	2.175	16.511	1.026	1.333	0.0313
258.2346	1203	9.388	0.794	12.892	1.074	0.728	0.0323
786.5967	1201	32.109	5.025	23.876	4.433	1.345	0.0328
881.7549	1203	71.086	7.371	55.715	3.669	1.276	0.0328
794.5419	1102	7.159	0.964	7.518	1.208	0.952	0.0336
338.2815	1203	63.145	7.870	40.778	4.091	1.549	0.034
781.497	1204	10.007	1.616	11.541	1.469	0.867	0.0347
184.1255	1203	5.899	0.492	5.451	0.380	1.082	0.0379
684.8037	1204	4.384	0.914	2.189	0.523	2.003	0.038
851.5686	1102	9.495	1.864	10.212	2.466	0.930	0.0392
880.7514	1203	127.222	15.132	97.678	7.882	1.302	0.0392
809.5934	1201	46.591	7.354	32.314	5.156	1.442	0.0408
148.1257	1203	7.366	0.769	6.753	0.637	1.091	0.0415
850.6899	1203	4.878	1.364	3.038	1.308	1.606	0.0417
161.1051	1101	4.609	0.742	3.270	0.598	1.409	0.0437
534.3186	1201	5.171	0.837	3.421	0.384	1.512	0.0444
852.5724	1102	5.230	0.908	5.566	1.313	0.940	0.0449
785.4799	1204	18.864	3.666	12.479	2.715	1.512	0.0455
207.8836	1101	6.333	0.187	5.541	0.328	1.143	0.0465
811.5718	1202	3.876	0.448	4.974	0.356	0.779	0.0466
793.5986	1201	6.760	0.935	5.649	0.925	1.197	0.0482
855.7361	1203	123.848	14.395	88.031	9.890	1.407	0.0482
720.4696	1204	6.438	0.985	4.636	0.451	1.389	0.0487
749.5762	1102	5.610	0.834	7.521	0.996	0.746	0.0489
283.903	1101	10.126	0.691	8.241	0.513	1.229	0.0491

[00196] Table 18: Accurate mass features differing between 10 clinically diagnosed SP- MULTIPLE SCLEROSIS patients and 10 controls ($p < 0.05$).

550.4602	1204	3.472	0.470	12.653	1.135	0.274	6.38E-07
551.4646	1204	1.329	0.135	4.604	0.425	0.289	8.16E-07
578.4923	1204	11.751	0.881	46.817	5.096	0.251	2.37E-06
579.4958	1204	4.875	0.338	18.473	2.036	0.264	3.46E-06
580.5089	1204	4.485	0.410	14.745	1.530	0.304	4.29E-06
577.4795	1204	4.976	0.560	19.699	2.206	0.253	4.39E-06
576.4757	1204	12.835	1.389	53.107	6.179	0.242	5.45E-06
597.5068	1202	1.268	0.146	6.421	0.806	0.197	6.27E-06
597.5062	1204	12.819	1.444	77.858	10.610	0.165	9.69E-06
594.4848	1204	24.497	3.105	133.680	17.759	0.183	1.00E-05
598.5107	1204	3.097	0.299	17.223	2.358	0.180	1.27E-05
596.5012	1204	34.571	4.087	218.437	31.246	0.158	1.58E-05
595.4883	1204	10.287	1.282	53.490	7.294	0.192	1.59E-05
596.5053	1202	3.474	0.443	17.217	2.345	0.202	1.85E-05
616.4675	1201	1.137	0.137	5.287	0.713	0.215	2.01E-05
548.4438	1204	2.840	0.301	8.439	0.949	0.337	2.45E-05
563.5013	1204	2.423	0.366	5.674	0.449	0.427	2.55E-05
595.4928	1202	1.000	0.000	4.713	0.664	0.212	2.61E-05
581.5126	1204	2.311	0.187	5.718	0.588	0.404	3.07E-05
568.4723	1204	3.894	0.425	16.023	2.182	0.249	3.85E-05
558.4649	1204	14.333	1.960	47.412	5.855	0.302	4.31E-05
552.4784	1204	3.168	0.473	8.828	0.955	0.359	4.76E-05
493.385	1204	1.314	0.162	3.646	0.454	0.360	0.0001
508.3782	1204	2.236	0.288	6.212	0.525	0.429	0.0001
510.3937	1204	2.540	0.206	6.379	0.706	0.398	0.0001
522.4313	1204	6.350	0.553	17.283	2.060	0.368	0.0001
523.4337	1204	2.286	0.315	5.980	0.630	0.382	0.0001
534.4645	1204	1.957	0.270	3.828	0.264	0.511	0.0001
559.4688	1204	5.792	0.817	17.767	2.149	0.326	0.0001
562.4989	1204	6.687	0.835	14.101	1.150	0.474	0.0001
566.454	1204	2.751	0.317	9.505	1.345	0.289	0.0001
576.4765	1202	1.498	0.174	4.719	0.648	0.317	0.0001
594.4875	1202	2.493	0.219	10.648	1.606	0.234	0.0001
440.3526	1204	1.971	0.293	4.720	0.499	0.418	0.0002
446.341	1204	5.988	1.031	14.861	1.583	0.403	0.0002
448.3562	1204	5.280	0.797	11.199	0.980	0.471	0.0002
462.3346	1204	1.743	0.330	3.910	0.341	0.446	0.0002
469.3863	1204	2.014	0.324	5.088	0.570	0.396	0.0002
480.3473	1204	1.744	0.295	3.644	0.278	0.479	0.0002
492.3832	1204	3.984	0.530	10.155	1.249	0.392	0.0002
494.3968	1204	5.707	0.583	15.601	2.032	0.366	0.0002
502.4054	1204	2.918	0.427	7.569	0.875	0.386	0.0002
524.4448	1204	2.616	0.262	5.893	0.649	0.444	0.0002
532.4503	1204	2.510	0.310	5.587	0.587	0.449	0.0002
560.4821	1204	11.764	1.633	33.134	4.258	0.355	0.0002
561.4863	1204	4.978	0.592	13.298	1.719	0.374	0.0002
569.4769	1204	1.777	0.290	6.363	0.957	0.279	0.0002
610.482	1204	3.412	0.252	9.034	1.167	0.378	0.0002
466.3656	1204	6.322	0.835	14.385	1.618	0.439	0.0003
496.4157	1204	4.466	0.397	10.900	1.400	0.410	0.0003
513.4116	1204	2.094	0.261	6.272	0.900	0.334	0.0003
520.4131	1204	6.912	0.762	17.327	2.230	0.399	0.0003
540.4387	1204	11.594	1.267	37.542	5.765	0.309	0.0003
558.4663	1202	1.658	0.286	5.209	0.754	0.318	0.0003
570.4903	1204	1.413	0.178	5.019	0.800	0.282	0.0003
574.4594	1204	9.005	1.272	39.939	6.897	0.225	0.0003
618.4834	1201	1.781	0.291	7.306	1.195	0.244	0.0003
464.3524	1204	4.094	0.584	9.801	1.187	0.418	0.0004

484.3788	1204	3.708	0.358	8.598	1.058	0.431	0.0004
495.4018	1204	2.006	0.239	5.043	0.657	0.398	0.0004
538.4267	1204	10.128	1.453	29.656	4.218	0.342	0.0004
541.4422	1204	4.227	0.390	12.395	1.858	0.341	0.0004
482.3604	1204	2.285	0.295	5.003	0.565	0.457	0.0005
490.3676	1204	2.790	0.429	6.947	0.872	0.402	0.0005
504.4188	1204	3.750	0.382	8.301	0.997	0.452	0.0005
512.4079	1204	5.922	0.643	18.085	2.823	0.327	0.0005
590.4585	1204	3.711	0.299	11.905	1.918	0.312	0.0005
530.4379	1204	2.841	0.360	6.143	0.710	0.462	0.0006
572.4455	1204	1.967	0.311	5.226	0.721	0.376	0.0006
575.4628	1204	3.809	0.596	15.394	2.708	0.247	0.0006
468.384	1204	8.098	1.139	18.856	2.365	0.429	0.0007
592.4717	1204	9.648	1.214	44.121	8.330	0.219	0.0007
450.3729	1204	5.207	0.618	10.174	1.075	0.512	0.0008
557.4527	1204	2.309	0.458	5.268	0.578	0.438	0.0008
447.3433	1204	2.229	0.361	4.668	0.498	0.478	0.0009
474.3731	1204	2.438	0.304	5.290	0.652	0.461	0.0009
621.4188	1204	2.734	0.263	6.323	0.869	0.432	0.0009
556.4497	1204	5.393	1.130	12.738	1.469	0.423	0.0009
593.4736	1204	4.158	0.455	18.437	3.549	0.226	0.0009
478.4044	1204	2.134	0.241	4.281	0.489	0.498	0.001
564.4396	1204	1.437	0.188	3.702	0.547	0.386	0.001
662.4267	1204	4.620	0.544	7.882	0.625	0.586	0.001
438.3354	1204	1.835	0.244	3.452	0.337	0.532	0.0011
462.3716	1204	1.755	0.262	3.104	0.228	0.565	0.0011
467.3711	1204	2.039	0.329	4.570	0.573	0.446	0.0012
537.4142	1204	1.601	0.279	4.530	0.714	0.353	0.0013
539.4274	1204	2.869	0.646	9.921	1.707	0.299	0.0013
546.4298	1204	1.760	0.276	4.749	0.733	0.371	0.0013
634.3951	1204	4.180	0.858	9.781	1.192	0.428	0.0013
327.0307	1204	5.519	0.363	7.555	0.397	0.731	0.0014
518.3989	1204	4.953	0.798	13.176	2.085	0.376	0.0016
564.513	1204	2.491	0.383	5.132	0.600	0.485	0.0016
591.4614	1204	1.560	0.188	4.512	0.780	0.346	0.0017
780.5303	1204	7.051	0.543	10.487	0.764	0.672	0.0018
536.41	1204	4.849	0.589	11.044	1.618	0.439	0.002
476.3869	1204	2.592	0.310	5.647	0.807	0.459	0.0024
452.3868	1204	2.151	0.258	3.957	0.448	0.544	0.0026
684.6037	1204	4.873	0.976	1.387	0.258	3.513	0.0028
786.51	1204	28.214	4.540	47.124	3.077	0.599	0.0029
702.4175	1204	6.417	0.863	10.462	0.810	0.613	0.0031
1227.091	1203	20.780	2.685	7.389	2.870	2.812	0.0031
574.4635	1202	1.186	0.124	3.815	0.767	0.311	0.0033
590.4964	1204	2.971	0.514	5.910	0.704	0.503	0.0034
872.6715	1204	3.959	0.500	1.847	0.382	2.143	0.0035
534.3912	1204	2.266	0.293	4.670	0.675	0.483	0.0042
519.3998	1204	1.795	0.343	4.507	0.764	0.398	0.0045
566.3431	1102	4.903	0.724	8.188	0.735	0.599	0.0052
1253.124	1203	8.027	1.087	3.768	0.790	2.130	0.0053
1227.109	1203	1.000	0.000	6.463	1.747	0.155	0.0058
325.0805	1203	6.643	0.259	4.437	0.662	1.497	0.0061
565.3391	1102	16.189	2.611	27.761	2.683	0.583	0.0063
612.4994	1204	3.321	0.458	5.636	0.594	0.589	0.0064
428.3653	1201	8.555	1.414	4.218	0.485	2.028	0.0095
477.3218	1201	5.669	0.639	3.353	0.484	1.691	0.0098
786.5408	1204	69.825	12.284	114.325	9.307	0.611	0.0098
516.3324	1204	5.629	0.723	8.488	0.678	0.653	0.0099
787.5452	1204	32.885	5.441	52.366	4.008	0.828	0.0099
542.3447	1102	4.631	0.544	7.432	0.812	0.623	0.0103
716.4323	1204	9.262	1.168	13.477	0.896	0.687	0.0103
700.4371	1204	5.142	0.606	8.366	0.953	0.615	0.0105
780.4907	1204	8.547	1.272	13.226	1.042	0.646	0.0107
738.5448	1102	2.946	0.273	4.660	0.546	0.632	0.0116
758.4785	1204	48.142	6.547	70.822	4.748	0.680	0.0117

541.3415	1102	15.508	1.952	25.108	2.840	0.618	0.0122
832.5211	1204	4.586	0.642	6.726	0.423	0.682	0.0122
860.7729	1203	15.620	1.641	9.251	1.594	1.688	0.0123
772.5265	1204	80.234	11.146	116.817	7.132	0.687	0.0128
503.3194	1204	6.562	0.772	9.407	0.685	0.698	0.013
531.312	1102	4.205	0.597	6.308	0.485	0.667	0.0137
1226.078	1203	14.056	3.132	4.162	1.855	3.377	0.0141
1251.104	1203	7.607	1.295	3.433	0.835	2.216	0.0144
264.9759	1204	5.745	0.400	7.002	0.237	0.820	0.0145
569.3687	1102	4.709	0.805	7.253	0.491	0.649	0.0147
136.1258	1203	5.487	0.328	4.195	0.349	1.308	0.0148
468.3577	1201	4.938	0.656	2.794	0.451	1.767	0.0149
150.1413	1203	5.000	0.275	3.627	0.433	1.379	0.0154
610.5204	1204	6.706	0.995	16.481	3.554	0.407	0.0163
730.4535	1204	19.286	2.509	27.501	1.823	0.701	0.0163
1019.384	1102	4.944	0.506	6.898	0.538	0.717	0.0165
808.5264	1204	14.792	2.579	22.126	1.047	0.669	0.0168
812.6122	1201	6.364	0.891	3.185	0.813	1.998	0.0168
723.6395	1204	8.368	0.522	6.496	0.484	1.288	0.017
808.5225	1204	28.835	5.040	42.966	1.954	0.671	0.0176
748.5722	1102	9.684	0.818	15.985	2.287	0.606	0.0183
722.5244	1204	6.342	0.668	9.102	0.844	0.697	0.0196
368.1655	1102	1.185	0.185	3.565	0.913	0.332	0.0199
749.5762	1102	4.605	0.360	7.438	1.050	0.619	0.0201
828.5477	1201	7.432	0.648	5.374	0.484	1.383	0.0203
861.5265	1102	3.988	0.528	6.978	1.050	0.572	0.0203
170.11	1203	4.805	0.180	3.587	0.444	1.340	0.0204
506.4338	1204	2.346	0.274	3.523	0.376	0.666	0.021
728.5605	1204	6.704	0.883	10.409	1.177	0.644	0.0215
897.5729	1102	3.972	0.484	7.232	1.201	0.549	0.0215
859.7662	1203	43.601	5.294	27.020	3.953	1.614	0.0219
794.5126	1204	37.888	4.676	51.378	2.752	0.737	0.023
754.5724	1204	6.394	0.637	8.581	0.615	0.745	0.0238
858.7607	1203	96.197	12.594	58.519	8.619	1.644	0.0238
602.5287	1203	442.732	42.991	298.697	39.820	1.482	0.0243
793.5381	1102	9.315	1.086	16.255	2.606	0.573	0.0243
997.3968	1102	4.647	0.582	6.752	0.635	0.688	0.0251
886.5582	1102	5.308	0.511	8.358	1.146	0.635	0.0258
759.5145	1204	82.727	18.578	135.933	11.699	0.609	0.0261
603.532	1203	185.229	18.521	124.877	16.890	1.483	0.027
899.5871	1102	5.197	0.456	9.083	1.551	0.572	0.0272
602.3165	1204	23.972	3.305	34.085	2.678	0.703	0.0287
667.3547	1102	8.637	1.277	12.512	1.016	0.690	0.0289
194.0803	1203	4.092	0.903	10.360	2.484	0.395	0.0291
590.5287	1203	14.508	1.311	9.764	1.535	1.486	0.0304
784.5238	1204	61.252	12.412	93.059	5.519	0.658	0.0309
770.5108	1204	66.669	9.999	91.610	3.726	0.728	0.0312
134.11	1203	12.200	0.544	9.502	1.019	1.284	0.0313
833.5931	1201	13.307	1.528	8.766	1.203	1.518	0.0313
148.1257	1203	8.105	0.361	6.336	0.670	1.279	0.0321
781.497	1204	9.800	2.175	15.526	1.165	0.631	0.0322
835.6094	1201	5.359	0.712	3.476	0.392	1.542	0.0326
555.3102	1102	4.391	0.655	6.728	0.771	0.653	0.0329
729.5727	1204	6.021	0.638	8.508	0.871	0.708	0.0334
617.0921	1204	207.196	23.604	268.008	11.876	0.773	0.0335
576.51	1203	705.744	95.161	426.474	76.051	1.655	0.0341
788.5549	1204	13.572	1.893	19.042	1.460	0.713	0.0344
162.1412	1203	6.114	0.257	5.035	0.396	1.214	0.0346
758.5089	1204	123.531	26.082	193.464	16.177	0.639	0.0351
766.4759	1204	8.570	0.885	11.487	0.925	0.746	0.0351
779.5828	1204	15.783	1.495	21.637	2.090	0.729	0.0351
821.5714	1102	5.465	0.646	9.535	1.667	0.573	0.0352
888.5121	1204	7.012	1.000	11.029	1.455	0.636	0.0354
872.5557	1102	4.145	0.546	7.038	1.156	0.589	0.0362
827.5446	1201	14.730	1.373	10.914	0.983	1.350	0.0364

742.4745	1204	7.998	1.054	10.693	0.572	0.748	0.0375
378.9906	1204	3.444	0.208	4.104	0.207	0.839	0.0378
541.3141	1101	4.561	0.406	7.761	1.369	0.588	0.0379
785.5287	1204	31.042	6.421	46.658	2.719	0.665	0.038
830.7332	1203	24.238	3.418	14.491	2.708	1.673	0.0383
1226.099	1203	2.984	1.394	8.082	1.804	0.369	0.0383
184.1255	1203	6.461	0.228	5.272	0.481	1.226	0.0384
830.5881	1102	7.945	0.851	13.643	2.417	0.582	0.0392
727.5554	1204	13.082	1.418	20.151	2.857	0.649	0.0398
858.6843	1102	1.895	0.602	5.452	1.489	0.348	0.0399
474.2846	1204	6.895	0.821	9.629	0.862	0.726	0.04
488.2996	1204	4.260	0.569	5.842	0.432	0.729	0.04
829.5851	1102	17.810	1.960	31.309	5.797	0.569	0.0406
780.5872	1204	8.415	0.637	11.092	1.041	0.759	0.0410
519.3322	1101	4.242	0.386	6.534	0.971	0.649	0.0417
832.6022	1102	10.376	0.935	17.404	3.081	0.596	0.0425
172.1255	1203	6.804	0.268	5.742	0.462	1.202	0.043
699.5206	1204	4.894	0.732	8.253	1.358	0.593	0.043
577.5134	1203	257.861	34.886	159.577	28.694	1.616	0.0431
720.5081	1204	3.975	0.701	6.114	0.689	0.650	0.0431
281.2447	1204	20.974	3.610	29.938	1.991	0.701	0.0433
760.5231	1204	58.934	11.826	88.317	6.590	0.667	0.0436
744.4942	1204	127.238	21.145	181.030	13.054	0.703	0.0441
379.2536	1204	1.941	0.435	3.035	0.280	0.640	0.0448
633.3232	1102	3.047	0.515	4.505	0.439	0.676	0.0451
804.5715	1102	19.957	2.117	36.054	7.169	0.554	0.0451
591.5321	1203	5.846	0.612	3.933	0.649	1.486	0.0458
832.7499	1203	13.911	2.292	7.895	1.626	1.782	0.0462
461.2707	1204	2.330	0.223	2.933	0.174	0.794	0.0467
302.2255	1204	3.093	0.424	4.225	0.320	0.732	0.0471
198.1411	1203	5.029	0.251	4.192	0.304	1.200	0.048
280.2413	1204	107.053	18.645	152.351	10.396	0.703	0.048
803.5683	1102	49.513	5.444	91.298	18.924	0.542	0.048
794.5419	1102	4.724	0.480	7.467	1.204	0.633	0.0486
558.3761	1204	4.723	0.502	6.382	0.605	0.740	0.049
777.5878	1204	13.748	1.328	18.871	2.037	0.729	0.0494
834.5963	1201	6.953	0.783	4.966	0.528	1.400	0.0497

[00197] Table 19: Accurate mass features differing between 10 clinically diagnosed RR- MULTIPLE SCLEROSIS patients and SP- MULTIPLE SCLEROSIS controls (p<0.05).

448.3562	1204	12.828	0.687	4.778	0.539	2.685	2.71E-08
467.3711	1204	5.915	0.391	1.851	0.225	3.196	4.32E-08
468.3656	1204	17.861	1.282	5.789	0.482	3.085	6.87E-08
484.3788	1204	10.785	0.702	3.694	0.316	2.920	7.22E-08
450.3729	1204	11.105	0.649	4.587	0.398	2.421	9.95E-08
580.5089	1204	18.694	1.434	4.911	0.628	3.807	1.10E-07
578.4923	1204	56.509	4.849	12.547	1.621	4.504	2.10E-07
579.4958	1204	22.270	1.913	5.359	0.654	4.156	3.04E-07
452.3868	1204	4.802	0.284	1.901	0.220	2.526	3.12E-07
469.3863	1204	6.668	0.510	2.064	0.276	3.231	3.69E-07
494.3968	1204	19.451	1.654	5.659	0.537	3.437	3.84E-07
468.384	1204	25.447	2.013	8.057	0.922	3.158	4.29E-07
581.5126	1204	7.479	0.640	2.018	0.219	3.706	4.72E-07
508.3782	1204	6.735	0.505	2.211	0.220	3.046	5.07E-07
618.4834	1201	10.132	1.055	2.065	0.335	4.907	6.11E-07
510.3937	1204	8.743	0.780	2.767	0.254	3.160	6.55E-07
495.4018	1204	6.117	0.522	1.938	0.194	3.156	6.94E-07
513.4116	1204	7.003	0.525	2.222	0.332	3.152	7.35E-07
596.5053	1202	23.532	2.544	4.292	0.670	5.483	7.38E-07
598.5107	1204	23.078	2.490	3.714	0.543	6.214	7.79E-07
597.5068	1202	8.982	0.941	1.674	0.317	5.366	7.99E-07
522.4313	1204	20.427	1.832	6.314	0.627	3.235	8.26E-07
568.4723	1204	19.787	2.282	4.028	0.450	4.912	1.01E-06
569.4769	1204	7.849	0.822	1.841	0.247	4.263	1.16E-06
597.5082	1204	108.868	12.859	15.045	2.165	7.103	1.20E-06
537.4142	1204	5.533	0.453	1.639	0.226	3.376	1.33E-06
610.482	1204	13.296	1.385	3.518	0.246	3.779	1.37E-06
551.4646	1204	5.117	0.491	1.286	0.215	3.979	1.43E-06
596.5012	1204	302.332	35.546	41.055	6.198	7.364	1.50E-06
512.4079	1204	20.539	1.811	6.306	0.840	3.257	1.90E-06
446.341	1204	15.895	1.428	5.248	0.644	3.029	2.20E-06
550.4602	1204	13.875	1.423	3.330	0.563	4.167	2.31E-06
464.3524	1204	11.766	1.096	3.915	0.357	3.005	2.60E-06
492.3832	1204	12.904	1.371	3.723	0.342	3.466	2.75E-06
595.4883	1204	74.421	9.813	11.967	1.663	6.219	3.22E-06
590.4585	1204	14.400	1.343	4.037	0.398	3.567	3.23E-06
577.4795	1204	24.484	2.935	5.137	0.663	4.766	3.38E-06
536.41	1204	14.183	1.310	4.859	0.376	2.919	3.48E-06
594.4848	1204	187.278	25.540	29.026	4.176	6.452	3.49E-06
523.4337	1204	6.812	0.627	2.209	0.276	3.084	3.70E-06
576.4757	1204	66.906	8.279	13.283	1.720	5.037	3.77E-06
524.4448	1204	6.702	0.649	2.252	0.265	2.976	3.94E-06
440.3526	1204	5.499	0.452	2.035	0.245	2.702	4.58E-06
482.3604	1204	6.620	0.663	2.324	0.131	2.849	4.58E-06
616.4675	1201	7.494	0.826	1.454	0.257	5.154	4.99E-06
594.4875	1202	14.552	1.781	3.055	0.390	4.763	5.09E-06
476.3869	1204	5.204	0.291	2.461	0.311	2.115	5.20E-06
534.3912	1204	6.113	0.560	2.085	0.205	2.932	5.42E-06
520.4131	1204	21.750	2.595	6.355	0.508	3.423	5.98E-06
566.454	1204	12.984	1.754	2.754	0.294	4.715	6.04E-06
570.4903	1204	5.738	0.678	1.315	0.165	4.363	6.16E-06
541.4422	1204	21.125	2.994	4.380	0.457	4.823	6.30E-06
496.4157	1204	16.461	1.840	4.566	0.456	3.605	6.77E-06
540.4387	1204	62.261	8.879	12.111	1.401	5.141	7.87E-06
538.4257	1204	44.839	5.751	10.268	1.135	4.367	7.94E-06
518.3969	1204	15.673	1.778	4.418	0.460	3.548	8.09E-06

462.3346	1204	4.382	0.382	1.547	0.236	2.833	8.53E-06
595.4928	1202	6.217	0.788	1.316	0.215	4.724	1.07E-05
519.3998	1204	5.809	0.641	1.523	0.270	3.814	1.10E-05
438.3354	1204	4.529	0.392	1.772	0.224	2.556	1.15E-05
591.4614	1204	5.115	0.467	1.890	0.269	2.706	1.28E-05
521.4188	1204	7.695	0.847	2.676	0.199	2.876	1.47E-05
552.4784	1204	10.058	1.199	2.886	0.350	3.485	1.70E-05
474.3731	1204	5.983	0.690	1.993	0.275	3.002	1.71E-05
548.4438	1204	10.315	1.447	2.758	0.308	3.740	1.83E-05
564.4396	1204	4.615	0.499	1.374	0.202	3.359	1.95E-05
447.3433	1204	4.979	0.485	1.984	0.273	2.535	2.69E-05
592.4717	1204	49.484	6.153	11.401	1.794	4.340	2.92E-05
480.3473	1204	4.690	0.393	2.037	0.277	2.302	3.09E-05
493.385	1204	4.379	0.513	1.374	0.155	3.187	3.09E-05
593.4736	1204	20.853	2.576	5.040	0.818	4.138	3.10E-05
578.4765	1202	5.839	0.693	1.624	0.253	3.595	4.93E-05
502.4054	1204	7.483	0.761	2.843	0.469	2.632	0.0001
504.4188	1204	7.357	0.554	3.472	0.503	2.119	0.0001
532.4503	1204	5.403	0.420	2.505	0.389	2.157	0.0001
534.4645	1204	4.466	0.372	1.859	0.315	2.402	0.0001
539.4274	1204	14.474	2.478	3.226	0.548	4.487	0.0001
563.5013	1204	5.244	0.403	2.149	0.406	2.440	0.0001
572.4455	1204	5.719	0.661	1.843	0.336	3.103	0.0001
327.0307	1204	8.810	0.616	5.574	0.325	1.581	0.0002
490.3676	1204	7.893	1.158	2.550	0.235	3.095	0.0002
574.4594	1204	41.888	6.536	9.914	1.605	4.225	0.0002
558.4649	1204	49.295	7.166	14.427	2.763	3.417	0.0003
559.4688	1204	18.364	2.624	5.551	0.957	3.308	0.0003
562.4989	1204	14.237	1.259	6.489	1.040	2.194	0.0003
575.4628	1204	16.105	2.497	4.318	0.711	3.730	0.0003
560.478	1203	9.998	0.752	6.482	0.540	1.542	0.0004
478.4044	1204	4.058	0.374	2.096	0.311	1.936	0.0006
530.4379	1204	6.139	0.698	2.640	0.429	2.325	0.0006
546.4298	1204	5.422	0.846	1.797	0.298	3.017	0.0006
557.4527	1204	4.841	0.521	2.044	0.378	2.271	0.0007
558.4663	1202	5.148	0.667	1.907	0.395	2.700	0.0007
612.4994	1204	8.097	0.941	3.294	0.473	2.458	0.0007
506.4338	1204	4.086	0.311	2.302	0.323	1.775	0.001
556.4497	1204	11.083	1.357	4.554	0.976	2.434	0.001
574.4635	1202	3.793	0.533	1.474	0.208	2.573	0.001
560.4821	1204	27.410	2.617	12.279	2.740	2.232	0.0011
561.4863	1204	11.089	1.073	5.196	1.045	2.134	0.0013
462.3716	1204	3.329	0.305	1.759	0.265	1.893	0.0031
856.6045	1102	7.793	1.319	6.545	0.691	1.191	0.0036
854.5884	1102	4.407	0.729	3.471	0.311	1.270	0.0042
634.3951	1204	11.113	2.036	4.104	0.836	2.708	0.0046
855.6009	1102	15.163	2.457	12.786	1.331	1.186	0.0047
519.3322	1101	5.790	0.368	4.417	0.460	1.311	0.0048
564.513	1204	5.141	0.573	2.558	0.458	2.010	0.005
611.3724	1201	5.966	0.607	3.351	0.968	1.780	0.0066
895.5575	1102	6.457	1.242	5.002	0.523	1.291	0.0109
853.5848	1102	8.298	1.332	6.301	0.585	1.317	0.0113
541.3141	1101	6.800	0.568	4.815	0.562	1.412	0.0121
610.5204	1204	18.173	4.912	7.931	1.596	2.291	0.0141
662.4267	1204	7.469	1.139	4.182	0.445	1.784	0.016
570.3766	1201	2.962	0.303	1.895	0.477	1.563	0.0171
827.5695	1102	26.452	5.262	20.283	2.505	1.304	0.0172
886.7804	1203	10.774	2.558	6.148	2.114	1.752	0.0192
546.3413	1204	4.855	0.708	2.749	0.410	1.766	0.0195
610.3691	1201	16.117	1.748	10.304	3.000	1.564	0.0223
378.9906	1204	4.292	0.313	3.160	0.309	1.358	0.0233
570.376	1203	3.876	0.750	1.590	0.318	2.438	0.0239
162.1412	1203	6.102	0.386	6.220	0.259	0.981	0.0244
810.5967	1201	17.404	2.577	28.950	3.444	0.581	0.0247
835.6094	1201	3.249	0.425	5.886	0.673	0.552	0.0251

785.4799	1204	10.252	2.193	21.576	4.246	0.475	0.026
606.4872	1204	7.177	1.330	4.058	0.558	1.769	0.0281
639.4037	1201	3.594	0.326	2.974	0.457	1.208	0.0281
797.5257	1204	47.197	3.685	69.761	11.106	0.677	0.0296
264.9759	1204	7.268	0.440	5.734	0.361	1.268	0.0311
809.5934	1201	34.411	5.094	58.757	6.906	0.586	0.0313
828.5732	1102	12.694	2.293	9.840	1.101	1.317	0.0337
744.55	1203	9.050	1.592	6.269	1.174	1.444	0.0347
831.5758	1101	6.461	0.619	7.739	0.919	0.835	0.035
590.4964	1204	5.436	1.167	3.151	0.719	1.725	0.0382
795.5083	1204	30.689	2.188	40.150	5.184	0.764	0.039
769.4929	1204	65.355	6.025	88.784	12.171	0.736	0.0397
743.5475	1203	19.547	3.139	14.226	2.394	1.374	0.041
181.9806	1102	5.873	0.891	2.752	0.581	2.134	0.0422
748.5722	1102	12.839	2.174	9.039	0.655	1.420	0.0437
200.1566	1203	8.329	0.174	7.692	0.276	1.083	0.044
729.5727	1204	9.886	1.145	5.989	0.859	1.651	0.044
638.4003	1201	9.386	1.059	7.767	1.259	1.208	0.0446
832.6027	1202	26.976	3.806	14.960	2.367	1.803	0.0484
160.1256	1203	7.941	0.470	8.253	0.387	0.962	0.0489
566.3433	1202	18.924	1.569	15.150	0.718	1.249	0.0497

[00198] Table 20: Accurate mass features differing between 10 RR- MULTIPLE SCLEROSIS patients transitioning to SP- MULTIPLE SCLEROSIS and 10 clinically diagnosed RR- MULTIPLE SCLEROSIS patients ($p < 0.05$).

580.5089	1204	4.522	0.282	18.548	1.208	0.244	7.78E-09
450.3729	1204	4.558	0.182	11.434	0.666	0.399	1.23E-08
578.4923	1204	12.054	1.113	56.683	4.013	0.213	1.66E-08
679.4958	1204	4.683	0.442	22.165	1.599	0.211	2.08E-08
448.3562	1204	5.180	0.347	13.087	0.778	0.396	3.80E-08
495.4018	1204	1.616	0.212	6.267	0.490	0.258	1.48E-07
581.5126	1204	1.907	0.240	7.358	0.552	0.259	1.68E-07
577.4795	1204	5.058	0.413	25.693	2.730	0.197	4.26E-07
550.4602	1204	3.541	0.391	14.589	1.385	0.243	4.63E-07
576.4757	1204	12.734	1.088	70.329	7.812	0.181	5.18E-07
484.3788	1204	3.788	0.327	10.428	0.737	0.363	5.84E-07
466.3656	1204	5.919	0.636	17.958	1.455	0.330	8.21E-07
551.4646	1204	1.432	0.191	5.465	0.509	0.262	9.00E-07
523.4337	1204	1.899	0.208	7.018	0.643	0.271	9.27E-07
561.4863	1204	5.053	0.475	12.026	0.798	0.420	9.30E-07
560.4821	1204	12.218	1.146	29.622	1.986	0.412	9.53E-07
494.3968	1204	5.859	0.504	19.493	1.724	0.301	9.75E-07
569.4769	1204	1.544	0.195	7.995	0.867	0.193	9.80E-07
512.4079	1204	5.855	0.647	19.615	1.654	0.288	1.05E-06
522.4313	1204	6.534	0.460	21.358	1.990	0.306	1.37E-06
598.5107	1204	3.116	0.474	22.331	2.356	0.140	1.53E-06
568.4723	1204	3.894	0.304	19.878	2.356	0.196	1.71E-06
504.4188	1204	3.450	0.155	7.928	0.540	0.435	1.91E-06
562.4989	1204	7.498	0.721	15.447	0.968	0.485	2.12E-06
521.4168	1204	2.171	0.302	8.244	0.832	0.263	2.22E-06
510.3937	1204	2.563	0.308	8.533	0.806	0.300	2.37E-06
595.4883	1204	11.511	1.238	75.048	9.692	0.153	2.53E-06
468.384	1204	8.388	0.798	25.639	2.296	0.327	2.54E-06
594.4848	1204	28.935	3.071	189.742	25.019	0.152	2.70E-06
478.3869	1204	2.471	0.222	5.512	0.344	0.448	2.71E-06
467.3711	1204	2.178	0.201	5.899	0.500	0.369	3.30E-06
559.4688	1204	5.645	0.677	20.172	2.260	0.280	3.64E-06
596.5012	1204	36.743	5.025	292.674	33.424	0.126	3.78E-06
548.4438	1204	2.435	0.275	10.928	1.436	0.223	4.42E-06
597.5062	1204	14.057	1.806	102.408	12.272	0.137	4.58E-06
452.3868	1204	1.961	0.284	4.714	0.297	0.416	4.71E-06
513.4116	1204	2.091	0.224	6.369	0.543	0.328	4.96E-06
469.3863	1204	2.527	0.300	6.732	0.561	0.375	5.15E-06
536.41	1204	4.313	0.401	14.234	1.396	0.303	5.27E-06
520.4131	1204	6.574	0.399	22.631	2.629	0.290	5.54E-06
496.4157	1204	4.416	0.420	16.630	1.813	0.266	5.67E-06
524.4448	1204	2.358	0.254	6.952	0.680	0.339	5.81E-06
558.4649	1204	16.025	1.572	53.745	6.329	0.280	6.09E-06
532.4503	1204	2.369	0.411	5.793	0.358	0.409	6.62E-06
610.482	1204	3.393	0.359	13.517	1.487	0.251	7.29E-06
596.5053	1202	3.960	0.749	21.759	2.655	0.182	7.39E-06
502.4054	1204	3.334	0.240	8.085	0.715	0.412	8.26E-06
492.3832	1204	4.090	0.389	13.443	1.453	0.304	9.14E-06
519.3998	1204	1.708	0.258	6.099	0.655	0.280	9.59E-06
597.5068	1202	1.593	0.313	8.239	1.011	0.193	1.17E-05
508.3782	1204	2.261	0.381	6.852	0.600	0.330	1.24E-05
534.3912	1204	1.750	0.274	6.155	0.599	0.284	1.41E-05
440.3526	1204	2.084	0.264	5.422	0.458	0.384	1.42E-05
541.4422	1204	3.411	0.379	20.925	3.127	0.163	1.44E-05
572.4465	1204	1.914	0.316	6.467	0.678	0.296	1.47E-05

590.4585	1204	4.211	0.410	15.828	1.796	0.266	1.60E-05
518.3969	1204	5.017	0.475	16.609	1.870	0.302	1.61E-05
557.4527	1204	2.731	0.271	5.357	0.471	0.510	1.62E-05
566.454	1204	3.340	0.340	13.368	1.792	0.250	1.63E-05
552.4784	1204	3.072	0.232	10.252	1.204	0.300	1.65E-05
482.3604	1204	2.387	0.293	6.828	0.684	0.360	1.72E-05
540.4387	1204	11.181	0.938	62.043	9.171	0.180	1.76E-05
594.4875	1202	2.688	0.532	13.976	1.893	0.192	1.78E-05
438.3354	1204	1.732	0.203	4.427	0.393	0.391	1.83E-05
464.3524	1204	4.499	0.474	12.320	1.229	0.365	2.20E-05
480.3473	1204	2.129	0.160	5.005	0.415	0.425	2.20E-05
537.4142	1204	1.910	0.208	5.646	0.543	0.338	2.22E-05
447.3433	1204	2.234	0.247	5.264	0.490	0.424	2.81E-05
563.5013	1204	3.247	0.289	5.586	0.321	0.581	3.14E-05
570.4903	1204	1.448	0.184	5.518	0.681	0.262	3.18E-05
618.4834	1201	1.560	0.468	8.422	1.166	0.185	3.74E-05
595.4928	1202	1.284	0.199	6.004	0.829	0.214	4.03E-05
591.4614	1204	1.801	0.306	5.617	0.656	0.321	4.23E-05
474.3731	1204	2.347	0.232	6.267	0.722	0.375	4.26E-05
478.4044	1204	2.307	0.091	4.386	0.373	0.526	4.46E-05
538.4257	1204	10.245	0.848	43.969	6.163	0.233	4.48E-05
446.341	1204	7.325	0.584	17.044	1.670	0.430	4.74E-05
462.3346	1204	2.077	0.306	4.729	0.424	0.439	0.0001
493.385	1204	1.231	0.155	4.527	0.575	0.272	0.0001
506.4338	1204	1.982	0.248	4.196	0.336	0.472	0.0001
534.4645	1204	2.424	0.211	4.518	0.357	0.537	0.0001
556.4497	1204	6.643	0.820	12.691	1.212	0.523	0.0001
574.4594	1204	9.293	1.092	47.951	7.271	0.194	0.0001
575.4628	1204	3.830	0.446	18.688	2.716	0.205	0.0001
576.4765	1202	1.606	0.237	5.861	0.728	0.274	0.0001
592.4717	1204	11.135	1.295	53.655	7.672	0.208	0.0001
593.4736	1204	4.494	0.518	22.565	3.209	0.199	0.0001
546.4298	1204	1.747	0.279	6.038	0.878	0.289	0.0002
558.4663	1202	1.899	0.327	5.605	0.822	0.339	0.0002
564.4396	1204	1.712	0.210	4.715	0.573	0.363	0.0002
616.4675	1201	1.703	0.290	6.509	0.964	0.262	0.0003
490.3676	1204	2.714	0.463	8.438	1.224	0.322	0.0005
327.0307	1204	6.681	0.286	9.139	0.560	0.731	0.0008
530.4379	1204	3.253	0.519	6.685	0.667	0.487	0.0009
612.4994	1204	3.768	0.404	8.583	0.926	0.439	0.0009
574.4635	1202	1.366	0.201	4.237	0.626	0.322	0.001
590.4964	1204	3.310	0.570	5.836	1.084	0.567	0.0011
564.513	1204	2.805	0.511	5.374	0.449	0.522	0.0012
606.4872	1204	5.701	0.704	7.360	1.168	0.775	0.0014
610.5204	1204	10.685	2.101	18.659	4.768	0.573	0.0021
539.4274	1204	3.329	0.516	13.055	2.878	0.255	0.0027
804.5715	1102	23.445	5.126	46.944	8.322	0.499	0.0042
871.5526	1102	10.664	2.211	20.033	3.436	0.532	0.0042
803.5683	1102	57.160	13.215	116.696	21.332	0.490	0.0068
733.6414	1204	22.549	2.616	15.527	1.847	1.452	0.0125
829.5851	1102	20.138	2.786	39.735	6.642	0.507	0.0131
872.5557	1102	5.121	1.092	9.357	1.681	0.547	0.0132
569.369	1202	17.157	1.105	13.563	0.716	1.265	0.0147
603.5305	1201	1.757	0.345	4.631	1.000	0.379	0.0154
899.5871	1102	6.673	1.249	11.918	1.901	0.560	0.0163
576.5115	1201	4.861	1.392	15.148	3.461	0.321	0.0165
804.5428	1201	1.324	0.133	2.470	0.401	0.536	0.0183
601.515	1201	1.735	0.300	3.771	0.676	0.460	0.0204
707.6248	1204	12.429	2.115	8.439	1.276	1.473	0.0207
859.7715	1201	2.313	0.462	6.009	1.199	0.385	0.0227
856.7627	1201	4.217	1.499	13.836	3.049	0.305	0.0229
602.5271	1201	3.644	0.789	9.846	2.283	0.370	0.0239
600.5115	1201	3.502	0.762	9.257	2.057	0.378	0.0245

577.5148	1201	2.335	0.557	6.227	1.401	0.375	0.0253
719.6222	1204	17.497	2.578	12.859	1.904	1.361	0.0258
734.6429	1204	12.854	1.562	9.336	1.071	1.377	0.0258
687.4916	1204	22.840	3.156	37.698	5.234	0.606	0.0267
757.5622	1101	39.148	3.393	48.113	5.667	0.814	0.0274
784.5809	1101	9.351	0.779	10.960	1.260	0.853	0.0284
296.2357	1204	9.614	0.724	11.283	0.721	0.852	0.0288
574.4958	1201	4.124	0.837	9.107	1.722	0.453	0.0288
634.3951	1204	6.136	1.289	10.597	1.635	0.579	0.0302
758.5656	1101	20.676	1.792	25.277	2.995	0.818	0.0306
830.5881	1102	9.353	1.250	17.793	2.882	0.626	0.0309
260.2137	1203	5.046	0.345	5.892	0.213	0.856	0.031
854.737	1201	3.193	0.952	9.115	1.923	0.350	0.031
462.3716	1204	2.409	0.361	3.400	0.244	0.709	0.0313
686.4879	1204	53.749	7.668	90.851	13.784	0.592	0.0338
239.939	1102	4.050	0.706	5.407	0.631	0.749	0.0347
611.3724	1201	3.100	0.337	4.945	0.611	0.627	0.0349
673.4765	1204	7.595	1.041	10.654	0.589	0.713	0.0355
721.6382	1204	21.966	3.380	15.879	2.227	1.383	0.0366
660.4958	1201	1.889	0.508	5.003	1.096	0.378	0.036
857.7557	1201	3.271	0.869	8.666	1.828	0.377	0.0362
897.5729	1102	5.114	0.690	9.533	1.628	0.536	0.0387
735.6554	1204	22.156	3.657	15.289	1.767	1.449	0.0388
712.5074	1204	21.951	1.994	34.704	4.558	0.633	0.042
438.2993	1204	1.840	0.322	1.128	0.128	1.631	0.043
830.5634	1201	6.463	0.700	4.807	0.558	1.344	0.044
834.5372	1204	16.169	1.535	13.200	1.532	1.225	0.044
705.6086	1204	7.811	1.337	5.498	0.728	1.421	0.0461
598.4959	1201	2.137	0.341	3.489	0.474	0.612	0.0492

[00199] Table 21: Accurate mass features differing between 10 RR- MULTIPLE SCLEROSIS patients transitioning to SP- MULTIPLE SCLEROSIS and 10 clinically diagnosed SP- MULTIPLE SCLEROSIS patients ($p < 0.05$).

761.529	1204	44.870	1.690	17.080	1.359	2.627	5.04E-10
760.5231	1204	109.100	5.078	36.556	3.639	2.984	2.17E-09
690.4843	1204	10.221	0.348	2.876	0.621	3.554	1.80E-08
758.5089	1204	220.433	11.793	77.867	7.929	2.831	1.86E-08
759.5145	1204	157.006	9.162	49.129	5.813	3.196	2.08E-08
784.5238	1204	110.088	6.582	37.629	3.285	2.926	2.20E-08
732.4929	1204	64.274	3.909	19.285	2.366	3.333	2.35E-08
742.4745	1204	11.943	0.592	5.829	0.243	2.049	3.27E-08
785.6287	1204	55.424	3.489	19.112	1.576	2.900	3.71E-08
812.5559	1204	27.252	1.631	10.799	0.817	2.524	7.82E-08
786.5408	1204	131.847	8.335	48.312	4.378	2.729	1.00E-07
787.5452	1204	60.356	3.628	23.335	2.087	2.586	1.08E-07
744.4942	1204	201.989	11.571	84.491	6.859	2.391	1.29E-07
733.501	1204	37.219	2.804	10.249	1.378	3.631	1.44E-07
731.4898	1204	139.168	10.350	37.560	5.674	3.705	1.55E-07
809.5264	1204	26.756	1.700	9.447	1.067	2.832	1.58E-07
770.5108	1204	106.431	6.784	45.101	2.511	2.360	1.74E-07
808.5225	1204	51.852	3.250	18.231	2.201	2.844	1.78E-07
734.508	1204	18.236	1.160	3.835	1.208	4.755	2.00E-07
788.5549	1204	22.912	1.296	10.186	0.813	2.249	2.60E-07
452.2536	1204	6.868	0.446	2.179	0.383	3.152	5.17E-07
780.4907	1204	16.058	1.092	6.003	0.661	2.675	5.41E-07
772.5265	1204	133.872	8.862	58.392	4.253	2.293	7.09E-07
757.5008	1204	61.528	5.970	15.151	1.849	4.061	1.03E-06
746.5118	1204	82.872	4.102	31.127	5.409	2.662	1.11E-06
810.5394	1204	89.171	6.119	33.859	4.327	2.634	1.34E-06
811.5436	1204	43.601	3.103	16.663	1.976	2.617	1.44E-06
836.5534	1204	10.314	0.657	4.772	0.395	2.161	1.68E-06
688.4658	1204	11.938	0.707	4.376	0.766	2.728	1.98E-06
794.5126	1204	60.457	4.186	26.308	2.449	2.298	2.32E-06
758.491	1204	31.271	2.881	9.848	1.070	3.175	2.38E-06
814.498	1204	13.596	1.071	3.697	0.945	3.677	3.28E-06
813.5617	1204	12.536	0.749	6.919	0.599	2.118	3.34E-06
781.497	1204	19.276	1.841	5.709	0.839	3.376	4.04E-06
779.4829	1204	12.474	1.000	3.574	0.849	3.490	4.20E-06
766.4759	1204	13.713	0.797	6.185	0.774	2.217	4.49E-06
783.5127	1204	85.018	7.493	29.020	3.903	2.930	4.83E-06
782.5084	1204	137.563	11.554	50.082	6.507	2.747	5.22E-06
718.4736	1204	7.136	0.524	2.585	0.481	2.782	8.51E-06
617.0921	1204	309.018	21.686	158.144	10.806	1.954	1.03E-05
712.4676	1204	10.933	0.771	4.225	0.746	2.588	1.20E-05
807.5103	1204	26.009	1.784	11.300	1.589	2.302	1.30E-05
716.4987	1204	25.836	1.753	13.642	1.011	1.894	1.58E-05
806.5068	1204	47.988	3.146	21.094	3.089	2.275	1.62E-05
755.4854	1204	48.306	5.655	12.351	2.199	3.911	1.76E-05
796.5278	1204	99.834	7.696	46.021	4.705	2.169	1.80E-05
816.5159	1204	9.394	0.509	5.772	0.339	1.628	1.98E-05
717.5011	1204	11.886	0.757	6.288	0.571	1.890	2.16E-05
379.2536	1204	3.809	0.419	1.217	0.160	3.130	2.37E-05
768.4944	1204	140.648	11.913	61.674	6.965	2.281	2.87E-05
835.5417	1204	10.463	0.567	5.827	0.576	1.796	3.32E-05
154.0035	1204	29.040	2.027	16.527	1.444	1.757	0.0001
306.2668	1204	9.934	0.611	5.961	0.492	1.667	0.0001
420.2651	1204	4.402	0.296	2.352	0.283	1.872	0.0001
712.5074	1204	27.892	2.096	14.835	1.259	1.880	0.0001
721.6382	1204	21.940	1.994	9.600	1.382	2.285	0.0001
815.5045	1204	10.228	0.990	3.882	0.751	2.635	0.0001

832.5211	1204	7.326	0.502	3.616	0.495	2.026	0.0001
834.5372	1204	18.398	1.170	9.245	1.179	1.990	0.0001
713.5097	1204	12.216	0.831	6.999	0.727	1.745	0.0002
740.4966	1204	21.462	1.753	12.268	0.652	1.749	0.0002
765.4894	1204	16.458	1.357	8.044	1.028	2.046	0.0002
780.5303	1204	10.368	0.680	5.616	0.727	1.846	0.0002
788.4794	1204	12.719	0.884	5.698	1.146	2.232	0.0002
872.6715	1204	1.743	0.360	4.635	0.476	0.376	0.0002
313.2702	1101	2.194	0.770	7.148	0.719	0.307	0.0003
714.5221	1204	35.796	2.655	19.874	2.112	1.801	0.0003
569.3687	1102	7.896	0.726	3.826	0.545	2.064	0.0004
690.5475	1204	6.802	0.612	3.252	0.496	2.092	0.0004
737.5045	1204	6.986	0.614	3.496	0.476	1.998	0.0004
789.5658	1204	11.139	0.523	6.630	0.841	1.680	0.0004
792.4954	1204	54.222	4.145	29.312	3.671	1.850	0.0004
686.4879	1204	64.699	6.152	34.789	3.217	1.860	0.0005
738.5185	1204	25.027	2.301	13.047	1.625	1.918	0.0008
757.5637	1204	13.400	0.954	7.554	0.966	1.774	0.0006
707.6248	1204	12.254	1.130	5.983	0.952	2.048	0.0007
736.5031	1204	14.979	1.418	7.958	0.879	1.882	0.0007
742.5142	1204	25.071	2.044	15.355	1.078	1.633	0.0007
886.5582	1102	9.158	1.042	4.660	0.365	1.965	0.0008
784.6228	1204	25.713	4.239	6.186	2.336	4.157	0.0009
820.5294	1204	23.887	1.921	13.544	1.608	1.764	0.0009
313.7724	1101	6.709	1.105	1.625	0.625	4.129	0.001
687.4916	1204	26.877	2.482	15.555	1.326	1.728	0.001
997.3968	1102	7.284	0.537	4.347	0.490	1.676	0.001
747.5121	1204	54.593	3.233	37.326	2.827	1.463	0.0011
735.6554	1204	20.623	2.274	10.165	1.429	2.029	0.0013
745.4938	1204	130.623	9.940	84.200	6.457	1.551	0.0013
495.3322	1201	3.373	0.314	5.878	0.549	0.574	0.0014
783.6174	1204	41.517	7.569	8.434	4.065	4.923	0.0014
688.5048	1204	26.182	2.869	14.168	1.295	1.848	0.0015
689.5083	1204	10.722	0.992	6.291	0.591	1.704	0.0015
771.5075	1204	73.947	6.085	47.202	3.498	1.567	0.0015
633.3232	1102	4.844	0.484	2.605	0.335	1.860	0.0016
773.5257	1204	93.493	6.310	60.586	5.795	1.543	0.0016
748.5722	1102	16.264	1.898	8.476	0.842	1.919	0.0017
770.569	1204	53.299	4.207	34.160	2.822	1.560	0.0017
812.6122	1201	3.134	0.525	6.820	0.801	0.460	0.0018
302.2255	1204	4.500	0.351	2.451	0.424	1.836	0.0022
1019.384	1102	7.026	0.536	4.416	0.461	1.591	0.0022
565.3391	1102	25.394	2.373	13.580	2.172	1.870	0.0023
715.4864	1204	18.244	1.938	10.740	0.788	1.699	0.0024
566.3431	1102	7.841	0.753	4.138	0.807	1.847	0.0025
794.5718	1204	28.883	2.201	18.210	1.958	1.586	0.0025
567.3547	1102	13.308	1.127	7.473	1.149	1.781	0.0026
833.5929	1101	4.412	0.340	3.011	0.199	1.465	0.0027
738.5448	1102	4.708	0.497	2.639	0.309	1.784	0.0028
719.6222	1204	16.138	1.622	9.274	1.107	1.740	0.0031
341.2443	1204	1.816	0.203	1.073	0.073	1.692	0.0032
795.5083	1204	45.741	3.585	29.214	3.139	1.586	0.0035
766.5153	1204	16.517	1.630	9.478	1.247	1.743	0.0037
714.4837	1204	38.123	4.503	22.200	1.520	1.717	0.0039
854.588	1202	12.123	1.922	5.379	0.618	2.254	0.004
722.5244	1204	8.482	0.797	4.920	0.683	1.724	0.0041
872.5557	1102	7.287	1.023	3.637	0.395	2.003	0.0041
541.3415	1102	24.790	2.463	14.323	1.902	1.731	0.0042
694.4953	1204	4.820	0.383	3.093	0.338	1.558	0.0042
749.5762	1102	7.246	0.867	4.089	0.392	1.772	0.0043
747.5781	1204	15.794	1.109	10.817	0.989	1.460	0.0045
854.6884	1102	6.310	0.977	3.038	0.217	2.077	0.0045
887.797	1203	7.782	1.662	1.939	0.632	4.014	0.0045
711.4947	1204	16.498	1.821	9.295	1.190	1.775	0.0046
769.5638	1204	114.290	10.000	74.462	6.683	1.535	0.0046

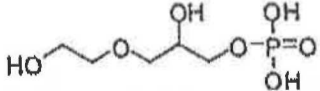
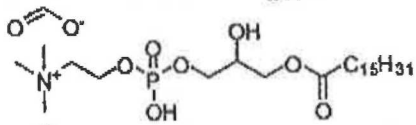
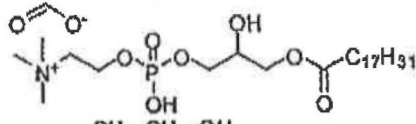
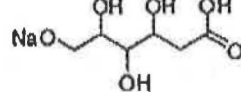
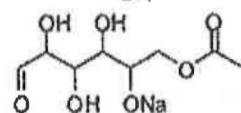
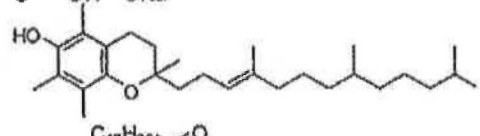
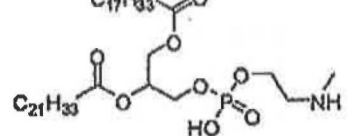
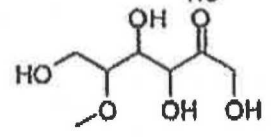
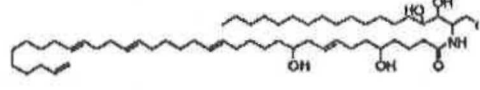
858.6843	1102	6.363	1.403	1.499	0.499	4.245	0.0047
861.5265	1102	6.739	0.942	3.423	0.386	1.989	0.0048
304.241	1204	35.178	2.794	23.043	2.384	1.527	0.0049
181.9806	1102	5.211	0.559	2.924	0.419	1.782	0.0051
280.2413	1204	153.282	19.301	81.661	11.092	1.877	0.0055
772.5842	1204	85.761	4.877	44.401	4.424	1.481	0.0056
830.5881	1102	14.539	2.244	7.110	0.677	2.045	0.0057
542.3447	1102	6.998	0.664	4.179	0.572	1.675	0.0059
281.2447	1204	29.936	3.733	16.187	2.186	1.849	0.006
696.4733	1204	10.289	2.321	2.124	1.124	4.844	0.006
744.5516	1204	132.030	11.098	87.905	8.188	1.502	0.006
699.4908	1204	7.332	0.636	4.685	0.530	1.565	0.0061
788.8128	1201	2.019	0.397	4.095	0.504	0.493	0.0062
853.5852	1202	23.738	3.949	10.596	1.444	2.240	0.0063
734.488	1204	15.884	2.046	8.871	0.912	1.791	0.0064
243.0719	1101	27.999	1.874	20.630	1.397	1.357	0.0066
256.24	1202	2.203	0.657	5.225	0.686	0.422	0.0066
345.8738	1101	6.834	0.780	3.405	0.747	2.007	0.0066
715.5228	1204	19.602	1.668	12.192	1.629	1.608	0.0066
746.5701	1204	50.317	4.437	33.781	2.855	1.489	0.0067
787.5995	1201	3.802	1.068	8.895	1.191	0.427	0.0067
897.5729	1102	7.234	1.101	3.570	0.427	2.026	0.0067
477.3218	1201	3.453	0.463	5.932	0.628	0.582	0.007
710.4916	1204	35.981	4.256	20.330	2.695	1.770	0.007
794.5419	1102	6.943	0.889	4.079	0.280	1.702	0.007
765.5704	1204	7.199	1.796	1.505	0.505	4.783	0.0073
829.5852	1202	51.097	8.889	23.637	2.394	2.162	0.0073
277.8861	1101	14.484	1.330	8.729	1.282	1.659	0.0075
803.5681	1202	107.752	18.098	50.203	5.542	2.146	0.0075
743.5461	1204	311.426	28.757	202.118	20.807	1.541	0.0076
830.5885	1202	22.923	3.777	11.038	1.041	2.077	0.0076
829.5851	1102	31.727	5.202	15.489	1.333	2.048	0.0077
825.5532	1202	4.799	0.761	2.209	0.389	2.172	0.008
827.5694	1202	62.787	10.158	29.412	4.314	2.135	0.008
773.5954	1204	28.422	2.040	19.815	1.943	1.434	0.0085
793.5381	1102	14.713	2.099	8.058	0.764	1.826	0.0086
828.5734	1202	29.459	4.653	14.227	2.071	2.071	0.0086
555.3102	1102	6.492	0.756	3.572	0.600	1.818	0.0087
804.5714	1202	45.461	7.448	22.273	2.328	2.041	0.0087
144.0944	1203	5.913	0.351	7.175	0.232	0.824	0.009
160.1256	1203	6.763	0.421	8.308	0.301	0.814	0.0093
809.5932	1101	15.522	1.647	10.274	0.674	1.511	0.0093
847.5315	1201	6.021	0.711	3.555	0.427	1.694	0.0093
1127.741	1204	3.217	0.760	1.000	0.000	3.217	0.0094
871.5526	1102	14.899	2.355	7.613	0.893	1.970	0.0096
531.312	1102	6.157	0.560	3.797	0.556	1.621	0.0097
824.5477	1201	2.172	0.368	1.078	0.078	2.015	0.0097
767.5827	1201	3.881	0.389	7.191	1.015	0.540	0.0099
634.4267	1201	1.649	0.225	1.000	0.000	1.649	0.01
757.5622	1201	98.678	7.301	147.948	14.540	0.667	0.01
797.5257	1204	75.239	6.508	49.628	5.701	1.516	0.0101
698.4885	1204	15.802	1.637	10.577	0.753	1.494	0.0104
305.2439	1204	7.384	0.624	5.028	0.505	1.469	0.0105
793.5986	1201	5.506	0.560	10.176	1.451	0.541	0.0108
805.5832	1202	35.953	5.189	20.118	1.873	1.787	0.0109
771.5792	1204	136.073	11.270	92.313	9.966	1.474	0.0113
798.6019	1204	14.243	1.131	9.780	1.037	1.456	0.0113
801.5543	1202	5.002	0.852	2.386	0.338	2.096	0.0113
856.6698	1102	6.227	0.961	3.318	0.356	1.877	0.0117
360.1467	1201	3.600	0.662	6.748	0.851	0.533	0.0118
739.4827	1204	11.114	1.363	6.723	0.728	1.653	0.012
742.5366	1204	14.637	1.420	9.248	1.224	1.583	0.0121
745.5643	1204	111.991	10.785	75.580	6.918	1.482	0.0123
806.5865	1202	15.199	2.072	8.845	0.893	1.718	0.0124
638.5138	1204	16.487	2.611	8.215	1.368	2.007	0.0129

821.5714	1102	9.269	1.328	5.241	0.559	1.769	0.0129
729.5727	1204	7.812	0.649	5.126	0.683	1.524	0.0131
852.5726	1202	9.104	1.529	4.422	0.890	2.059	0.0131
260.2507	1204	16.191	1.538	10.701	1.205	1.513	0.0135
593.3416	1204	1.607	0.248	3.779	0.717	0.425	0.0144
793.5663	1204	51.823	4.721	34.525	4.078	1.501	0.0148
758.5855	1201	49.966	3.650	73.507	7.476	0.680	0.015
796.5864	1204	43.191	4.094	28.514	3.400	1.515	0.0152
1225.093	1203	11.951	2.621	3.439	1.681	3.475	0.0154
501.3217	1201	3.393	0.263	5.242	0.602	0.647	0.0156
428.295	1204	5.120	0.550	2.795	0.634	1.832	0.0157
766.5372	1204	14.283	1.111	9.541	1.307	1.497	0.0159
817.5374	1102	6.843	1.095	3.764	0.350	1.818	0.0161
448.3194	1204	3.018	0.318	4.428	0.400	0.682	0.0162
832.6022	1102	19.054	3.035	10.450	1.065	1.823	0.0163
666.5449	1204	30.717	5.437	14.525	2.657	2.115	0.0167
278.2254	1204	14.617	3.007	5.989	1.213	2.441	0.017
134.11	1203	10.616	0.728	12.993	0.503	0.817	0.0171
736.4951	1201	1.974	0.296	3.164	0.318	0.626	0.0173
759.5779	1201	36.709	2.675	55.828	6.376	0.658	0.0173
787.6095	1201	3.839	1.154	8.481	1.266	0.453	0.0176
279.2284	1204	3.218	0.622	1.424	0.270	2.260	0.0177
338.2461	1204	2.353	0.284	3.961	0.516	0.594	0.018
851.6886	1202	17.271	3.006	8.654	1.290	1.996	0.018
853.5848	1102	11.637	2.132	5.890	0.529	1.976	0.0181
589.3398	1102	9.802	1.051	5.292	1.297	1.852	0.0182
440.308	1201	2.501	0.752	6.005	1.056	0.417	0.0184
769.4929	1204	100.639	10.775	65.749	7.683	1.531	0.019
759.5779	1101	20.424	2.309	13.845	1.041	1.475	0.0195
454.2969	1201	2.600	0.647	6.918	1.462	0.376	0.0196
283.2602	1204	39.989	6.317	20.305	4.126	1.969	0.0199
786.5967	1201	23.862	1.962	35.935	4.045	0.664	0.02
282.2572	1204	208.013	33.546	104.163	21.783	1.997	0.0204
785.5934	1201	49.962	3.937	74.679	8.353	0.669	0.0204
194.0803	1203	9.310	2.283	3.242	0.734	2.872	0.0219
612.4994	1204	4.837	0.644	3.004	0.327	1.610	0.0223
810.5986	1101	7.551	0.886	5.184	0.311	1.457	0.0223
146.11	1203	5.303	0.322	6.310	0.229	0.840	0.0227
722.486	1204	9.551	0.914	6.518	0.757	1.465	0.0229
741.5302	1204	27.862	3.508	17.261	2.288	1.615	0.0231
678.4528	1201	2.508	0.327	3.923	0.439	0.639	0.0232
765.5316	1204	26.722	2.313	17.611	2.680	1.517	0.0232
279.9312	1102	5.532	0.590	3.566	0.499	1.551	0.0234
681.5631	1204	5.859	1.282	2.128	0.748	2.753	0.0237
158.1101	1203	6.179	0.353	7.354	0.299	0.840	0.0238
831.5992	1102	42.798	7.551	23.095	2.383	1.853	0.0238
899.5871	1102	9.096	1.561	5.053	0.463	1.800	0.024
780.5811	1201	17.113	1.223	25.455	2.973	0.672	0.0242
799.5401	1204	22.934	2.034	16.153	1.745	1.420	0.0242
804.5715	1102	33.353	5.961	18.084	1.583	1.844	0.0242
856.6048	1202	18.380	2.965	10.122	1.484	1.816	0.0245
150.1413	1203	4.270	0.313	5.372	0.303	0.795	0.0247
640.5285	1204	19.392	3.594	9.651	1.655	2.009	0.0257
806.5863	1102	12.671	2.471	6.480	0.601	1.955	0.0262
462.3716	1204	2.921	0.398	1.744	0.262	1.675	0.0264
678.5469	1204	12.126	3.131	3.709	1.422	3.269	0.0265
797.5973	1204	28.011	2.531	19.814	2.126	1.414	0.0267
760.581	1101	9.590	1.054	6.767	0.476	1.417	0.0268
304.1111	1202	3.207	0.871	5.725	0.541	0.560	0.0269
738.4806	1204	20.803	3.350	11.861	1.518	1.754	0.0274
664.5313	1204	27.632	5.581	12.747	2.590	2.168	0.0281
803.5683	1102	83.168	15.663	44.392	3.974	1.873	0.0282
855.6011	1202	36.927	6.113	20.484	2.974	1.803	0.0282
674.4902	1204	8.970	0.728	6.706	0.575	1.338	0.0286
255.2283	1204	4.079	0.561	2.202	0.524	1.852	0.0291

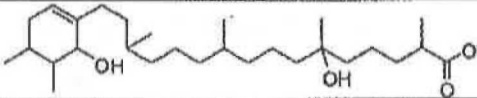
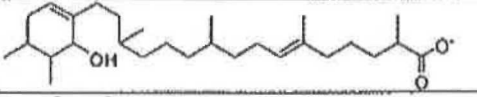
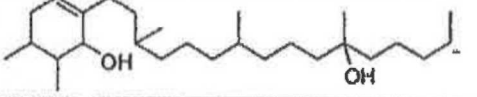
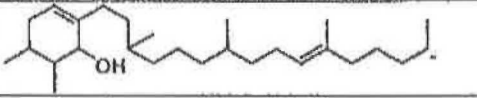
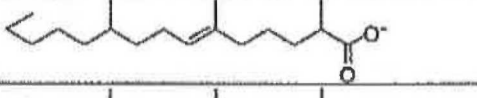
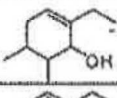
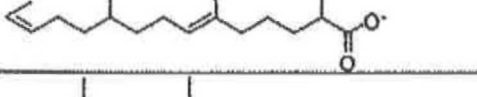
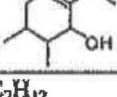
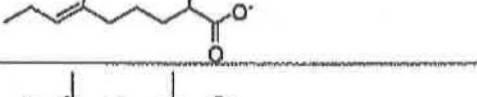
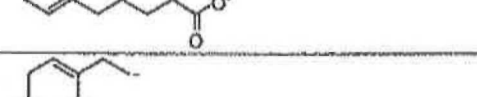
781.5619	1204	9.821	0.767	6.454	1.145	1.522	0.031
828.5732	1102	17.547	3.618	8.840	0.788	1.985	0.031
446.219	1201	2.849	0.439	4.644	0.593	0.813	0.0313
172.1255	1203	6.279	0.331	7.346	0.297	0.855	0.0314
330.2569	1204	3.124	0.365	2.136	0.203	1.463	0.0316
827.5695	1102	37.092	7.880	18.168	1.787	2.042	0.0316
794.802	1201	2.964	0.305	4.814	0.690	0.616	0.0318
700.6037	1204	17.186	1.533	12.592	1.177	1.365	0.0323
720.5081	1204	5.498	0.665	3.262	0.656	1.685	0.0324
832.6027	1202	26.701	4.628	15.120	1.758	1.766	0.0325
768.5525	1204	84.034	8.929	56.831	7.192	1.479	0.0328
246.1468	1201	3.198	0.597	6.521	1.233	0.491	0.0331
252.2096	1204	2.857	0.418	1.603	0.327	1.782	0.0332
832.5793	1101	4.657	0.579	3.138	0.294	1.484	0.0333
831.5995	1202	58.297	10.567	32.015	4.088	1.821	0.0338
664.4374	1201	2.661	0.255	3.734	0.370	0.713	0.0344
701.5064	1204	7.572	0.687	5.461	0.588	1.387	0.0355
752.4902	1201	2.178	0.299	3.211	0.322	0.678	0.0356
856.6045	1102	11.620	2.428	5.891	0.628	1.973	0.0356
609.324	1102	4.496	0.485	3.054	0.386	1.472	0.036
702.5676	1101	4.257	0.546	2.826	0.297	1.506	0.0361
340.2621	1204	3.898	0.427	5.495	0.529	0.709	0.0364
1227.091	1203	10.143	2.889	19.624	2.874	0.517	0.0369
782.565	1101	16.038	1.944	11.211	0.883	1.430	0.0385
452.244	1201	5.926	0.441	7.858	0.702	0.754	0.0388
805.5833	1102	28.734	6.015	15.001	1.230	1.916	0.0388
884.709	1204	6.673	1.627	2.508	0.862	2.661	0.0388
512.3347	1201	3.117	0.364	6.248	1.279	0.499	0.039
156.0943	1203	3.520	0.264	4.172	0.120	0.844	0.0396
662.5164	1204	15.920	3.832	6.974	1.209	2.283	0.0402
783.5778	1201	50.681	3.561	70.630	7.811	0.718	0.0404
795.5814	1204	81.007	8.721	55.971	6.822	1.447	0.0406
868.7141	1204	6.709	1.697	2.424	0.889	2.768	0.0408
734.6429	1204	12.043	1.235	8.376	1.054	1.438	0.0412
148.1257	1203	7.248	0.541	8.624	0.296	0.840	0.0414
303.1081	1102	4.733	0.404	2.761	0.755	1.714	0.0414
438.2924	1201	1.630	0.290	3.698	0.851	0.441	0.0429
786.5965	1101	13.950	1.980	9.293	0.754	1.501	0.0429
781.5617	1101	32.132	3.891	22.886	1.609	1.404	0.0433
833.5931	1201	8.927	1.272	14.084	1.894	0.634	0.0439
811.6093	1201	8.938	1.013	13.019	1.497	0.687	0.0442
807.5754	1101	23.270	3.292	15.705	1.123	1.482	0.0447
257.1709	1201	4.487	0.776	9.661	2.145	0.464	0.0453
228.1362	1201	23.118	2.605	44.117	8.870	0.524	0.0455
767.5473	1204	179.498	20.657	121.949	16.310	1.472	0.0469
823.5427	1204	10.161	0.841	7.863	0.633	1.292	0.0469
821.4768	1204	7.624	0.798	11.838	1.704	0.644	0.0471
702.5676	1201	2.900	0.357	3.945	0.316	0.735	0.0473
855.6009	1102	23.110	5.234	11.619	1.220	1.989	0.0473
718.5348	1204	6.442	0.685	4.178	0.783	1.542	0.0475
162.1412	1203	5.557	0.362	6.454	0.204	0.861	0.0476
606.4872	1204	7.849	1.273	4.484	0.885	1.750	0.0477
784.5811	1201	24.262	1.796	33.594	3.775	0.722	0.0477
895.5575	1102	8.641	1.842	4.628	0.422	1.867	0.0487
705.6086	1204	7.166	0.897	4.887	0.562	1.466	0.0489
242.1519	1201	4.098	0.605	6.579	0.951	0.623	0.0494
264.9759	1204	6.810	0.414	5.630	0.354	1.210	0.0494

[00200] Table 22: Accurate masses, mode of ionization, putative molecular formulae and proposed structures for multiple sclerosis biomarkers detected in aqueous and organic extracts of human serum.

	Detected Mass	Exact Mass	Mode	Formula	Proposed Structure
1	452.3868	452.3866	1204	$C_{28}H_{52}O_4$	
2	496.4157	496.4128	1204	$C_{30}H_{56}O_5$	
3	524.4448	524.4441	1204	$C_{32}H_{60}O_5$	
4	540.4387	540.4390	1204	$C_{32}H_{60}O_6$	
5	576.4757	576.4754	1204	$C_{36}H_{64}O_5$	
6	578.4923	578.4910	1204	$C_{36}H_{66}O_5$	
7	580.5089	580.5067	1204	$C_{36}H_{68}O_5$	
8	594.4848	594.4859	1204	$C_{36}H_{66}O_6$	
9	596.5012	596.5016	1204	$C_{36}H_{68}O_6$	
10	786.5408	786.5411	1204	$C_{43}H_{79}O_{10}P$	

11	216.04	216.0399	1102	$C_5H_{13}O_7P$	
12	541.3415	541.3379	1102	$C_{25}H_{52}NO_9P$	
13	565.3391	565.3380	1102	$C_{27}H_{52}NO_9P$	
14	202.0453	202.0453	1101	$C_6H_{11}O_6Na$	
15	244.0559	244.0559	1101	$C_8H_{13}O_7Na$	
16	428.3653	428.3654	1201	$C_{29}H_{48}O_2$	
17	805.5609	805.5621	1201	$C_{46}H_{80}NO_8P$	
18	194.0803	194.0790	1203	$C_7H_{14}O_6$	
19	857.7516	857.7472	1203	$C_{54}H_{99}NO_6$	

[00201] Table 23: MS/MS fragmentation of multiple sclerosis biomarker I,
452.3868 ($C_{28}H_{52}O_4$)

m/z	Formula	Molecular fragment	Fragment loss
451	$C_{28}H_{51}O_4$		$-H^+$
433	$C_{28}H_{49}O_3$		$-H_2O$
407	$C_{27}H_{51}O_2$		$-CO_2$
389	$C_{27}H_{49}O$		433 - CO_2
281	$C_{18}H_{33}O_2$		433 - 
279	$C_{18}H_{31}O_2$		433 - 
183	$C_{11}H_{19}O_2$		279 - C_7H_{12}
169	$C_{10}H_{17}O_2$		279 - C_8H_{14}
153	$C_{10}H_{17}O$		- phytol chain
139	$C_9H_{15}O$		153 - CH_4
125	$C_8H_{13}O$		139 - CH_4
111	$C_7H_{11}O$		125 - CH_4
97	C_6H_9O		111 - CH_4

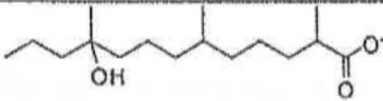
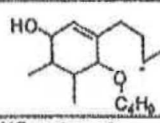
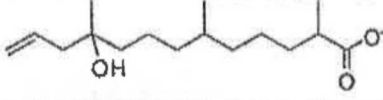
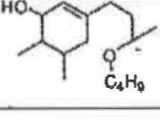
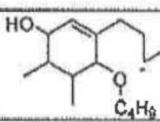
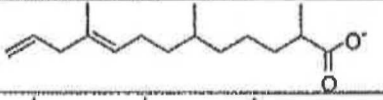
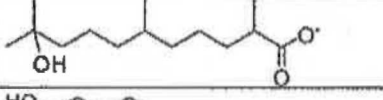
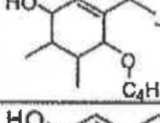
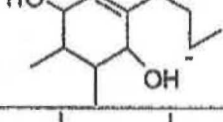
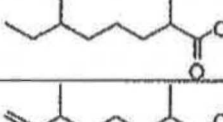
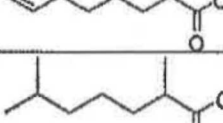
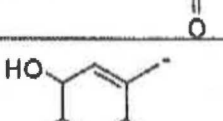
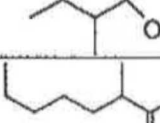
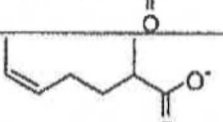
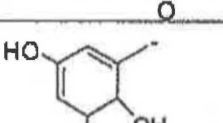
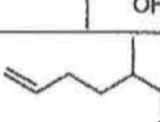
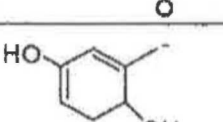
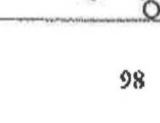
[00202] Table 24: MS/MS fragmentation of multiple sclerosis biomarker 2, 496.4157 (C₃₀H₅₆O₅)

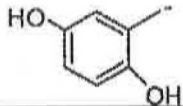
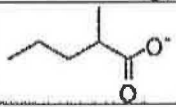
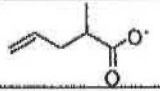
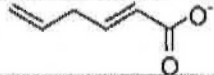
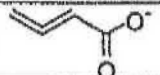
m/z	Formula	Molecular fragment	Fragment loss
495	C ₃₀ H ₅₅ O ₅		-H ⁺
477	C ₃₀ H ₅₃ O ₃		-H ₂ O
451	C ₂₉ H ₅₅ O ₃		-CO ₂
433	C ₂₉ H ₅₃ O ₂		-(CO ₂ + H ₂ O)
415	C ₂₉ H ₅₁ O		-(CO ₂ + 2H ₂ O)
307	C ₂₀ H ₃₅ O ₂		415 -
297	C ₁₈ H ₃₃ O ₃		-
279	C ₁₈ H ₃₁ O ₂		297 - H ₂ O
235	C ₁₆ H ₂₇ O		279 -
223	C ₁₄ H ₂₃ O ₂		477 -
215	C ₁₂ H ₂₃ O ₂		Fragmentation at C13-C14 and loss of CH ₃
197	C ₁₂ H ₂₁ O ₂		- phytol chain

179	$C_{12}H_{19}O$		197 - H_2O
181	$C_{13}H_{25}$		415 - 235
169	$C_{10}H_{17}O_2$		179 - C_2H_4
157	$C_8H_{13}O_3$		215 - C_4H_{10}
155	$C_9H_{15}O_2$		223 -
153	$C_{10}H_{17}O$		197 - C_2H_4O
141	$C_9H_{17}O$		451 -
139	$C_9H_{15}O$		153 - CH_4
127	$C_8H_{15}O$		141 - CH_2
125	$C_8H_{13}O$		184 - C_4H_8
113	$C_6H_9O_2$		157 - C_2H_4O

[00203] Table 25: MS/MS fragmentation of multiple sclerosis biomarker 3, 524.4448 ($C_{32}H_{60}O_5$)

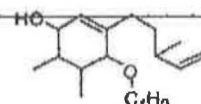
m/z	Formula	Molecular fragment	Fragment loss
523	$C_{32}H_{59}O_5$		- H^+
505	$C_{32}H_{57}O_4$		- H_2O
487	$C_{32}H_{55}O_3$		- $2 \times H_2O$
479	$C_{31}H_{59}O_3$		- CO_2
463	$C_{30}H_{55}O_3$		479 - CH_4
461	$C_{31}H_{57}O_2$		-($CO_2 + H_2O$)
443	$C_{31}H_{55}O$		-($CO_2 + 2H_2O$)
365	$C_{23}H_{41}O_3$		463 - C_7H_{13}
337	$C_{21}H_{37}O_3$		365 - C_2H_4
299	$C_{18}H_{33}O_3$		
297	$C_{18}H_{33}O_3$		
281	$C_{18}H_{33}O_2$		
279	$C_{18}H_{31}O_2$		297 - H_2O

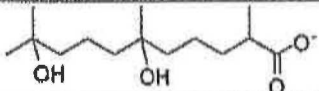
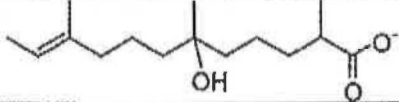
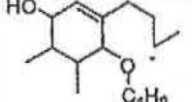
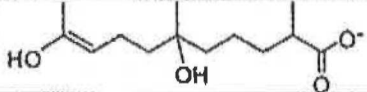
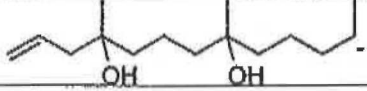
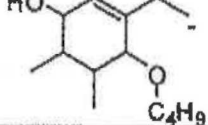
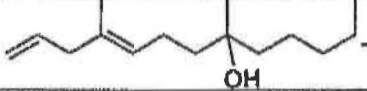
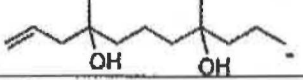
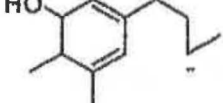
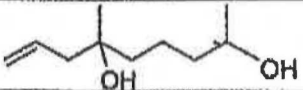
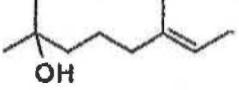
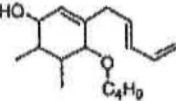
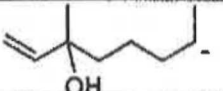
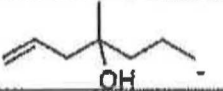
271	$C_{16}H_{31}O_3$		- 
269	$C_{16}H_{29}O_3$		- 
253	$C_{16}H_{29}O_2$		- 271
251	$C_{16}H_{27}O_2$		269 - H_2O
243	$C_{14}H_{27}O_3$		- 281
225	$C_{14}H_{25}O_2$		- phytol chain
197	$C_{12}H_{21}O_2$		253 - C_4H_8
171	$C_{10}H_{19}O_2$		251 - C_6H_8
169	$C_{10}H_{17}O_2$		251 - C_6H_{10}
157	$C_9H_{17}O_2$		271 - CH_2
155	$C_9H_{15}O_2$		197 - C_7H_6
143	$C_8H_{15}O_2$		157 - CH_2
141	$C_9H_{17}O$		157 - CH_4
139	$C_8H_{11}O_2$		155 - CH_4
127	$C_7H_{11}O_2$		143 - CH_4
125	$C_8H_{13}O_2$		139 - CH_3

123	$C_7H_7O_2$		139 - CH_4
115	$C_6H_{11}O_2$		141 - C_3H_6
113	$C_6H_9O_2$		141 - C_3H_4
111	$C_6H_7O_2$		127 - CH_4
83	$C_4H_5O_2$		111 - C_2H_4

[00204] Table 26: MS/MS fragmentation of multiple sclerosis biomarker 4, 540.4390 (C₃₂H₆₀O₆)

m/z	Formula	Molecular fragment	Fragment loss
539	C ₃₂ H ₅₉ O ₆		- H ⁺
521	C ₃₂ H ₅₇ O ₅		- H ₂ O
503	C ₃₂ H ₅₅ O ₄		- 2 x H ₂ O
495	C ₃₁ H ₅₉ O ₄		- CO ₂
477	C ₃₁ H ₅₇ O ₃		- (CO ₂ + H ₂ O)
461	C ₃₀ H ₅₃ O ₃		477 - CH ₄
459	C ₃₁ H ₅₅ O ₂		- (CO ₂ + 2 x H ₂ O)
419	C ₂₇ H ₄₇ O ₃		461 - C ₃ H ₆
335	C ₂₂ H ₃₉ O ₂		459 - C ₉ H ₁₆
315	C ₁₈ H ₃₅ O ₄		
313	C ₁₈ H ₃₃ O ₄		
297	C ₁₈ H ₃₃ O ₃		315 - H ₂ O
279	C ₁₈ H ₃₁ O ₂		297 - H ₂ O



259	$C_{14}H_{27}O_4$		-
255	$C_{15}H_{27}O_3$		297 - C_3H_6
253	$C_{16}H_{29}O_2$		503 - phytol chain
243	$C_{14}H_{27}O_3$		259 - CH_4
241	$C_{15}H_{29}O_2$		495 - 253
225	$C_{14}H_{25}O_2$		- phytol chain
223	$C_{14}H_{23}O_2$		241 - H_2O
213	$C_{13}H_{25}O_2$		241 - C_2H_4
179	$C_{12}H_{19}O$		253 - C_4H_9OH
171	$C_{10}H_{19}O_2$		213 - C_3H_6
155	$C_9H_{15}O_2$		419 - 
141	$C_8H_{13}O_2$		223 - C_6H_{10}
127	$C_8H_{15}O$		171 - C_2H_4O

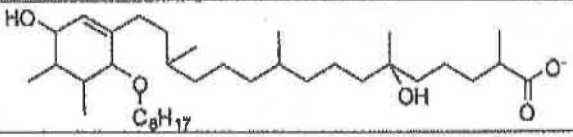
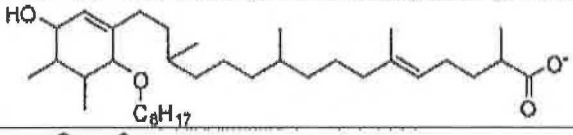
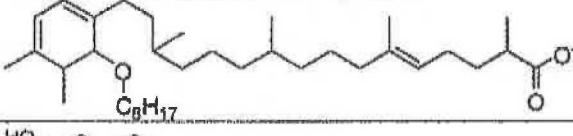
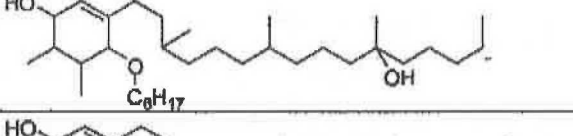
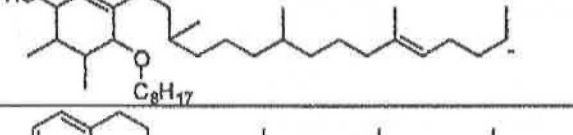
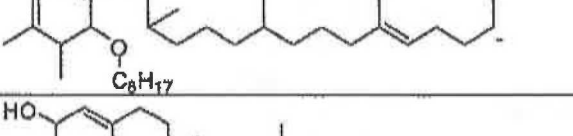
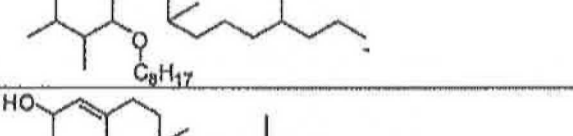
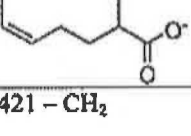
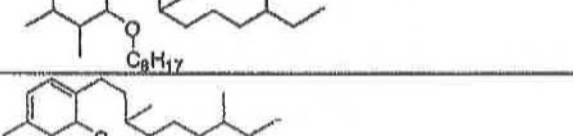
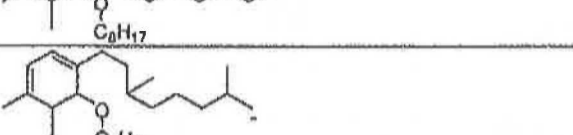
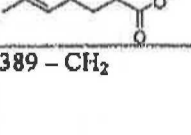
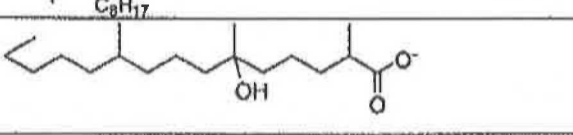
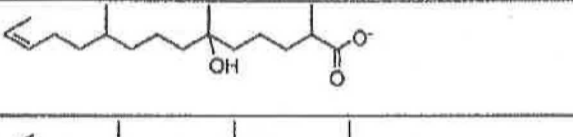
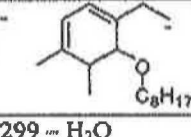
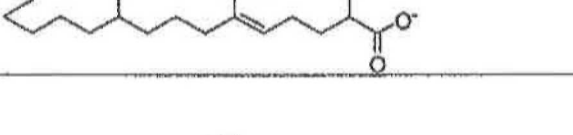
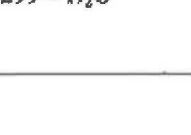
[00205] Table 27: MS/MS fragmentation of multiple sclerosis biomarker 5, 576.4757 (C₃₆H₆₄O₅)

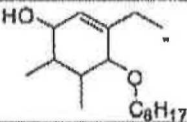
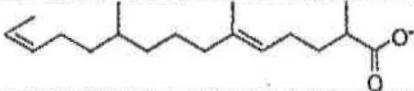
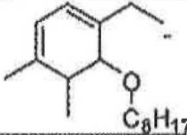
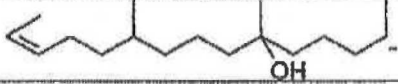
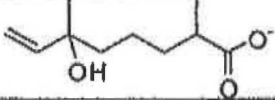
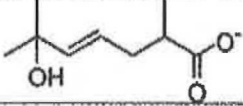
m/z	Formula	Molecular fragment	Fragment loss
575	C ₃₆ H ₆₃ O ₅		-H ⁺
557	C ₃₆ H ₆₁ O ₄		-H ₂ O
539	C ₃₆ H ₅₉ O ₃		-2xH ₂ O
531	C ₃₅ H ₆₃ O ₃		-CO ₂
513	C ₃₅ H ₆₁ O ₂		557 - CO ₂
495	C ₃₅ H ₅₉ O		531 - CO ₂
417	C ₂₈ H ₄₉ O ₂		-
403	C ₂₈ H ₄₇ O ₂		417 - CH ₂
371	C ₂₆ H ₄₃ O		387 - CH ₂
297	C ₁₈ H ₃₃ O ₃		-
279	C ₁₈ H ₃₃ O ₂		-
279	C ₁₈ H ₃₁ O ₂		- phytol chain
251	C ₁₆ H ₂₇ O ₂		557 -
183	C ₁₁ H ₁₉ O ₂		557 -

[00206] Table 28: MS/MS fragmentation of multiple sclerosis biomarker 6, 578.4848 ($C_{36}H_{66}O_5$)

m/z	Formula	Molecular fragment	Fragment loss
577	$C_{36}H_{65}O_5$		$-H^+$
559	$C_{36}H_{63}O_4$		$-H_2O$
541	$C_{36}H_{61}O_3$		$-2xH_2O$
533	$C_{35}H_{65}O_3$		$-CO_2$
515	$C_{35}H_{63}O_2$		$559 - CO_2$
497	$C_{35}H_{61}O$		$533 - CO_2$
419	$C_{29}H_{51}O_2$		$-CH_2OHCOO^-$
405	$C_{28}H_{49}O_2$		$419 - CH_2$
387	$C_{27}H_{47}O$		$405 - H_2O$
373	$C_{26}H_{45}O$		$387 - CH_2$
297	$C_{18}H_{33}O_3$		$-3,4,5\text{-trimethylphenyl-C}_8\text{H}_{17}\text{-O-}$
281	$C_{18}H_{33}O_2$		
279	$C_{18}H_{31}O_2$		$297 - H_2O$
279	$C_{18}H_{31}O_2$		$- \text{phytol chain}$

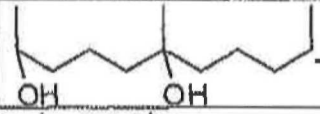
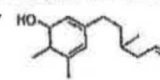
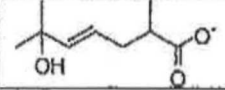
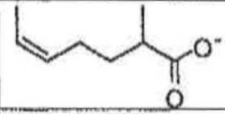
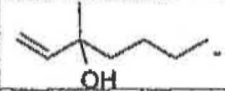
[00207] Table 29: MS/MS fragmentation of multiple sclerosis biomarker 7, 580.5089 ($C_{36}H_{68}O_5$)

m/z	Formula	Molecular fragment	Fragment loss
579	$C_{36}H_{67}O_5$		$-H^+$
561	$C_{36}H_{65}O_4$		$-H_2O$
543	$C_{35}H_{65}O_3$		$-2 \times H_2O$
535	$C_{35}H_{67}O_3$		$-CO_2$
517	$C_{35}H_{65}O_2$		$561 - CO_2$
499	$C_{35}H_{63}O$		$535 - CO_2$
421	$C_{28}H_{53}O_2$		$561 -$ 
407	$C_{27}H_{51}O_2$		$421 - CH_2$
389	$C_{27}H_{49}O$		$543 -$ 
375	$C_{26}H_{47}O$		$389 - CH_2$
299	$C_{18}H_{35}O_3$		$-$ 
297	$C_{18}H_{33}O_3$		$-$ 
281	$C_{18}H_{33}O_2$		$299 - H_2O$

281	$C_{18}H_{33}O_2$		• phytol chain
279	$C_{18}H_{31}O_2$		297 - H_2O
263	$C_{18}H_{31}O$		543 - phytol chain
253	$C_{17}H_{33}O$		535 - 263
185	$C_{10}H_{17}O_3$		299 - C_8H_{19}
171	$C_9H_{15}O_3$		

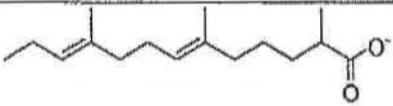
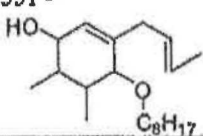
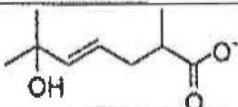
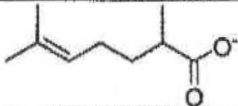
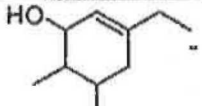
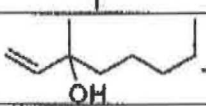
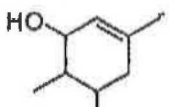
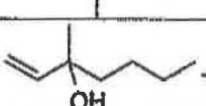
[00208] Table 30: MS/MS fragmentation of multiple sclerosis biomarker 8, 594.4848 (C₃₆H₆₆O₆)

m/z	Formula	Molecular fragment	Fragment loss
593	C ₃₆ H ₆₅ O ₆		-H ⁺
575	C ₃₆ H ₆₅ O ₅		-H ₂ O
557	C ₃₆ H ₆₅ O ₄		-2xH ₂ O
549	C ₃₅ H ₆₅ O ₄		-CO ₂
531	C ₃₅ H ₆₅ O ₃		575 - CO ₂
513	C ₃₅ H ₆₅ O ₂		549 - CO ₂
495	C ₃₅ H ₆₁ O		495 - H ₂ O
421	C ₂₇ H ₄₉ O ₃		531 - C ₈ H ₁₆ O
371	C ₂₆ H ₄₃ O		513 -
315	C ₁₈ H ₃₅ O ₄		-
297	C ₁₈ H ₃₃ O ₃		495 - H ₂ O
279	C ₁₈ H ₃₁ O ₂		421 - H ₂ O
279	C ₁₈ H ₃₁ O ₂		- phytol chain

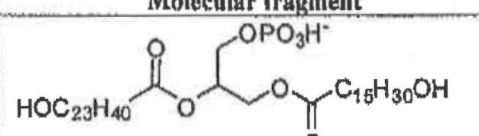
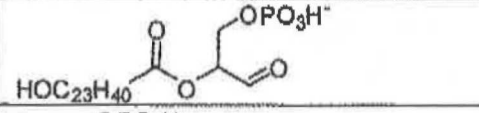
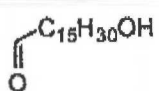
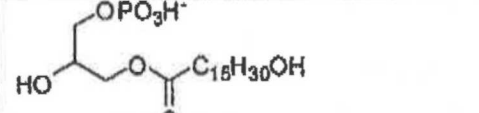
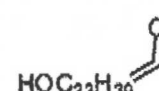
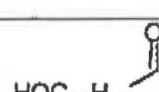
201	$C_{12}H_{23}O_2$		421 - 
171	$C_9H_{15}O_3$		
141	$C_8H_{13}O_2$		
127	$C_8H_{15}O$		

[00209] Table 31: MS/MS fragmentation of multiple sclerosis biomarker 9, 596.5012 (C₃₆H₆₈O₆)

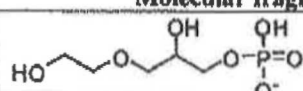
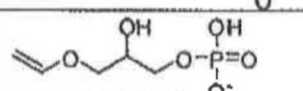
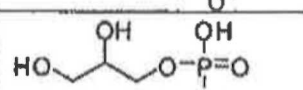
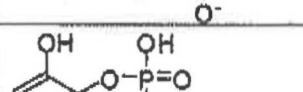
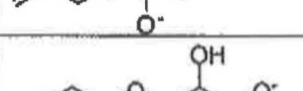
m/z	Formula	Molecular fragment	Fragment loss
595	C ₃₆ H ₆₇ O ₆		-H ⁺
577	C ₃₆ H ₆₅ O ₅		-H ₂ O
559	C ₃₆ H ₆₃ O ₄		-2xH ₂ O
551	C ₃₅ H ₆₇ O ₂		-CO ₂
515	C ₃₅ H ₆₃ O ₂		559 - CO ₂
497	C ₃₃ H ₆₁ O		515 - H ₂ O
423	C ₂₇ H ₅₁ O ₃		515 - C ₈ H ₁₆ O
373	C ₂₆ H ₄₅ O		497 -
315	C ₁₈ H ₃₅ O ₄		-
297	C ₁₈ H ₃₃ O ₃		315 - H ₂ O
281	C ₁₈ H ₃₂ O ₂		- phytol chain
279	C ₁₈ H ₃₁ O ₂		297 - H ₂ O
269	C ₁₆ H ₂₉ O ₃		577 -

251	$C_{16}H_{27}O_2$		551 - 
171	$C_9H_{15}O_3$		
155	$C_9H_{15}O_2$		
153	$C_{10}H_{17}O$		
141	$C_9H_{17}O$		
139	$C_9H_{15}O$		
127	$C_8H_{15}O$		

[00210] Table 32: MS/MS fragmentation of multiple sclerosis biomarker 10, 786.5408 (C₄₃H₇₉O₁₀P)

m/z	Formula	Molecular fragment	Fragment loss
785	C ₄₃ H ₇₈ O ₁₀ P		- H ⁺
529	C ₂₇ H ₄₆ O ₈ P		- 
425	C ₁₉ H ₃₈ O ₈ P		- 
169	C ₃ H ₆ O ₆ P		529 - 
97	H ₂ PO ₄		

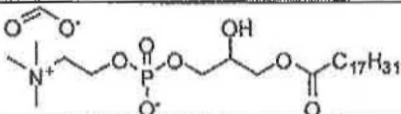
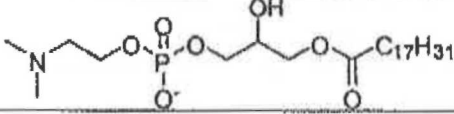
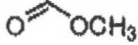
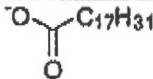
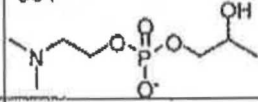
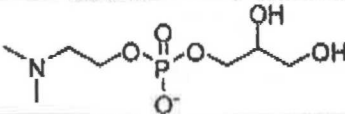
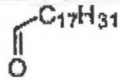
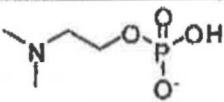
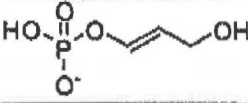
[00211] Table 33: MS/MS fragmentation of multiple sclerosis biomarker 11, 216.04 (C₅H₁₃O₇P)

m/z	Formula	Molecular fragment	Fragment loss
215	C ₅ H ₁₂ O ₇ P		- H ⁺
197	C ₅ H ₁₀ O ₆ P		- H ₂ O
171	C ₃ H ₈ O ₆ P		197 - C ₂ H ₂
153	C ₃ H ₆ O ₅ P		171 - H ₂ O
135	C ₃ H ₁₁ O ₄		- 

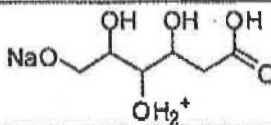
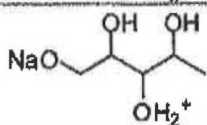
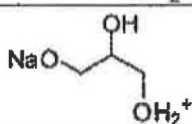
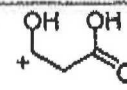
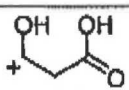
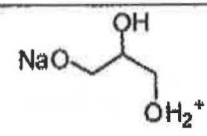
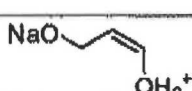
[00212] Table 34: MS/MS fragmentation of multiple sclerosis biomarker 12, 541.3415 (C₂₅H₅₂NO₉P)

m/z	Formula	Molecular fragment	Fragment loss
540	C ₂₅ H ₅₁ NO ₉ P		-H ⁺
480	C ₂₃ H ₄₇ NO ₇ P		-
255	C ₁₆ H ₃₁ O ₂		480 -
242	C ₇ H ₁₇ NO ₆ P		480 -
224	C ₇ H ₁₅ NO ₅ P		242 - H ₂ O
168	C ₄ H ₁₁ NO ₄ P		224 -
153	C ₃ H ₆ O ₅ P		224 -
79	PO ₃		

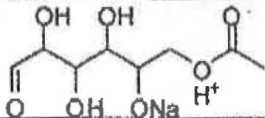
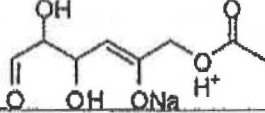
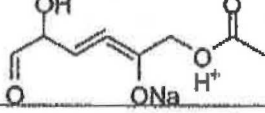
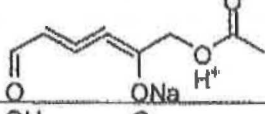
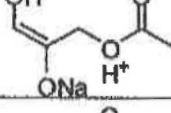
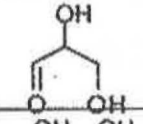
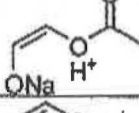
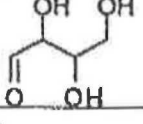
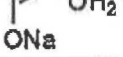
[00213] Table 35: MS/MS fragmentation of multiple sclerosis biomarker 13, 565.3391 (C₄₇H₈₃NO₁₃P)

m/z	Formula	Molecular fragment	Fragment loss
564	C ₂₇ H ₅₁ NO ₉ P		-H ⁺
504	C ₂₅ H ₄₅ NO ₈ P		- 
279	C ₁₈ H ₃₁ O ₂		504 - 
242	C ₇ H ₁₇ NO ₆ P		504 - 
224	C ₇ H ₁₅ NO ₅ P		242 - H ₂ O
168	C ₄ H ₁₁ NO ₄ P		224 - 
153	C ₃ H ₆ O ₅ P		224 - 
79	PO ₃		

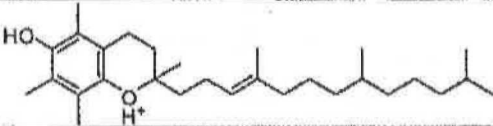
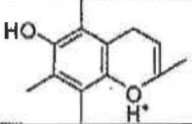
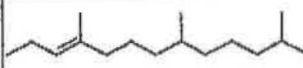
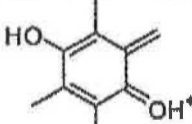
[00214] Table 36: MS/MS fragmentation of multiple sclerosis biomarker 14, 202.0453 ($C_6H_{11}O_6Na$)

m/z	Formula	Molecular fragment	Fragment loss
203	$C_6H_{12}O_6Na$		+ H^+
159	$C_5H_{12}O_4Na$		- CO_2
115	$C_3H_8O_3Na$		- 
89	$C_3H_5O_3$		- 
97	$C_3H_6O_2Na$		115 - H_2O

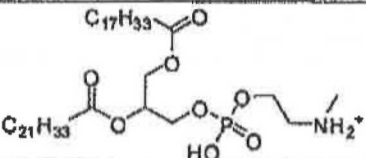
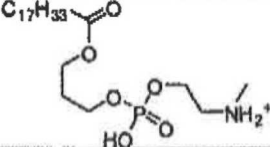
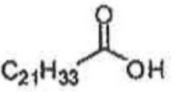
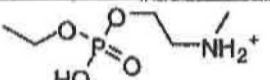
[00215] Table 37: MS/MS fragmentation of multiple sclerosis biomarker 15, 244.0559 ($C_8H_{13}O_7Na$)

m/z	Formula	Molecular fragment	Fragment loss
245	$C_8H_{14}O_7Na$		+ H^+
227	$C_8H_{12}O_6Na$		- H_2O
209	$C_8H_{10}O_5Na$		227 - H_2O
191	$C_8H_8O_4Na$		209 - H_2O
155	$C_5H_8O_4Na$		- 
125	$C_4H_6O_3Na$		- 
83	$C_2H_4O_2Na$		125 - 

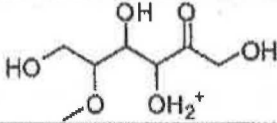
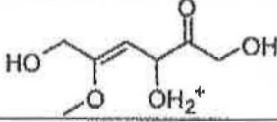
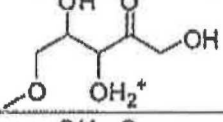
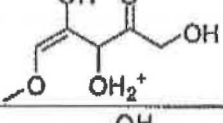
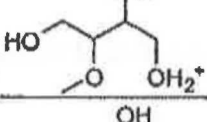
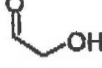
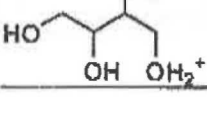
[00216] Table 38: MS/MS fragmentation of multiple sclerosis biomarker 16, 428.3653 ($C_{29}H_{48}O_2$)

m/z	Formula	Molecular fragment	Fragment loss
429	$C_{29}H_{49}O_2$		$+H^+$
205	$C_{13}H_{17}O_2$		
165	$C_{10}H_{13}O_2$		

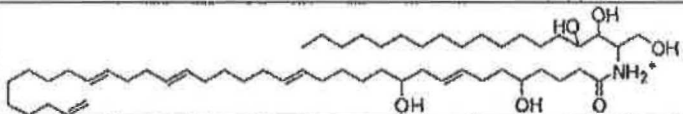
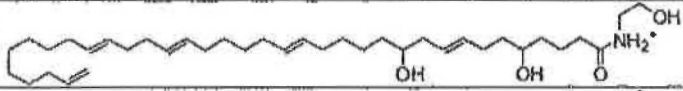
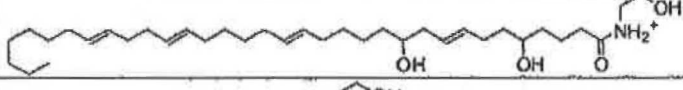
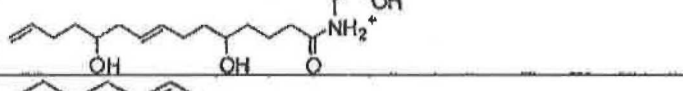
[00217] Table 39: MS/MS fragmentation of multiple sclerosis biomarker 17, 805.5609 ($C_{46}H_{80}NO_8P$)

m/z	Formula	Molecular fragment	Fragment loss
806	$C_{46}H_{81}NO_8P$		$+H^+$
478	$C_{24}H_{49}NO_6P$		
237	$C_{17}H_{33}$		
184	$C_5H_{15}NO_4P$		

[00218] Table 40: MS/MS fragmentation of multiple sclerosis biomarker 18, 194.0803 ($C_7H_{14}O_6$)

m/z	Formula	Molecular fragment	Fragment loss
195	$C_7H_{15}O_6$		$+H^+$
177	$C_7H_{13}O_5$		$-H_2O$
165	$C_6H_{13}O_5$		$-CH_2O$
163	$C_6H_{11}O_5$		$-CH_3OH$
137 Observed 138	$C_5H_{13}O_4$		
123	$C_4H_{11}O_4$		$137 - CH_2$

[00219] Table 41: MS/MS fragmentation of multiple sclerosis biomarker 19, 857.7516 ($C_{54}H_{99}NO_6$)

m/z	Formula	Molecular fragment	Fragment loss
858	$C_{54}H_{100}NO_6$		+ H^+
602	$C_{38}H_{68}NO_4$		- $C_{16}H_{34}O_2$
576	$C_{36}H_{66}NO_4$		602 - C_2H_2
314	$C_{17}H_{32}NO_4$		576 - $C_{19}H_{34}$
165	$C_{12}H_{21}$		
151	$C_{11}H_{19}$		
95	C_7H_{11}		

REFERENCES

- [00220] 1. MacDonald, B.K., et al. *The incidence and lifetime prevalence of neurological disorders in a prospective community-based study in the UK*. Brain, 2000. **123**: p. 665-76.
- 5 [00221] 2. Martyn, C. *Epidemiology*. McAlpine's Multiple Sclerosis. New York: Churchill Livingstone, 1992.
- [00222] 3. Sayetta, R.B. *Theories of the etiology of multiple sclerosis: a critical review*. J Clin Lab Immunol, 1986. **21**: p. 55-70.
- [00223] 4. Matthews, W.B. *Clinical aspects*. McAlpine's Multiple Sclerosis.
10 New York: Churchill Livingstone, 1992.
- [00224] 5. Poser, C.M. *Trauma and multiple sclerosis: an hypothesis*. J Neuro, 1987. **234**: p. 155-9.
- [00225] 6. Chelmicka-Schorr, E., and Arnason, B.G. *Nervous stem-immune system interactions and their role in multiple sclerosis*. Ann Neurol, 1994. **36**: p.29-
15 s32.
- [00226] 7. Noseworthy, J.H., et al. *Multiple sclerosis*. N Engl J Med, 2000. **343**: p.938-52.
- [00227] 8. ffrench-Constant, C. *Pathogenesis of multiple sclerosis*. Lancet, 1994. **343**: p.271-5.
- 20 [00228] 9. Prineas, J.W., and McDonald, W.I. *Demyelinating diseases*. Arnold: London, 1997.
- [00229] 10. Kidd, D., et al. *Spinal cord MRI using multi-array coils and fast spin echo. II. Findings in multiple sclerosis*. Neurology, 1993. **43**: p. 2632-7.
- [00230] 11. Lucchinetti, C.W., et al. *Heterogeneity of multiple sclerosis lesions: implication for the pathogenesis of demyelination*. Ann Neurol, 2000. **47**: p. 707-17.
25

- [00231] 12. Keegan, B.M., and Noseworthy, J.H. *Multiple Sclerosis*. Annu Rev Med, 2002. **52**: 285-302.
- [00232] 13. McDonald, W.I, et al. *Recommended diagnostic criteria for multiple sclerosis: guidelines from the international panel on the diagnosis of multiple sclerosis*. Ann Neurol, 2001. **50**: 121-7.
- [00233] 14. Reo, N.V., NMR-based metabolomics. Drug Chem Toxicol, 2002. **25**(4): p. 375-82.
- [00234] 15. Fiehn, O., et al. *Metabolite profiling for plant functional genomics*. Nat Biotechnol, 2000. **18**(11): p. 1157-61.
- [00235] 16. Hirai, M.Y., et al., *Integration of transcriptomics and metabolomics for understanding of global responses to nutritional stresses in Arabidopsis thaliana*. Proc Natl Acad Sci U S A, 2004. **101**(27): p. 10205-10.
- [00236] 17. Roessner, U., et al., *Metabolic profiling allows comprehensive phenotyping of genetically or environmentally modified plant systems*. Plant Cell, 2001. **13**(1): p. 11-29.
- [00237] 18. Castrillo, J.I., et al., *An optimized protocol for metabolome analysis in yeast using direct infusion electrospray mass spectrometry*. Phytochemistry, 2003. **62**(6): p. 929-37.
- [00238] 19. Fiehn, O., *Metabolomics---the link between genotypes and phenotypes*. Plant Mol Biol, 2002. **48**(1-2): p. 155-71.
- [00239] 20. Aharoni, A., et al., *Nontargeted metabolome analysis by use of Fourier Transform Ion Cyclotron Mass Spectrometry*. Omics, 2002. **6**(3): p. 217-34.
- [00240] 21. Hirai, M.Y., et al., *Elucidation of gene-to-gene and metabolite-to-gene networks in arabidopsis by integration of metabolomics and transcriptomics*. J Biol Chem, 2005. **280**(27): p. 25590-5.

- [00241] 22. Murch, S.J., et al., *A metabolomic analysis of medicinal diversity in Huang-qin (Scutellaria baicalensis Georgi) genotypes: discovery of novel compounds*. Plant Cell Rep, 2004. 23(6): p. 419-25.
- [00242] 23. Tohge, T., et al., *Functional genomics by integrated analysis of metabolome and transcriptome of Arabidopsis plants over-expressing an MYB transcription factor*. Plant J, 2005. 42(2): p. 218-35.
- [00243] 24. Tibshirani, R., et al., *Diagnosis of multiple cancer types by shrunken centroids of gene expression*. Proc Natl Acad Sci U S A, 2002. 99(10): p. 6567-72.
- 10 [00244] 25. 15. Wu, B., et al., *Comparison of statistical methods for classification of ovarian cancer using mass spectrometry data*. Bioinformatics, 2003. 19(13): p. 1636-43.
- [00245] 26. Hager, J.W., et al., *High-performance liquid chromatography-tandem mass spectrometry with a new quadrupole/linear ion trap instrument*. J Chromatogr A, 2003. 1020(1): p. 3-9.
- 15 [00246] 27. Hopfgartner, G., et al., *Triple quadrupole linear ion trap mass spectrometer for the analysis of small molecules and macromolecules*. J Mass Spectrom, 2004. 39(8): p. 845-55.
- [00247] 28. Xia, Y.Q., et al., *Use of a quadrupole linear ion trap mass spectrometer in metabolite identification and bioanalysis*. Rapid Commun Mass Spectrom, 2003. 17(11): p. 1137-45.
- 20 [00248] 29. Zhang, M.Y., et al., *Hybrid triple quadrupole-linear ion trap mass spectrometry in fragmentation mechanism studies: application to structure elucidation of buspirone and one of its metabolites*. J Mass Spectrom, 2005. 40(8): p. 1017-1029.
- 25

WHAT IS CLAIMED:

- I. A method for diagnosing a patient's multiple sclerosis health state, or change in multiple sclerosis health state, or for diagnosing another demyelinating disorder, or for diagnosing the risk of multiple sclerosis or another demyelinating disorder in a patient, the method comprising the steps of:
 - a) analyzing a blood sample from said patient to obtain quantifying data for one or more than one metabolite marker;
 - b) comparing the quantifying data for said one or more than one metabolite marker to corresponding data obtained from one or more than one reference sample to identify an increase or decrease in the level of said one or more than one metabolite marker in said sample; and
 - c) using said increase or decrease in the level of said one or more than one metabolite marker in said sample to diagnose said patient's multiple sclerosis health state, or change in multiple sclerosis health state, or for diagnosing another demyelinating disorder, or for diagnosing the risk of multiple sclerosis or another demyelinating disorder,

wherein the one or more than one metabolite marker comprises one or more molecule selected from the molecular formulae consisting of $C_{28}H_{52}O_4$, $C_{30}H_{56}O_5$, $C_{32}H_{60}O_5$, $C_{32}H_{60}O_6$, $C_{36}H_{64}O_5$, $C_{36}H_{66}O_5$, $C_{36}H_{68}O_5$, $C_{36}H_{66}O_6$, $C_{36}H_{68}O_6$, $C_{43}H_{79}O_{10}P$, $C_5H_{13}O_7P$, $C_{25}H_{52}NO_9P$, $C_{27}H_{52}NO_9P$, $C_6H_{11}O_6Na$, $C_8H_{13}O_7Na$, $C_{29}H_{48}O_2$, $C_{46}H_{80}NO_8P$, $C_7H_{14}O_6$ and $C_{54}H_{99}NO_6$,

the metabolite having the molecular formula of $C_{28}H_{52}O_4$ being characterized by a collision induced dissociation (CID) MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 451, 433, 407, 389, 281, 279, 183, 169, 153, 139, 125, 111 and 97,

the metabolite having the molecular formula of $C_{30}H_{56}O_5$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 495, 477, 451, 433, 415, 307, 297, 279, 235, 223, 215, 197, 179, 181, 169, 157, 155, 153, 141, 139, 127, 125 and 113,

the metabolite having the molecular formula of $C_{32}H_{60}O_5$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 523, 505, 487, 479, 463, 461, 443, 365, 337, 299, 297,

281, 279, 271, 269, 253, 251, 243, 225, 197, 171, 169, 157, 155, 143, 141, 139, 127, 125, 123, 115, 113, 111 and 83,

the metabolite having the molecular formula of $C_{32}H_{60}O_6$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 539, 521, 503, 495, 477, 461, 459, 419, 335, 315, 313, 297, 279, 259, 255, 253, 243, 241, 225, 223, 213, 179, 171, 155, 141 and 127,

the metabolite having the molecular formula of $C_{36}H_{64}O_5$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 575, 557, 539, 531, 513, 495, 417, 403, 371, 297, 279, 251 and 183,

the metabolite having the molecular formula of $C_{36}H_{66}O_5$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 577, 559, 541, 533, 515, 497, 419, 405, 387, 373, 297, 281 and 279,

the metabolite having the molecular formula of $C_{36}H_{68}O_5$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 579, 561, 543, 535, 517, 499, 421, 407, 389, 375, 299, 297, 281, 279, 263, 253, 185 and 171,

the metabolite having the molecular formula of $C_{36}H_{66}O_6$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 593, 575, 557, 549, 531, 513, 495, 421, 371, 315, 297, 279, 201, 171, 141 and 127,

the metabolite having the molecular formula of $C_{36}H_{68}O_6$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 595, 577, 559, 551, 515, 497, 423, 373, 315, 297, 281, 279, 269, 251, 171, 155, 153, 141, 139 and 127,

the metabolite having the molecular formula of $C_{43}H_{79}O_{10}P$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 785, 529, 425, 169 and 97,

the metabolite having the molecular formula of $C_5H_{13}O_7P$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 215, 197, 171, 153 and 135,

the metabolite having the molecular formula of $C_{25}H_{52}NO_9P$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 540, 480, 255, 242, 224, 168, 153 and 79,

the metabolite having the molecular formula of $C_{27}H_{52}NO_9P$ being characterized by a CID MS/MS fragmentation pattern under negative ionization comprising the following daughter ions: 564, 504, 279, 242, 224, 168, 153 and 79,

the metabolite having the molecular formula of $C_6H_{11}O_6Na$ being characterized by a CID MS/MS fragmentation pattern under positive ionization comprising the following daughter ions: 203, 159, 115, 89 and 97,

the metabolite having the molecular formula of $C_8H_{13}O_7Na$ being characterized by a CID MS/MS fragmentation pattern under positive ionization comprising the following daughter ions: 245, 227, 209, 191, 155, 125 and 83,

the metabolite having the molecular formula of $C_{29}H_{48}O_2$ being characterized by a CID MS/MS fragmentation pattern under positive ionization comprising the following daughter ions: 429, 205 and 165,

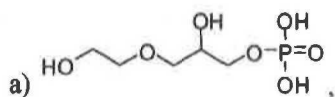
the metabolite having the molecular formula of $C_{46}H_{80}NO_8P$ being characterized by a CID MS/MS fragmentation pattern under positive ionization comprising the following daughter ions: 806, 478, 237 and 184,

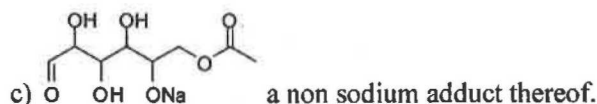
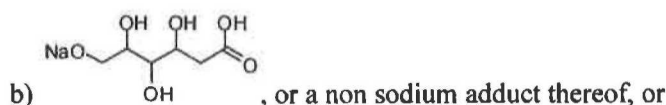
the metabolite having the molecular formula of $C_7H_{14}O_6$ being characterized by a CID MS/MS fragmentation pattern under positive ionization comprising the following daughter ions: 195, 177, 165, 163, 138 and 123, and

the metabolite having the molecular formula of $C_{54}H_{99}NO_6$ being characterized by a CID MS/MS fragmentation pattern under positive ionization comprising the following daughter ions: 858, 602, 576, 314, 165, 151 and 95.

2. The method of claim 1, wherein the multiple sclerosis metabolite markers are relapsing-remitting as compared to a normal reference sample and are selected from the molecular formulae consisting of a) $C_{30}H_{56}O_5$, b) $C_{32}H_{60}O_5$, c) $C_{32}H_{60}O_6$, d) $C_{36}H_{68}O_5$, e) $C_{36}H_{66}O_6$, f) $C_{36}H_{68}O_6$, and g) $C_{36}H_{66}O_5$, including combinations thereof.
3. The method of claim 2, wherein the samples are analyzed by MS/MS transition and one or more daughter ions are selected for:

- a) $C_{30}H_{56}O_5$ selected from 495, 477, 451, 433, 415, 307, 297, 279, 235, 223, 215, 197, 179, 181, 169, 157, 155, 153, 141, 139, 127, 125 and 113;
- b) $C_{32}H_{60}O_5$ selected from 523, 505, 487, 479, 463, 461, 443, 365, 337, 299, 297, 281, 279, 271, 269, 253, 251, 243, 225, 197, 171, 169, 157, 155, 143, 141, 139, 127, 125, 123, 115, 113, 111 and 83;
- c) $C_{32}H_{60}O_6$ selected from 539, 521, 503, 495, 477, 461, 459, 419, 335, 315, 313, 297, 279, 259, 255, 253, 243, 241, 225, 223, 213, 179, 171, 155, 141 and 127;
- d) $C_{36}H_{68}O_5$ selected from: 579, 561, 543, 535, 517, 499, 421, 407, 389, 375, 299, 297, 281, 279, 263, 253, 185 and 171;
- e) $C_{36}H_{66}O_6$ selected from 593, 575, 557, 549, 531, 513, 495, 421, 371, 315, 297, 279, 201, 171, 141 and 127;
- f) $C_{36}H_{68}O_6$ selected from 595, 577, 559, 551, 515, 497, 423, 373, 315, 297, 281, 279, 269, 251, 171, 155, 153, 141, 139 and 127; or
- g) $C_{36}H_{66}O_5$ selected from 577, 559, 541, 533, 515, 497, 419, 405, 387, 373, 297, 281 and 279.
4. The method of claim 1, wherein the multiple sclerosis metabolite markers are primary-progressive as compared to a normal reference sample and are selected from the molecular formulae consisting of a) $C_5H_{13}O_7P$, b) $C_6H_{11}O_6Na$, c) $C_8H_{13}O_7Na$, and d) $C_{54}H_{99}NO_6$, including combinations thereof.
5. The method of claim 4, wherein the samples are analyzed by MS/MS transition and one or more daughter ions are selected for:
- a) $C_5H_{13}O_7P$ selected from 215, 197, 171, 153 and 135;
- b) $C_6H_{11}O_6Na$ selected from 203, 159, 115, 89 and 97;
- c) $C_8H_{13}O_7Na$ selected from 245, 227, 209, 191, 155, 125 and 83; or
- d) $C_{54}H_{99}NO_6$ selected from 858, 602, 576, 314, 165, 151 and 95.
6. The method of claim 4 or 5, wherein the one or more than one metabolite is further characterized by the structure

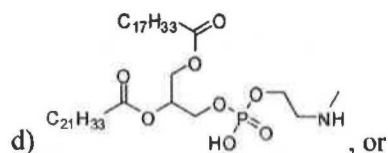
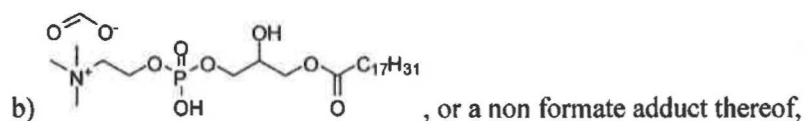
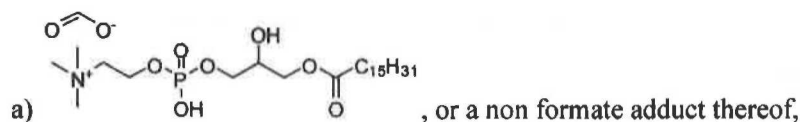


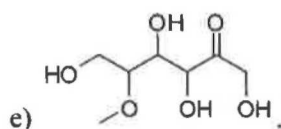


7. The method of claim 1, wherein the multiple sclerosis metabolite markers are secondary-progressive as compared to a normal reference sample and are selected from the molecular formulae consisting of a) $C_{25}H_{52}NO_9P$, b) $C_{27}H_{52}NO_9P$, c) $C_{29}H_{48}O_2$, d) $C_{48}H_{80}NO_8P$, e) $C_7H_{14}O_6$, and f) $C_{36}H_{66}O_5$.

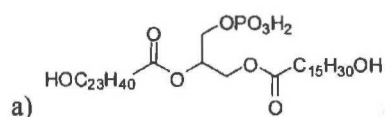
8. The method of claim 7, wherein the samples are analyzed by MS/MS transition and one or more daughter ions are selected for:
- a) $C_{25}H_{52}NO_9P$ selected from 540, 480, 255, 242, 224, 168, 153 and 79;
 - b) $C_{27}H_{52}NO_9P$ selected from 564, 504, 279, 242, 224, 168, 153 and 79;
 - c) $C_{29}H_{48}O_2$ selected from 429, 205 and 165;
 - d) $C_{46}H_{80}NO_8P$ selected from 806, 478, 237 and 184;
 - e) $C_7H_{14}O_6$ selected from 195, 177, 165, 163, 138 and 123; or
 - f) $C_{36}H_{66}O_5$ selected from 577, 559, 541, 533, 515, 497, 419, 405, 387, 373, 297, 281 and 279.

9. The method of claim 7 or 8, wherein the one or more than one metabolite is further characterized by the structure:





10. The method of claim 1, wherein the multiple sclerosis metabolite markers are relapsing-remitting as compared to a secondary-progressive reference sample and are selected from the molecular formulae consisting of a) $C_{32}H_{60}O_6$, and b) $C_{36}H_{64}O_5$, including combinations thereof.
11. The method of claim 2, wherein the samples are analyzed by MS/MS transition and one or more daughter ions are selected for:
- a) $C_{32}H_{60}O_6$ selected from 539, 521, 503, 495, 477, 461, 459, 419, 335, 315, 313, 297, 279, 259, 255, 253, 243, 241, 225, 223, 213, 179, 171, 155, 141 and 127; or
- b) $C_{36}H_{64}O_5$ selected from 575, 557, 539, 531, 513, 495, 417, 403, 371, 297, 279, 251 and 183.
12. The method of claim 1, wherein the multiple sclerosis metabolite markers are relapsing-remitting transitioning to secondary-progressive as compared to a secondary-progressive reference sample and comprise the molecular formula of a) $C_{43}H_{79}O_{10}P$.
13. The method of claim 12, wherein the samples are analyzed by MS/MS transition and one or more daughter ions are selected for the metabolite having the molecular formula of a) $C_{43}H_{79}O_{10}P$ selected from 785, 529, 425, 169 and 97.
14. The method of claim 12 or 13, wherein the metabolite is further characterized by the structure



15. The method of claim 1, wherein the multiple sclerosis metabolite markers are relapsing-remitting transitioning to secondary-progressive as compared to a relapsing-remitting reference sample and are selected from the molecular formulae consisting of a) $C_{36}H_{64}O_5$, and b) $C_{36}H_{66}O_5$, including combinations thereof.

16. The method of claim 15, wherein the samples are analyzed by MS/MS transition and one or more daughter ions are selected for:
 - a) $C_{36}H_{64}O_5$ selected from 575, 557, 539, 531, 513, 495, 417, 403, 371, 297, 279, 251 and 183; or
 - b) $C_{36}H_{66}O_5$ selected from 577, 559, 541, 533, 515, 497, 419, 405, 387, 373, 297, 281 and 279.
17. The method of claim 1, wherein the increase or decrease in the level of said one or more than one metabolite marker in said sample is identified in the comparing step (b) using a colorimetric chemical assay, an antibody-based enzyme-linked immunosorbant assay (ELISA), a dipstick chemical assay, or a mass spectrometry-based system.
18. The method of any one of claims 1 to 17, wherein the sample is whole blood, plasma, serum, or a subfraction of whole blood.
19. The method of any one of claims 1 to 18, wherein said one or more than one reference sample is a plurality of samples obtained from control individuals; one or more than one baseline sample obtained from the patient at an earlier date; or a combination thereof.
20. A method for diagnosing a patient's multiple sclerosis health state, or change in multiple sclerosis health state, or for diagnosing another demyelinating disorder, or for diagnosing the risk of multiple sclerosis or another demyelinating disorder in a patient, the method comprising the steps of:
 - a) analyzing a blood sample from said patient by high resolution mass spectrometry to obtain accurate mass intensity data;
 - b) comparing the accurate mass intensity data to corresponding data obtained from one or more than one reference sample to identify an increase or decrease in accurate mass intensity; and
 - c) using said increase or decrease in accurate mass intensity to diagnose said patient's multiple sclerosis health state, or change in multiple sclerosis health state, or for diagnosing another demyelinating disorder or for diagnosing the risk of multiple sclerosis or another demyelinating disorder,

wherein the accurate mass intensity is measured at or ± 5 ppm of a hydrogen and electron adjusted accurate mass or neutral accurate mass in Daltons selected from the group consisting of: 161.1051, 174.1411, 181.9806, 188.0143, 188.1567, 194.0803, 200.1566, 202.0453, 202.1721, 203.1157, 205.8867, 207.8836, 209.8525, 211.8495, 214.1721, 216.04, 216.1877, 217.9124, 218.0192, 218.0371, 226.0688, 239.939, 243.0719, 244.0559, 244.2189, 246.1467, 246.1468, 246.1472, 246.2345, 249.9677, 253.8165, 255.8135, 256.24, 257.8106, 260.004, 260.2137, 263.8453, 264.9759, 265.8424, 267.811, 268.1287, 269.132, 269.8081, 271.8051, 273.8743, 273.9985, 274.1777, 274.1778, 275.1811, 275.8713, 277.8861, 278.1494, 280.2413, 281.2447, 282.1444, 282.2572, 283.2602, 290.1737, 296.1601, 296.2357, 297.1634, 303.108, 303.1081, 304.2398, 305.8792, 306.2568, 310.1758, 311.7754, 312.231, 314.2461, 318.1421, 326.1708, 327.0307, 328.2393, 328.2415, 331.957, 333.8302, 333.9539, 340.2407, 345.8738, 364.3123, 366.3284, 367.3325, 368.1656, 374.2613, 376.2759, 377.2801, 378.2921, 378.9906, 379.2536, 379.2957, 380.3079, 381.311, 382.1084, 384.3399, 385.3428, 392.294, 428.3653, 432.2365, 432.3252, 438.3354, 440.3526, 446.2525, 446.341, 447.3433, 448.3194, 448.3562, 450.3729, 452.3868, 460.2681, 462.3346, 462.3716, 464.3524, 466.3656, 467.3711, 468.3577, 468.384, 469.3863, 474.2846, 474.3731, 476.3869, 477.3218, 478.4044, 480.3473, 482.3604, 484.3788, 488.2996, 490.3676, 492.3832, 493.385, 494.3968, 495.3321, 495.3322, 495.4018, 496.4157, 501.3217, 502.3165, 502.4054, 503.3194, 504.4188, 506.2853, 506.4338, 508.3782, 510.3937, 512.4079, 513.4116, 517.3141, 518.3969, 519.3322, 519.3998, 520.4131, 521.3478, 521.4188, 522.4313, 523.3637, 523.4337, 524.296, 524.4448, 530.4379, 532.1851, 532.4503, 534.3912, 534.4645, 536.41, 537.4142, 538.4257, 539.4274, 539.4286, 540.4387, 541.3141, 541.3415, 541.4422, 542.3447, 544.3636, 546.3413, 546.4298, 548.4438, 550.4602, 550.4954, 551.4646, 551.4976, 551.4991, 552.3273, 552.4784, 552.5008, 555.3102, 556.4497, 557.4527, 558.3761, 558.4649, 558.4663, 559.4688, 560.478, 560.4821, 561.4863, 562.4989, 563.5013, 564.4396, 564.513, 565.3391, 566.3431, 566.454, 567.3547, 568.4723, 569.3687, 569.4769, 570.4903, 572.4455, 574.4594, 574.4635, 574.4958, 575.4628, 576.4757, 576.4765, 577.4795, 578.4923, 579.4958, 580.5089, 581.5126, 584.2641, 590.4585, 590.4964, 591.4614, 592.4717, 593.4736, 594.4848, 594.4875, 595.4883, 595.4928, 596.5012, 596.5053, 597.5062, 597.5068, 598.5107, 606.4872, 610.3691, 610.482, 610.5204, 611.3724, 612.4994, 613.3404, 616.4675, 617.0921, 618.4834, 622.4973, 623.5003,

632.5038, 633.3232, 634.3951, 646.5702, 648.5861, 649.5895, 662.4267, 684.6037, 686.4879, 687.4916, 688.5048, 689.5083, 693.631, 694.6323, 702.4175, 708.4632, 712.5074, 713.5097, 715.5228, 716.4323, 716.4987, 718.4736, 720.4696, 722.479, 723.6395, 726.5438, 729.5727, 730.4535, 731.5464, 733.501, 734.488, 734.508, 736.5031, 737.5045, 739.5146, 741.5302, 742.4745, 742.5366, 743.5461, 743.5475, 744.55, 744.5516, 745.5643, 746.5118, 746.5701, 747.5121, 747.5761, 757.5008, 758.4785, 760.5231, 761.529, 763.5153, 764.5196, 765.5316, 766.5051, 766.5372, 767.5473, 768.5468, 768.5525, 769.5638, 770.5108, 770.569, 772.5265, 775.5514, 779.5828, 780.5303, 780.5872, 781.5619, 781.6001, 783.6174, 784.5238, 784.6228, 785.5287, 786.51, 786.5408, 787.5452, 788.5549, 789.5658, 791.5488, 792.555, 793.4936, 793.5663, 794.5718, 804.5476, 805.5609, 806.5643, 808.5225, 809.5264, 813.5871, 816.5159, 819.5831, 820.589, 821.5288, 823.5427, 827.5445, 827.5446, 828.5477, 828.5479, 828.7213, 830.5634, 831.5758, 831.5992, 832.5211, 832.6022, 833.5931, 833.7541, 834.5372, 834.5963, 835.5417, 835.6094, 839.6019, 851.5686, 852.5724, 854.737, 855.6009, 855.6798, 856.6045, 856.7527, 872.6715, 886.5582, 886.7896, 899.5871, 1016.9279, 1017.9316, 1018.9399, 1228.1101, 1251.1042, 1253.1236 and 1254.1311, or any combination thereof.

21. The method of claim 20, wherein the blood sample is whole blood, plasma, serum, or a subfraction of whole blood.
22. The method claim 20 or 21, wherein a decrease in accurate mass intensity is identified in the comparing step (b).
23. The method claim 20 or 21, wherein an increase in accurate mass intensity is identified in the comparing step (b).
24. The method of any one of claims 20 to 23, wherein a liquid/liquid extraction is performed on the blood sample whereby non-polar metabolites are dissolved into an organic solvent and polar metabolites are dissolved into an aqueous solvent.
25. The method of claim 24, wherein the accurate mass intensity data is obtained using a Fourier Transform Ion Cyclotron Resonance, time of flight, magnetic sector, quadrupole or triple quadrupole mass spectrometer.

26. The method of claim 24, wherein accurate mass intensity data is obtained from the ionization of the extracted samples using an ionization method selected from the group consisting of: positive electrospray ionization, negative electrospray ionization, positive atmospheric pressure chemical ionization, and negative atmospheric pressure chemical ionization, or a combination thereof.
27. The method of any one of claims 20 to 26, wherein said one or more than one reference sample is a plurality of samples obtained from control individuals; one or more than one baseline sample obtained from the patient at an earlier date; or a combination thereof.
28. The method of any one of claims 20 to 27, wherein the increase or decrease in accurate mass intensity is used to diagnose:
- relapsing-remitting as compared to a normal reference sample wherein the accurate mass intensity is measured at or ± 5 ppm of a hydrogen and electron adjusted accurate mass or neutral accurate mass in Daltons selected from the group consisting of: 188.0143, 202.0453, 205.8867, 207.8836, 209.8525, 211.8495, 218.0192, 243.0719, 244.0559, 246.1467, 249.9677, 253.8165, 255.8135, 256.24, 257.8106, 260.004, 263.8453, 265.8424, 267.811, 268.1287, 269.8081, 271.8051, 273.8743, 274.1778, 275.1811, 277.8861, 280.2413, 281.2447, 282.2572, 283.2602, 290.1737, 296.1601, 296.2357, 304.2398, 306.2568, 311.7754, 312.231, 314.2461, 318.1421, 327.0307, 328.2393, 333.8302, 333.9539, 345.8738, 364.3123, 366.3284, 367.3325, 374.2613, 376.2759, 377.2801, 378.2921, 378.9906, 379.2957, 380.3079, 381.311, 382.1084, 384.3399, 385.3428, 392.294, 428.3653, 432.3252, 438.3354, 440.3526, 446.2525, 446.341, 447.3433, 448.3194, 448.3562, 450.3729, 452.3868, 460.2681, 462.3346, 462.3716, 464.3524, 466.3656, 467.3711, 468.3577, 468.384, 469.3863, 474.2846, 474.3731, 476.3869, 478.4044, 480.3473, 482.3604, 484.3788, 488.2996, 490.3676, 492.3832, 493.385, 494.3968, 495.3322, 495.4018, 496.4157, 502.3165, 502.4054, 503.3194, 504.4188, 506.2853, 508.3782, 510.3937, 512.4079, 513.4116, 518.3969, 519.3998, 520.4131, 521.4188, 522.4313, 523.4337, 524.296, 524.4448, 530.4379, 532.4503, 534.3912, 534.4645, 536.41, 537.4142, 538.4257, 539.4274, 540.4387, 541.4422, 542.3447, 544.3636, 546.3413, 546.4298, 548.4438, 550.4602, 550.4954, 551.4646, 551.4976, 552.3273, 552.4784, 552.5008, 556.4497, 557.4527, 558.3761, 558.4649, 558.4663, 559.4688, 560.478, 562.4989, 563.5013, 564.4396,

564.513, 566.454, 568.4723, 569.4769, 570.4903, 572.4455, 574.4594, 574.4635, 575.4628, 576.4757, 576.4765, 577.4795, 578.4923, 579.4958, 580.5089, 581.5126, 590.4585, 591.4614, 592.4717, 593.4736, 594.4848, 594.4875, 595.4883, 595.4928, 596.5012, 596.5053, 597.5062, 597.5068, 598.5107, 610.3691, 610.482, 611.3724, 612.4994, 616.4675, 618.4834, 632.5038, 634.3951, 646.5702, 648.5861, 649.5895, 686.4879, 687.4916, 688.5048, 689.5083, 693.631, 694.6323, 702.4175, 708.4632, 712.5074, 713.5097, 715.5228, 716.4323, 722.479, 730.4535, 734.488, 736.5031, 737.5045, 739.5146, 741.5302, 742.5366, 743.5461, 743.5475, 744.55, 744.5516, 745.5643, 746.5701, 747.5761, 758.4785, 763.5153, 765.5316, 766.5051, 766.5372, 767.5473, 768.5525, 769.5638, 770.569, 775.5514, 781.5619, 781.6001, 789.5658, 791.5488, 792.555, 793.5663, 794.5718, 831.5992, 832.6022, 833.7541, 855.6009, 855.6798, 856.6045, 886.5582, 1016.9279, 1017.9316, 1018.9399, 1251.1042, 1253.1236 and 1254.1311 or a combination thereof,

primary-progressive as compared to a normal reference sample wherein the accurate mass intensity is measured at or ± 5 ppm of a hydrogen and electron adjusted accurate mass or neutral accurate mass in Daltons selected from the group consisting of: 161.1051, 188.0143, 202.0453, 203.1157, 216.04, 217.9124, 218.0192, 218.0371, 226.0688, 243.0719, 244.0559, 246.1467, 246.1472, 253.8165, 257.8106, 260.004, 263.8453, 268.1287, 269.132, 273.8743, 273.9985, 274.1778, 275.8713, 278.1494, 282.1444, 290.1737, 296.1601, 297.1634, 310.1758, 318.1421, 326.1708, 328.2393, 328.2415, 333.9539, 432.2365, 446.2525, 460.2681, 474.2846, 488.2996, 503.3194, 506.2853, 506.4338, 517.3141, 546.3413, 610.3691, 611.3724, 613.3404, 622.4973, 623.5003, 634.3951, 791.5488, 792.555, 793.4936, 819.5831, 827.5445, 828.5479, 831.5992, 832.6022, 899.5871, 1228.1101 or a combination thereof,

secondary-progressive as compared to a normal reference sample wherein the accurate mass intensity is measured at or ± 5 ppm of a hydrogen and electron adjusted accurate mass or neutral accurate mass in Daltons selected from the group consisting of: 194.0803, 202.0453, 243.0719, 244.0559, 246.1468, 249.9677, 253.8165, 255.8135, 311.7754, 312.231, 328.2393, 328.2415, 382.1084, 384.3399, 385.3428, 428.3653, 438.3354, 440.3526, 446.341, 448.3194, 448.3562, 464.3524, 466.3656, 467.3711, 468.3577, 468.384, 469.3863, 474.3731, 480.3473, 482.3604, 484.3788, 490.3676, 492.3832, 493.385, 494.3968, 495.3321, 495.4018, 508.3782, 510.3937, 512.4079, 513.4116, 518.3969, 519.3322, 519.3998, 520.4131, 521.3478, 521.4188,

522.4313, 523.3637, 523.4337, 534.3912, 536.41, 537.4142, 538.4257, 539.4274, 540.4387, 541.3141, 541.3415, 541.4422, 542.3447, 546.4298, 548.4438, 550.4602, 551.4646, 552.4784, 558.4663, 560.4821, 564.4396, 565.3391, 566.3431, 566.454, 568.4723, 569.3687, 569.4769, 574.4594, 574.4635, 575.4628, 576.4757, 576.4765, 577.4795, 578.4923, 579.4958, 580.5089, 590.4585, 590.4964, 592.4717, 593.4736, 594.4848, 594.4875, 595.4883, 595.4928, 596.5012, 596.5053, 597.5062, 597.5068, 598.5107, 606.4872, 610.482, 610.5204, 613.3404, 616.4675, 618.4834, 622.4973, 684.6037, 694.6323, 720.4696, 723.6395, 763.5153, 764.5196, 765.5316, 766.5372, 768.5468, 791.5488, 792.555, 793.4936, 805.5609, 806.5643, 819.5831, 820.589, 821.5288, 827.5446, 828.5477, 851.5686, 852.5724, 855.6798, 872.6715, 886.5582 and 1228.1101, or a combination thereof,

relapsing-remitting as compared to secondary-progressive wherein the accurate mass intensity is measured at or ± 5 ppm of a hydrogen and electron adjusted accurate mass or neutral accurate mass in Daltons selected from the group consisting of: 194.0803, 264.9759, 274.1777, 327.0307, 340.2407, 366.3284, 368.1656, 378.9906, 428.3653, 438.3354, 440.3526, 446.341, 447.3433, 448.3562, 450.3729, 452.3868, 462.3346, 462.3716, 464.3524, 466.3656, 467.3711, 468.384, 469.3863, 474.3731, 476.3869, 477.3218, 478.4044, 480.3473, 482.3604, 484.3788, 490.3676, 492.3832, 493.385, 494.3968, 495.4018, 496.4157, 501.3217, 502.4054, 504.4188, 506.4338, 508.3782, 510.3937, 512.4079, 513.4116, 518.3969, 519.3998, 520.4131, 521.4188, 522.4313, 523.4337, 524.4448, 530.4379, 532.1851, 532.4503, 534.3912, 534.4645, 536.41, 537.4142, 538.4257, 539.4274, 539.4286, 540.4387, 541.3141, 541.4422, 546.4298, 548.4438, 550.4602, 551.4646, 551.4991, 552.4784, 556.4497, 557.4527, 558.4649, 558.4663, 559.4688, 560.4821, 561.4863, 562.4989, 563.5013, 564.4396, 564.513, 566.454, 568.4723, 569.4769, 570.4903, 572.4455, 574.4594, 574.4635, 575.4628, 576.4757, 576.4765, 577.4795, 578.4923, 579.4958, 580.5089, 581.5126, 590.4585, 590.4964, 591.4614, 592.4717, 593.4736, 594.4848, 594.4875, 595.4883, 595.4928, 596.5012, 596.5053, 597.5062, 597.5068, 598.5107, 606.4872, 610.482, 610.5204, 612.4994, 616.4675, 618.4834, 632.5038, 634.3951, 646.5702, 662.4267, 729.5727, 780.5303, 786.51, 793.4936, 805.5609, 806.5643, 821.5288, 833.5931, 834.5963, 835.6094, 851.5686, 852.5724, 856.6045, 872.6715 and 886.7896, or a combination thereof,

relapsing-remitting transiting to secondary-progressive as compared to relapsing-remitting wherein the accurate mass intensity is measured at or ± 5 ppm of a hydrogen and electron adjusted accurate mass or neutral accurate mass in Daltons selected from the group consisting of: 174.1411, 181.9806, 188.1567, 200.1566, 202.1721, 214.1721, 216.1877, 239.939, 244.2189, 246.2345, 260.2137, 277.8861, 303.108, 327.0307, 331.957, 333.9539, 345.8738, 378.9906, 384.3399, 438.3354, 440.3526, 446.341, 447.3433, 448.3562, 450.3729, 452.3868, 462.3346, 462.3716, 464.3524, 466.3656, 467.3711, 468.384, 469.3863, 474.3731, 476.3869, 478.4044, 480.3473, 482.3604, 484.3788, 490.3676, 492.3832, 493.385, 494.3968, 495.4018, 496.4157, 502.4054, 504.4188, 506.4338, 508.3782, 510.3937, 512.4079, 513.4116, 518.3969, 519.3998, 520.4131, 521.4188, 522.4313, 523.4337, 524.4448, 530.4379, 532.4503, 534.3912, 534.4645, 536.41, 537.4142, 538.4257, 539.4274, 540.4387, 541.3141, 541.4422, 546.3413, 546.4298, 548.4438, 550.4602, 551.4646, 552.4784, 556.4497, 557.4527, 558.4649, 558.4663, 559.4688, 560.478, 560.4821, 561.4863, 562.4989, 563.5013, 564.4396, 564.513, 566.454, 568.4723, 569.4769, 570.4903, 572.4455, 574.4594, 574.4635, 575.4628, 576.4757, 576.4765, 577.4795, 578.4923, 579.4958, 580.5089, 581.5126, 584.2641, 590.4585, 590.4964, 591.4614, 592.4717, 593.4736, 594.4848, 594.4875, 595.4883, 595.4928, 596.5012, 596.5053, 597.5062, 597.5068, 598.5107, 610.482, 610.5204, 612.4994, 616.4675, 618.4834, 718.4736, 726.5438, 731.5464, 743.5475, 744.55, 747.5121, 779.5828, 780.5872, 781.6001, 783.6174, 784.6228, 793.4936, 804.5476, 813.5871, 816.5159, 821.5288, 823.5427, 828.5477, 830.5634, 831.5758, 832.5211, 834.5372, 834.5963, 835.5417 and 839.6019, or a combination thereof, and

relapsing-remitting transiting to secondary-progressive as compared to secondary-progressive wherein the accurate mass intensity is measured at or ± 5 ppm of a hydrogen and electron adjusted accurate mass or neutral accurate mass in Daltons selected from the group consisting of: 239.939, 246.1468, 249.9677, 303.1081, 305.8792, 331.957, 333.9539, 379.2536, 541.3141, 555.3102, 565.3391, 566.3431, 567.3547, 574.4958, 617.0921, 633.3232, 684.6037, 716.4987, 718.4736, 733.501, 734.508, 742.4745, 746.5118, 747.5121, 757.5008, 760.5231, 761.529, 770.5108, 772.5265, 783.6174, 784.5238, 784.6228, 785.5287, 786.5408, 787.5452, 788.5549, 808.5225, 809.5264, 828.7213, 854.737, 856.7527 and 872.6715, or a combination thereof.

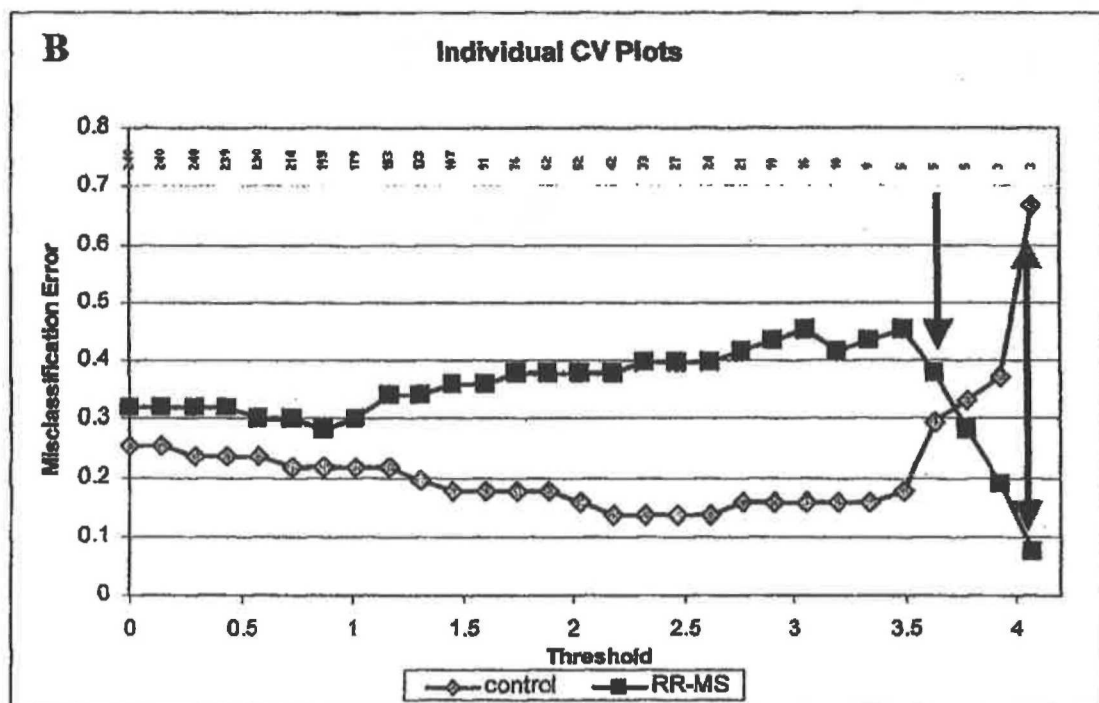
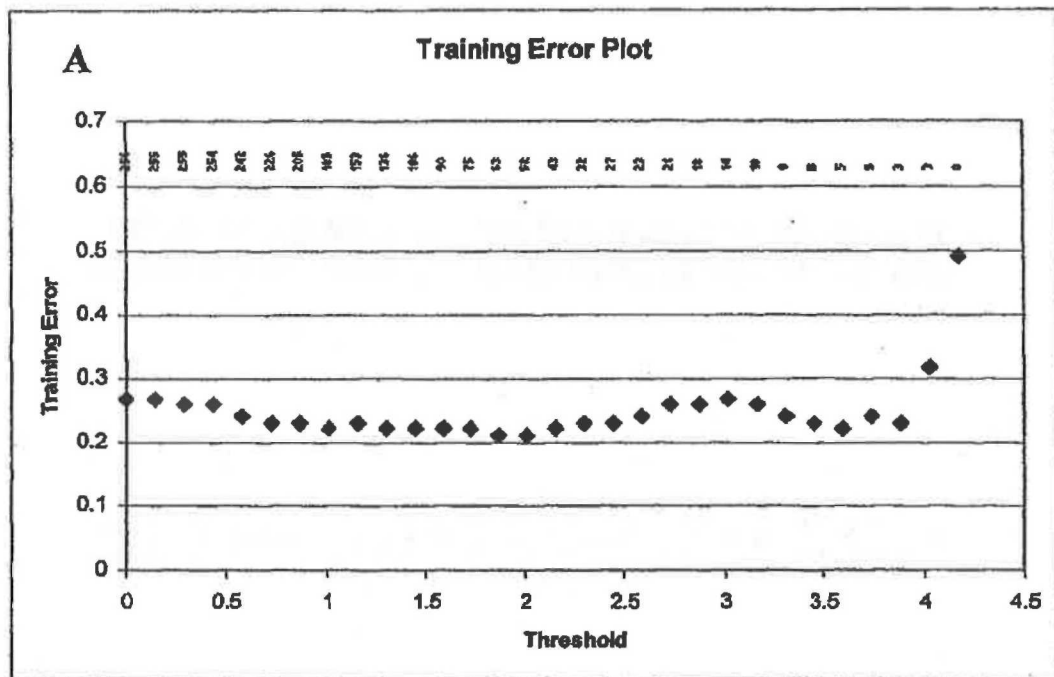


Fig. 1

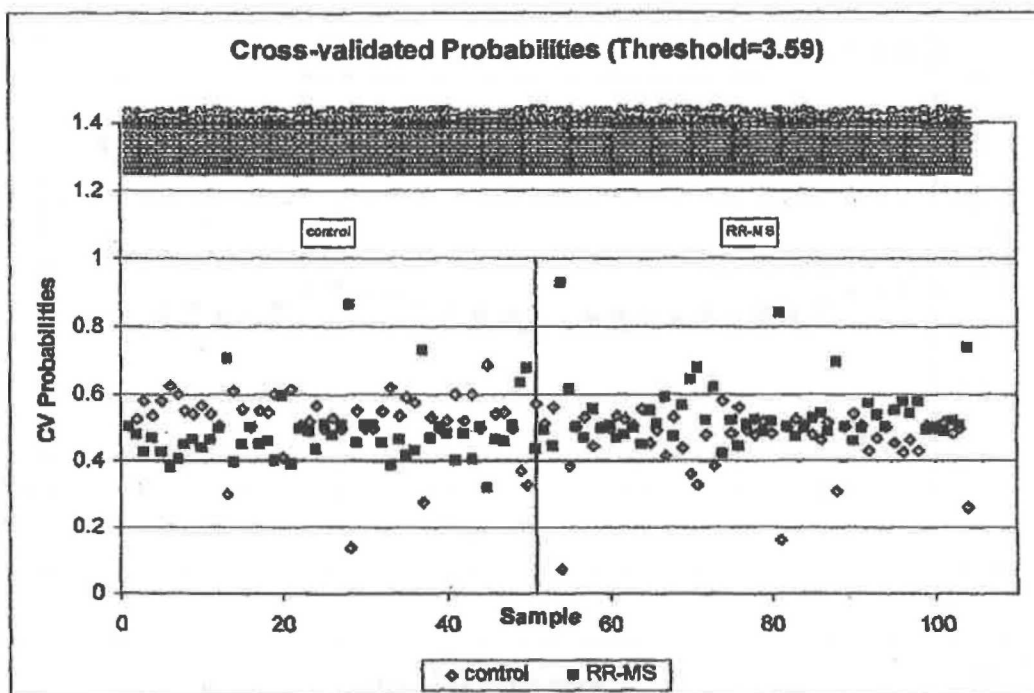


Fig. 2

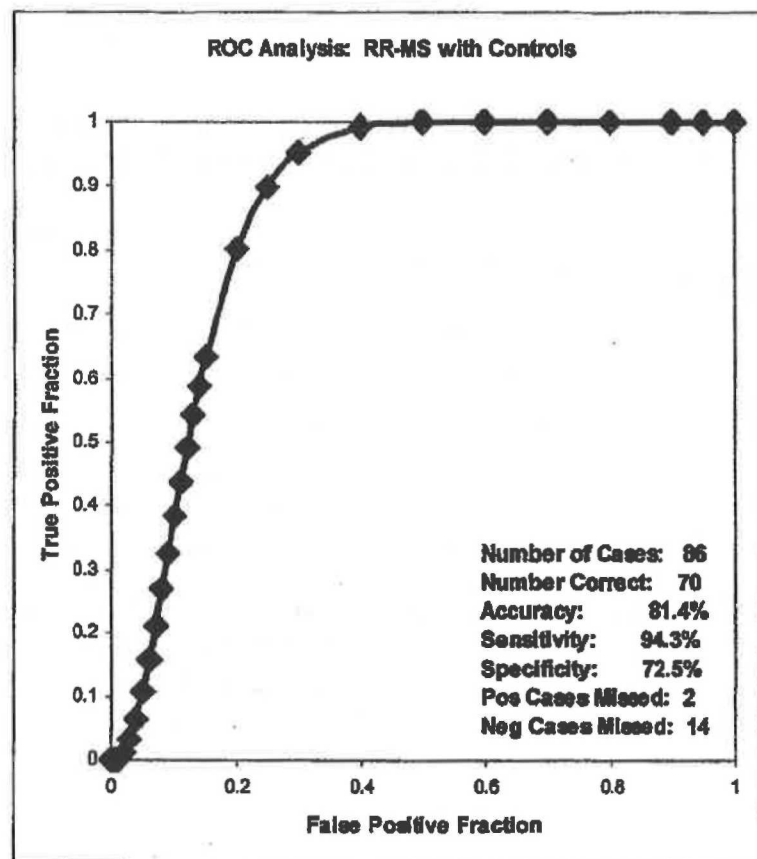


Fig. 3

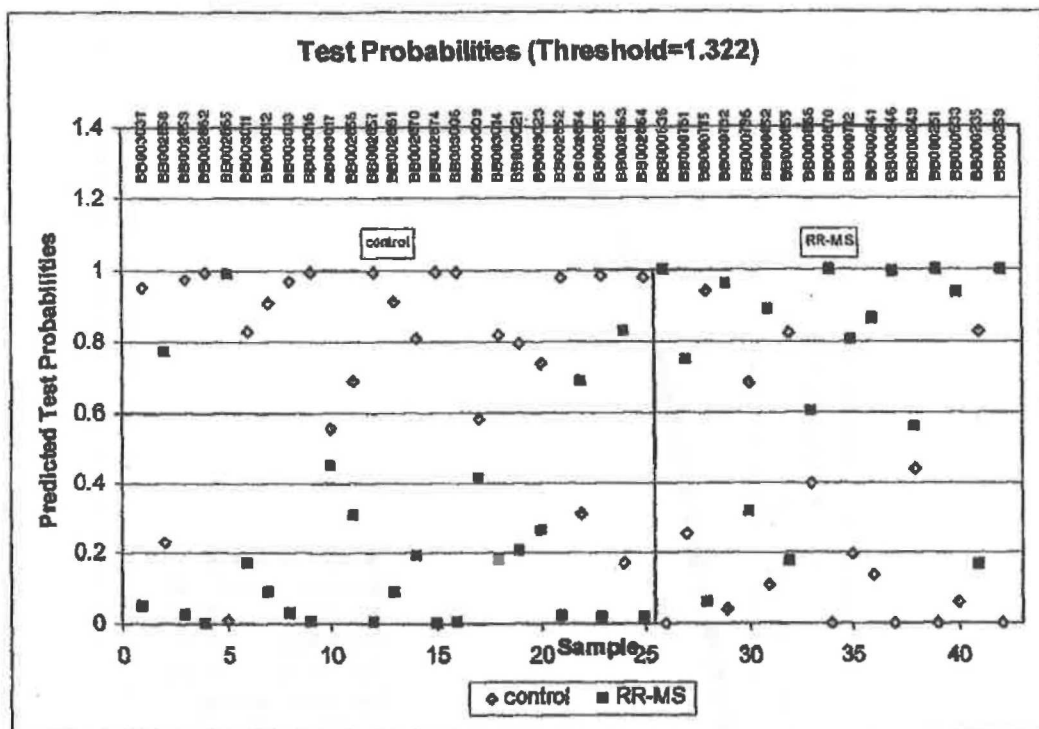


Fig. 4

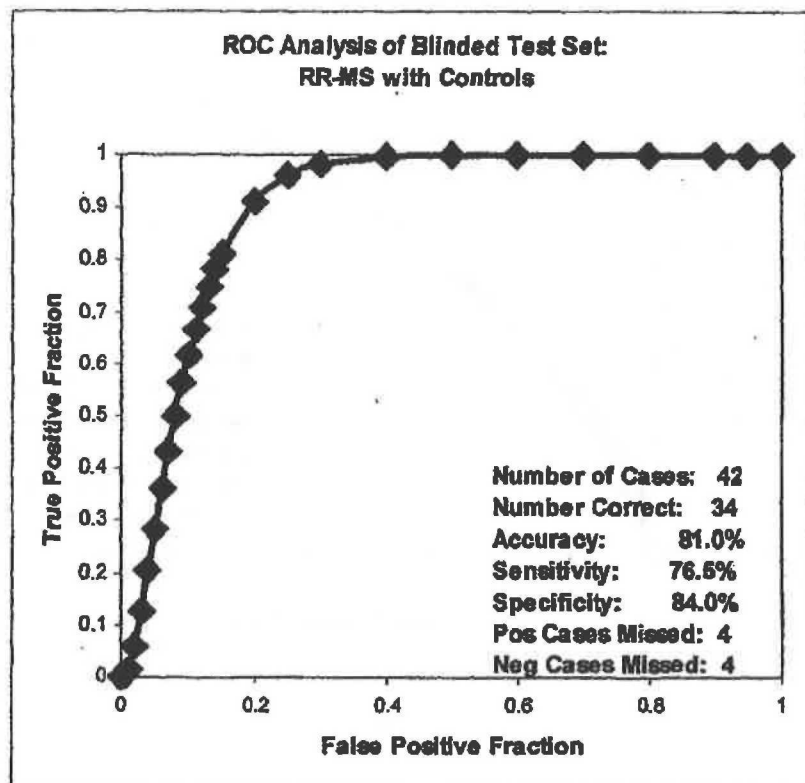


Fig. 5

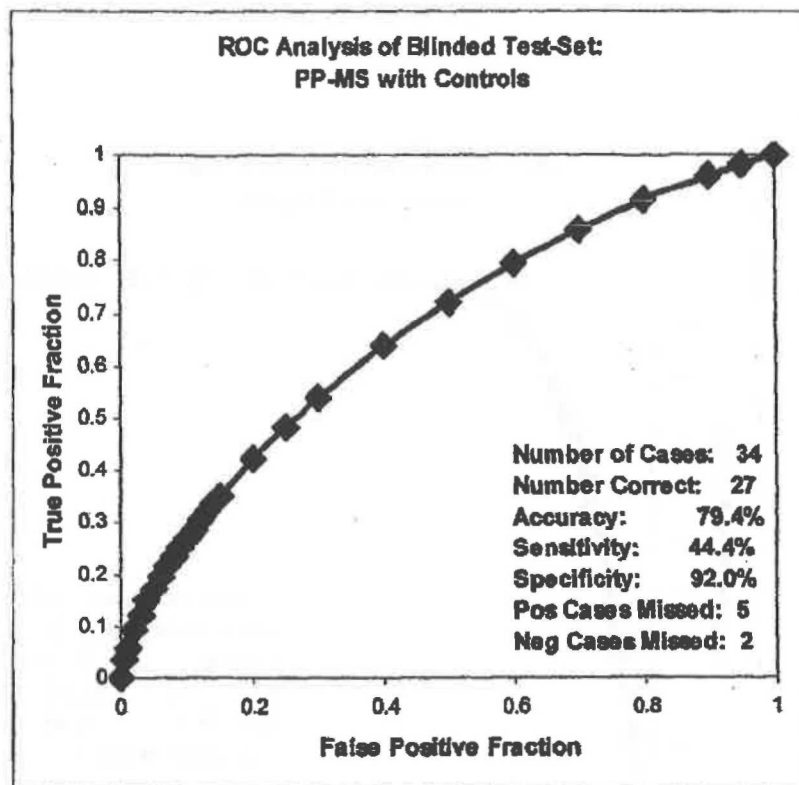


Fig. 6

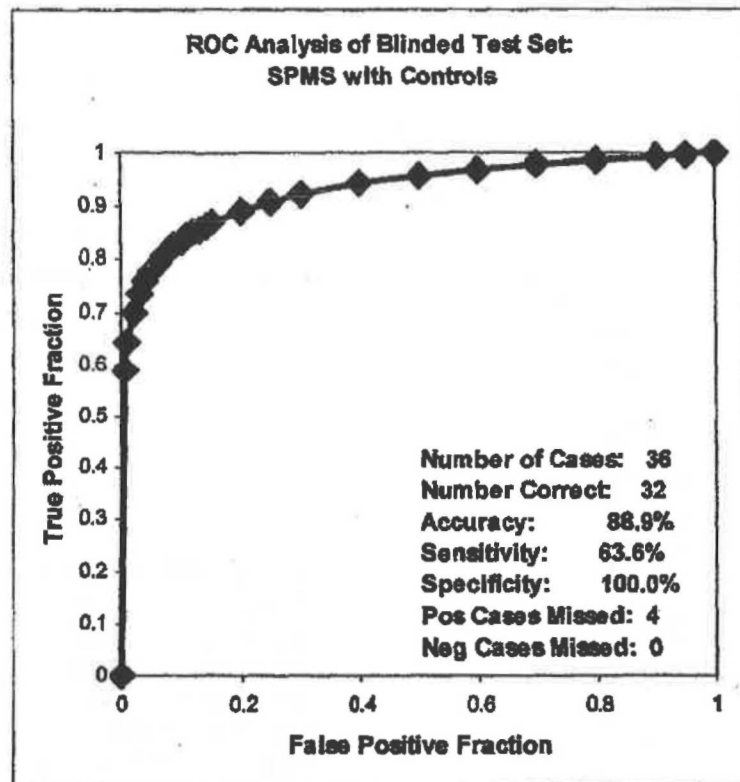


Fig. 7

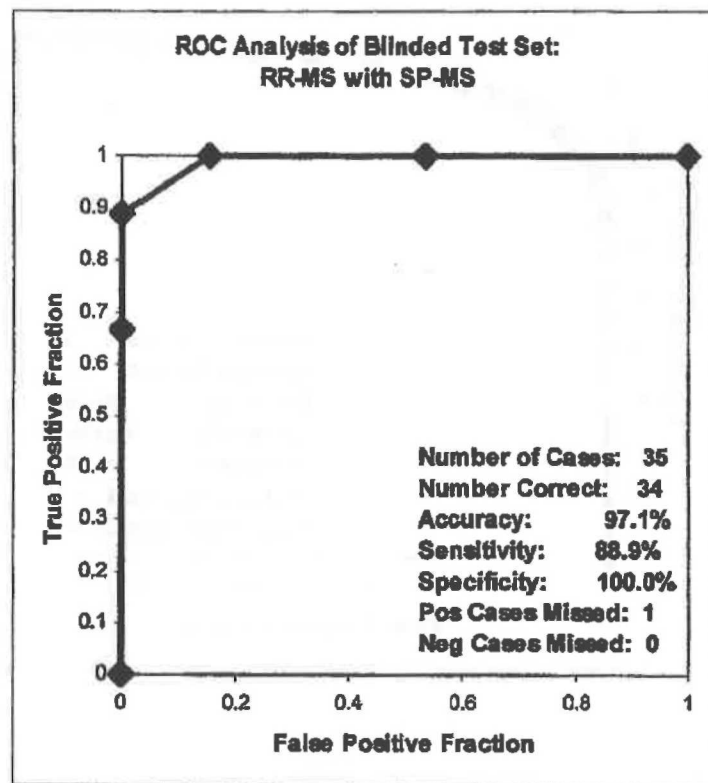


Fig. 8

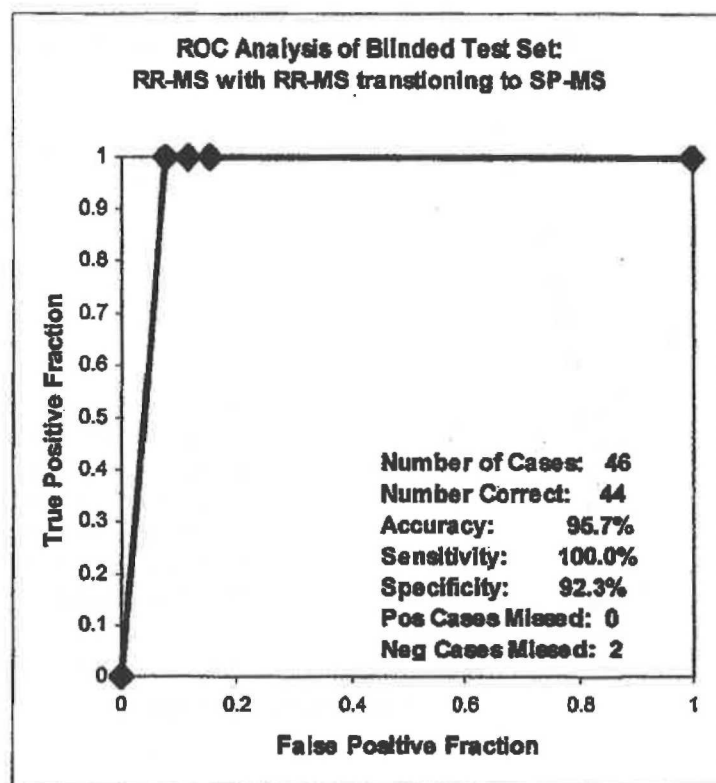


Fig. 9

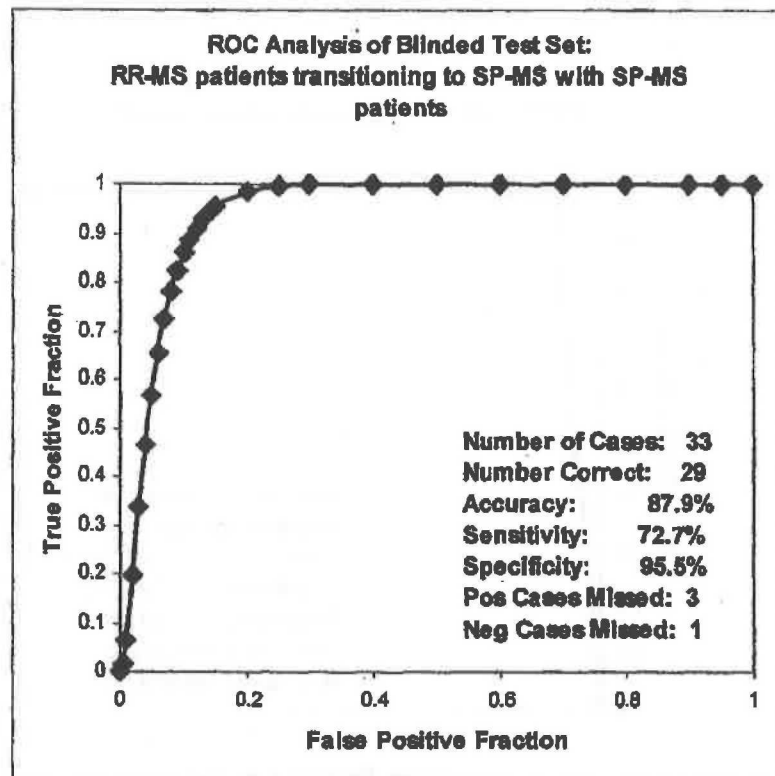


Fig. 10

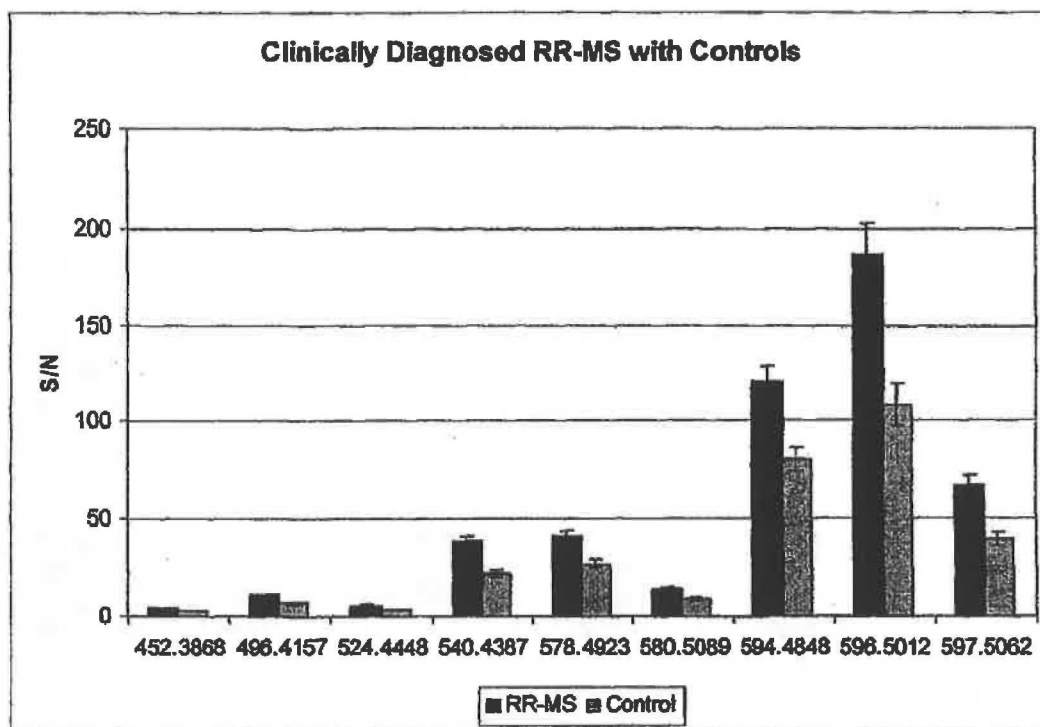


Fig. 11

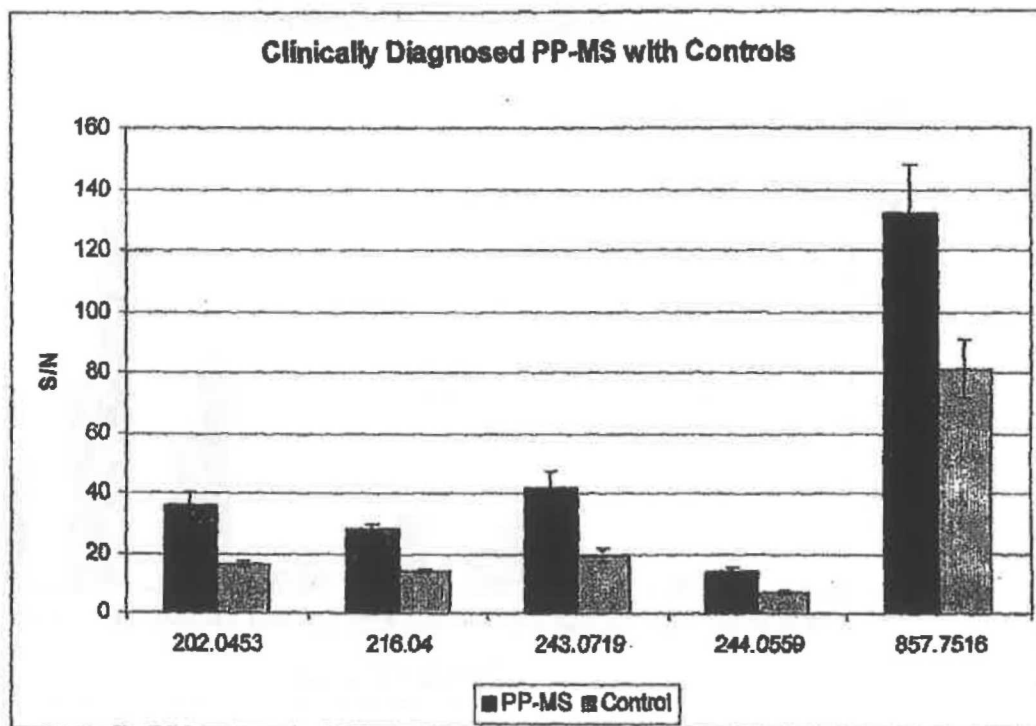


Fig. 12

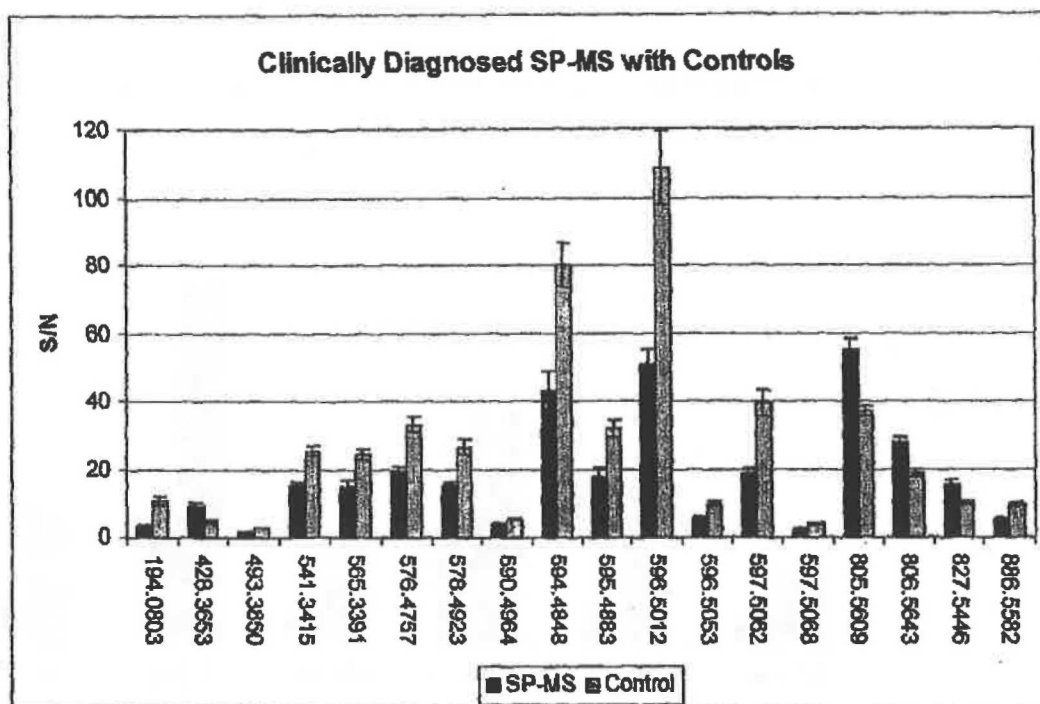


Fig. 13

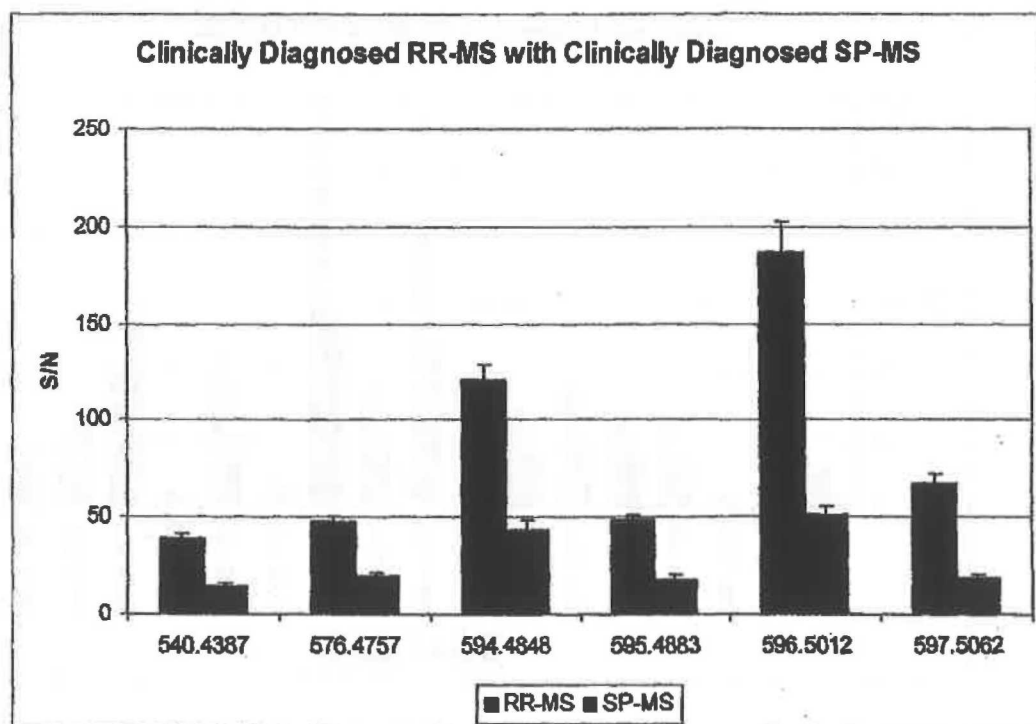


Fig. 14

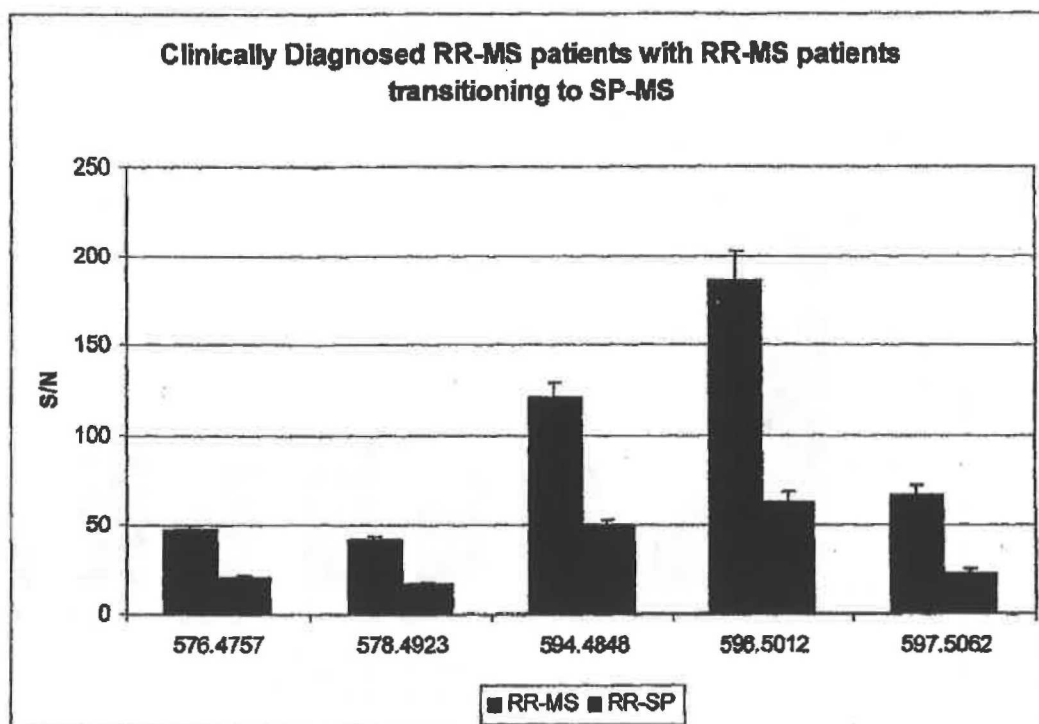


Fig. 15

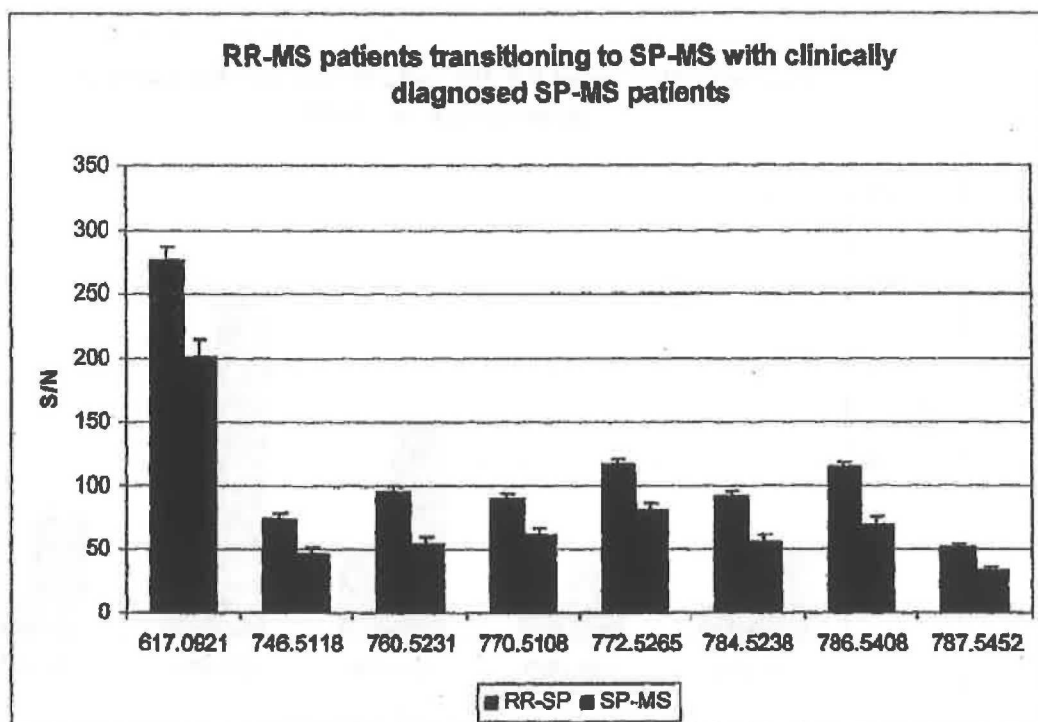


Fig. 16